

**Medtronic**

# MiniMed™ 780G

with Guardian™ 4 Sensor

System User Guide



Includes technology developed by **dreamed**  
diabetes ai



# **MiniMed™ 780G**













## System User Guide












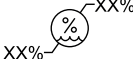








Medtronic and Medtronic logo are trademarks of Medtronic.™\*\* Third-party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.



## Icon table

	Bluetooth® wireless technology or Bluetooth® enabled
	Follow instructions for use or electronic instructions for use
	Manufacturer
	Date of manufacture (DoM)
	Country of manufacture (and Date of manufacture when a date appears beside)
	Importer
	Do not re-sterilize
	Do not re-use
	Do not use if package is damaged and consult instructions for use
FCC ID	Complies with United States regulations for RF devices
	Fragile, handle with care
IPX8	Protected against the effects of continuous immersion in water
	Keep dry
	Batch code

	Medical device
	Magnetic Resonance (MR) Unsafe
	Non-pyrogenic
	Type BF applied part
(1x)	One per container/package
	Open here
	Recyclable, contains recycled content
	Catalogue number
RF	Identification number for global radio frequency certification
R <sub>x</sub> Only	Requires prescription in the USA
	Single patient multiple use
	Single sterile barrier system
	Serial number
	Sterilized using irradiation
	Humidity limits

	Temperature limits
	Use-by date
	Do not dispose of this product in unsorted municipal waste stream
	Unique Device Identifier
	Non-ionizing electromagnetic radiation
	Caution: consult instructions for use for important warnings or precautions not found on the label



# Contents

## ■ Safety and indications

23	Important system information
24	Using this guide
24	Conventions
25	Emergency kit
26	User safety
26	Intended use
29	Risks and side effects
33	General warnings
33	Pump
38	Reservoir and infusion sets
38	Sensor and serter
41	Transmitter
41	Meter
41	Exposure to magnetic fields and radiation
42	General precautions
43	Waterproof capabilities
43	Electrostatic discharge
44	Extreme temperatures
44	Skin care products
44	Infusion sets and sites, transmitter, sensor, and meter

44	Adverse reactions
44	Security precautions
46	Insulin guidelines
47	Consumables
48	Other MiniMed 780G system devices
48	Optional items

## ■ System overview

53	What are the components of the MiniMed 780G system?
54	Modes
55	Manual mode without CGM
55	Manual mode with CGM
56	SmartGuard
57	Delivery settings

## ■ Pump basics

62	Using the buttons
62	Sleep mode
63	Pump delivery system
64	Infusion set
65	Reservoir
65	Pump
65	Battery cap
65	Identifying the correct battery cap for your pump
66	Inserting the battery
68	Removing the battery
69	Startup settings
71	Home screen in Manual mode
72	Shortcuts from the Home screen
72	Status icons
74	Menu screen

- 75 Menu map
- 75 Sound & Vibration screen
- 76 Status screen

## ■ Basal and bolus insulin delivery

- 81 Setting up basal insulin
- 81 Basal rate
- 82 Max basal rate
- 84 Basal patterns
- 87 Settings covering a 24-hour period
- 91 Viewing basal delivery information
- 91 Suspending all insulin delivery and resuming basal insulin delivery
- 94 Temp basal rates
- 96 Entering a blood glucose (BG) meter reading
- 97 Setting up bolus delivery
- 97 Manual mode bolus delivery options
- 97 Max bolus
- 99 Bolus Wizard feature in Manual mode
- 105 Delivering a Normal bolus
- 108 Stopping a bolus delivery

## ■ Reservoir and infusion set

- 114 Setting up the reservoir and infusion set
- 114 Removing the reservoir and rewinding the pump
- 117 Filling the reservoir and connecting it to the infusion set tubing
- 121 Placing the reservoir into the pump and filling the tubing with insulin
- 125 Inserting the infusion set into the body
- 126 Filling the cannula
- 129 Disconnecting the infusion set
- 129 Reconnecting the infusion set

## ■ Paired devices

- 133 Setting up the Accu-Chek Guide Link meter
- 133 Pairing the pump and meter
- 135 Unpair and delete the meter
- 135 Unpairing a meter from the pump
- 136 Deleting the pump from a meter
- 137 Pairing the pump and transmitter
- 139 Transmitter and sensor care
- 139 Unpairing the transmitter from the pump
- 140 Disconnecting the transmitter from the sensor
- 140 Removing the sensor
- 140 Cleaning the transmitter
- 140 Storing the transmitter
- 141 MiniMed Mobile app
- 141 Updating the pump software
- 141 Downloading the pump software update
- 141 Preparing to install the pump software update
- 142 Installing the pump software update
- 144 Completing the pump software update
- 144 Uploading device data to CareLink software
- 145 Sharing device data with the CareLink Connect app

## ■ Continuous glucose monitoring

- 149 CGM overview
- 149 What is CGM
- 151 What is blood glucose (BG) and sensor glucose (SG)?
- 153 Calibrating the sensor
- 153 CGM settings
- 153 Sensor alert settings



168	Manually resuming basal insulin delivery during a Suspend before low or Suspend on low event
168	Silencing sensor alerts
170	Guardian 4 sensor
170	Inserting the sensor
172	Connecting the transmitter to the sensor
172	Starting the sensor
172	Reconnecting the sensor
173	Deactivating the Sensor feature
173	Turning on the Sensor feature
174	Removing the sensor
174	Using CGM
174	Home screen with CGM in Manual mode
175	Sensor graph

## ■ SmartGuard

179	Introduction
181	Auto Basal
181	Auto Correction
181	Giving a bolus when the SmartGuard feature is active
181	Preparing to set up the SmartGuard feature
183	Setting up the SmartGuard feature
183	Conditions to activate the SmartGuard feature
184	Suspending manually while using the SmartGuard feature
184	Suspend before low and Suspend on low features while using the SmartGuard feature
184	SmartGuard Checklist
187	Home screen with the SmartGuard feature
187	Using the SmartGuard feature
187	The sensor graph with the SmartGuard feature
189	Entering a BG value in the SmartGuard feature

190	Delivering a bolus in the SmartGuard feature
196	Setting a temp target
197	Staying in the SmartGuard feature
198	Exiting the SmartGuard feature
199	Returning to the SmartGuard feature after an exit
200	Using Block mode with the SmartGuard feature
200	Alert silence feature

## ■ General settings

203	Time and date
203	Display options
204	Block mode
206	Self Test
207	Manage Settings
207	Saving the settings
207	Restoring the settings
208	Clearing the settings
209	Clearing the active insulin
210	Viewing the pump setting history
210	Max basal rate
211	Auto suspend
211	Language

## ■ History and graph

215	Introduction
215	History & Graph menu
215	History
215	Summary screen
216	Understanding the Summary screen
219	Daily History screen
219	Alarm History screen

- 220 Paired Sensors screen
- 220 SG Review screen
- 222 Graph screen
- 222 Time in Range screen

## ■ Notifications and reminders

- 227 Notifications in the MiniMed Mobile app
- 227 Reminders
- 228 Personal reminders
- 229 Bolus BG Check reminder
- 229 Missed Meal Bolus reminder
- 230 Low Reservoir reminder
- 231 Set Change reminder
- 232 Alarms, alerts, and messages
- 233 Alarms
- 234 Alerts
- 235 Messages
- 236 Pump alarms, alerts, and messages

## ■ Additional basal features

- 239 Preset temp basal rates
- 241 Starting a preset temp basal delivery
- 242 Canceling a temp basal or preset temp basal
- 242 Additional basal patterns
- 242 Adding an additional basal pattern
- 243 Editing, copying, or deleting a basal pattern
- 244 Changing from one basal pattern to another

## ■ Additional bolus features

- 249 Bolus types
- 250 Bolus type example

250	Bolus settings
250	Bolus increment
251	Bolus speed
252	Changing the Bolus Wizard settings in Manual mode
252	Changing the carb ratio
252	Changing the insulin sensitivity factor
253	Changing the BG target
253	Changing the active insulin time
254	Square Wave bolus
255	Turning the Square Wave bolus feature on or off
255	Delivering a Square Wave bolus using the Bolus Wizard feature
257	Delivering a Square Wave bolus using the Manual bolus feature
258	Dual Wave bolus
259	Turning the Dual Wave bolus feature on or off
259	Delivering a Dual Wave bolus using the Bolus Wizard feature
261	Delivering a Dual Wave bolus using the Manual bolus feature
262	Easy bolus
263	Setting up the Easy bolus feature
264	Delivering a bolus using the Easy bolus feature
265	Preset bolus
265	Setting up and managing preset bolus deliveries
268	Delivering a preset bolus
268	Stopping a Square Wave or Dual Wave bolus delivery

## ■ Troubleshooting

276	Pump issues
279	Sensor issues

## ■ Pump Care

283	Clean, store, and dispose of the pump
283	Cleaning the pump

284 Storing the pump

286 Pump disposal

## ■ Appendix A: List of alarms, alerts, and messages

291 Pump alarms, alerts, and messages

306 CGM (sensor) alarms, alerts, and messages

315 SmartGuard feature alerts and messages

320 CareLink software alert and message

## ■ Appendix B: Product specifications

325 Specifications and default settings

325 Alarm and alert escalation

326 Altitude range

326 Backlight

326 Basal delivery

327 BG meter reading

327 Bolus delivery

327 Bolus Wizard feature in Manual mode default settings

328 Bolus Wizard feature in Manual mode specifications

331 Carb ratios

331 Easy bolus feature

332 Environmental conditions

332 Essential performance

332 Expected service life

333 Filling the infusion set and cannula

333 Insulin delivery default settings

334 Low Reservoir reminder

334 Max bolus

334 Normal bolus

334 Percent temp basal

335 Program safety checks

335	Pump dimensions
335	Pump memory
335	Pump performance characteristics
338	Pump weight
338	Sensor default settings
339	IEC 60601-1-2
339	IEC 60601-1-2, Special EMC Precautions for Medical Electrical Equipment
340	IEC 60601-1-2
340	IEC 60601-1-10: PCLCS
341	Guidance and manufacturer's declaration
344	Wireless communication
344	FCC notice
346	Open Source Software disclosure

## ■ **Appendix C: Menu map**

349	Menu map
-----	----------

## ■ **Appendix D: Performance data**

355	MiniMed 780G system device performance
356	The SmartGuard feature
366	The Suspend before low feature
367	Guardian 4 System Performance
367	Clinical study overview
368	Sensor accuracy
369	Agreement when CGM reads "Below 50 mg/dL" or "Above 400 mg/dL"
370	Concurrence of SG and YSI values
371	Trend accuracy
372	Accuracy over time
372	Precision
372	Sensor life
373	Safety

373	Alert performance
374	Glucose TRUE Alert Rate
375	Glucose FALSE Alert Rate
377	Glucose Correct Detection Rate
378	Glucose Missed Detection Rate

## ■ **Glossary**

## ■ **Index**







1



# 1

## Safety and indications

This user guide describes the operation of the MiniMed 780G system with smart device connectivity and SmartGuard technology. SmartGuard technology adjusts insulin delivery based on sensor glucose (SG) values without the need to enter a blood glucose (BG) meter reading for confirmation. The MiniMed 780G insulin pump operates in Manual mode when the SmartGuard feature is not active.

The MiniMed 780G system with the Guardian 4 transmitter is designed to operate without requiring calibration.

Consult a healthcare professional before starting insulin pump therapy.

### Important system information

Only use rapid-acting U-100 insulin with the MiniMed 780G system. For more information see *Insulin guidelines*, page 46.

The MiniMed 780G system uses the Guardian 4 sensor and Guardian 4 transmitter for continuous glucose monitoring. For more information see *Continuous glucose monitoring*, page 149.

The Guardian 4 sensor does not require calibration. However, the system is designed to use every blood glucose (BG) meter reading either entered manually or received from a linked glucose meter to calibrate the sensor. For more information see *Calibrating the sensor*, page 153.

The Guardian 4 sensor is indicated for arm insertion only. For more information see *Inserting the sensor*, page 170.



**CAUTION:** The Guardian 4 sensor is indicated for arm use only. Do not use the Guardian 4 sensor in the abdomen or other body sites including the buttocks, due to unknown or different performance that could result in hypoglycemia or hyperglycemia.





Only use MiniMed or Medtronic reservoirs and infusion sets that are specifically designed for use with the MiniMed 780G system. For more information on compatible reservoirs and infusion sets see *Consumables, page 47*.

## Using this guide

Use the table of contents at the beginning of the user guide and the index at the end of the user guide to locate specific information.

Refer to the glossary for definitions of terms and acronyms used.

## Conventions

Convention	Definition
Select	Press  to activate a screen item, accept a value, or initiate an action.
Select and hold	Press and hold  to perform an action.
Press	Press and release a button.
Press and hold	Press and hold a button.
Bold text	Indicates screen items and buttons, such as “Select <b>Next</b> to continue.”
X	Indicates a value that might appear differently on the pump screen.
Note	<b>Note:</b> A note provides helpful information.
Caution	 <b>CAUTION:</b> A caution informs of a potential hazard which, if not avoided, might result in minor or moderate injury, or damage to the equipment.
WARNING	 <b>WARNING:</b> A warning informs of a potential safety hazard which, if not avoided, may result in serious injury or death. It may also describe potential serious adverse reactions.

For instructions about setting up devices on the MiniMed 780G system, such as a sensor or infusion set, refer to the user guide for the related device.

## Emergency kit

Keep an emergency kit available at all times and confirm that necessary supplies are available and not expired. Tell a family member or friend where to find the emergency kit.

When traveling, check blood glucose more frequently to accommodate for changes in activity levels and meal times.

Consult your healthcare professional on which of the following items to include in your emergency kit:

- Rapid-acting glucose
- Blood glucose (BG) testing supplies
- Urine or blood ketone monitoring supplies
- Extra infusion set and reservoir
- Extra new AA lithium or alkaline batteries, or fully charged NiMH batteries
- Insulin syringe
- Short-acting insulin, long-acting insulin, or both (with dosage instructions from a healthcare professional)
- Adhesive dressing
- Glucagon



**WARNING:** Do not use the Bolus Wizard feature to calculate a bolus for a period of time after giving a manual injection of insulin by syringe or pen. Manual injections are not accounted for in the active insulin amount. Using the Bolus Wizard feature too soon after a manual injection may result in over-delivery of insulin and may cause hypoglycemia. Consult a healthcare professional for how long to wait after a manual injection before using the Bolus Wizard feature.



**WARNING:** Do not use the SmartGuard feature for a period of time after giving a manual injection of insulin by syringe or pen. Manual injections are not accounted for in the active insulin amount. Using the SmartGuard feature too soon after a manual injection may result in over-delivery of insulin and may cause hypoglycemia. Consult a healthcare professional for how long to wait after a manual injection before using the SmartGuard feature.

## User safety



**WARNING:** Do not use the MiniMed 780G system until appropriate training has been received from a healthcare professional. Training is essential to ensure the safe use of the MiniMed 780G system.

## Intended use

The MiniMed 780G system is intended for the continuous delivery of basal insulin at selectable rates, and the administration of insulin boluses at selectable amounts for the management of type 1 diabetes mellitus in persons 7 years of age and older requiring insulin. The system is also intended to continuously monitor glucose values in the fluid under the skin. The MiniMed 780G system includes SmartGuard technology, which can be programmed to automatically adjust insulin delivery based on continuous glucose monitoring (CGM) sensor glucose values and can suspend delivery of insulin when the SG value falls below or is predicted to fall below predefined threshold values.

The MiniMed 780G system consists of the following devices:

- MiniMed 780G insulin pump
- Guardian 4 transmitter
- Guardian 4 sensor
- One-press server
- Accu-Chek™\* Guide Link blood glucose meter
- Accu-Chek Guide Test Strips

The system requires a prescription from a healthcare professional.



**WARNING:** Do not use the Suspend before low or Suspend on low features to prevent or treat low glucose. Always follow the instructions of a healthcare professional to treat low glucose. Using Suspend before low or Suspend on low features to prevent or treat low BG may result in prolonged hypoglycemia.

Guardian 4 sensor

The Guardian 4 sensor is intended for use with the MiniMed 780G system and the Guardian 4 transmitter to monitor glucose levels for the management of diabetes.

The sensor is intended for single use and requires a prescription. The Guardian 4 sensor is indicated for **up to** 7 days of continuous use.

The Guardian 4 sensor is not intended to be used directly to make therapy adjustments while the MiniMed 780G is operating in manual mode. All therapy adjustments in Manual mode should be based on measurements obtained using a blood glucose meter and not on values provided by the Guardian 4 sensor.

The Guardian 4 sensor has been studied and is approved for use in the systems, insertion sites, and ages listed in the following table.

System	Approved Age	Sensor Insertion Site
MiniMed 780G system	7 years and older	Arm

One-pressserter

Theserter is used as an aid for inserting the sensor. It is indicated for single-patient use and is not intended for multiple patient use.

Guardian 4 transmitter

The Guardian 4 transmitter is intended for use with the MiniMed 780G system and Guardian 4 sensor to monitor glucose levels for the management of diabetes.

## Accu-Chek Guide Link Blood Glucose Monitoring System

The Accu-Chek Guide Link Blood Glucose Monitoring System is comprised of the Accu-Chek Guide Link meter and the Accu-Chek Guide test strips. The Accu-Chek Guide Link Blood Glucose Monitoring System is intended to quantitatively measure glucose in fresh capillary whole blood from the fingertip, palm, and upper arm as an aid in monitoring the effectiveness of glucose control.

The Accu-Chek Guide Link Blood Glucose Monitoring System is intended for in vitro diagnostic single-patient use by people with diabetes. The Accu-Chek Guide Link Blood Glucose Monitoring System is intended to be used by a single person and should not be shared.

This system is not for use in diagnosing or screening for diabetes mellitus and not for neonatal use. Alternative site testing should be done only during steady-state times (when glucose is not changing rapidly). The Accu-Chek Guide Link Blood Glucose Monitoring System is intended to be used to wirelessly transmit glucose values to the MiniMed 780G system and MiniMed 770G system with Bluetooth™\* wireless technology through the use of Bluetooth low energy communication.



### WARNINGS:

- Do not use Alternative Site Testing to calibrate a continuous glucose monitoring system.
- Do not use Alternative Site Testing to make insulin dosing calculations.

## Contraindications

**The MiniMed 780G system is contraindicated for use in persons under age seven.**

Do not use the inserter to insert sensors other than the Guardian 4 sensor. Medtronic cannot guarantee the safety or efficacy of this product if used with other sensors.

The reservoir is contraindicated for the infusion of blood or blood products.

Infusion sets are indicated for subcutaneous use only and not for intravenous (IV) infusion.

Infusion sets are not indicated for the infusion of blood or blood products.



Pump therapy is not recommended for persons who are unwilling to or unable to perform BG meter readings.



**WARNING:** Do not use the SmartGuard feature for people who require less than 8 units or more than 250 units of total daily insulin per day. A total daily dose of at least 8 units, but no more than 250 units, is required to use the SmartGuard feature.

Pump therapy is not recommended for people who are unwilling or unable to maintain contact with their healthcare professional.

### Intended target population

The intended target population for the MiniMed 780G insulin pump includes children and adolescents ages 7-17 years and adults ages 18 years and older who are responsive to insulin delivered subcutaneously.

### Risks and side effects

#### Risks related to insulin administration and pump use

Risks related to insulin infusion and potential interruptions of insulin delivery include:

- Hypoglycemia
- Hyperglycemia
- Diabetic ketoacidosis
- Seizure
- Coma
- Death

#### Risks related to insulin pump infusion set

Risks related to insulin pump infusion set use include:

- Localized infection
- Skin irritation or redness

- Bruising
- Discomfort or pain
- Bleeding
- Irritation
- Rash
- Occlusions that may interrupt insulin delivery and lead to hyperglycemia and diabetic ketoacidosis

Follow the instructions in the provided user guides for the insertion and care of infusion sets. If an infusion site becomes irritated or inflamed, dispose of the infusion set in a sharps container, and select a different location to insert a new infusion set.

### **Risks related to sensor use**

Risks related to sensor use include:

- Skin irritation or other reactions
- Bruising
- Discomfort
- Redness
- Bleeding
- Pain
- Rash
- Infection
- Raised bump
- Appearance of a small freckle-like dot where needle was inserted
- Fainting secondary to anxiety or fear of needle insertion
- Allergic reaction
- Soreness or tenderness
- Swelling at insertion site

- Sensor filament fracture, breakage, or damage
- Minimal blood splatter associated with sensor needle removal
- Residual redness associated with adhesive or tapes or both
- Scarring

### Specific risks related to sensor use

Do not use continuous glucose monitoring if hydroxyurea, also known as hydroxycarbamide, is taken. Hydroxyurea is used to treat certain diseases, such as cancer and sickle cell anemia. Hydroxyurea use results in higher sensor glucose readings compared to blood glucose readings. Taking hydroxyurea while using continuous glucose monitoring can result in hypoglycemia caused by over-delivery of insulin, inaccurate or missed alarms and alerts, delay or loss of sensor-enabled insulin suspension, and substantially higher sensor glucose readings in reports than actual blood glucose readings.

Always check the label of any medication being taken to confirm if hydroxyurea or hydroxycarbamide is an active ingredient. If hydroxyurea is taken, consult a healthcare professional. Turn the Sensor feature off to disable continuous glucose monitoring. For more information, see *Deactivating the Sensor feature*, page 173. Use additional blood glucose meter readings to verify glucose levels.

Consult a healthcare professional if a medication that contains acetaminophen or paracetamol is taken while wearing the sensor. Medications that contain acetaminophen or paracetamol can falsely raise sensor glucose readings. The level of inaccuracy depends on the amount of acetaminophen or paracetamol active in the body and can differ for each person. Falsely elevated sensor readings can result in over-delivery of insulin, which can cause hypoglycemia. Medications that contain acetaminophen or paracetamol include, but are not limited to, cold medicines and fever reducers. Check the label of any medications being taken to see if acetaminophen or paracetamol is an active ingredient. Use additional blood glucose meter readings to confirm blood glucose levels.

While the SmartGuard feature is active, if acetaminophen or paracetamol is taken, program a temp target for up to 8 hours, or the amount of time recommended by a healthcare provider. For more information, see *Setting a temp target*, page 196. Use blood

glucose values instead of sensor glucose readings to calculate a meal bolus or correction bolus for up to 8 hours, or the duration recommended by a healthcare provider, after taking acetaminophen or paracetamol.

Do not use SG values to make treatment decisions, including delivering a bolus, while the pump is in Manual mode. When the SmartGuard feature is active and you are no longer in Manual mode, the pump uses an SG value, when available, to calculate a bolus amount. However, if your symptoms do not match the SG value, use a BG meter to confirm the SG value. Failure to confirm glucose levels when your symptoms do not match the SG value can result in the infusion of too much or too little insulin, which may cause hypoglycemia or hyperglycemia. For more information on using CGM, see *Continuous glucose monitoring, page 149*. For more information on using the SmartGuard feature, see *SmartGuard, page 179*.

For persons 7 years of age and older, use the upper arm insertion site. Do not insert the sensor into any other location.

### Risks related to meter use

- For the most current risks, see the User's Manual that came with the device.
- A list of warnings for the meter are provided in the meter section, see *Meter, page 41*.



#### **WARNINGS:**

- Do not use Alternative Site Testing to calibrate a continuous glucose monitoring system.
- Do not use Alternative Site Testing to make insulin dosing calculations.

### Risks related toserter use

General risks withserter use may include skin infection around the area where theserter is used.

## Risks related to the MiniMed 780G system

- Hypoglycemia
- Hyperglycemia
- Diabetic ketoacidosis
- Seizure
- Coma
- Death

## Removing the pump for temporary storage

If there is a need or desire to remove the pump, use the following guidelines:

- Write down the current basal rates and use the Save Settings feature. For more information, see *Saving the settings*, page 207.
- Remove the battery. For more information, see *Storing the pump*, page 284.
- If the pump is disconnected for less than one hour, an insulin adjustment may not be required. If the pump is disconnected for more than one hour, consult a healthcare professional to determine an alternate method of insulin delivery.

## General warnings

### Pump

- Do not use the pump in the presence of anesthetic mixtures that include oxidizing agents such as oxygen or nitrous oxide. Exposure to these conditions may damage the pump and result in serious injury.
- Always use the fingertip for blood samples when entering a BG meter reading into the pump. All BG values are used for calibration. Do not use blood samples from the palm for BG values entered into the pump. The palm has not been studied for use with the SmartGuard feature and the performance of the system using such blood samples is not known.
- When the SmartGuard feature is active, SG readings are used to calculate basal insulin delivery and correction boluses. Do not use SG readings to make treatment decisions while the pump is in Manual mode. SG and BG values may differ. Sensor

performance may occasionally vary from sensor to sensor and in different situations for a sensor, such as on the first day of use.

A BG meter reading is required in the following situations:

- Before a correction bolus is given in Manual mode.
- The SG reading is lower than expected.
- The SG reading is higher than expected.
- Suspected hypoglycemia or symptoms of hypoglycemia.
- Suspected hyperglycemia or symptoms of hyperglycemia.
- Suspected diabetic ketoacidosis or symptoms of diabetic ketoacidosis.

Do not use SG readings to make treatment decisions while the pump is in Manual mode.

- The low SG alert functionality is distinct from the automated insulin dosing function of the MiniMed 780G system. When using the SmartGuard feature, the MiniMed 780G system has been shown to be safe and effective for its intended use in this population. However, do not rely solely on the use of the Low SG alarm, or the use of “Alert on Low” and “Alert before Low” when those alerts are set at or below 60 mg/dL. At these BG levels, a low SG alarm or alert may not reflect the user’s true BG, and you may not be notified. Do not ignore symptoms of low glucose. Always confirm SG readings with a BG meter, and treat according to the recommendation of a healthcare professional. Solely relying on these SG alerts and readings for treatment decisions could result in missing severe hypoglycemia (low BG) events.
- Do not rely on the pump tones or vibrations to navigate the pump screens or menus. Relying on pump tones or vibrations may result in incorrect menu or setting selection. Always view the pump screen when selecting menus and entering information into the system.
- Only use rapid-acting U-100 insulin prescribed by a healthcare professional for use with an infusion pump. For a list of compatible insulins, see *Insulin guidelines, page 46*. Use of any other drug or medication in the reservoir can cause serious injury.

- Confirm that the infusion set is disconnected from the body before rewinding the pump or filling the infusion set tubing. Never insert the reservoir into the pump while the tubing is connected to the body. Doing so may result in an accidental infusion of insulin, which may cause hypoglycemia.
- Do not insert the reservoir before rewinding the pump. Doing so may result in an accidental infusion of insulin, and may result in hypoglycemia.
- Do not use the MiniMed 780G insulin pump or additional system devices next to other electrical equipment, which may cause interference. This includes mobile communication devices such as cell phones that are not paired with the MiniMed 780G system, GPS navigation systems, anti-theft systems, and any electrical equipment that has an output transmitter power greater than 1 W. The recommended separation distance between the insulin pump and common RF emitters is 12 in (30 cm). For more information about recommended separation distance guidelines between the insulin pump and common RF emitters, see *Guidance and manufacturer's declaration, page 341*. Other electrical equipment that may compromise normal system operation has been contraindicated. For more information, see *Exposure to magnetic fields and radiation, page 41*.
- Do not unscrew or retighten the tubing connector on the reservoir while the infusion set is connected to the body. Doing so may result in an accidental infusion of insulin, and may cause hypoglycemia.
- Do not use Luer sets with the MiniMed 780G system. Only use MiniMed or Medtronic reservoirs and infusion sets that are specifically designed for use with the MiniMed 780G system.
- Do not change or modify the MiniMed or Medtronic reservoir and infusion set. Modification of these components may cause serious injury, interfere with device operation, and void the warranty.
- Do not rely on pump alerts, alarms, or reminders alone to check BG levels. Set additional reminders on other devices, such as a cell phone.
- Do not change or modify the internal RF transmitter or antenna. Doing so may interfere with the safe operation of the equipment.

- The MiniMed 780G system is approved for use with the Guardian 4 transmitter with Bluetooth wireless technology (MMT-7841). Use of a transmitter not approved for communication with the pump may cause damage to system components and may result in inaccurate SG readings.
- If other devices that employ radio frequencies are in use, such as cell phones that are not paired with the MiniMed 780G system, cordless phones, walkie-talkies, and wireless networks, they may prevent communication between the sensor and the insulin pump. This interference does not cause any incorrect data to be sent and does not cause any harm to devices. Moving away from, or turning off, these other devices may enable communication. Contact 24-Hour Technical support if RF interference continues.
- Special Precautions regarding Electromagnetic Compatibility (EMC): This body-worn device is intended to be operated within a residential, domestic, public or work environment, where common levels of radiated "E" (V/m) or "H" fields (A/m) exist. Technologies that emit these fields include: cellular phones that are not paired with the MiniMed 780G system, wireless technology, electric can openers, microwaves, and induction ovens. The MiniMed 780G system generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the provided instructions, may cause harmful interference to radio communications.
- Portable and mobile RF communications equipment can affect the operation of the MiniMed 780G system. If interference occurs, move away from the RF transmitter.
- The MiniMed 780G insulin pump can generate, use, and radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. If the MiniMed 780G insulin



pump does cause interference to radio or television reception, try to correct the interference by one or more of the following measures:

- Decrease the distance between the sensor and the insulin pump to 6 feet (1.8 meters) or less.
  - Decrease the distance between the meter and the insulin pump to 6 feet (1.8 meters) or less.
  - Increase the separation between the sensor and the device that is receiving/emitting interference.
- The safety of the MiniMed 780G system has not been studied in persons with impaired kidney function. Persons with kidney disease should consult a healthcare professional to determine if the potential benefits of pump therapy outweigh the risks.
  - Monitor for diabetic retinopathy. During the beginning of insulin pump therapy, rapid improvement in glucose control and reduction in A1c may result in worsening of existing diabetic retinopathy. Use of the MiniMed 780G system has been associated with rapid improvement in glucose control. Monitor for diabetic retinopathy with retinal eye examinations and if necessary adequate treatment must be performed by a healthcare professional before beginning a treatment with the MiniMed 780G insulin pump.
  - The safety of the MiniMed 780G system has not been studied in pregnant women, persons with type 2 diabetes, or in persons using other anti-hyperglycemic therapies that do not include insulin. Persons in these situations should consult a healthcare professional to determine if the potential benefits of pump therapy outweigh the risks.
  - The safety of using the Suspend before low and Suspend on low features in patients who have no pump experience is not known. The Suspend before low and Suspend on low features should not be used if insulin pump settings have not been previously established. Insulin pump settings include basal rates, insulin to carb ratio, and insulin sensitivity factors. Consult a healthcare professional before using the Suspend before low or Suspend on low features.

- If a serious incident related to the device occurs, immediately report the incident to a healthcare professional. Serious incidents may include death, temporary or permanent serious decline in health, or a serious public health threat. Immediately report any serious incident to the manufacturer and local competent authority.
- The user must have adequate vision and hearing to recognize all functions of the pump, including alerts, alarms, and reminders. Not recognizing an alert, alarm, or reminder could result in a hypoglycemic or hyperglycemic event.

## Reservoir and infusion sets

See the user guides that came with the device for the most current warnings related to the reservoir and infusion set.

- If insulin, or any other liquid, gets inside the tubing connector, it can temporarily block the vents that allow the pump to properly fill the infusion set. This may result in the infusion of too little or too much insulin, and may result in hyperglycemia or hypoglycemia. If this occurs, start over with a new reservoir and infusion set.
- If a BG reading is unexpectedly high during the infusion of insulin or if an occlusion alarm occurs, check the infusion set for clogs and leaks.

If in doubt, change the infusion set in case the soft cannula is dislodged, crimped, or partially clogged. Consult a healthcare professional to create a plan for rapid insulin replacement in the event this occurs. Check BG to confirm that the appropriate amount of insulin has been administered.

- Only use reservoir and infusion sets manufactured or distributed by Medtronic Diabetes. The pump has been tested to operate when used with compatible reservoirs and infusion sets. Medtronic Diabetes cannot guarantee appropriate operation if the pump is used with reservoirs or infusion sets offered by third parties. Medtronic Diabetes is not responsible for any injury or pump malfunction that may occur in association with the use of incompatible components.

## Sensor andserter

For the most current warnings, see the user guide that came with the device.

- Do not attempt to connect a transmitter that is not compatible with the sensor. The sensor is designed to work with approved transmitters only. Connecting the

sensor to a transmitter that is not approved for use with the sensor may damage the components.

- Do not use continuous glucose monitoring if hydroxyurea, also known as hydroxycarbamide, is taken. Hydroxyurea is used to treat certain diseases, such as cancer and sickle cell anemia. Hydroxyurea use results in higher sensor glucose readings compared to blood glucose readings. Taking hydroxyurea while using continuous glucose monitoring can result in hypoglycemia caused by over-delivery of insulin, inaccurate or missed alarms and alerts, delay or loss of sensor-enabled insulin suspension, and substantially higher sensor glucose readings in reports than actual blood glucose readings.

Always check the label of any medication being taken to confirm if hydroxyurea or hydroxycarbamide is an active ingredient. If hydroxyurea is taken, consult a healthcare professional. Turn the Sensor feature off to disable continuous glucose monitoring. For more information, see *Deactivating the Sensor feature*, page 173. Use additional blood glucose meter readings to verify glucose levels.

- Consult a healthcare professional if a medication that contains acetaminophen or paracetamol is taken while wearing the sensor. Medications that contain acetaminophen or paracetamol can falsely raise sensor glucose readings. The level of inaccuracy depends on the amount of acetaminophen or paracetamol active in the body and can differ for each person. Falsely elevated sensor readings can result in over-delivery of insulin, which can cause hypoglycemia. Medications that contain acetaminophen or paracetamol include, but are not limited to, cold medicines and fever reducers. Check the label of any medications being taken to see if acetaminophen or paracetamol is an active ingredient. Use additional blood glucose meter readings to confirm blood glucose levels.

While the SmartGuard feature is active, if acetaminophen or paracetamol is taken, program a temp target for up to eight hours, or the amount of time recommended by a healthcare provider. For more information, see *Setting a temp target*, page 196. Use blood glucose values instead of sensor glucose readings to calculate a meal bolus or correction bolus for up to eight hours, or the duration recommended by a healthcare provider, after taking acetaminophen or paracetamol.

- Always inspect the packaging for damage prior to use. Sensors are sterile and non-pyrogenic, unless the package has been opened or damaged. If the sensor

packaging is open or damaged, discard the sensor directly into a sharps container. Use of a non-sterile sensor may result in infection at the insertion site.

- Instructions for using the One-pressserter (MMT-7512) are different from other Medtronic insertion devices. Failure to follow directions, or using a differentserter, may result in improper insertion, pain, or injury.
- Confirm that the sensor is securely placed in theserter to avoid improper insertion, pain, or minor injury.
- Do not allow children to put small parts in their mouth. This product poses a choking hazard for young children.
- Healthcare professionals and caregivers:
  - Always wear gloves to insert the sensor. A retractable needle is attached to the sensor. Minimal bleeding may occur.
  - Cover the sensor with sterile gauze to remove the needle housing from the sensor.
- Place the needle housing directly into a sharps container after sensor insertion to prevent accidental needlestick injury.
- Watch for bleeding at the insertion site (under, around, or on top of the sensor). If bleeding occurs, do the following:
  1. Apply steady pressure, using sterile gauze or a clean cloth placed on top of the sensor, for up to three minutes. The use of unsterile gauze can cause site infection.
  2. If bleeding stops, connect the transmitter to the sensor. If bleeding does not stop, do not connect the transmitter to the sensor because blood can get into the transmitter connector, and may damage the device.
- If bleeding continues, causes excessive pain or discomfort, or blood is significantly visible in the plastic base of the sensor, do the following:
  1. Remove the sensor and continue to apply steady pressure until the bleeding stops. Discard the sensor in a sharps container.

2. Check the site for redness, bleeding, irritation, pain, tenderness, or inflammation. Treat based on instructions from a healthcare professional.
  3. Insert a new sensor in a different location.
- Wear gloves when inserting the sensor into someone other than yourself to avoid contact with patient blood. Minimal bleeding may occur. Contact with patient blood can cause infection.
  - Never point a loadedserter toward any body part where insertion is not desired. An accidental button-push may cause the needle to inject the sensor in an undesired location and cause minor injury.
  - The safety of sensor use in critically ill patients is not known. Sensor use in critically ill patients is not recommended.

## Transmitter

See the user guide included with the device for the most current warnings related to transmitter use.

Do not allow children to put small parts in their mouth. This product poses a choking hazard for young children.

## Meter

For the most current warnings, see the User's Manual that came with the device.

Always use the fingertip for blood samples when entering a BG meter reading into the pump. All BG values are used for calibration. Do not use blood sample from the palm for BG values entered into the pump. The palm has not been studied for use with the SmartGuard feature and the performance of the system using such blood samples is not known.

## Exposure to magnetic fields and radiation

- Do not expose the pump, transmitter, or sensor to MRI equipment, diathermy devices, or other devices that generate strong magnetic fields (for example, x-ray, CT scan, or other types of radiation). Strong magnetic fields can cause the system to malfunction, and result in serious injury. If the pump is exposed to a strong

magnetic field, discontinue use and contact 24-Hour Technical Support for further assistance.

Magnetic fields, and direct contact with magnets, may affect the accurate functioning of the system which may lead to health risks such as hypoglycemia or hyperglycemia.

- Remove the pump, transmitter, sensor, and meter before entering a room with x-ray, MRI, diathermy, or CT scan equipment. The magnetic fields and radiation in the immediate vicinity of this equipment can make the devices nonfunctional or damage the part of the pump that regulates insulin delivery, possibly resulting in over-delivery and severe hypoglycemia.
- Do not expose the pump to a magnet, such as pump cases that have a magnetic clasp. Exposure to a magnet may interfere with the motor inside the pump. Damage to the motor can cause the device to malfunction, and result in serious injury.
- Do not send the pump, transmitter, or sensor through an x-ray scanning machine. The radiation can damage the pump components that regulate insulin delivery, and may result in over-delivery of insulin and hypoglycemia.

All system components, including the pump, transmitter, and sensor, must be removed before being screened with a full-body scanner. To avoid system removal, request an alternative screening method, if necessary.

- Carry the Medical emergency card provided with the device when traveling. The Medical emergency card provides critical information about airport security systems and pump use on an airplane. Not following the guidance on the Medical emergency card may result in serious injury.

## General precautions

The pump does not notify the user of leaks in the infusion set or degradation of insulin. If BG is too high, check the pump and the infusion set to confirm that the necessary amount of insulin is being delivered.

Check for adverse reactions where the pump comes into contact with skin. These reactions include redness, swelling, irritation, sensitization, rash, and other allergic

reactions. Do not allow the pump to come into contact with skin wounds, as the pump materials have only been evaluated for safe contact with intact skin.



**Note:** If you drop your pump, be sure to monitor your glucose levels for the next four hours.

### Waterproof capabilities

- The pump is waterproof at the time of manufacture and when the reservoir and tubing are properly inserted. It is protected against the effects of being underwater to a depth of up to 8 feet (2.4 meters) for up to 30 minutes.
- If the pump is dropped, hit against a hard object, or otherwise damaged, the waterproof characteristics of the outer casing of the pump may be compromised. If the pump is dropped or might be damaged, carefully inspect it to confirm that there are no cracks before exposing the pump to water.
- This waterproof capability rating applies only to the pump.
- If water may have entered the pump or other pump malfunction is observed, check BG and treat high BG as necessary using an alternative source of insulin. Contact 24-Hour Technical Support for further assistance, and consult a healthcare professional about high or low BG levels or with any other questions about care.

### Electrostatic discharge

- Very high levels of electrostatic discharge (ESD) can result in a reset of the pump's software and a pump error alarm. After clearing the alarm, confirm that the pump is set to the correct date and time, and that all other settings are programmed to the desired values. Following a pump reset, the SmartGuard feature may be unavailable for 5 hours to allow active insulin to be updated.
- For more information on pump alarms, see *Pump alarms, alerts, and messages*, page 291. Contact 24-Hour Technical Support with any problems entering pump settings.

## Extreme temperatures

Exposure to extreme temperatures can damage the device. Avoid the following conditions:

- Pump storage temperature above 122 °F (50 °C) or below -4 °F (-20 °C).
- Pump operating temperature above 98.6 °F (37 °C) or below 41 °F (5 °C).
- Insulin solutions freeze near 32 °F (0 °C) and degrade at temperatures higher than 98.6 °F (37 °C). In cold weather, wear the pump close to the body and cover it with warm clothing. In a warm environment, take measures to keep the pump and insulin cool.
- Do not steam, sterilize, autoclave, or otherwise heat the pump.

## Skin care products

Some skin care products, such as lotion, sunscreen, and insect repellents, can damage the plastic in the pump case. After using skin care products, wash hands prior to handling the pump. If a skin care product comes into contact with the pump, wipe it off as soon as possible with a damp cloth and mild soap. For instructions on cleaning the pump, see *Cleaning the pump*, page 283.

## Infusion sets and sites, transmitter, sensor, and meter

Refer to the corresponding device user guide for all warnings, precautions, and instructions relating to the device. Failure to reference the corresponding device user guide can result in minor injury, or damage to the device.

## Adverse reactions

Refer to the sensor user guide for adverse reactions related to sensor use. Failure to reference the sensor user guide may result in minor injury, or damage to the sensor.

## Security precautions

The MiniMed 780G insulin pump system is designed with security features to help keep the system and the data secure. These security features in the insulin pump system are set in the factory and ready to use when the insulin pump is received. For example, when the pump communicates with other devices in the system, such as a compatible BG meter, transmitter, or compatible mobile device, the data that it sends and receives



is encrypted and protected by cyclic redundancy checks. This helps prevent other people from being able to see system data, or to interfere with insulin pump therapy.

To help keep the system secure, follow these instructions:

- Do not leave the insulin pump or paired devices unattended.
- Do not share the pump, transmitter, sensor, or BG meter serial number.
- Do not connect the pump to any third-party devices not authorized by Medtronic.
- Do not use any software not authorized by Medtronic to control the system.
- Be attentive to pump notifications, alarms, and alerts because they may indicate that someone else is trying to connect to or interfere with the device.
- Disconnect the Blue Adapter from the computer whenever it is not being used.
- Use good cyber security practices; use anti-virus software and keep computer software up to date.
- Refer to the MiniMed Mobile App User Guide for information on how to keep the compatible mobile device safe for use with the Medtronic devices.

The pump only communicates with paired devices. The short time that it takes to pair the pump with other devices is a sensitive time for security. During this time, it is possible for an unintended device to pair with the pump. While Medtronic has designed security features into the system to prevent this, always follow these instructions to keep the system safe during pairing:

- Pair the transmitter, BG meter, or the compatible mobile device with the pump away from other persons and devices.
- When the transmitter successfully pairs with the pump, the green LED on the transmitter stops blinking. If the green LED on the transmitter continues to blink for several minutes or more after it is successfully paired, it may have been paired with an unintended device. See *Unpairing the transmitter from the pump*, page 139 to delete the transmitter from the pump and then follow the steps to pair it again.
- After pairing the BG meter or the compatible mobile device with the pump, make sure that the BG meter or compatible mobile device indicates that pairing was successful.

If there are symptoms of severe hypoglycemia or diabetic ketoacidosis or if unexpected changes of insulin pump settings or insulin delivery are suspected, consult a healthcare professional.

If there is a concern that someone else is trying to connect to or interfere with the device, stop using it and contact 24-Hour Technical Support immediately.

## Insulin guidelines



**WARNING:** Do not insert an insulin-filled reservoir into the pump, or connect an insulin-filled infusion set into the body while training with the system. Doing so may result in the unintentional infusion of insulin, which may result in hypoglycemia. Start insulin therapy only when directed by a healthcare professional.

The MiniMed 780G system has been studied with, and is intended for use with, the following rapid-acting U-100 insulins:

- U-100 Admelog<sup>TM\*</sup>
- U-100 Humalog<sup>TM\*</sup>
- U-100 NovoLog<sup>TM\*</sup>

Some insulin products are labeled for use in any pump that is compatible with the insulins listed in this section. To see if another insulin not listed in this section can be used, refer to section 2.2 of the prescribing information for that insulin product.



**WARNING:** Only use rapid-acting U-100 insulin from those listed in this section, as prescribed by a healthcare professional, in the MiniMed 780G system. Use of the incorrect type of insulin, or insulin with a greater or lesser concentration may result in over-delivery or under-delivery of insulin, which may result in hypoglycemia or hyperglycemia. Consult a healthcare professional with any questions about the type of insulin that is compatible with the pump.

Consumables

The pump uses disposable, single-use MiniMed and Medtronic reservoirs and infusion sets for insulin delivery.



**WARNING:** Only use reservoir and infusion sets manufactured or distributed by Medtronic Diabetes. The pump has undergone extensive testing to confirm appropriate operation when used with compatible reservoirs and infusion sets manufactured or distributed by Medtronic Diabetes. Medtronic Diabetes cannot guarantee appropriate operation if the pump is used with reservoirs or infusion sets offered by third parties and therefore Medtronic Diabetes is not responsible for any injury or malfunctioning of the pump that may occur in association with such use.

- **Reservoirs**–If using a Medtronic Extended infusion set, use the Medtronic Extended reservoir MMT-342, 3.0 mL (300-unit). Otherwise, use the MiniMed reservoir MMT-332A, 3.0 mL (300-unit).
- **Infusion sets**–Contact a healthcare professional for help in choosing a Medtronic Diabetes infusion set. Change the infusion set per the duration of use in the infusion set user guide.

The following table lists the compatible infusion sets. The MMT numbers may change if other compatible infusion sets become available.



**Note:** For infusion sets, MMT numbers that include an A (such as MMT-396A, MMT-396-AT) are compatible with the pump system. MMT numbers that do not include an A are no longer compatible with the pump system.

Type	MMT number
MiniMed Quick-set infusion set	MMT-386A, MMT-387A, MMT-394A, MMT-396A, MMT-397A, MMT-398A, MMT-399A
MiniMed Silhouette infusion set	MMT-368A, MMT-377A, MMT-378A, MMT-381A, MMT-382A, MMT-383A, MMT-384A

Type	MMT number
MiniMed Sure-T infusion set	MMT-862A, MMT-864A, MMT-866A, MMT-874A, MMT-876A, MMT-884A, MMT-886A
MiniMed Mio 30 infusion set	MMT-905A, MMT-906A
MiniMed Mio infusion set	MMT-921A, MMT-923A, MMT-925A, MMT-941A, MMT-943A, MMT-945A, MMT-965A, MMT-975A
MiniMed Mio Advance infusion set	MMT-213A, MMT-242A, MMT-243A, MMT-244A
Medtronic Extended infusion set	MMT-430A, MMT-431A, MMT-432A, MMT-433A, MMT-440A, MMT-441A, MMT-442A, MMT-443A

## Other MiniMed 780G system devices

- **Accu-Chek Guide Link meter**—The meter sends BG meter readings to the pump.
- **Guardian 4 transmitter (MMT-7841)**—The transmitter pairs with the pump, collects data measured by the sensor, and wirelessly sends this data to monitoring devices. This device is required for CGM.
- **Guardian 4 sensor (MMT-7040)**—The sensor is a disposable, single-use device inserted just below the skin to measure glucose levels in interstitial fluid. This device is required for CGM. The Guardian 4 glucose sensor is the only sensor compatible with the Guardian 4 transmitter.

## Optional items

The following items may be used with the MiniMed 780G system.

- **Pump clip**—The pump clip attaches to a belt and can be used to open the battery compartment.
- **Activity guard**—The activity guard helps to prevent the reservoir from being rotated or removed from the pump during physical activities.
- **MiniMed Mobile app (MMT-6101 for Android or MMT-6102 for iOS)**—The app provides a secondary display of insulin pump data and uploads system data to CareLink software. The app can check for eligible and available software updates for the pump. The Update pump feature in the app allows you to update the pump software remotely. The app can be installed on multiple mobile devices, but only one mobile device can be paired with the pump at a time.

- **CareLink Connect app (MMT-6111 for Android or MMT-6112 for iOS)**—The app can be downloaded onto compatible mobile devices from the app store. Refer to the app user guide for setup and operation within the app. This optional app is available to care partners to view patient therapy data and to be notified of selected patient alerts. This app does not replace the real-time display of insulin pump data on the primary display device. All therapy decisions should be based on the primary display device. Refer to the local Medtronic Diabetes website for information about supported devices and operating systems.
- **Blue Adapter**—The Blue Adapter uploads system data to CareLink software through a USB port on a computer. Refer to the CareLink software user guide for setup and operation of the Blue Adapter.



# 2







# 2

## System overview


In this chapter, you will learn about the components of the system and some important concepts and terminology that you will need to understand when using the system.

### What are the components of the MiniMed 780G system?

The following items are the main system components:

- **MiniMed 780G pump**—The pump delivers insulin into your body through the infusion set, based on the settings provided by your healthcare professional.
- **Infusion sets**—An infusion set connects to both the pump and your body. It carries the insulin as it is pushed out of the pump and delivers it.
- **Reservoirs**—The reservoir is filled with insulin and placed in the pump so that insulin can be delivered into your body through the infusion set.
- **Sensors and transmitter**—The sensor measures glucose in the fluid under your skin. The transmitter communicates with the pump through a wireless connection. The sensor and transmitter make up the continuous glucose monitoring (CGM) system.
- **Accu-Chek Guide Link meter**—Use this meter to measure the glucose in your blood. The meter sends this blood glucose (BG) information to your pump through a wireless connection.

The following diagram shows what the pump, meter, sensor, and transmitter look like and how you may wear them on your body. A diagram later in chapter 3 will show you more details about the infusion set and reservoir.

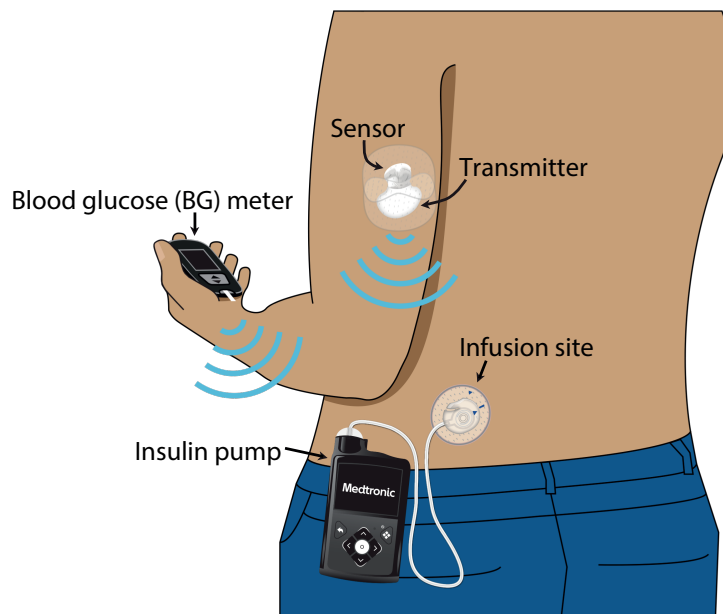

## Modes

Your pump operates in two different modes: Manual mode and SmartGuard mode (also called the SmartGuard feature).

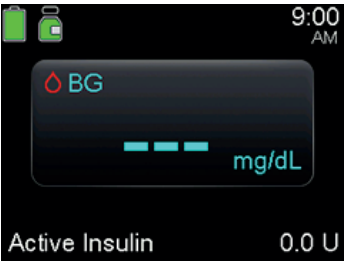
When you first use your MiniMed 780G insulin pump, it is in Manual mode.

Manual mode refers to a group of features that requires your input to deliver boluses for meals and to correct glucose levels. You may use Manual mode with or without CGM. When using CGM in Manual mode, you can see sensor glucose trends, receive low and high sensor glucose alerts, and suspend insulin delivery according to your settings.

After a few days of use in Manual mode, and at the direction of your healthcare provider, you can use SmartGuard. When in SmartGuard, the pump automatically adjusts and delivers basal insulin and can also deliver automatic correction boluses to regulate glucose levels to a target SG value. You will still need to enter carbs that you eat to deliver a food bolus.

The following tables show the main features of Manual mode and SmartGuard. There are details on each of these topics throughout this guide.

## Manual mode without CGM



Bolus delivery	Basal delivery	Suspend
<ul style="list-style-type: none"><li>• Bolus Wizard calculates a bolus based on your settings<ul style="list-style-type: none"><li>– A blood glucose (BG) meter reading is needed for a correction bolus</li><li>– A carb entry is needed for a food bolus</li></ul></li><li>• Manual bolus<ul style="list-style-type: none"><li>– You enter the number of units of insulin to cover food, high BG, or both</li></ul></li></ul>	<ul style="list-style-type: none"><li>• Programmed basal delivery settings</li><li>• A Temp basal rate can be used to temporarily increase or decrease basal insulin delivery</li></ul>	<ul style="list-style-type: none"><li>• Suspend all delivery<ul style="list-style-type: none"><li>– Choose this option to stop all delivery of insulin</li></ul></li></ul>

## Manual mode with CGM



Bolus delivery	Basal delivery	Suspend
<ul style="list-style-type: none"> <li>Same as Manual mode without CGM</li> </ul>	<ul style="list-style-type: none"> <li>Same as Manual mode without CGM</li> </ul>	<ul style="list-style-type: none"> <li>Suspend all delivery <ul style="list-style-type: none"> <li>Same as Manual mode without CGM</li> </ul> </li> <li>Suspend before low <ul style="list-style-type: none"> <li>Suspends insulin delivery and alerts based on your settings</li> </ul> </li> <li>Suspend on low <ul style="list-style-type: none"> <li>Suspends insulin delivery and alerts based on your settings</li> </ul> </li> </ul>

SmartGuard



Bolus delivery	Basal delivery	Suspend
<ul style="list-style-type: none"> <li>SmartGuard bolus feature delivers bolus insulin based on sensor glucose (SG) values and carb entries</li> <li>A blood glucose (BG) meter reading may be required when a sensor glucose (SG) value does not appear on the Bolus screen</li> </ul>	<ul style="list-style-type: none"> <li>The pump automatically delivers basal insulin based on recent insulin delivery needs, Sensor glucose (SG) values, and your glucose target</li> <li>A Temp target can be set when less insulin is needed, such as for exercise</li> </ul>	<ul style="list-style-type: none"> <li>Suspend all delivery <ul style="list-style-type: none"> <li>Same as Manual mode without CGM</li> </ul> </li> </ul>

Bolus delivery	Basal delivery	Suspend
<ul style="list-style-type: none"> <li>The bolus amount cannot be adjusted</li> <li>The pump may automatically deliver an Auto Correction bolus to maximize the time in range.</li> </ul>		

### Delivery settings

The delivery settings table describes whether or not each setting applies to the SmartGuard feature and Manual mode. Consult your healthcare professional before changing delivery settings.

Delivery setting	Impacts SmartGuard therapy	Impacts Manual mode
Active insulin time	Yes	Yes
Basal pattern and basal rates	No	Yes
BG target in Bolus Wizard	No	Yes
Bolus increment	No	Yes
Bolus speed	Yes	Yes
Carb ratio	Yes	Yes
Dual/Square Wave bolus	No	Yes
Insulin sensitivity factor	No	Yes
Max basal	No	Yes
Max bolus	No	Yes
Preset bolus	No	Yes
Preset temp	No	Yes



# 3







# 3

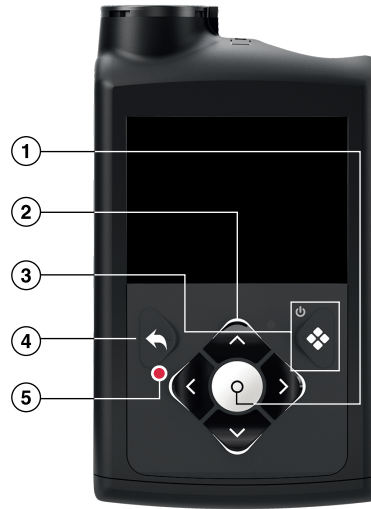
## Pump basics

This chapter provides information about the basic features, buttons, and screens of the MiniMed 780G insulin pump.













**CAUTION:** Do not use sharp objects to press the pump buttons. The use of sharp objects can damage the pump.

## Using the buttons

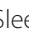


The following table describes the notification light and how to use the pump buttons.

Item	Description
①	Press  to go to the Menu screen from the Home screen and to select a highlighted menu option.
②	Press  or  to scroll up or down, highlight an item on a screen, and to increase or decrease the value of a setting. Press  or  to move left or right on certain screens and to highlight the icons on the Menu screen.
③	Press  to access the Graph screen. Press and hold  to put the pump in Sleep mode.
④	Press  to go back to the previous screen. Press and hold  to return to the Home screen.
⑤	The notification light  flashes when the pump has an alarm or alert. The notification light is not visible unless it flashes.


## Sleep mode

The pump screen goes dark if there is no activity for the duration specified in the Backlight setting. Two minutes after the screen goes dark, the pump enters Sleep mode to conserve battery power. Sleep mode does not affect insulin delivery. Press any

button to wake up the pump. Press and hold  for two seconds to manually enter Sleep mode.

## Pump delivery system

The following diagram shows the parts of the pump delivery system, including the infusion set\*, reservoir, and pump.

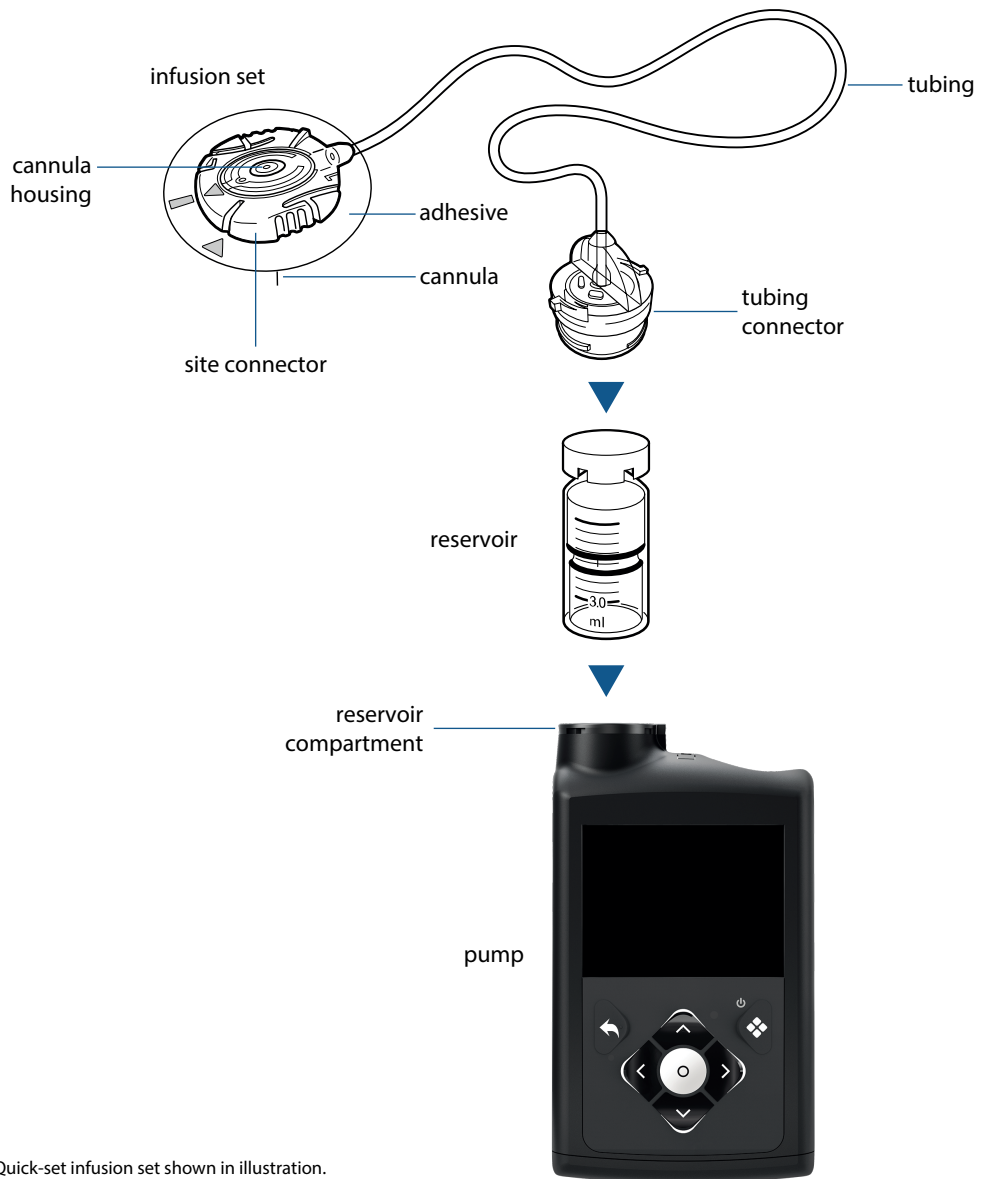

\*Quick-set infusion set shown in illustration.

## Infusion set

The infusion set consists of the following components:

- The tubing carries insulin from the reservoir into the body.
- The tubing connector attaches to the reservoir.
- The insertion piece attaches to the body.

- The cannula is a small, flexible tube inserted into the body. Some infusion sets use a small needle instead of a cannula.
- The adhesive holds the infusion set in place.

Change the infusion set according to the user guide provided with the infusion set.

## Reservoir

The reservoir stores insulin for delivery and is inserted into the pump reservoir compartment.

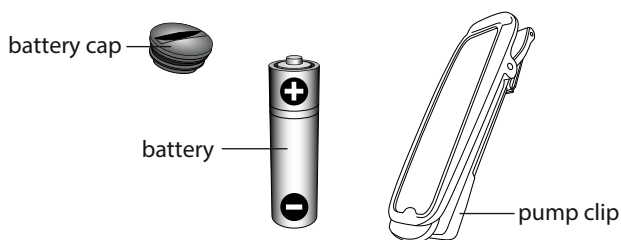
## Pump

Underneath the reservoir compartment, a piston pushes up on the bottom of the reservoir to move insulin into the tubing, through the cannula, and into the body.

The pump delivers small doses of insulin. The smallest dose of insulin is 0.025 units. The piston inside the pump must be rewound each time a newly filled reservoir is inserted into the reservoir compartment.

## Battery cap

The battery cap is located in the pump box with the accessories. It seals the battery compartment and ensures the pump is connected to the battery.

















## Identifying the correct battery cap for your pump

Only use the correct battery cap that came with your pump. Using an incorrect battery cap may cause your pump to not seal properly or may result in an incomplete battery connection.

The correct battery cap depends on your pump model. Identify your pump model by looking on the back of your pump to see if the serial number is laser-etched or printed

on a silver sticker. Use the diagram to identify the correct battery cap for your pump model.

			
Laser-etched serial number		Silver sticker with printed serial number	
Pump model (ACC-1528)		Pump model (ACC-1529)	
<b>✓Correct battery cap</b>		<b>✗Incorrect battery cap</b>	<b>✓Correct battery cap</b>
	 ✓No frame around "AA"		 ✗No frame around "AA"
	 ✓Large metal battery contact		 ✗Small metal battery contact
	 ✓Frame around "AA"		 ✓Large metal battery contact

If you need a spare battery cap, please contact your local Medtronic support representative for assistance, questions, or concerns.

Inserting the battery

The pump requires one new AA (1.5 V) battery. For best results, use a new AA lithium (FR6) battery. The pump also accepts an AA alkaline (LR6) or a fully charged AA NiMH (HR6) nickel-metal hydride rechargeable battery.



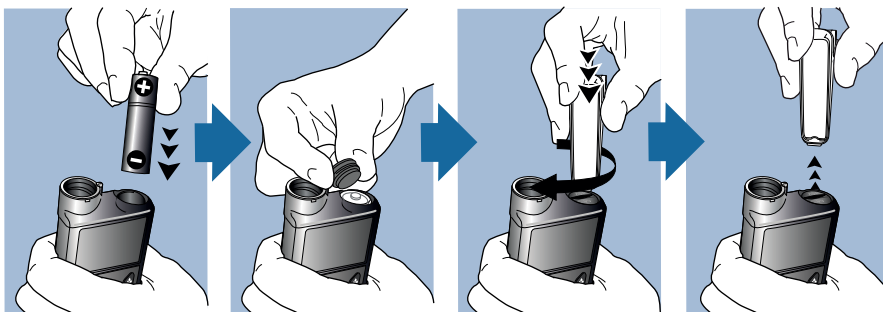
**CAUTION:** Do not use a carbon zinc battery in the pump. Carbon zinc batteries are not compatible with the pump and can cause the pump to report inaccurate battery levels.



**Note:** Do not use cold batteries because the battery life may incorrectly appear low. Allow cold batteries to reach room temperature before they are inserted into the pump.

## To insert the battery:


1. Insert a new or fully charged AA battery. Make sure to insert the negative end (–) first.



2. Place the battery cap onto the pump. Use the bottom edge of the pump clip or a coin to tighten the cap.



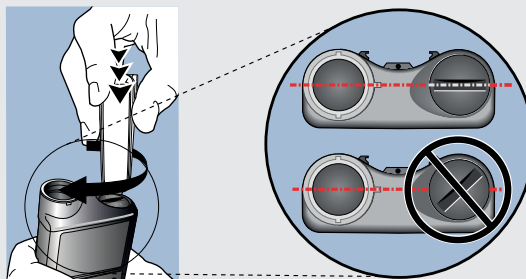**CAUTION:** Do not overtighten or undertighten the battery cap. A battery cap that is too tight can cause damage to the pump case. A battery cap that is too loose can prevent detection of the new battery. Turn the battery cap clockwise until the cap slot is aligned horizontally with the pump case, as shown in the following example.



The first time a battery is inserted into the pump, the Startup Wizard begins. Any other time a battery is inserted into the pump, the Home screen appears and the pump resumes basal delivery.


## Removing the battery



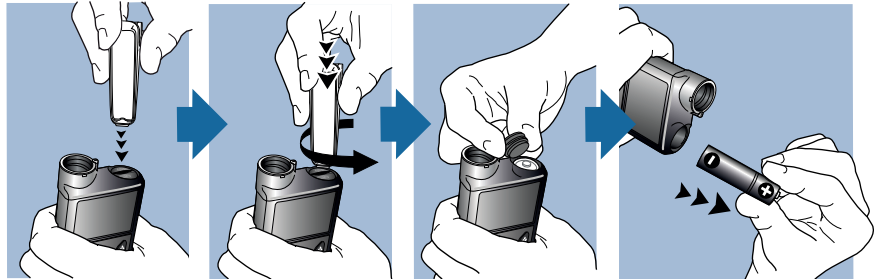**CAUTION:** Do not remove the battery unless a new battery needs to be inserted or to store the pump. The pump cannot deliver insulin while the battery is removed. After an old battery is removed, make sure to replace it with a new battery within 10 minutes to clear the Insert battery alarm and avoid a Power loss alarm. If power loss occurs for an extended period of time, the time and date settings must be re-entered.

### To remove the battery:

1. Before a battery is removed from the pump, clear any active alarms or alerts.
2. Use the pump clip or a coin to loosen and remove the battery cap.
3. Remove the battery.



4. Dispose of old batteries in an appropriate container and in accordance with local laws for battery disposal.
5. After a battery is removed, wait until the Insert Battery screen appears before inserting a new battery.

If a battery is removed to place the pump in storage, see *Storing the pump*, page 284 for more information.



## Startup settings

The Startup Wizard appears after a battery is inserted for the first time. Use the Startup Wizard to set the language, time format, current time and date, and to rewind the pump. To re-enter these settings later, see *Pump issues*, page 276.

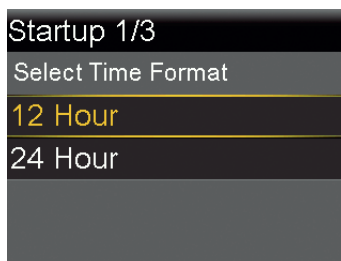
### To use the Startup Wizard:

1. On the Select Language screen, select a language.




The Select Time Format screen appears.

2. Select a time format.



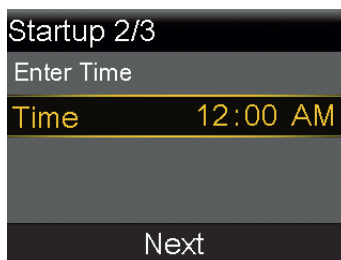3. Enter the current time, and then select **Next**.



The Enter Date screen appears.

4. Enter the current date, and then select **Next**.

Startup 3/3	
Enter Date	
Year	2021
Month	Jan
Day	1, Fri
Next	

A "Rewinding" message appears. The piston returns to its start position in the reservoir compartment. This may take several seconds.



When rewinding is complete, a message appears to confirm the startup is complete.

5. Select **OK** to go to the Home screen.


Startup
Congratulations! Startup Wizard complete
✓
OK

## Home screen in Manual mode

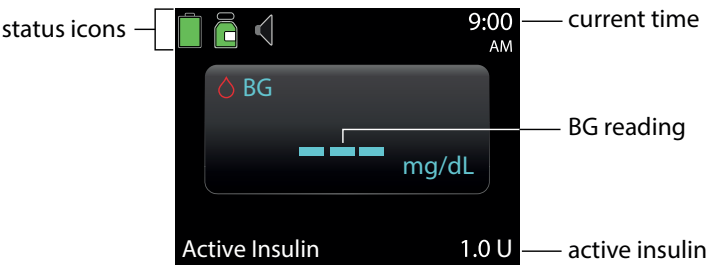The Home screen appears after the battery is changed, when the pump wakes from Sleep mode, and when another screen is not actively being used.



**Note:** This example shows the Home screen in Manual mode when the Sensor feature is turned off. For information about the Home screen when the Sensor feature is turned on, see *Home screen with CGM in Manual mode, page 174*. For information about the Home screen with the SmartGuard feature, see *Home screen with the SmartGuard feature, page 187*.






The following items appear on the Home screen:

Item	Description
Status icons	The status icons show a quick status of the pump system. For more information, see <i>Status icons, page 72</i> .
Current time	For details on setting the time, see <i>Time and date, page 203</i> .
BG readings	The current blood glucose (BG) meter reading is shown. The BG is either entered manually, or it is sent from a paired Accu-Chek Guide Link meter.
Active insulin	Active insulin is bolus insulin delivered by the insulin pump that continues to lower BG levels. Active insulin is not necessarily reflective of the pharmacokinetics and pharmacodynamics of rapid acting insulins. For more details on active insulin, see the description of Active Insulin Time in <i>Bolus Wizard feature in Manual mode , page 99</i> .








## Shortcuts from the Home screen














The following table describes shortcuts that can be used to quickly access certain pump functions. These shortcuts only work on the Home screen.









Shortcut	Description
	Press this button to access the Status screen.
	Press this button to access the Time in Range screen when the Sensor feature is turned on.
	Press this button to access the Bolus screen. The Bolus screen that appears varies depending on the bolus feature that is currently active.

## Status icons

The status icons provide the current status of the pump system. For information on viewing detailed status screens, see *Status screen, page 76*.

Icon name	Description
Active Insulin re-set to zero	After the Active Insulin reset to zero alarm occurs,  appears on the Home screen and Bolus screens until the time shown in the alarm. For more information, see <i>Pump issues, page 276</i> .
Alert silence	The Alert silence icon  indicates that the Alert Silence feature is turned on and some alerts will not make a sound or vibration. Sensor alerts can be silenced for a specific duration using the Alert silence feature. For more information, see <i>Silencing sensor alerts, page 168</i> .
Battery	<p>The color and fill level of the icon indicate the charge level of the pump battery. As the battery is used, the icon changes from solid green in the following order: </p> <ul style="list-style-type: none"><li> The battery is full.</li><li> The battery is low.</li><li> The battery can be used for less than 30 minutes and needs to be replaced.</li></ul>
Block mode	The Block mode icon  shows that the pump is locked. For more information about Block mode, see <i>Block mode, page 204</i> .
Connection	The connection icon shows the following information:


Icon name	Description
	<ul style="list-style-type: none"> <li>•  The Sensor feature is on and communicating.</li> <li>•  The Sensor feature is on, but the sensor is not communicating with the pump.</li> </ul>
Reservoir	<p>The reservoir icon shows the fill status of the MiniMed or Medtronic 3.0 mL (300-unit) reservoir.</p> <ul style="list-style-type: none"> <li>•  Approximately 85%–100% of the insulin remains in the reservoir.</li> <li>•  Approximately 71%–84% of the insulin remains in the reservoir.</li> <li>•  Approximately 57%–70% of the insulin remains in the reservoir.</li> <li>•  Approximately 43%–56% of the insulin remains in the reservoir.</li> <li>•  Approximately 29%–42% of insulin remains in the reservoir.</li> <li>•  Approximately 15%–28% of the insulin remains in the reservoir.</li> <li>•  Approximately 1%–14% of the insulin remains in the reservoir.</li> <li>•  The amount of insulin remaining in the reservoir is unknown.</li> </ul>
Sensor life	<p>The number on the sensor life icon indicates the number of days that remain in the life of the sensor. When the Sensor feature is turned on, the icon appears on the Status screen. After a new sensor is inserted, the icon is solid green. When one day remains in the life of the sensor, the icon turns red and appears on the Home screen. When the sensor expires, the icon turns solid black with an X.</p> <div data-bbox="639 1067 1130 1133">  </div> <p>If the number of days that remain in the life of the sensor is not yet available, such as when the sensor is warming up, the sensor life icon appears with three dots. </p> <p>If the number of days that remain in the life of the sensor is unknown, the sensor life icon appears with a question mark. </p>
Sensor status	<p>The sensor status icon shows whether the sensor is warming up, is monitoring sensor glucose (SG) values, a BG is required, or the sensor status is unavailable. The icon appears only when the Sensor feature is turned on.</p>

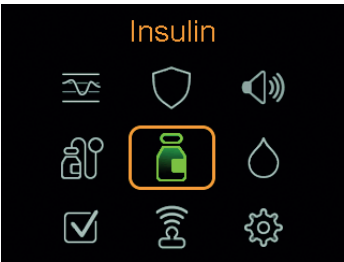
Icon name	Description
	<ul style="list-style-type: none"> <li> An icon with a green circle around it means the sensor is working and no action is required.</li> <li> A red icon means a BG reading is required.</li> <li> A question mark icon with a blue circle around it means that sensor information is unavailable.</li> <li> An icon with three white dots on a black background means the pump is waiting for the sensor status to update.</li> </ul>
Suspend by sensor	<p>When the current low alert time segment has either the Suspend before low or Suspend on low feature turned on, the Suspend by sensor icon appears on the Home screen.</p> <p></p> <p>When Suspend before low or Suspend on low feature suspends insulin delivery, the icon flashes.</p> <p></p> <p>If the Suspend before low or Suspend on low feature is turned on but unavailable, the icon has a red X.</p> <p>This can be due to a recent suspend by sensor event or when no SG values are available.</p> <p></p> <p>For more information, see <i>The Suspend before low feature, page 160</i> and <i>The Suspend on low feature, page 162</i>.</p>
Temporary network connection	<p> The temporary network connection icon shows when the pump is temporarily connected to a remote upload device.</p>












**Note:** Status icons provide limited information. For example, the reservoir icon may indicate the reservoir is low on insulin. The Status screen shows more detailed information about how many units are left. For more information about the status screens, see *Status screen, page 76*.

## Menu screen

Use the menu to go to screens that show various features and functions of the system. Press  from the Home screen to go to the menu. The highlighted menu option appears in color. All other menu options appear in black and gray.



Use the menu to go to the following screens:

Menu selection	Menu icon	Description
Insulin		Deliver a bolus, set up and deliver basal insulin, suspend insulin delivery, and stop bolus during bolus delivery.
History & Graph		View history, SG review, graph, and time in range.
SmartGuard		Set up the SmartGuard feature.
Sound & Vibration		Set sound, vibrate, and volume options for notifications.
Reservoir & Set		Set up a new reservoir and infusion set, and fill a cannula.
Blood Glucose		Enter a BG value.
Status		View the status of the pump and other system features.
Paired Devices		Pair devices or CareLink software.
Settings		Set up device settings, delivery settings, and alert settings.










Menu map

Refer to *Menu map*, page 349 to see the menu map diagrams.

Sound & Vibration screen

The sound and vibration options are set on the Sound & Vibration screen. Sensor alerts can also be temporarily silenced. For information about silencing alerts, see *Silencing sensor alerts*, page 168. The alert silence icon on the Home screen indicates when alerts are silenced. For more information, see *Status icons*, page 72.

**To adjust the sound and vibration settings:**

- 1. From the Home screen, press , and then select .
- 2. Adjust the volume:
  - a. Select **Volume**.
  - b. Press , , , or , and then press .
- 3. Select **Sound**, and then press  to turn the sound on or off.
- 4. Select **Vibration**, and then press  to turn the vibration on or off.

**Status screen**

The Status screen provides access to information about the pump and information about the sensor, if applicable. The Status screen also provides the option to suspend all insulin delivery or resume basal insulin delivery.

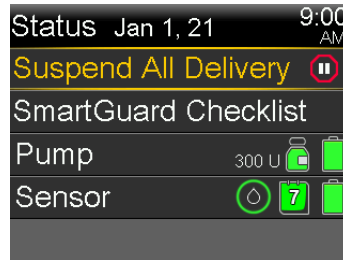
Use the Status screen to access the following screens or options:

Screen or option	Description
Stop Bolus	This option appears when a bolus delivery is in progress. Select <b>Stop Bolus</b> to stop the active bolus.
Suspend All Delivery or Resume Basal	This option indicates whether insulin delivery is currently suspended. Select <b>Suspend All Delivery</b> to suspend insulin delivery. Select <b>Resume Basal</b> to resume basal insulin delivery. For more information see <i>Suspending all insulin delivery and resuming basal insulin delivery, page 91</i> .
SmartGuard Checklist screen	This screen shows a list of the required conditions before the pump can use the SmartGuard feature. For more information, see <i>SmartGuard Checklist, page 184</i> .
Pump status screen	This screen shows a detailed view of the pump status, the reservoir and infusion set status, battery status, pump serial number, pump name, model number, and other pump details.
Sensor status screen	This screen appears when the Sensor feature is turned on. The Sensor status screen includes sensor life, transmitter battery life, and shows the serial number and version number of the transmitter.



### To view the status screens:

1. From the Home screen, press , and then select .



2. Press  or  to highlight a status screen, and then press .



# 4



Basal and bolus insulin delivery



# 4

## Basal and bolus insulin delivery

This chapter explains how to use different types of insulin delivery.

### Setting up basal insulin

Basal insulin is the “background” insulin that the body needs throughout the day and night to maintain target blood glucose (BG) meter readings when food is not eaten. Basal insulin accounts for approximately one half of daily insulin requirements. The MiniMed 780G insulin pump simulates a pancreas by delivering insulin continuously over 24 hours.



**WARNING:** The pump is intended to be used with a basal pattern. The basal pattern must be manually entered and saved into the pump. The pump will operate with a basal rate of 0.0 U/hr until a basal pattern is entered and saved. There is no reminder message to program basal rates. Consult a healthcare professional to determine what basal pattern is needed. For more information about basal patterns, see *Basal patterns, page 84*.

### Basal rate

Basal rate is the specific amount of basal insulin that the pump continuously delivers each hour. While some people use one basal rate all day, others require different rates at different times of the day.

Basal rates are set in one or more basal patterns. Each basal pattern covers 24 hours. For specific information about basal patterns, see *Basal patterns, page 84*.

## Max basal rate

The Max basal rate setting applies to Manual mode only and does not affect how much insulin SmartGuard delivers. When using the SmartGuard feature, delivery limits are determined automatically.

The Max basal rate setting limits the maximum amount of basal insulin per hour that can be programmed. Consult with your healthcare professional to personalize your Max basal rate setting.

It is not possible to set a basal rate, a temp basal rate, or a preset temp basal rate that would exceed the Max basal rate limit. After the basal patterns or preset temp basal rates are set, the Max basal rate setting cannot be set lower than any existing basal rate.



If you require a Max basal rate higher than 6 units per hour, a message appears indicating that a Max basal rate higher than 6 units per hour is outside of the typical range. A Max basal rate of 6 U/hr can allow basal delivery of up to 144 units per day in Manual mode. This amount of insulin may not be safe. Consult with your healthcare professional to ensure your basal setting is appropriate for you.

If you require a Max basal rate higher than 10 units per hour, the system will not let you increase the Max basal rate setting from **Delivery Settings > Max Basal Rate**. The only time you will see a rate higher than 10 U/hr through **Delivery Settings > Max Basal Rate** is when you have either:

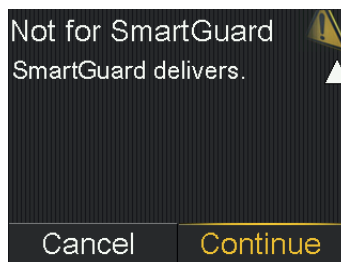
- Updated your pump software using the **MiniMed Mobile app**, and your previous Max basal rate was larger than 10 U/hr.
- Increased your Max basal rate as directed by your healthcare professional.

For more information on how to increase your Max basal rate above 10 U/hr, see *Max basal rate*, page 210.

### To set the Max basal rate:


1. From the Home screen, press , and then select .
2. Select **Delivery Settings > Max Basal/Bolus**.

The Not for SmartGuard screen appears.




3. Select **Continue**.

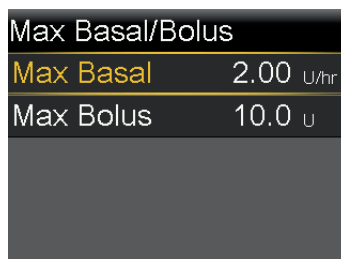The Max Basal/Bolus screen appears.



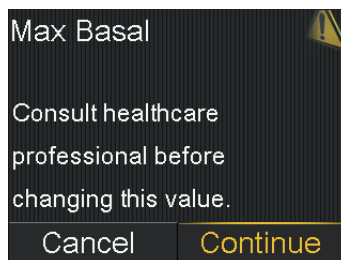4. Select **Max Basal**.



5. To continue to the Max Basal Rate screen, select **Continue**.
6. Select **Max Basal**, and then set the maximum number of basal insulin units per hour.

Max Basal Rate	
Max Basal	2.00 U/hr
This applies to Manual mode. SmartGuard delivery not affected.	
Save	

- If you lower the Max basal rate setting below an existing basal rate, a Max Basal Error will appear.
- If you raise the Max basal rate setting above 6 U/hr, you are warned that the amount of insulin may not be safe.
- The system will stop you from raising the Max basal rate setting above 10 U/hr. If your healthcare professional has increased your rate above 10 U/hr, see *Max basal rate*, page 210.

7. Select **Save**.


## Basal patterns

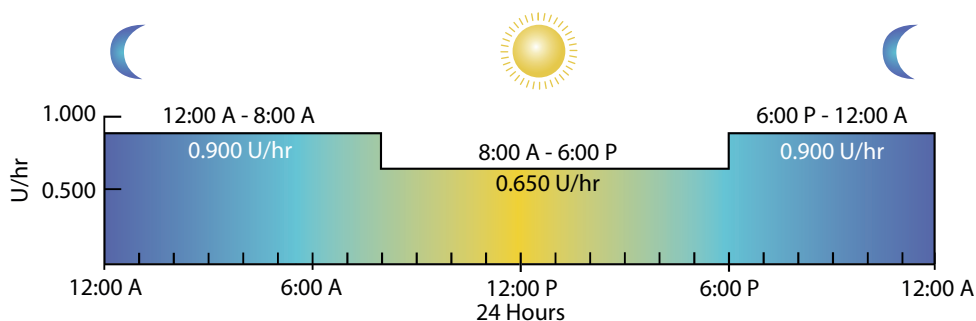Basal pattern settings apply to Manual mode only and do not affect how much insulin SmartGuard delivers.

The basal pattern determines the amount of basal insulin delivered throughout the day and night while in Manual mode. A basal pattern is made up of one to 48 basal rates that are set to cover a full 24-hour period. Because basal insulin needs can vary, up to eight basal patterns can be set.

The following example represents one basal pattern with three basal rates set for three different time periods.





Consult a healthcare professional to determine the basal pattern. The basal pattern must be manually entered into the pump. There will be no reminder message to program basal rates.





**WARNING:** Confirm a basal pattern is entered. If a basal pattern is needed but not entered and saved, this could result in an under-delivery of basal insulin. Under-delivery of insulin can potentially cause severe hyperglycemia, which may lead to diabetic ketoacidosis.

## Setting up a basal pattern

This procedure shows how to set up a basal pattern for the first time. To add an additional basal pattern, see *Adding an additional basal pattern*, page 242.


### To set up a basal pattern:

1. From the Home screen, press , and then select .
2. In Manual mode, select **Delivery Settings** > **Basal Pattern Setup**.
3. The system reminds you that basal patterns only apply to Manual mode. Select **Continue**.

Basal Pattern Setup		
Basal 1	0.0	U ✓
Add New		

4. Select **Basal 1**.
5. Select **Options**, and then select **Edit**.

Edit Basal 1		
Start	End	U/hr
12:00 A	12:00 A	0.025
Review		

6. For one basal rate, the End time does not need to change. Press  on the 12:00 A.




**Note:** For instructions on setting up multiple basal rates over a 24-hour period, see *Settings covering a 24-hour period*, page 87.


7. Enter the unit value for the time period.
8. Select **Review**.

If your basal pattern rate will deliver significantly more insulin than you typically need, the system displays a message. You can edit the basal pattern or, after consulting with a healthcare professional, you can continue.

Basal 1		
24 hr Total: 24 U		
Start	End	U/hr
12:00 A	12:00 A	1.00
Save		

Review the basal pattern. Press  to return to the previous screen to make changes.



**Note:** If  is pressed and **Save** is not selected, the changes are not saved.

9. Select **Save**. If you do not select Save, your changes are not saved.  
If this is an added basal pattern and you want to activate it, see *Changing from one basal pattern to another*, page 244.



**CAUTION:** If you have not pressed Save after settings are entered and the screen goes dark, the entered settings will not be saved.



**Note:** Programming a basal pattern is an important part of setting up the insulin pump for use. Please review the settings to confirm that these are programmed accurately based on settings provided from a healthcare professional.

## Settings covering a 24-hour period

Some pump functions allow settings to change over a 24-hour period. Basal rates are one of those settings.

Setting up multiple values over a 24-hour period applies to the following settings:

- Basal patterns  
See *Setting up a basal pattern*, page 85

- High SG settings  
See *Setting up the high SG settings, page 155*
- Low SG settings  
See *Setting up the low SG settings, page 165*
- Carb ratios, insulin sensitivities, and BG targets in the Bolus Wizard feature in Manual mode  
See *Setting up the Bolus Wizard feature in Manual mode, page 101*

This screen is an example of a basal pattern with different rates of basal insulin for specific times of the day:

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	8:00 <sub>A</sub>	0.900
8:00 <sub>A</sub>	6:00 <sub>P</sub>	0.650
6:00 <sub>P</sub>	12:00 <sub>A</sub>	0.900
Review		

### To set up values over a 24-hour period:

1. On the appropriate settings screen, select the End time and enter the end time for the first time period. In this example, the first desired time period is 8 hours. The start time always begins at 12:00 A. To set an 8-hour period, enter an end time of 8:00 A.

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	12:00 <sub>A</sub>	0.025
Review		

2. Enter the unit value for the first time period.

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	8:00 <sub>A</sub>	0.900
Review		

3. Press .

The start time for the next time period appears.

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	8:00 <sub>A</sub>	0.900
8:00 <sub>A</sub>	8:30 <sub>A</sub>	---
Review		

4. Enter the end time for the next time period.

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	8:00 <sub>A</sub>	0.900
8:00 <sub>A</sub>	6:00 <sub>P</sub>	---
Review		

5. Enter the unit value for the next time period.

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	8:00 <sub>A</sub>	0.900
8:00 <sub>A</sub>	6:00 <sub>P</sub>	0.650
Review		

6. Press .

The start time for the next time period appears.

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	8:00 <sub>A</sub>	0.900
8:00 <sub>A</sub>	6:00 <sub>P</sub>	0.650
6:00 <sub>P</sub>	12:00 <sub>A</sub>	---
Review		


7. Repeat steps 3-5 for every desired time period until the end time of 12:00 A is reached. This completes the 24-hour duration.

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	8:00 <sub>A</sub>	0.900
8:00 <sub>A</sub>	6:00 <sub>P</sub>	0.650
6:00 <sub>P</sub>	12:00 <sub>A</sub>	0.900
Review		


8. Select **Review**.

If your basal pattern rate will deliver significantly more insulin than you typically need, the system displays a message. You can edit the basal pattern or, after consulting with a healthcare professional, you can continue.

Edit Basal 1		
Start	End	U/hr
12:00 <sub>A</sub>	8:00 <sub>A</sub>	0.900
8:00 <sub>A</sub>	6:00 <sub>P</sub>	0.650
6:00 <sub>P</sub>	12:00 <sub>A</sub>	0.900
Review		

Review the basal pattern. Press  to return to the previous screen to make changes.





**Note:** If  is pressed and **Save** is not selected, the changes are not saved.

9. Select **Save**.



## Viewing basal delivery information

### To view the current basal rate:

1. From the Home screen, press , and then select .
2. Select **Basal**.

The current basal rate appears at the top of the screen.

### To view basal patterns:

1. From the Home screen, press , and then select .
2. Select **Basal**.
3. Select **Basal Patterns**.

The Basal Patterns screen shows a list of configured basal patterns and the 24-hour insulin total for each basal pattern. A check mark appears next to the active basal pattern.

4. To view details for a basal pattern, select the basal pattern.

For more information about basal patterns, see *Basal patterns*, page 84.

## Suspending all insulin delivery and resuming basal insulin delivery

Use this feature to suspend all active basal and bolus insulin deliveries. A reminder that insulin is not being delivered occurs every 15 minutes while this feature is active. The pump beeps, vibrates, or both every 15 minutes as a reminder that insulin is not being delivered.



**Note:** The first reminder occurs 15 minutes after the pump display times out. The pump beeps, vibrates, or both 15 minutes after the pump display times out. If a button is pressed to wake up the pump, the pump beeps, vibrates, or both 15 minutes after the pump display times out again. To adjust the timeout setting, see *Display options*, page 203.

To restart basal insulin delivery, use the Resume Basal feature. The pump starts the programmed basal pattern but does not start any previously programmed bolus deliveries.



**Note:** To stop a bolus delivery without stopping the basal delivery, see *Stopping a bolus delivery*, page 108.




**WARNING:** If insulin delivery is suspended during a bolus, check the pump daily history to determine the amount of insulin that was delivered before resuming insulin delivery. Bolus delivery and fill cannula do not restart when insulin delivery is resumed. If needed, program a new bolus or fill the cannula. Failure to resume basal insulin delivery can result in hyperglycemia and diabetic ketoacidosis.



**WARNING:** Do not rely solely on the sound or vibration notifications when using the sound or vibrate options. These notifications may not occur as expected if the speaker or vibrator in the pump malfunctions. A missed notification may result in the delivery of too much or too little insulin. This is most common when using the Easy bolus feature or when the pump is in manual suspend. Contact 24-Hour Technical Support with any concerns.

### To suspend all insulin delivery:

1. From the Home screen, press , and then select .
2. Select **Suspend All Delivery**.

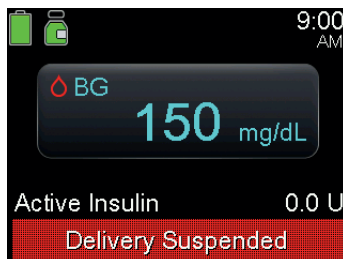


A confirmation message appears.



3. Select **Yes** to suspend all insulin delivery.

The pump functions are limited until insulin delivery is resumed.


The Delivery Suspended banner appears on the Home screen while insulin is suspended.



#### To resume basal insulin delivery:

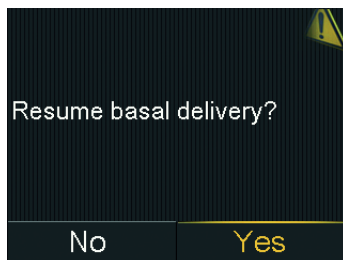1. While insulin delivery is suspended, from the Home screen press , and then select .
2. Select **Resume Basal**.

A confirmation message appears.



3. To resume basal insulin delivery, select **Yes**.

If a temp basal was active when the pump was suspended, it resumes, provided the time is still within the duration set.



**Note:** If a bolus delivery that was in progress before delivery was suspended is needed, check the Daily History screen for the actual bolus units delivered and the intended bolus amount. Then set up a new bolus amount as needed. For details about using the Daily History screen, see *Daily History screen, page 219*.

## Temp basal rates

The temp basal feature is used to set and start a temporary basal rate that can be used immediately to manage blood glucose (BG) during short-term activities or conditions, such as exercise or meals.

The duration of the temp basal rate can range from 30 minutes to 24 hours. After the temp basal rate delivery is completed or canceled, the programmed basal pattern resumes. The temp basal rates and preset temp basal rates can be defined using either a percentage of the current basal pattern or by setting a specific rate, as described in the table:



Temp basal rate type	Description
Percent	Percent delivers a percentage of the basal rates programmed in the active basal pattern for the duration of the temp basal rate. The temp basal amount is rounded down to the next 0.025 units if the basal rate is set at less than 1 unit per hour, or to the next 0.05 units if the basal rate is set at more than 1 unit per hour. Temp basal rates can be set to deliver from 0% to 200% of the scheduled basal rate. The percentage used is based on the largest basal rate scheduled during the temp basal rate duration and is limited by the Max basal rate.
Rate	Rate delivers a fixed basal insulin rate in units per hour for the duration of the temp basal rate. The amount set is limited by the Max basal rate.

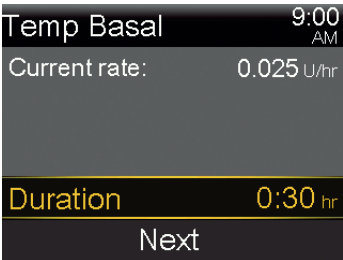
Preset temp basal rates can be set for recurring short term situations. For more information on Preset temp basal rates, see *Preset temp basal rates, page 239*.

Starting a temp basal rate


When a temp basal rate starts, basal delivery changes to the temp basal rate for the set duration. When the duration completes, the basal insulin returns to the active basal pattern.

To start a temp basal rate:

- 1. From the Home screen, press , and then select .
- 2. Select **Basal** > **Temp Basal**.
- 3. Set the **Duration**.



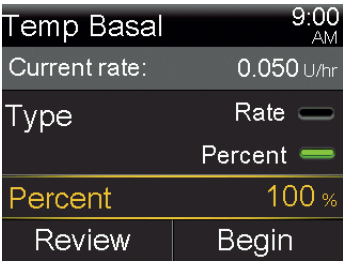- 4. Select **Next**.
- 5. Select **Type** to select Rate or Percent.



- 6. Depending on the type selected, do one of the following:
  - Enter a percentage.
  - Enter a basal rate.

Select **Review** to review the temp basal setting.

- 7. Select **Begin** to start the temp basal rate.

The Temp Basal banner appears on the Home screen during delivery.





## Entering a blood glucose (BG) meter reading

The system may request a blood glucose (BG) meter reading to continue use of sensor. Additionally, a blood glucose (BG) meter reading can be entered at any time, if desired.

The BG screen allows manual entry of a blood glucose (BG) meter reading. Previously entered manual or meter BG readings do not appear on the BG screen. A blood glucose (BG) meter reading received from a linked meter appears in a separate BG Meter screen that requires confirmation.

### To manually enter blood glucose (BG) meter readings:

1. From the Home screen, press , and then select .
2. Enter a blood glucose (BG) meter reading. Do not enter a sensor glucose (SG) value in place of a blood glucose (BG) meter reading. A blood glucose (BG) meter reading must always come from a blood glucose meter. The entered glucose value is used to calibrate the sensor.
3. Select **Save**.

### To manually enter blood glucose (BG) meter readings on the Bolus Wizard screen in Manual mode:

From the Bolus Wizard screen in Manual mode, select **BG**.

### To confirm a blood glucose (BG) meter reading from a blood glucose meter:

When the BG Meter screen with the message Confirm BG? shows, select **Yes** to confirm the blood glucose (BG) meter reading.

The BG received message shows.

## Setting up bolus delivery

A bolus is given for two reasons: to cover food that contains carbohydrates or to correct glucose levels that are above the target range.

### Manual mode bolus delivery options

Three types of bolus are available using the Bolus Wizard or Manual bolus feature, including normal bolus, Square Wave bolus, and Dual Wave bolus. Discuss these options with a healthcare professional to determine what is best. For more information, see *Bolus types*, page 249.



**Note:** Do not use a blood glucose (BG) meter reading if more than 12 minutes have passed since the last BG meter reading was taken. That BG meter reading and the calculated bolus amount may no longer be accurate.

The following table describes how to deliver a bolus using the Bolus Wizard feature or Manual bolus feature. These bolus options are only available in Manual mode.

Feature	Description
Bolus Wizard feature in Manual mode	Enter the BG meter value or the amount of carbs expected from a meal, or both. Then the Bolus Wizard feature in Manual mode calculates an estimated bolus amount based on the individual settings. For details about using the Bolus Wizard feature, see <i>Bolus Wizard feature in Manual mode</i> , page 99.
Manual bolus feature	Calculate and manually enter the bolus amount. For details about using the Manual bolus feature, see <i>Delivering a normal bolus using the Manual bolus feature</i> , page 107.

### Max bolus

The Max bolus setting limits the amount of insulin that can be programmed by the user for a single bolus in Manual mode. The pump prevents single bolus insulin deliveries


that exceed the Max bolus amount. The Max bolus can be set from 0 to 25 units. Set the Max bolus value as indicated by a healthcare professional.

If the Max bolus is set up after the preset bolus deliveries are set, the Max bolus cannot be set lower than any of the existing preset bolus amounts.

The Max bolus setting applies to boluses programmed by the user in Manual mode.

When the SmartGuard feature is active, SmartGuard determines the limits for each auto correction bolus.

### To set the Max bolus:

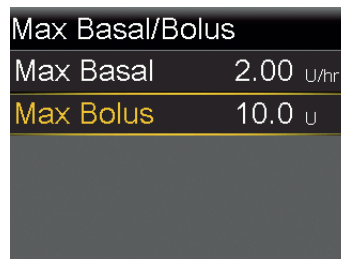
1. From the Home screen, press , and then select .
2. Select **Delivery Settings** > **Max Basal/Bolus**.

The Not for SmartGuard screen appears.




3. Select **Continue**.


The Max Basal/Bolus screen appears.



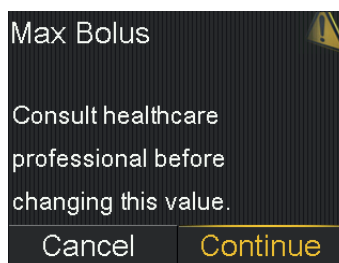4. Select **Max Bolus**.



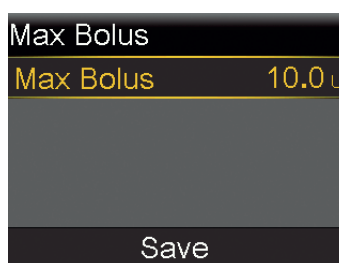5. To continue to the Max Bolus screen, select **Continue**.
6. Select **Max Bolus**, and then set the maximum number of insulin units the pump can deliver in one bolus.



7. Select **Save**.

## Bolus Wizard feature in Manual mode

In Manual mode, the Bolus Wizard feature uses Bolus Wizard settings to calculate an estimated bolus amount based on the BG readings and carbs that are entered.

After the Bolus Wizard feature is set up, use a normal bolus to deliver a food bolus, a correction bolus, or a food plus correction bolus. For more information, see *Delivering a normal bolus with the Bolus Wizard feature*, page 105.

The Bolus Wizard feature can also be used to deliver a Dual Wave bolus or a Square Wave bolus. For more information, see *Bolus types*, page 249.

## Bolus Wizard settings in Manual mode

To use the Bolus Wizard feature, consult a healthcare professional to determine the personal settings that should be used. The carb ratio, insulin sensitivity factor, BG target, and the active insulin time are needed to complete the setup. Always consult a healthcare professional before changes are made to the Bolus Wizard settings. The

setup procedure begins on *Setting up the Bolus Wizard feature in Manual mode*, page 101.


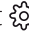
Manual mode Setting	Description
Active Insulin Time	<p>Active insulin is the bolus insulin that has been delivered by the pump and is still working to lower glucose levels. In the Bolus Wizard and SmartGuard Bolus feature, the Active Insulin Time setting is used to calculate a correction bolus by subtracting the estimated active insulin from each bolus. In SmartGuard, auto correction boluses are delivered up to every 5 minutes. A shorter Active Insulin Time setting may result in more insulin being delivered in correction boluses.</p> <p>A healthcare professional provides the personalized active insulin time based on historic glycemic control data for the individual user. When using SmartGuard, the recommended initial setting is an Active Insulin Time of 2-3 hours. The Active Insulin Time setting in the MiniMed 780G system is not necessarily reflective of the physiological insulin metabolism. Adjustments are not based on the pharmacokinetics and pharmacodynamics of the rapid-acting insulin. The current active insulin amount appears on the Home screen and includes only the bolus insulin received.</p>
BG Target	<p>In Manual mode, the Bolus Wizard feature calculates the estimated bolus based on the BG target range. The high and low values set are the values to which the BG is corrected. To use a single target value rather than a range, set the same value for the high and low value of the BG target.</p> <p>If the BG reading is above the high target value, a correction dose is calculated. If the BG reading is below the low target value, a negative correction is calculated and subtracted from the food bolus.</p>
Carb Ratio	<p>The carb ratio setting is used for food bolus calculations. The number of carb grams that are covered by 1 unit of insulin.</p>
Insulin Sensitivity Factor	<p>The insulin sensitivity factor setting is used to calculate correction bolus amounts.</p> <p>The insulin sensitivity factor is the amount that BG is reduced by 1 unit of insulin.</p>




## Setting up the Bolus Wizard feature in Manual mode

In Manual mode, to use the Bolus Wizard feature to calculate a bolus, first turn on the Bolus Wizard feature and enter the Bolus Wizard settings. There are four settings needed to set up the Bolus Wizard. Each setting is shown using 1/4, 2/4, 3/4, and 4/4 on the screens.

### To set up the Bolus Wizard feature:


1. From the Home screen, press , and then select .
2. Select **Delivery Settings > Bolus Wizard Setup**.

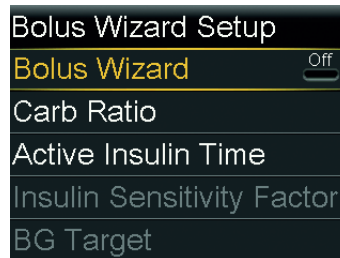This feature is only available in Manual mode. The Bolus Wizard Setup screen appears.




3. Select **Bolus Wizard** to turn on the feature.

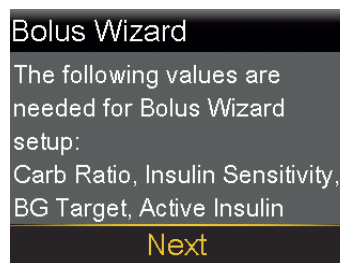If this is the first time the Bolus Wizard feature has been turned on, the following screen appears.




4. Confirm the values needed are ready to be entered, then select **Next**.

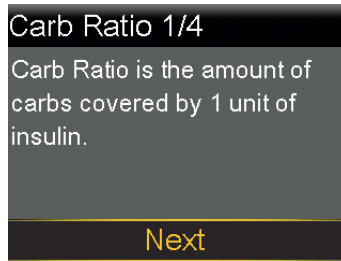The Carb Ratio 1/4 screen appears.



5. Select **Next**.

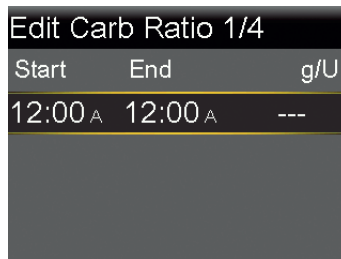The Edit Carb Ratio 1/4 screen appears.



6. To enter one carb ratio, enter the g/U, and then press .



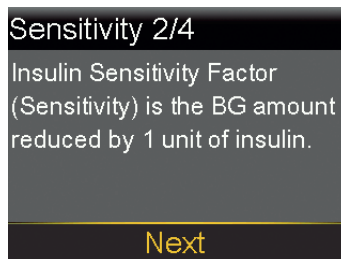
**Note:** For instructions on setting up more than one carb ratio over a 24-hour period, see *Settings covering a 24-hour period*, page 87.

7. Select **Next**.




**Note:** If the values are outside of the value range, a message asks to confirm the settings.

The Sensitivity 2/4 screen appears.



8. Select **Next**.

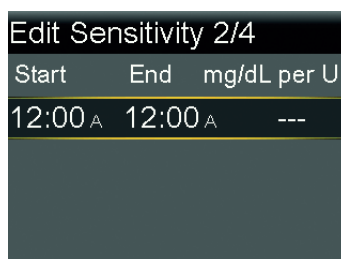The Edit Sensitivity 2/4 screen appears.




9. For one sensitivity factor, enter the mg/dL per U, and then press .



**Note:** For instructions on setting up more than one sensitivity factor over a 24-hour period, see *Settings covering a 24-hour period*, page 87.

10. Select **Next**.

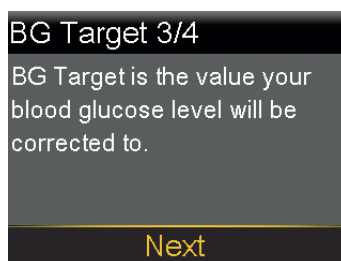The BG Target 3/4 screen appears.



11. Select **Next**.

The Edit BG Target 3/4 screen appears.

Edit BG Target 3/4		
Start	End	Lo-Hi (mg/dL)
12:00 A	12:00 A	--- - ---

12. For one BG target range, enter the Lo and Hi target, and then press .



**Note:** For instructions on setting up more than one BG target range over a 24-hour period, see *Settings covering a 24-hour period*, page 87.

13. Select **Next**.

The Active Insulin Time 4/4 screen appears.

Active Insulin Time 4/4	
Active Insulin Time is the length of time bolus insulin lowers glucose levels.	
<b>Next</b>	

14. Select **Next**.

The Active Insulin Time 4/4 screen appears.

Active Insulin Time 4/4	
<b>Duration</b>	4:00 hr
<b>Save</b>	

15. Enter the **Duration** of the active insulin time, and then press .



16. Select **Save**.

The Bolus Wizard feature setup is now complete.

Turning the Bolus Wizard feature off

The Bolus Wizard feature can be turned off at any time. The Bolus Wizard settings remain in the pump. When the Bolus Wizard feature is turned off, the Bolus Wizard menu selection does not appear on the Bolus screen, and the insulin sensitivity factor or BG target settings can not be edited on the Bolus Wizard Setup screen.

To turn the Bolus Wizard feature off:

- 1. From the Home screen, press , and then select .
- 2. Select **Delivery Settings > Bolus Wizard Setup**.
- 3. Select **Bolus Wizard** to turn the feature off.

Delivering a Normal bolus



In Manual mode, a normal bolus provides a single immediate dose of insulin. Use a normal bolus to cover food intake or to correct a high BG meter reading or both.




**Note:** The pump can deliver a normal bolus while a Square Wave bolus or the Square portion of a Dual Wave bolus is being delivered.

Delivering a normal bolus with the Bolus Wizard feature

In Manual mode, the Bolus Wizard screen shows the most recent BG reading, if available. The table indicates the different ways that the Bolus Wizard screen shows the BG reading.

Bolus Wizard screen	Glucose reading information
	The  icon indicates that a recent blood glucose (BG) meter reading is used by the Bolus Wizard feature to calculate a correction bolus. DO NOT enter a sensor glucose (SG) value in place of a blood glucose (BG) meter reading.

## Bolus Wizard screen

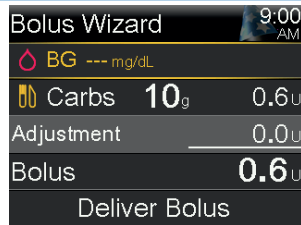

Bolus Wizard 9:00 AM

BG --- mg/dL

Carbs 10g 0.6U

Adjustment 0.0U



Bolus 0.6U

Deliver Bolus


## Glucose reading information

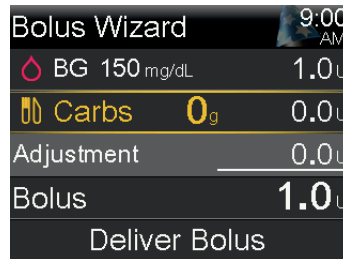The BG appears as dashes when no BG is available for the Bolus Wizard feature to calculate a correction bolus.

### To deliver a normal bolus using the Bolus Wizard feature:

1. From the Home screen, press , and then select .
2. Select **Bolus** > **Bolus Wizard**.

This feature is only available in Manual mode. The Bolus Wizard screen appears.



Bolus Wizard 9:00 AM

BG 150 mg/dL 1.0U

Carbs 0g 0.0U

Adjustment 0.0U

Bolus 1.0U

Deliver Bolus

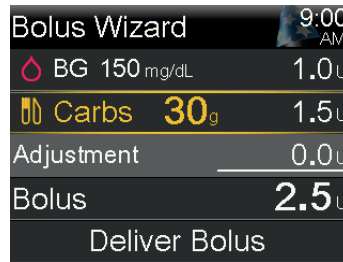
3. For a correction bolus or a food bolus with a correction, use a blood glucose (BG) meter for a blood glucose (BG) meter reading. Do not enter a sensor glucose (SG) value in place of a blood glucose (BG) meter reading. A blood glucose (BG) meter reading must always come from a blood glucose meter. The entered glucose value is used to calibrate the sensor.



**Note:** A blood glucose (BG) meter reading can be entered on the Bolus Wizard screen. On the Bolus Wizard screen, select **BG**.

4. For a food bolus, select **Carbs** to enter the carb count of the meal. For a correction bolus where no food was eaten, leave the carbs value at 0.

The calculated bolus appears in the Bolus field.



- If a change to the bolus amount is needed, select **Bolus** and modify the bolus amount.





- Select **Deliver Bolus** to start the bolus.

The pump beeps or vibrates and a message appears when the bolus starts. The Home screen shows the bolus amount as it is being delivered. The pump beeps or vibrates when bolus delivery is complete.


## Delivering a normal bolus using the Manual bolus feature

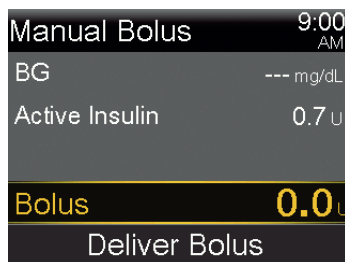The following procedure describes how to deliver a normal bolus using the Manual bolus feature.

### To deliver a normal bolus using the Manual bolus feature:

- From the Home screen, press , and then select .
- In Manual mode, do one of the following:
  - Select **Bolus** if the Bolus Wizard feature is turned off.
  - Select **Bolus** > **Manual Bolus** if the Bolus Wizard feature is turned on.

The Manual Bolus screen appears.





3. Select **Bolus** to set the bolus delivery amount in units.
4. Select **Deliver Bolus** to start the bolus.



## Stopping a bolus delivery

These procedures describe how to stop a bolus.





**WARNING:** Always press , select , and then select **Stop Bolus** to stop bolus insulin delivery. Do not use the Suspend All Delivery feature to stop bolus insulin. The Suspend All Delivery feature stops both basal insulin and bolus insulin delivery. Failure to resume basal insulin delivery could result in too little insulin, which may cause high BG.

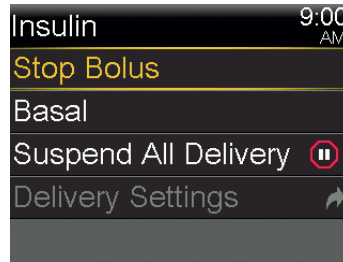


**Note:** To stop all insulin delivery, use the Suspend All Delivery feature (press , select , and then select **Suspend All Delivery**). For more information on using the Suspend All Delivery feature, see *Suspending all insulin delivery and resuming basal insulin delivery*, page 91.

### To stop a bolus delivery:


1. While the pump delivers a bolus, press  and then select .
- The Insulin menu appears.





2. Select **Stop Bolus**.

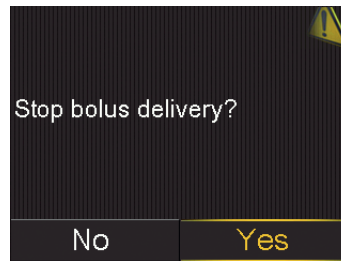A message appears confirming if bolus delivery should be stopped.



3. Select **Yes** to confirm.

The Bolus Stopped screen appears and shows the amount of bolus delivered, and the original bolus amount set up.



4. Select **Done**.



**Note:** The delivered amount can be viewed in the insulin delivery history screen after the procedure is closed. For more information, see *Daily History screen*, page 219.



5



# 5

## Reservoir and infusion set

The pump has options to change the reservoir and infusion set, reservoir only, or infusion set only. This chapter provides information about setting up the reservoir and infusion set with the Reservoir & Set option.

If the reservoir runs out of insulin and the infusion set has not been used for the duration of use indicated for the infusion set, the New Reservoir Only option may be used to change the reservoir. If only the infusion set needs to be changed, the New Set Only option may be used to change the infusion set.

Refer to the infusion set user guide for the duration of use indicated for the infusion set. Refer to the reservoir user guide for the duration of use indicated for the reservoir.

Do not begin the steps to replace the reservoir and infusion set until training has been received.



**WARNING:** Always confirm that the infusion set tubing is disconnected from the body before doing the following steps:

- placing the reservoir into the pump
- rewinding the pump
- loading the reservoir
- filling the infusion set tubing

Failing to disconnect the infusion set tubing from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

## Setting up the reservoir and infusion set

Confirm that the time and date on the pump are correct before insulin is used with the pump for the first time. For information about how to change the time and date on the pump, see *Time and date*, page 203. Consult a healthcare professional to determine the appropriate pump settings before insulin is used with the pump.

The following items are needed:

- MiniMed 780G insulin pump
- vial of rapid-acting U-100 insulin
- MiniMed or Medtronic reservoir
- MiniMed or Medtronic infusion set and its user guide



**WARNING:** Do not use the pump to deliver insulin for the first time until the active insulin has been cleared. If the pump has been used for training with bolus delivery before insulin is used, the active insulin value may be inaccurate. This may result in inaccurate insulin delivery, and serious injury. For details, see *Clearing the active insulin*, page 209.



**Note:** Different infusion sets may have different instructions for insertion into the body. All the procedures in the sections within this chapter must be followed to change the reservoir and infusion set.

## Removing the reservoir and rewinding the pump

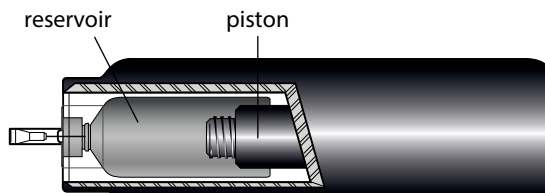
If this is the first time a reservoir is inserted into the pump, proceed to the pump rewind instructions. For more information about the reservoir see the reservoir user guide.




**WARNING:** Always confirm that the infusion set is disconnected from the body before rewinding the pump or filling the infusion set tubing. Never insert the reservoir into the pump while the tubing is connected to the body. Doing so may result in an unintentional infusion of insulin, and may cause hypoglycemia.

When the pump rewinds, the piston in the reservoir compartment returns to its starting position and allows a new reservoir to be placed into the pump.

The piston is located in the reservoir compartment of the pump. It engages the reservoir and pushes insulin through the tubing.




### Start here:


1. Wash hands with soap and water. On the pump, press  to go to the Menu screen.



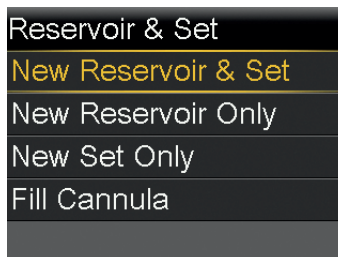2. Select , and then select **New Reservoir & Set**.




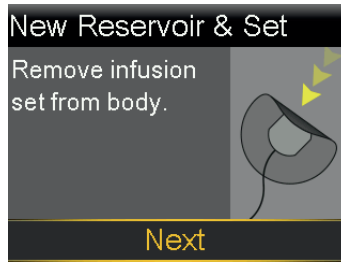3. Remove the infusion set by loosening the adhesive and pulling the set away from the body. Select **Next**.




**Note:** For instructions on how to remove the infusion set from the body refer to the user guide that came with the infusion set.

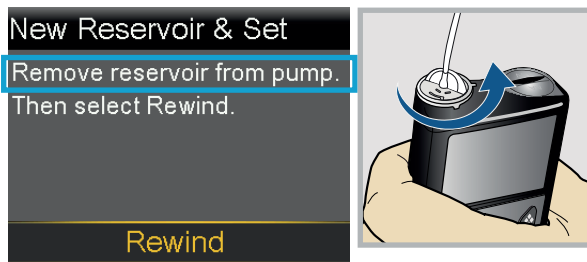4. If the optional activity guard is attached to the reservoir compartment on the pump, remove it now.
5. Remove the used reservoir from the pump.



6. Dispose of the used reservoir and infusion set per the disposal information in the corresponding user guide.
7. Select **Rewind**.

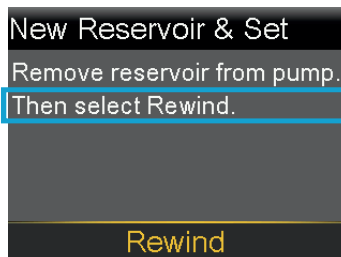**Do not** connect the infusion set to the body.








**WARNING:** Always confirm that the infusion set is disconnected from the body before rewinding the pump. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.



8. Follow the next steps to fill the new reservoir with insulin and to connect the infusion set tubing.

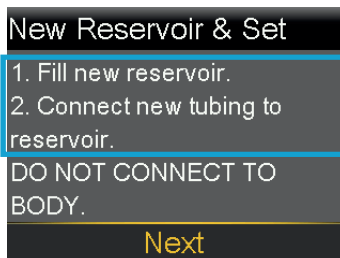Do not select **Next**.



### Filling the reservoir and connecting it to the infusion set tubing

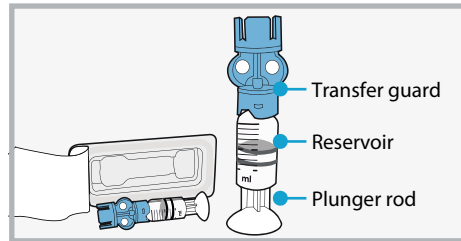


**WARNING:** Always allow the insulin to reach room temperature before use. Cold insulin may cause air bubbles in the reservoir and tubing, which may result in inaccurate insulin delivery.


The following procedures must be performed in the order presented.

**To fill the reservoir and connect it to the infusion set tubing:**


1. Remove the reservoir from the package. Make sure the insulin vial is at room temperature to reduce the risk of air bubbles.



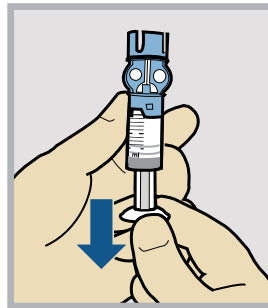
2. Pull the plunger down based on the planned insulin fill amount for the duration of use indicated for the reservoir.



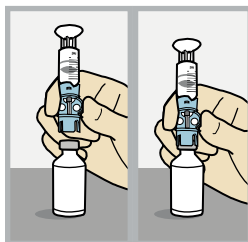
3. Wipe the top of the vial with alcohol. Place the vial on a sturdy flat surface. Firmly press the transfer guard onto the vial.



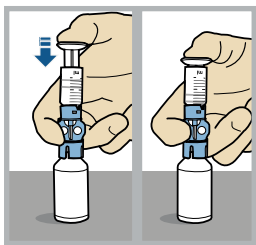
4. Push and hold the plunger down.




**WARNING:** Do not use the reservoir or infusion set if insulin or any liquid gets on the top of the reservoir or inside the tubing connector, as shown in the image. Insulin or any liquid may temporarily block the vents. This may result in the delivery of too little or too much insulin, which may cause hyperglycemia or hypoglycemia. If insulin or any liquid gets on the top of the reservoir or inside the tubing connector, start over with a new reservoir and infusion set.

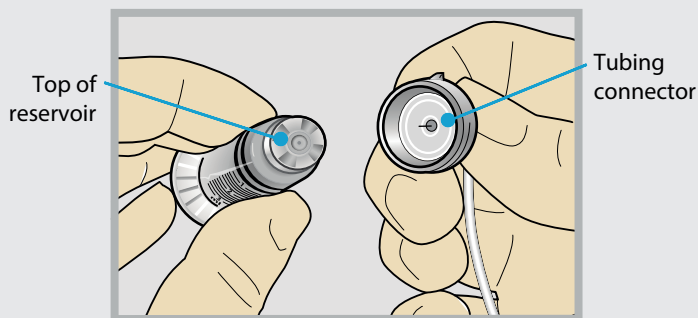

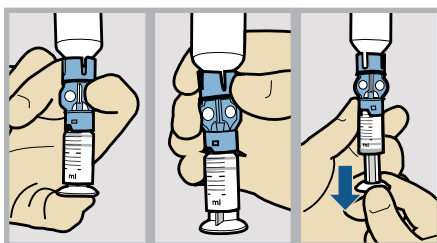
5. Keeping a thumb on the plunger, flip the vial over so the vial is on top. Release the thumb and pull the plunger down to fill the reservoir with insulin.



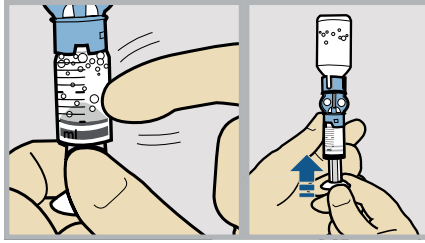
6. Tap the reservoir to move air bubbles to top of reservoir. Push the plunger up to move air into vial.



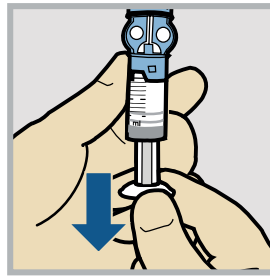
7. Pull the plunger back down to allow the reservoir to fill with the amount of insulin needed for the duration of use indicated for the reservoir.



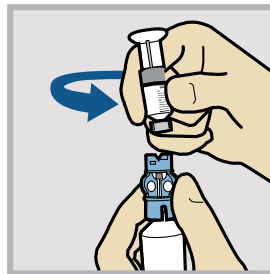
8. To avoid getting insulin on the top of the reservoir, **flip the vial over again so the reservoir is on top**. Hold the transfer guard and turn the reservoir counter-clockwise and remove the reservoir from the transfer guard.



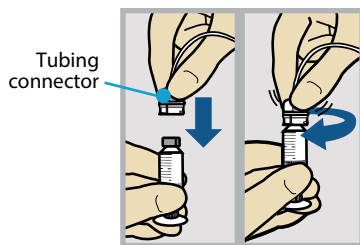
9. Follow the instructions in the infusion set user guide to access the infusion set tubing.
10. Gently push the tubing connector onto the reservoir. Turn the connector clockwise until it is locked into place.



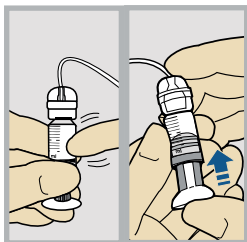11. Tap the reservoir to move any air bubbles to the top. Push the plunger slightly to move the bubbles into the tubing.



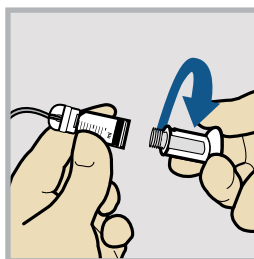12. Twist the plunger counter-clockwise to loosen it and to remove it.




### Placing the reservoir into the pump and filling the tubing with insulin



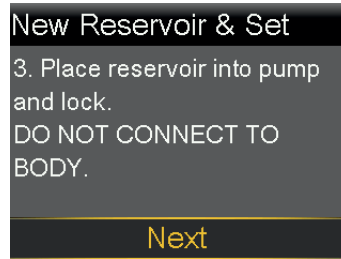
**WARNING:** Always rewind the pump before placing a new reservoir. Failing to rewind the pump may result in an unintentional infusion of insulin, which may cause hypoglycemia.

## To place the reservoir into the pump and fill the tubing with insulin:




**Note:** The backlight may have turned off. Press any button to turn the screen back on. Press  to go to the Menu screen, and then select .

1. Select **Next**.




2. Place the reservoir into the pump.

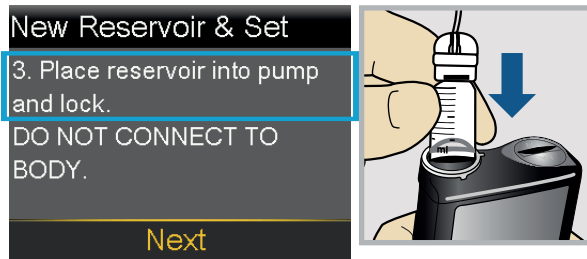**Do not** connect the infusion set to the body.



**WARNING:** Always confirm that the infusion set is disconnected from the body before placing the reservoir into the pump. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

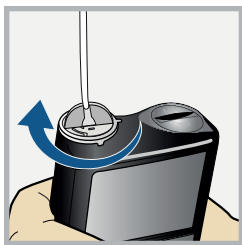
3. Turn the reservoir clockwise until the reservoir locks into place, and select **Next**.



### New Reservoir & Set

3. Place reservoir into pump and lock.  
DO NOT CONNECT TO BODY.

Next


4. Select **Load** and hold  until the checkmark appears on the screen.  
**Do not** connect the infusion set to the body.



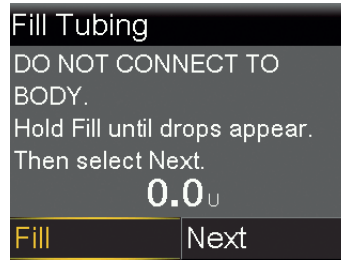
5. When the checkmark appears, select **Next**.



**WARNING:** Always confirm that the infusion set is disconnected from the body before loading the reservoir and filling the tubing. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

6. Select **Fill** and keep holding  until there are no air bubbles visible in the tubing, and there are drops at the end of the tubing.

**Do not** connect the infusion set to the body.




**WARNING:** Always check the tubing for air bubbles. Continue to press Fill until no bubbles remain in the tubing. Air bubbles may result in inaccurate insulin delivery.

7. After drops appear, press > and select **Next**.




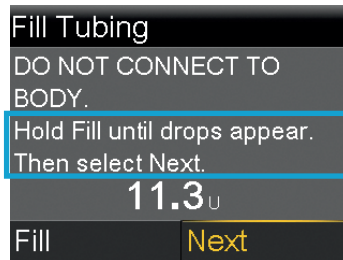**Note:** The location of the infusion set needle may be different depending on the type of infusion set being used.



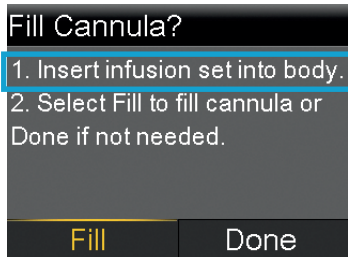8. Follow the steps in the infusion set user guide to insert the infusion set into the body before proceeding with the steps on the pump screen.





**Note:** If an infusion set with a steel cannula is used, the cannula does not need to be filled, and **Done** may be selected.


### Inserting the infusion set into the body

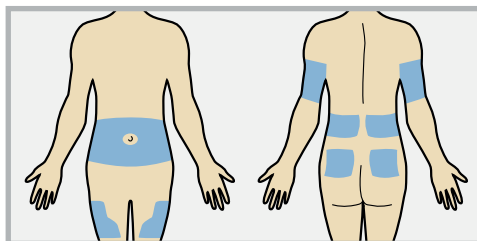Always refer to the infusion set user guide and the serter user guide, if needed, for instructions about how to insert an infusion set into the body.



**WARNING:** Do not remove the reservoir from the pump while the infusion set is connected to the body. Doing so may result in the delivery of too little or too much insulin, which may cause hyperglycemia or hypoglycemia.

Choose an insertion site from the shaded areas. Clean the insertion site with alcohol or other antiseptic as directed by a healthcare professional.



**CAUTION:** Do not use the same infusion set insertion site for an extended period of time. This may cause the site to become overused. Rotate the infusion set insertion sites regularly.



**CAUTION:** Always change the infusion set as indicated by the infusion set user guide. Using the same infusion set for an extended period of time beyond its product labeling can cause infusion set occlusion or site infection.

After the infusion set is inserted into the body follow the steps in the following section to fill the cannula.

### Filling the cannula

Filling the soft cannula with insulin is required after the infusion set is inserted into the body and the introducer needle is pulled out. The insulin amount required to fill the cannula depends on the type of infusion set used. Refer to the user guide that came with the infusion set for more information.



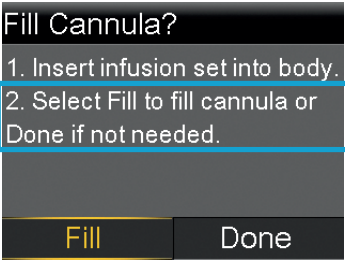
**Note:** The Fill Cannula action is not required during a reservoir only change. If performing a reservoir only change, select **Done** on the **Fill Cannula?** screen.



**WARNING:** Never leave the pump on the Fill Cannula? screen. Insulin delivery is suspended while on the Fill Cannula? screen. Always finish filling the cannula or return to the Home screen, to avoid continued insulin delivery suspension. Prolonged suspension of insulin delivery may cause hyperglycemia.

### To fill the cannula:


1. After the infusion set is inserted into the body, select **Fill**.




**Note:** Always verify that the amount shown in the **Fill amount** field is correct. The pump will remember the fill amount last used. Change the **Fill amount** if needed.

- If the Fill amount is correct, press **✓** to select **Fill Now** and then press **⊙**.
- If the Fill amount is incorrect, press **⊙**. Change to the correct amount and press **⊙**. Then select **Fill Now**.

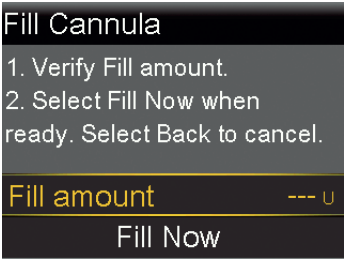2. Select **Fill amount** and enter the amount per the infusion set user guide.  
After entering the cannula size, press **⊙**.



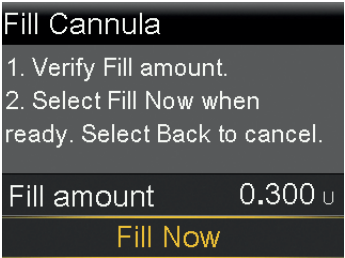3. Select **Fill Now**.



The screen displays the insulin amount as insulin fills the cannula.

The reservoir and infusion set change is now complete.


Always check blood glucose (BG) using a blood glucose meter one to three hours after changing the infusion set or reservoir.



**Note:** Use the following procedure only when it is necessary to stop filling the cannula.

#### To stop filling the cannula:

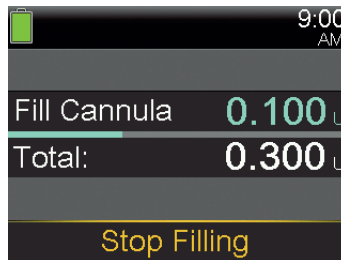1. Select **Stop Filling** to stop filling the cannula.



2. Select **Yes**.

The Fill Stopped screen appears.



3. Select **Done**.

### Disconnecting the infusion set

Refer to the infusion set user guide for instructions on how to disconnect the infusion set.

### Reconnecting the infusion set

Refer to the infusion set user guide for instructions on how to reconnect the infusion set.



# 6







# 6

## Paired devices

This chapter explains how to pair the MiniMed 780G insulin pump with compatible devices.

### Setting up the Accu-Chek Guide Link meter

The MiniMed 780G insulin pump with smart device connectivity can pair only with an Accu-Chek Guide Link meter to automatically receive blood glucose (BG) meter readings. If the Accu-Chek Guide Link meter is not paired with the pump, enter BG readings manually. The pump beeps, vibrates, or simultaneously beeps and vibrates when the pump receives a BG reading. Confirm the BG reading and deliver a bolus, if necessary. If a BG reading is not confirmed within 12 minutes, the BG will not be stored. If the BG reading is outside the range of 70 mg/dL to 250 mg/dL, an alert appears. Follow instructions from a healthcare professional to treat low BG or high BG.

To pair the pump and meter, use the following items:

- MiniMed 780G insulin pump with smart device connectivity
- Accu-Chek Guide Link meter

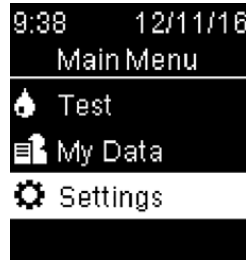
### Pairing the pump and meter

The MiniMed 780G insulin pump with smart device connectivity can pair with up to four Accu-Chek Guide Link meters.

#### To prepare the meter to pair with the pump:

1. Press the **OK** button on the meter to turn on the meter.


2. Select **Settings**.



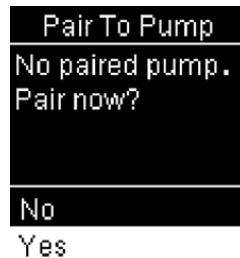3. Select **Wireless**.






4. Select **Yes** if the confirmation screen appears on the meter screen. Or, if the confirmation screen does not appear, select **Pairing**.



The meter serial number appears on the meter screen. The meter is now ready to pair with the pump.

#### To prepare the pump to pair with the meter:

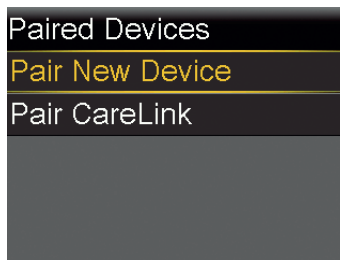1. From the Home screen, press , and then select .
2. Select **Pair New Device**.



The Searching... screen appears. After the pump is done searching, the Select Device screen appears.

3. Select the meter that matches the serial number that displays on the meter screen.

If the correct serial number does not appear, select **Search Again**.



If the connection is successful, a "Pairing successful!" message appears on the pump. A "Paired with pump" message with the serial number of the pump appears on the meter screen. If a Device not found alert appears, see *Pump alarms, alerts, and messages*, page 291 for more information.

## Unpair and delete the meter

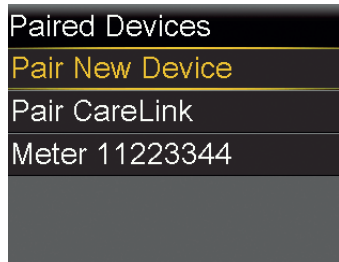
### Unpairing a meter from the pump

Follow this procedure to unpair the Accu-Chek Guide Link meter from the pump.

#### To unpair the meter from the pump:


1. From the Home screen, press , and then select .

The Paired Devices screen appears.




2. Select the serial number of the meter to unpair the device. The Accu-Chek Guide Link meter serial number is located on the back of the meter.

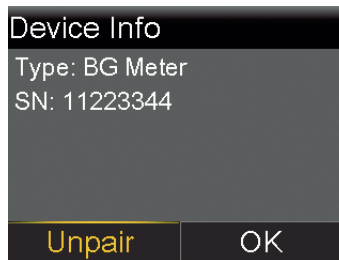The Device Info screen appears.



3. Select **Unpair**.

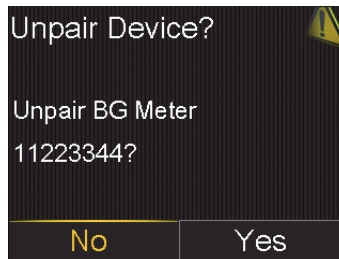The Unpair Device? screen appears.



4. Select **Yes** to confirm. Select **No** to cancel.

## Deleting the pump from a meter

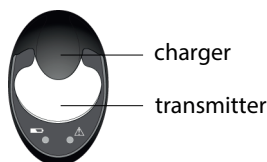
For steps to delete the pump from a meter, see the Accu-Chek Guide Link User's Manual.

## Pairing the pump and transmitter



The pump and transmitter must be paired to use the sensor. When paired, the pump and transmitter communicate with each other through a wireless connection. Only one transmitter can be paired with the pump. If a transmitter is already paired with the pump, delete the transmitter, and then continue. For instructions on how to delete a transmitter from the pump, see *Unpairing the transmitter from the pump*, page 139.


### To pair the pump and transmitter:

1. Attach the transmitter to the charger. Fully charge the transmitter. Keep the transmitter attached to the charger.




**Note:** Both lights on the charger are off when the transmitter is fully charged. For more information, see the transmitter user guide.

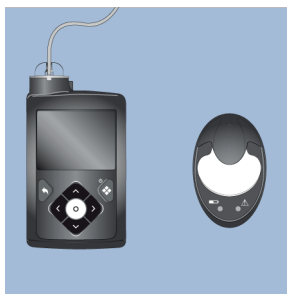2. From the Home screen, press , and then select .
3. Place the transmitter (still attached to the charger) next to the pump.



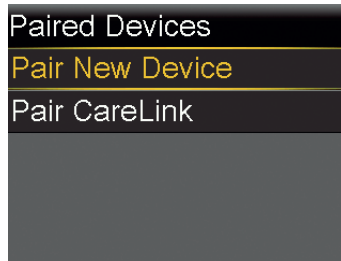4. Select **Pair New Device**.




The Searching... screen appears.



**Note:** The search process can take up to 20 seconds.

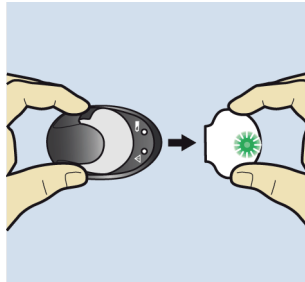
5. Remove the transmitter from the charger. The transmitter light flashes 10 times and turns off.




The Select Device screen appears with a list of available devices.

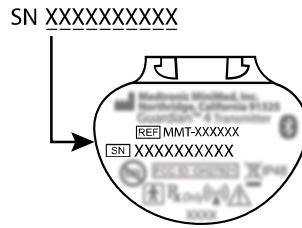6. Select the CGM device that matches the serial number indicated on the back of the transmitter.





If the correct serial number does not appear, select **Search Again**.



If the connection is successful, a “Pairing successful!” message appears on the pump. When the transmitter is communicating with the pump, the Sensor feature is turned on and  appears on the Home screen. For information on using the sensor with the transmitter, see *Connecting the transmitter to the sensor*, page 172. If a Device not found alert appears, see *Pump alarms, alerts, and messages*, page 291 for more information.


## Transmitter and sensor care
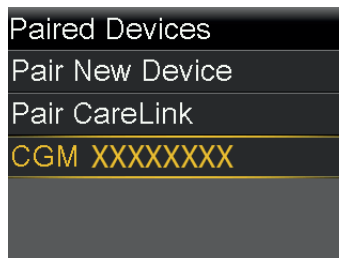
### Unpairing the transmitter from the pump


Follow this procedure to unpair the transmitter from the pump, including when the transmitter needs to be replaced.

#### To unpair the transmitter from the pump:

1. From the Home screen, press , and then select .  
The Paired Devices screen appears.




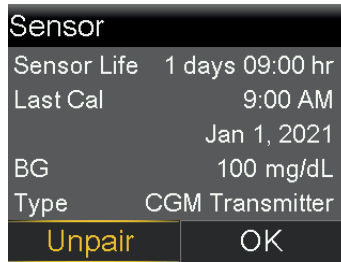2. Select CGM with the correct serial number.  
The Sensor screen appears.



3. Select **Unpair**.

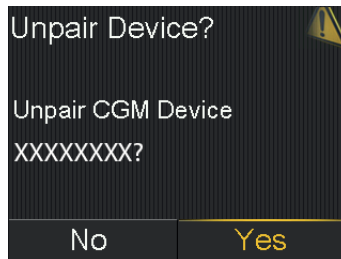The Unpair Device? screen appears.



4. Select **Yes** to confirm. Select **No** to cancel.

When the transmitter is unpaired from the pump, a No Paired CGM banner appears on the Home screen.

## Disconnecting the transmitter from the sensor

Refer to the transmitter user guide for instructions on how to disconnect the transmitter from the sensor.

## Removing the sensor

Refer to the sensor user guide for instructions on how to remove the sensor.

## Cleaning the transmitter

Refer to the transmitter user guide for instructions on how to clean the transmitter.

## Storing the transmitter

Refer to the transmitter user guide for instructions on how to store the transmitter.



## MiniMed Mobile app

The MiniMed Mobile app is an optional accessory that is compatible with the MiniMed 780G system. The app provides a secondary display that allows the user to view pump data. A compatible smartphone is required for the app to function. The app is available for both iOS and Android platforms. Consult the MiniMed Mobile app user guide for installation instructions.

## Updating the pump software



**Note:** This feature may not be available in all geographies.

After an eligibility message for a pump software update is received, use the MiniMed Mobile app to perform the pump software update. Training materials will be provided to guide you on which app to use. The app provides instructions for each step of the process. Follow the instructions provided on the app screens to perform the update.



**CAUTION:** A stable internet connection is required throughout the entire update process. Avoid the use of unsecure Wi-Fi networks or public Wi-Fi hotspots.

## Downloading the pump software update

To check if an update is available, ensure you are logged into the MiniMed Mobile app. The Software is Ready screen appears on the app when the download is complete.

## Preparing to install the pump software update

**To prepare to install the pump software update:**



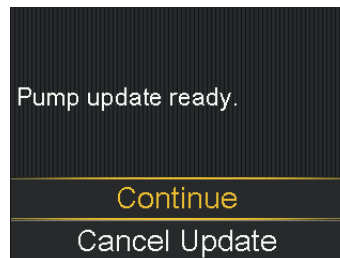
**Note:** After installation is complete, the SmartGuard feature requires a 5-hour warm-up period before it is active.

- Ensure glucose is within target before starting the update.
- Clear active alerts or alarms.


- If the pump is Suspended on low or Suspended before low, wait until insulin delivery resumes and BG recovers before starting the update.
- If a bolus delivery is in progress, wait until the bolus delivery completes before installing the pump software update.
- If the battery is low, the pump software update will not install. If the battery icon is not green, replace the battery before installing the pump software update.
- Insulin is not delivered and sensor glucose (SG) values are not shown for up to 20 minutes during the pump software installation. Manual injections are not accounted for in the active insulin amount. If an injection is needed during the software update, consult a healthcare professional for how long to wait after a manual injection before using the Bolus Wizard feature. Refer to *Emergency kit*, page 25 for necessary supplies to use for backup insulin delivery if needed.

## Installing the pump software update

1. When instructed by the app, go to the Home screen on the pump. On the pump, a screen appears when the pump is ready for the software update.
2. Select **Continue**.



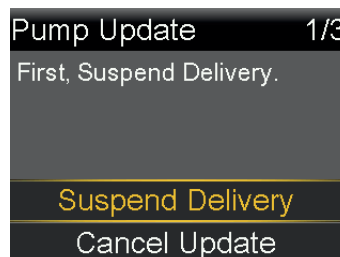
3. Select **Suspend Delivery** to suspend bolus and basal insulin delivery.




4. Disconnect the infusion set from the body, and then select **Confirm**.



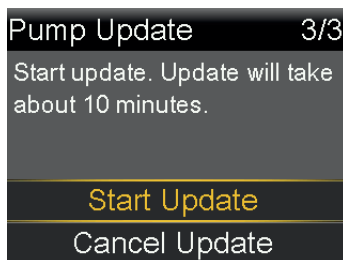
5. Select **Start Update**.



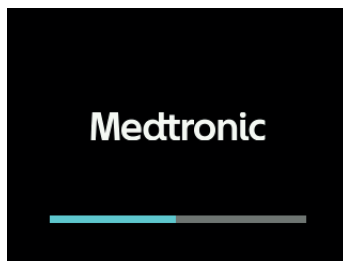
While the pump updates, a screen shows the progress.



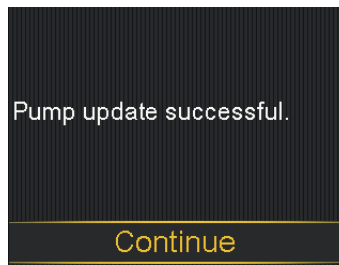6. Select **Continue**.



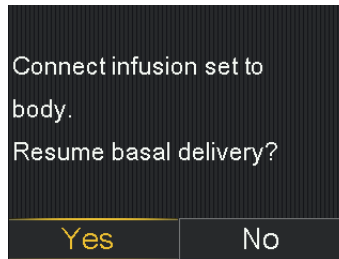7. Reconnect the infusion set to the body.
8. Select **Yes** to resume basal insulin delivery.



**Note:** The previous version of the software is retained if the update is not successful.



## Completing the pump software update

Follow the instructions on the app to complete the pump software update.

## Uploading device data to CareLink software

Upload system data to CareLink software with the MiniMed Mobile app or the Blue Adapter. Follow the instructions found on the CareLink software to upload system data with the Blue Adapter. Refer to the MiniMed Mobile app user guide for instructions to upload MiniMed 780G system data to CareLink software with the app.

### To prepare the pump to upload to CareLink software:

1. From the Home screen, press , and then select .
2. Select **Pair CareLink**.

Follow instructions on the CareLink uploader to complete steps.

## Sharing device data with the CareLink Connect app

The CareLink Connect app works with CareLink software. Through the CareLink Connect app, care partners can see information sent from a connected MiniMed Mobile app. A compatible smartphone is required for the app to function. The app is available for both iOS and Android platforms.

For more information about sharing data with the CareLink Connect app, see the MiniMed Mobile app user guide and the CareLink Connect app user guide.



7







# 7 Continuous glucose monitoring

This chapter explains how to enter sensor settings and set up continuous glucose monitoring (CGM). CGM requires these items:

- MiniMed 780G insulin pump
- Sensor alert settings provided by a healthcare professional
- Guardian 4 sensor
- Guardian 4 transmitter

## CGM overview

### What is CGM

CGM is a tool that uses a sensor to continuously measure the amount of glucose in interstitial fluid. CGM consists of the following:

- Sensor Glucose (SG) readings that are displayed every 5 minutes.
- Alerts based on current and predicted high and low glucose levels.
- Graphs that show glucose trends over time.
- Trend arrows that show the rate at which the most recent SG readings are (has been) rising or falling.

The sensor does not require calibration for use with the system. However, every blood glucose (BG) meter reading either entered manually or received from a paired meter is used to calibrate the sensor.



**WARNING:** Do not use SG values to make treatment decisions, including delivering a bolus, while the pump is in Manual mode. When the SmartGuard feature is active and you are no longer in Manual mode, the pump uses an SG value, when available, to calculate a bolus amount. However, if your symptoms do not match the SG value, use a BG meter to confirm the SG value. Failure to confirm glucose levels when your symptoms do not match the SG value can result in the infusion of too much or too little insulin, which may cause hypoglycemia or hyperglycemia. For more information on using the SmartGuard feature, see *SmartGuard*, page 179.



**WARNING:** Do not use continuous glucose monitoring if hydroxyurea, also known as hydroxycarbamide, is taken. Hydroxyurea is used to treat certain diseases, such as cancer and sickle cell anemia. Hydroxyurea use results in higher sensor glucose readings compared to blood glucose readings. Taking hydroxyurea while using continuous glucose monitoring can result in hypoglycemia caused by over-delivery of insulin, inaccurate or missed alarms and alerts, delay or loss of sensor-enabled insulin suspension, and substantially higher sensor glucose readings in reports than actual blood glucose readings.

Always check the label of any medication being taken to confirm if hydroxyurea or hydroxycarbamide is an active ingredient. If hydroxyurea is taken, consult a healthcare professional. Turn the Sensor feature off to disable continuous glucose monitoring. For more information, see *Deactivating the Sensor feature*, page 173. Use additional blood glucose meter readings to verify glucose levels.



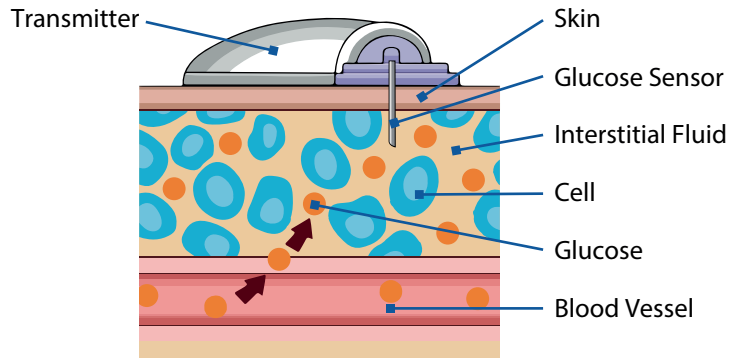
**WARNING:** Consult a healthcare professional if a medication that contains acetaminophen or paracetamol is taken while wearing the sensor. Medications that contain acetaminophen or paracetamol can falsely raise sensor glucose readings. The level of inaccuracy depends on the amount of acetaminophen or paracetamol active in the body and can differ for each person. Falsely elevated sensor readings can result in over-delivery of insulin, which can cause hypoglycemia. Medications that contain acetaminophen or paracetamol include, but are not limited to, cold medicines and fever reducers. Check the label of any medications being taken to see if acetaminophen or paracetamol is an active ingredient. Use additional blood glucose meter readings to confirm blood glucose levels.

While the SmartGuard feature is active, if acetaminophen or paracetamol is taken, program a temp target for up to eight hours, or the amount of time recommended by a healthcare provider. For more information, see *Setting a temp target*, page 196. Use blood glucose values instead of sensor glucose readings to calculate a meal bolus or correction bolus for up to eight hours, or the duration recommended by a healthcare provider, after taking acetaminophen or paracetamol.

## What is blood glucose (BG) and sensor glucose (SG)?

Blood glucose and sensor glucose are measured in different places. It is important to understand the differences between the two, as there are times when the system requires you to enter a blood glucose and there are other times when the system will use a sensor glucose.

Glucose travels between the blood and interstitial fluid. The blood glucose (BG) meter measures glucose levels in your blood. The glucose sensor measures glucose in the interstitial fluid. Blood glucose (BG) meter readings and sensor glucose (SG) readings will be close but will rarely exactly match. This difference is normal and should be expected.



**IMPORTANT:** When a glucose value is entered into the pump, it must be from a blood glucose (BG) meter.

The system automatically uses the entered glucose value to calibrate the sensor, unless the system gives you the option to calibrate the sensor.

The following table shows when to use a blood glucose (BG) meter reading:

When to use a BG	Examples	
Anytime glucose is entered into the pump, it needs to be a blood glucose (BG) meter reading, not a sensor glucose (SG) value.		
Anytime you deliver a bolus in Manual Mode and you want to use a glucose for a correction.		
Anytime the system requests a blood glucose (BG) meter reading.		



**Note:** See when to use a blood glucose (BG) meter reading when SmartGuard is active in *Entering a BG value in the SmartGuard feature*, page 189.

## Calibrating the sensor

Calibration is the process of using a blood glucose (BG) meter reading to help the sensor glucose (SG) readings more closely match the glucose measured in your blood. This blood glucose (BG) meter reading is not used until you confirm it on the screen. For more information, see *Entering a blood glucose (BG) meter reading*, page 96.

When you are using the MiniMed 780G system with the Guardian 4 CGM, you do not need to calibrate. However, the system is designed to use every blood glucose (BG) meter reading either entered manually or received from glucose meter to calibrate the sensor.

## CGM settings

### Sensor alert settings

An SG alert occurs when an SG reading changes at a particular rate, reaches a specified high or low limit, or before a high or low limit is reached. The pump can also be set to suspend insulin delivery before or when a low limit is reached.

### High SG settings




High SG settings provide alerts under the following conditions:

- When SG rises rapidly (Rise Alert).
- When SG approaches the high limit (Alert before high).
- When SG reaches the high limit (Alert on high).

The following graph shows the types of high SG settings.



 High SG alert settings



High glucose setting	Description
High limit	The high limit is used as a basis for some high SG settings. The high limit can be set from 100 to 400 mg/dL, for up to eight different time segments.
Alert before high	This setting provides an alert when SG is predicted to reach the high limit, raising awareness of potential high SG.
Time before high	This setting determines how long an Alert before high occurs before the high limit may be reached. It can be set between 5 and 30 minutes.
Alert on high	This setting provides an alert when SG reaches or exceeds the high limit.
High SG alert	This setting provides an alert when SG is at 250 mg/dL or higher for 3 hours. This is a fixed setting and cannot be changed.
Rise Alert	<p>This setting provides an alert when glucose is rising rapidly, such as after a meal or if a bolus is missed. Set the rise rates to match the trend arrows, as shown below, or to a custom rise rate.</p> <ul style="list-style-type: none"> <li>•  - SG is rising at a rate of 1 mg/dL per minute or more.</li> <li>•  - SG is rising at a rate of 2 mg/dL per minute or more.</li> <li>•  - SG is rising at a rate of 3 mg/dL per minute or more.</li> <li>• <b>Custom</b> - SG is rising at a custom rate, set from 1.0 mg/dL to 5.0 mg/dL per minute.</li> </ul>
Rise Limit	This setting determines when a Rise Alert occurs.

## Setting up the high SG settings

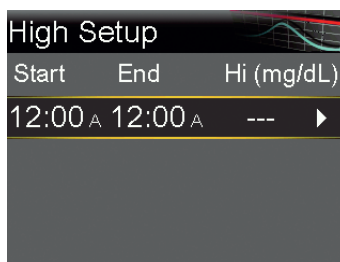
The Sensor feature must be turned on to set up the sensor settings. For more information, see *Turning on the Sensor feature*, page 173.

For details about high SG settings, see *High SG settings*, page 153.

### To set up the high SG settings:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **High Alert**.

The High Setup screen appears.




3. Select the time segment. The end time flashes.  
The start time of the first time segment is always 12:00 A. Up to eight time segments can be set, each with a different high limit. All the time segments must add up to a 24-hour period.
4. Set the End time.
5. Set your High limit. You can enter a value from 100 to 400 mg/dL, in increments of 5 mg/dL.
6. Select the arrow beside the End time to select the high alerts for the time segment.

A screen appears and shows the high alerts for the selected time segment.



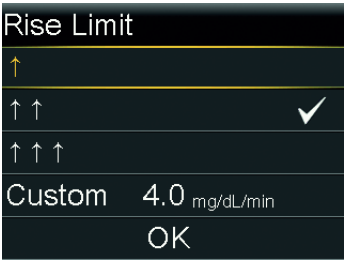
7. Set the following alerts, as desired:
- a. Select **Alert before high** to receive an alert before the high limit is reached.
  - b. Set the **Time before high** option between 5 to 30 minutes to receive an alert before the high limit is reached.
  - c. Select **Alert on high** to receive an alert when the high limit is reached.
  - d. Select **Rise Alert** to receive an alert when SG is rising quickly.
8. If Rise Alert is on, perform the following steps to set up the Rise Limit. Otherwise, proceed to step 9.
- a. Scroll down and select **Rise Limit**.
- The Rise Limit screen appears.



- b. Select one, two, or three arrows for the rise rate, or enter a custom rate.

Arrow selection	Minimum rate that SG is rising when an alert occurs.
↑	SG is rising at a rate of 1 mg/dL per minute or more.
↑↑	SG is rising at a rate of 2 mg/dL per minute or more.



Arrow selection	Minimum rate that SG is rising when an alert occurs.
	SG is rising at a rate of 3 mg/dL per minute or more.



**Note:** These arrows appear on the Home screen to indicate the rate at which SG is rising.



- c. To enter a custom rate, select **Custom**, enter the Rise Limit on the Custom Limit screen, and then select **OK**.
- d. Select **OK** again to confirm the Rise Limit settings.
9. Select **Next**.
10. If necessary, enter the remaining time segments to complete the 24-hour period.



**Note:** For instructions on setting up more than one high limit over a 24-hour period, see *Settings covering a 24-hour period*, page 87.

11. Select **Review**.
12. Review the high SG settings and select **Save**.



### To change the high SG settings:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **High Alert**.  
The High Setup screen appears.
3. Select **Edit**.
4. Select and adjust the time segment.
5. Select any alert setting to make adjustments, or to turn the setting on or off.
6. Select **Next**.
7. Select **Review**.
8. Review the high SG settings and select **Save**.

## High Snooze

The High Snooze feature sets the amount of time before a high alert repeats. The pump shows the high alert again if the high alert condition still exists after the specified snooze time.

### To set the High Snooze:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Snooze High & Low**.  
The Snooze screen appears.
3. Select **High Snooze** and enter a time in 5-minute increments from 5 minutes to 3 hours.
4. Select **Save**.


## Low SG settings

Low SG settings alert or suspend insulin delivery when SG either approaches or reaches the low limit.



**Note:** The MiniMed Mobile app may be used to view the sensor graph on a mobile device. Always read and acknowledge all alarms and alerts on the pump. If the pump simultaneously generates more than one alarm or alert, only one of the alarms or alerts appears on the mobile device.

The following graph shows the available low SG settings.

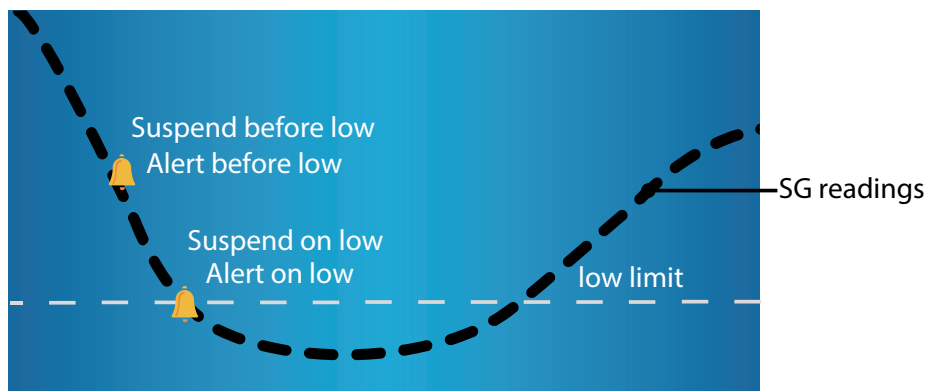

Low SG alert and suspend settings



**WARNING:** The Suspend before low and Suspend on low features are not intended to treat low BG. Suspending insulin delivery when SG is low may not bring BG back to the target range for several hours, which may cause hypoglycemia. If symptoms do not match the SG value, use a BG meter to confirm glucose levels before making treatment decisions.

For information about how to program low SG settings in Manual mode, see *Setting up the low SG settings*, page 165. The sensor must be turned on before low SG settings can be programmed.

## Low limit

The low limit is used as a basis for optional low SG alerts and Suspend by sensor features. The low limit can be set from 50 mg/dL to 90 mg/dL, for up to eight different time segments.

The Low SG alarm appears when SG readings fall below 64 mg/dL. This threshold may be higher than other sensors. Glucose detection and alert rates are less reliable below this threshold. This is a fixed setting and cannot be changed. When the alarm appears, it shows the SG reading next to the Low SG alarm. This alarm does not stop insulin delivery.

## The Suspend before low feature

The Suspend before low feature stops insulin delivery when SG is approaching the low limit. This feature can help minimize the amount of time spent with low glucose.



**WARNING:** Do not use the Suspend before low feature without first reading the information in this user guide and receiving training from a healthcare professional. The Suspend before low feature temporarily suspends insulin delivery for a maximum of two hours. Under some conditions of use, the pump can suspend insulin delivery again, resulting in under-delivery. Prolonged under-delivery of insulin may increase the risk of hyperglycemia and diabetic ketoacidosis. Always be aware of symptoms. If symptoms don't match SG readings, confirm SG with a BG meter reading.

The Suspend before low feature is turned off by default. Consult a healthcare professional before the Suspend before low feature is used.

If the Suspend before low feature is turned on, Alert on low is automatically turned on. Enabling Alert before low is optional.


- If Alert before low is off, a Suspend before low alert occurs, but the pump does not beep or vibrate when insulin delivery is suspended.
- The Suspend before low and Suspend on low features cannot be on at the same time. When either feature is on, the Resume basal alert can be activated.

## Suspend before low conditions

When a Suspend before low event occurs, insulin delivery is suspended. A Suspend before low event occurs if both of the following conditions are met:

- SG reading is at the low limit or is within 70 mg/dL above the low limit.
- SG is predicted to reach or fall below a level that is 20 mg/dL above the low limit within approximately 30 minutes.

The following image is an example of what can happen during a Suspend before low event.

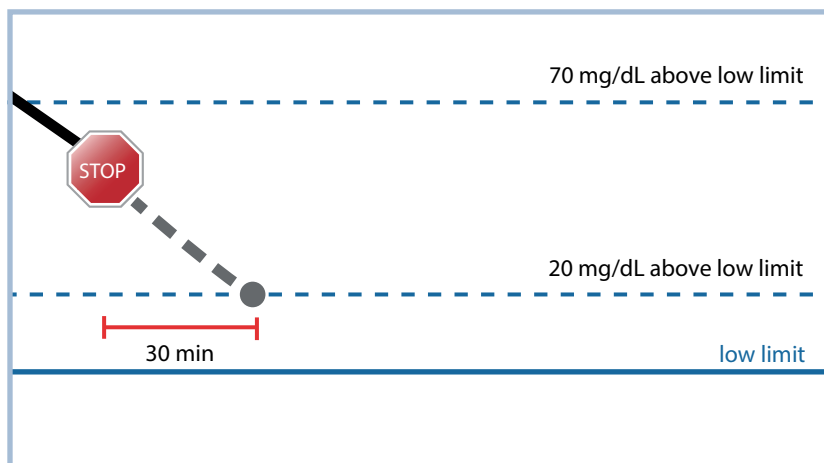

### Responding to a Suspend before low event

When the Suspend before low feature suspends insulin delivery, the icon flashes. If SG reaches the low limit, an Alert on low occurs.

When a Suspend before low event occurs, insulin delivery can be suspended for a minimum of 30 minutes or up to a maximum of two hours. Basal insulin delivery can be manually resumed at any time. For details, see *Manually resuming basal insulin delivery during a Suspend before low or Suspend on low event*, page 168. After 30 minutes, basal insulin delivery resumes if both of the following conditions are met:

- SG is at least 20 mg/dL above the low limit.
- SG is predicted to be more than 40 mg/dL above the low limit within 30 minutes.

If the Suspend before low alert is not cleared within two hours, the pump resumes insulin delivery and displays a Basal delivery resumed alert.

### Alert before low

Alert before low provides an alert when SG is predicted to reach the low limit, and increases awareness of potential low SG.

The Alert before low feature works as follows:

- If Alert before low is on, and both suspend features are off, Alert before low occurs 30 minutes before the low limit is reached.
- If the Suspend on low feature is on and Alert before low is on, Alert before low occurs 30 minutes before the low limit is reached.
- If the Suspend before low feature is on and Alert before low is on, a Suspend before low alert occurs when insulin delivery is suspended. For details, see *The Suspend before low feature*, page 160.

## The Suspend on low feature

The Suspend on low feature stops insulin delivery when SG readings reach or fall below the low limit. When a Suspend on low event occurs, insulin delivery is suspended. This feature is for situations when a person cannot respond to a low glucose condition and can help minimize the amount of time spent with low glucose.




**WARNING:** Do not use the Suspend on low feature without first reading the information in this user guide and receiving training from a healthcare professional. The Suspend on low feature temporarily suspends insulin delivery for a maximum of two hours. Under some conditions of use, the pump can suspend insulin delivery again, resulting in under-delivery. Prolonged suspension of insulin delivery may increase the risk of serious hyperglycemia, ketosis, and ketoacidosis.

The Suspend on low feature is off by default. Consult a healthcare professional for guidance before the Suspend on low feature is used.

When the Suspend on low feature is on, Alert on low is activated automatically. For more information, see *Alert on low*, page 164.

The following image is an example of what can happen during a Suspend on low event.

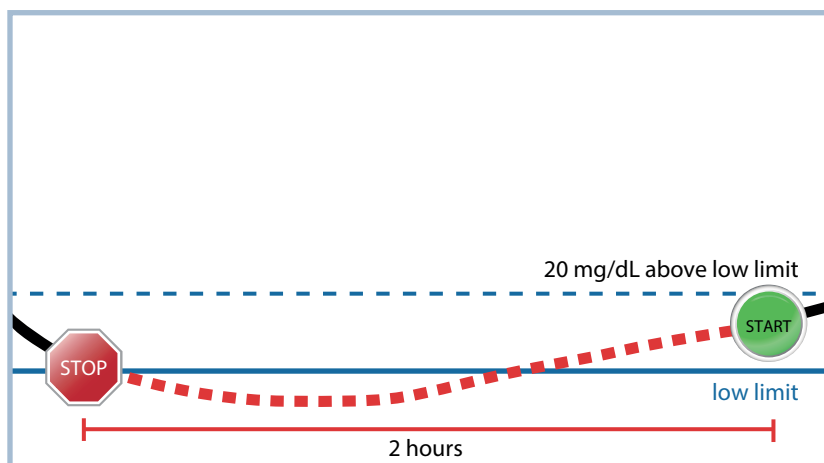

### Responding to a Suspend on low event

When the Suspend on low feature suspends insulin delivery, the icon flashes.


When a Suspend on low event occurs, a pump alarm occurs and insulin delivery remains suspended for a minimum of 30 minutes, up to a maximum of two hours. Insulin delivery can be resumed manually at any time. For details, see *Manually resuming basal insulin delivery during a Suspend before low or Suspend on low event*, page 168. After 30 minutes, basal insulin delivery resumes under the following conditions:

- SG is at least 20 mg/dL above the low limit.
- SG is predicted to be more than 40 mg/dL above the low limit within 30 minutes.

If the Suspend on low alarm is not cleared within two hours, the pump resumes insulin delivery and displays an emergency message.

### When the Suspend before low or Suspend on low features are unavailable

After a Suspend before low or Suspend on low event, both features are not active for a period of time to help prevent prolonged suspension of insulin delivery. Insulin delivery is suspended for a maximum of two hours. Insulin delivery can be manually suspended at any time. For details, see *Suspending all insulin delivery and resuming basal insulin delivery*, page 91.

When the Suspend before low and the Suspend on low features are unavailable, the suspend by sensor icon on the Home screen appears with a red X .

<b>Response to Suspend before low or Suspend on low events</b>	<b>Duration that the Suspend before low or Suspend on low feature is unavailable</b>
The alert is cleared within two hours and the pump stays suspended for the maximum two-hour suspend time.	The feature is unavailable for 30 minutes after basal insulin delivery resumes.
The alert is cleared within two hours and insulin delivery automatically resumes due to rising SG levels.	The feature is unavailable for 30 minutes after basal insulin delivery resumes.
The alert is cleared within two hours and basal insulin delivery is manually resumed.	The feature is unavailable for 30 minutes after basal insulin delivery resumes.
The alert is not cleared within 2 hours.	Basal insulin delivery automatically resumes and the feature is available.
The alert is cleared within 30 minutes after basal insulin delivery is automatically resumed.	The feature is unavailable for the remaining time left in the 30 minutes after basal insulin delivery resumed.
The alert is cleared between 30 minutes and four hours after basal insulin delivery is resumed.	The feature is available.
The alert is not cleared.	The feature is unavailable for four hours after basal delivery automatically resumes.

## Alert on low

The Suspend before low and the Suspend on low features automatically activate Alert on low. When Alert on low is on, the pump displays an alert when SG reaches or falls below the low limit. If insulin delivery is suspended by one of these features and the alert is not cleared, an emergency message appears.

## Automatically resuming basal insulin delivery after a Suspend before low or Suspend on low event

If insulin delivery is suspended by either the Suspend before low or the Suspend on low feature, basal insulin delivery automatically resumes under one of the following conditions:

- If insulin delivery is suspended for a minimum of 30 minutes and SG readings are at least 20 mg/dL above the low limit and expected to be more than 40 mg/dL above the low limit within 30 minutes
- After a maximum of two hours



## Resume basal alert

The Resume basal alert indicates when basal insulin is resumed automatically. When basal insulin delivery resumes and the Resume basal alert is off, a message appears indicating that basal insulin delivery has resumed.



If basal insulin delivery resumes after the maximum suspend time of two hours, an alert appears even if the Resume basal alert is off.

## Setting up the low SG settings

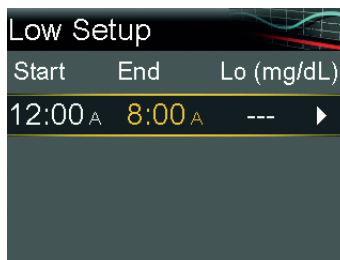
The Sensor feature must be turned on to set up the sensor settings. For more information, see *Turning on the Sensor feature*, page 173.

For information about the low SG settings, see *Low SG settings*, page 158.

### To set up the low SG settings:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Low Alert**.

The Low Setup screen appears.




3. Select the time segment. The end time flashes.

The start time of the first time segment is always 12:00 A. Up to eight time segments can be set, each with a different low limit. All the time segments must add up to a 24-hour period.

4. Set the End time.
5. Set your low limit. You can enter a value from 50 to 90 mg/dL, in increments of 5 mg/dL.

6. Select the arrow beside the End time to select the low SG settings for the time segment.

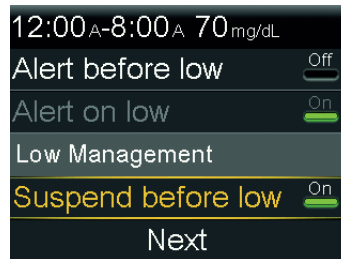A screen appears and shows the available settings for the selected time period.



7. Set the following alerts, as desired:
  - a. Select **Suspend before low** to set the pump to suspend insulin delivery before the low limit is reached.
  - b. Select **Alert before low** to receive an alert before the low limit is reached.
  - c. Select **Suspend on low** to set the pump to suspend insulin delivery when SG reaches or falls below the low limit.
  - d. Select **Alert on low** to receive an alert when SG reaches or falls below the low limit.
  - e. Select **Resume basal alert** to receive an alert when basal insulin delivery resumes during a suspend event. When this alert is off, the Basal delivery resumed message still appears.



**Note:** The Suspend before low and the Suspend on low features cannot both be on during the same time segment.



8. Select **Next**.
9. If necessary, enter the remaining time segments to complete the 24-hour period.



**Note:** For instructions on setting up more than one low limit over a 24-hour period, see *Settings covering a 24-hour period*, page 87.

10. Select **Review**.
11. Review the low SG settings, and select **Save**.


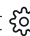
### To change the low SG settings:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Low Alert**.  
The Low Setup screen appears.
3. Select **Edit**.
4. Select and adjust the time segment.
5. Select any alert setting to make adjustments, or to turn the setting on or off.
6. Select **Next**.
7. Select **Review**.
8. Review the low SG settings, and select **Save**.

### Low Snooze

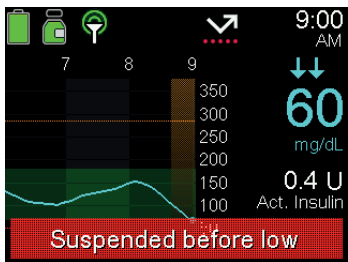
The Low Snooze feature sets the amount of time before a low alert repeats. The pump shows the low alert again if the low alert condition still exists after the specified snooze time.

### To set the Low Snooze:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Snooze High & Low**.  
The Snooze screen appears.
3. Select **Low Snooze** and enter a time in 5-minute increments from 5 minutes to 1 hour.
4. Select **Save**.



## Manually resuming basal insulin delivery during a Suspend before low or Suspend on low event

When the pump suspends insulin due to a Suspend before low or Suspend on low event, the Home screen shows which event is active in the red banner.




Basal insulin delivery automatically resumes when certain conditions are met. Basal delivery can be manually resumed at any time.

### To manually resume basal delivery:

1. From the Home screen, press , and then select .
2. Select **Resume Basal**.
3. Select **Yes** to resume basal insulin delivery.

## Silencing sensor alerts

The Alert silence feature silences certain sensor alerts for a set period of time. When using this option, the Alert silence icon  appears on the Home screen. The system still displays any alerts that occur, but there is no sound or vibration if they are silenced. This information can be reviewed in the Alarm History screen.



The Alert silence feature does not silence:

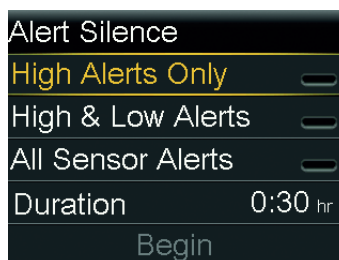
- **High SG alert**—When your sensor glucose (SG) value is above 250 mg/dL for more than three hours
- **Low SG alarm**—When your sensor glucose (SG) value falls below 64 mg/dL
- **SmartGuard exit alert**—When the pump exits the SmartGuard feature

The following table describes the sensor alerts that are silenced with each option.

Option	Silences these alerts
High Alerts Only	Alert on high, Alert before high, and Rise Alert
High & Low Alerts	Alert on high, Alert before high, Rise alert, Alert on low, Alert before low, Suspend before low, and Resume basal alert  <b>Note:</b> Alert on low cannot be silenced if the Suspend before low or Suspend on low features are turned on
All Sensor Alerts	All alerts listed previously for High & Low Alerts, as well as the following: <ul style="list-style-type: none"> <li>• All calibration alerts, reminders, or error messages that may result from entering a BG reading</li> <li>• All alerts related to sensor insertion, including alerts about sensor, changing the sensor, sensor expiration, sensor updating, and connection issues</li> <li>• All alerts related to the transmitter, including transmitter battery alerts and connection issues</li> </ul>

### To silence sensor alerts:

1. From the Home screen, press , and then select .
2. Select **Silence Sensor Alerts**.





3. Select **High Alerts Only**, **High & Low Alerts**, or **All Sensor Alerts**. Refer to the previous table for details about the alerts silenced with each selection.




**Note:** Silencing **All Sensor Alerts**, prevents the sound and vibration of most alerts related to SG readings, and the sensor. Silencing **All Sensor Alerts** does not silence the SmartGuard exit alert, the High SG alert, or the Low SG alarm for when SG is below 64 mg/dL.

4. Set the **Duration**. The duration can be set in 15-minute increments from 30 minutes to 24 hours.
5. Select **Begin**.

### To cancel Alert Silence:

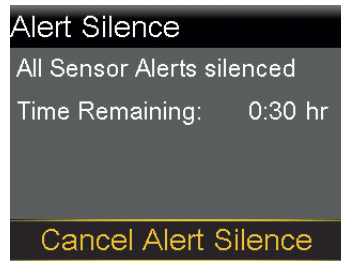1. From the Home screen, press , and then select .
2. Select **Alert Silence**.



3. Select **Cancel Alert Silence**.


## Guardian 4 sensor

### Inserting the sensor


The images that follow show insertion sites for ages 7 years and older.

Target the shaded areas shown in the image, and make sure that the insertion site has a sufficient amount of fat.

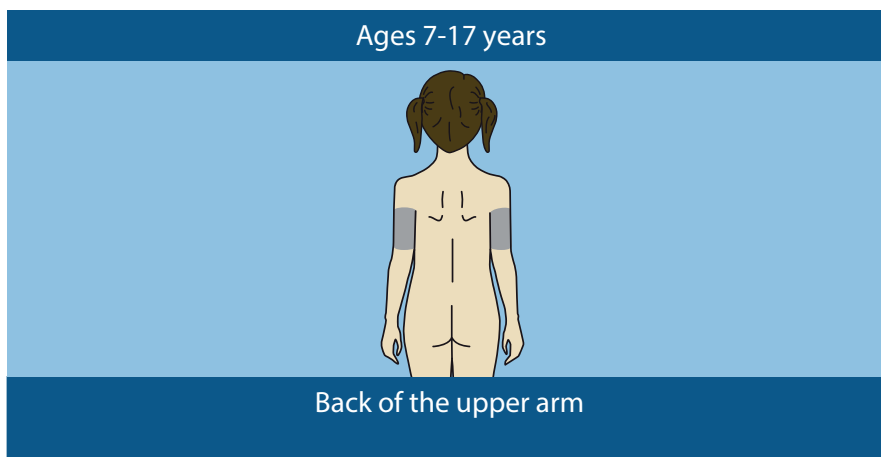Approved Age	Sensor Insertion Site
7 and older	Arm



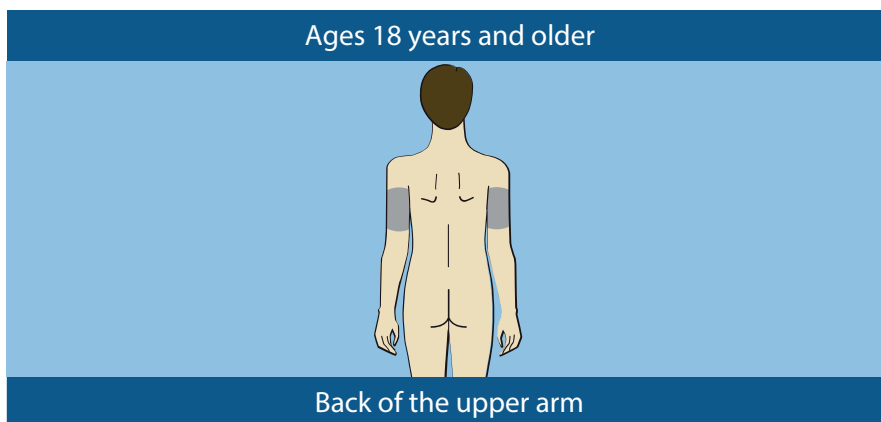Insertion on abdomen for ages 7-17 years has not been evaluated for accuracy.



**Note:** Assistance will likely be needed for sensor insertion into the back of the upper arm. Some users find it difficult to insert the sensor into their arm by themselves.

## Connecting the transmitter to the sensor

Refer to the transmitter user guide for instructions on how to connect the transmitter to the sensor.

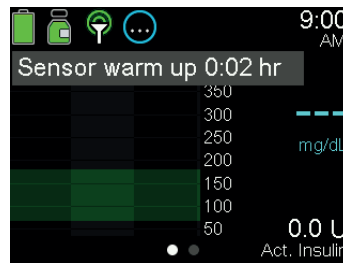
## Starting the sensor

After the sensor is inserted and paired with the transmitter, the pump will display a Start New Sensor screen.

### To start a new sensor:

1. Select **Start New Sensor** when it appears on the pump screen.

The “Sensor warm up X:XX hr” message appears.



**Note:** It may take up to five minutes for the “Sensor warm up X:XX hr” message to appear. The warm up period lasts two hours.

2. Select **OK**.

The “Sensor warm up X:XX hr” message appears on the Home screen until the sensor warm up is complete.

After the warm up is complete, the pump begins receiving SG readings.

## Reconnecting the sensor

If the transmitter is removed from a sensor while the sensor is inserted in the body, the pump detects when the transmitter is reconnected to the sensor and a “Sensor connected” message appears.



### To reconnect a sensor:

1. Select **Reconnect Sensor**.

The “Sensor warm up X:XX hr” message appears.



**Note:** It may take up to five minutes for the “Sensor warm up X:XX hr” message to appear. The warm up period lasts two hours.

2. Select **OK**.



The “Sensor warm up X:XX hr” message appears on the Home screen until the sensor warm up is complete.

After the warm up is complete, the pump begins receiving SG readings.

### Deactivating the Sensor feature

The Sensor feature may be turned off at any time. When the sensor is not used, turn off the Sensor feature to avoid a sensor alert. The Sensor feature must be turned on again before settings can be changed. The Sensor feature is turned on when the sensor is paired to the pump.



### To deactivate the Sensor feature:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Sensor**.
3. Select **Sensor**.
4. Select **Yes** to turn off the Sensor feature.

### Turning on the Sensor feature

The Sensor feature automatically turns on when a Guardian 4 sensor is paired to the pump. These instructions can help you if you turn the Sensor feature off during sensor use and need to turn the feature back on. The Sensor feature must be on before sensor alerts can be set up and SG levels can be monitored.

## To turn on the Sensor feature:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Sensor**.
3. Select **Sensor** to turn the feature on or off.

## Removing the sensor

### To remove the sensor:

1. Gently peel the sensor adhesive away from the body.
2. Dispose of the Guardian 4 sensor in accordance with all local laws and regulations. See the Disposal section in the sensor user guide.

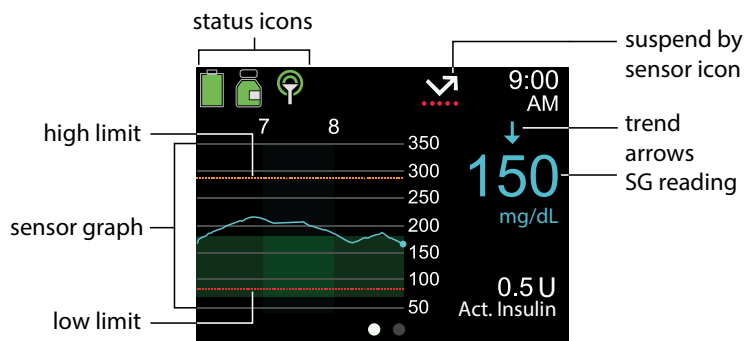
## Using CGM

### Home screen with CGM in Manual mode

When the Sensor feature is active, the Home screen displays a real-time graph that shows CGM information.




**Note:** To see the Home screen while the SmartGuard feature is active, see *Home screen with the SmartGuard feature, page 187*.



For more information about the icons that appear on the Home screen with CGM in Manual mode, see *Status icons, page 72*.

## Sensor graph

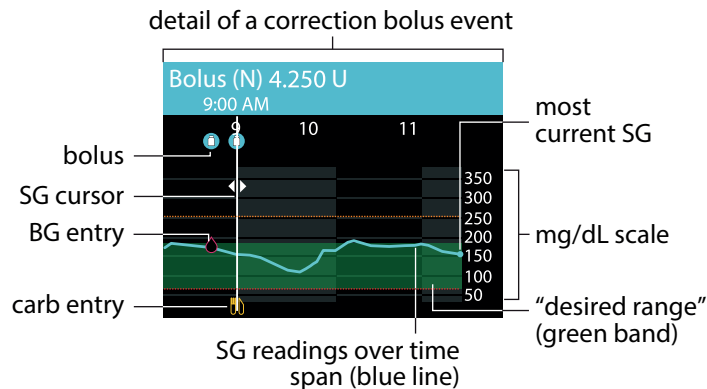

The sensor graph includes the following information:

- The most recent SG reading.
- Historical SG readings for the last 3-hour, 6-hour, 12-hour, or 24-hour periods.
- High and low SG limits.
- Carb entries.
- Delivered boluses.
- Suspend events caused by Suspend before low or Suspend on low.
- BG entries.





There are several reasons why an SG reading may not appear on the graph:

- A recently inserted sensor is warming up.
- A recently connected or reconnected transmitter is not ready.
- An error condition or a sensor-related alert has occurred or is occurring. For a list of sensor alerts, see *CGM (sensor) alarms, alerts, and messages*, page 306.

### To view the sensor graph:







1. From the Home screen, press the  button.

A full-screen view of the 3-hour graph appears.

2. Press  to navigate to the 6-hour, 12-hour, and 24-hour graphs.
3. Press  to view SG readings and event details.
4. To exit the full-screen view, press , or press the  button again.

## Trend arrows

The trend graph indicates how sensor glucose (SG) may have recently changed. The trend arrows indicate the rate at which the most recent SG readings are rising or falling. SG readings may trend up or down during certain activities, such as eating, giving a bolus, or when exercising. These icons appear only when the sensor feature is turned on.

- SG with no trend arrows: Rate of change of SG is less than 1 mg/dL per minute.
-  or : SG has been rising or falling at a rate of 20-40 mg/dL over the last 20 minutes, or 1-2 mg/dL per minute.
-  or : SG has been rising or falling at a rate of 40-60 mg/dL over the last 20 minutes, or 2-3 mg/dL per minute.
-  or : SG has been rising or falling at a rate of more than 60 mg/dL over the last 20 minutes, or more than 3 mg/dL per minute.

Consider any active insulin that is available. Active insulin may cause SG to decrease and can affect treatment decisions. For more information about active insulin, see *Bolus Wizard settings in Manual mode*, page 99.

# 8





# 8

## SmartGuard

This chapter provides information about how to set up and start using the SmartGuard feature. The SmartGuard feature uses sensor glucose (SG) values provided by the Guardian 4 sensor to automatically adjust insulin delivery.



**WARNING:** Do not use the SmartGuard feature for people who require less than 8 units or more than 250 units of total daily insulin per day. A total daily dose of at least 8 units, but no more than 250 units, is required to use the SmartGuard feature.

### Introduction

The SmartGuard feature uses carb information, SG, and SmartGuard target values to control insulin delivery. It also can automatically deliver a correction bolus to help correct a high SG reading. The MiniMed 780G insulin pump requires a minimum of eight units and a maximum of 250 units per day to operate using the SmartGuard feature.



**Note:** The Auto correction feature uses SG values to determine bolus insulin doses. Auto correction boluses are delivered without user acknowledgment. The accuracy of SG values can be lower than the accuracy of blood glucose (BG) meter readings, which are checked with a blood glucose (BG) meter.

The SmartGuard feature is designed to maximize the amount of time that glucose levels stay in the range of 70 mg/dL to 180 mg/dL. The following table describes features that the system uses to maximize time in range.

Feature name	Description
SmartGuard target: 100 mg/dL, 110 mg/dL, or 120 mg/dL	Consult a healthcare provider to determine which SmartGuard target to use to maximize time in range. The default setting is 100 mg/dL.
Auto Basal	When using the SmartGuard feature, basal insulin is automatically delivered based on SG readings and recent insulin delivery needs.
Target for Auto correction bolus based on SG: 120 mg/dL	The MiniMed 780G system may deliver a bolus automatically, as frequently as every 5 minutes, if the SmartGuard feature determines that a correction bolus is necessary. The default setting for Auto correction is set to On.
Temp Target: 150 mg/dL	A temp target can be set for events such as exercise or other times when less insulin is needed. If a temp target is used for exercise, consider starting it one to two hours before beginning the exercise. Auto correction boluses are not delivered while a temp target is active.



**Note:** When using the SmartGuard feature, meal boluses are still required.

The SmartGuard feature requires accurate sensor measurements and carb information to deliver insulin for meals. This insulin therapy requires the use of the Bolus feature to deliver boluses to cover meals.

When using the SmartGuard feature:

- If an Enter BG now alert occurs, enter a BG meter reading.
- Do not enter an SG reading when the system requests a BG reading.
- Bolus amount cannot be modified when delivering a bolus in the SmartGuard feature. If SG readings do not match with symptoms, enter a BG value from a BG meter.



## Auto Basal

When the SmartGuard feature is active, the basal insulin dose is calculated using SG values from the sensor. The automatic delivery of insulin is called Auto Basal.



**Note:** In Manual Mode, basal insulin will not be delivered if basal settings are not entered and saved. There will be no message that the basal rates are not programmed.

## Auto Correction

The pump may deliver a bolus automatically when the SmartGuard feature determines it is needed for correction, to maximize the time in range, between 70 mg/dL and 180 mg/dL. Because this is an automated bolus, no action is required. The Home screen shows when an Auto Correction bolus occurs.

## Giving a bolus when the SmartGuard feature is active

A meal bolus can be delivered while using the SmartGuard feature. For more information, see *Delivering a bolus in the SmartGuard feature*, page 190.



**WARNING:** Always confirm an SG value that does not match your symptoms. When the SmartGuard feature is active and you are no longer in Manual mode, the pump uses an available SG value to calculate a bolus amount. However, if your symptoms do not match, the SG value can result in the infusion of too much or too little insulin, which may cause hypoglycemia or hyperglycemia.

## Preparing to set up the SmartGuard feature

The SmartGuard feature requires a 48-hour warm-up period before activation. This warm-up period begins at midnight after the pump starts delivering insulin and it does not require sensor use. During the warm-up period, the pump collects and processes data for use by the SmartGuard feature.



**Note:** A basal pattern must be programmed for use during the warm-up period and for instances when the pump is in manual mode. During the warm-up period the pump should also be used to give boluses.

### To prepare the pump for the SmartGuard feature:

1. Cancel any active Temp Basal rates. See *Canceling a temp basal or preset temp basal*, page 242.
2. Confirm that insulin delivery is not suspended. See *Suspending all insulin delivery and resuming basal insulin delivery*, page 91.
3. Set the carb ratio. See *Changing the carb ratio*, page 252.
4. Review the high and low limit settings. High and low limit settings apply when in Manual mode and when using the SmartGuard feature. See *Sensor alert settings*, page 153 for details.
5. Enter a new BG reading.





**WARNING:** If the pump has been used in the last 21 days to practice button pressing, or if insulin that was programmed into the pump was not the user's actual insulin delivery, clear active insulin and the total daily doses tracked by the SmartGuard feature before using the SmartGuard feature. Failure to do so may result in the delivery of too little or too much insulin, which can cause hyperglycemia or hypoglycemia. The SmartGuard feature uses the recent delivery history on the pump to determine the insulin delivery amount.

Consult with your healthcare professional about using the Clear Active Insulin feature in the Manage Settings menu to clear both active insulin and the total daily dose for the SmartGuard feature.

## Setting up the SmartGuard feature

The SmartGuard feature requires 48-hours of insulin delivery before the feature can be used. This warm up period begins at the first midnight after delivery has started. For more information, see *Preparing to set up the SmartGuard feature, page 181*.

### To set up the SmartGuard feature:

1. From the Home screen, press , and then select .
2. Select **SmartGuard** to turn the feature on or off.



**Note:** Certain additional requirements must be met before the SmartGuard feature activates. For more information, see *SmartGuard Checklist, page 184*.

3. Select **SmartGuard Settings** and enter the following information:
  - Select the SmartGuard target: 100 mg/dL, 110 mg/dL, or 120 mg/dL.
  - Confirm that **Auto Correction** is on to activate automatic correction boluses.



**Note:** The Auto correction feature is turned on by default. When this setting is on, the pump automatically delivers correction boluses to help correct a high SG reading. For information, see *Delivering a bolus in the SmartGuard feature, page 190*.

4. Select **Save**.

### Conditions to activate the SmartGuard feature

If the pump is turned off for more than 2 weeks and is turned back on, the pump requires a 48 hour warm-up period before the SmartGuard feature activates.

If the pump has been off for 2 weeks or less and is turned back on, a 5 hour warm-up period is required before the SmartGuard feature activates.

If the SmartGuard feature is on but not active, the SmartGuard Checklist screen indicates the requirements needed to activate the SmartGuard feature. See *SmartGuard Checklist*, page 184.

The system requires 5 hours for the SmartGuard active insulin amount to update. This update time begins under the following conditions:

- A complete pump reset caused by a loss of power or a software error.
- When the insulin is resumed after being manually suspended for 4 hours or longer.
- When you update the pump software.

SmartGuard active insulin information is valid until one of the conditions listed above occurs, which restarts the 5 hour update time. The SmartGuard feature is unavailable during this time.

### Suspending manually while using the SmartGuard feature




For information about manually suspending insulin delivery, see *Suspending all insulin delivery and resuming basal insulin delivery*, page 91.


### Suspend before low and Suspend on low features while using the SmartGuard feature

When the SmartGuard feature is active, the Suspend before low and the Suspend on low features are unavailable and automatically turn off. If the system exits the SmartGuard feature, the Suspend before low and the Suspend on low features return to the state they were in before using the SmartGuard feature. For information about turning on the Suspend before low or the Suspend on low feature, see *Low SG settings*, page 158.

### SmartGuard Checklist

The SmartGuard Checklist screen indicates the requirements necessary to start or continue using the SmartGuard feature. For more information, see *Staying in the SmartGuard feature*, page 197.

The following table shows what to do when the wait icon  or the question icon  appear by items on the SmartGuard Checklist screen. In general, the wait icon 

indicates the system is processing or waiting for information. In general, the question icon  indicates the system needs an action from you.

1

2

3

4

5

6

7

8

9

SmartGuard Checklist

Enter BG

SmartGuard turned off

Sensor not ready

Bolus in progress

Delivery suspended










Carb ratio not set



Basal rate for Manual mode not set

SmartGuard updating

SmartGuard warming up

Line	Item	Instructions
1	Calibrating...	The system is using the recent BG meter reading to calibrate the sensor.
	Enter BG	Enter a new BG meter reading.
	Wait to calibrate	The system requires a BG reading and will ask when it is ready.
2	SmartGuard turned off	Turn on the SmartGuard feature.

Line	Item	Instructions
3	Sensor not ready 	<ul style="list-style-type: none"> <li>Confirm the pump shows a sensor serial number on the Paired Devices screen. Example: CGM XXXXXXXX Make sure the pump is paired with a sensor. For more information, see <i>Pairing the pump and transmitter</i>, page 137.</li> <li>Check the Home screen. If  displays, move the pump and sensor closer together. It may take 15 minutes to find the sensor signal. If after 30 minutes the pump and sensor are still not communicating, a Lost sensor signal alert appears. Check that the sensor is still inserted in the skin. Move the pump closer to the sensor.</li> </ul>
	Sensor off 	Turn on the Sensor feature in Settings > Device Settings.
	No paired CGM 	Pair the pump and sensor. For more information, see <i>Pairing the pump and transmitter</i> , page 137.
4	Bolus in progress 	Wait until the bolus is complete or stop the bolus before the SmartGuard feature can be used.
5	Delivery suspended 	If insulin delivery is suspended, the SmartGuard feature cannot be used. Treat low BG as instructed by a healthcare professional.
6	Carb ratio not set 	Enter a carb ratio in the Bolus Wizard feature or in the Bolus Wizard Setup screen.
7	Basal rate for Manual mode not set 	When no basal pattern is set, the SmartGuard checklist will display “Basal rate for manual mode not set”. You must program, confirm, and save a basal pattern before the pump will enter SmartGuard.
	Temp Basal rate 	When a basal pattern is set and the pump is currently running a temp basal, the SmartGuard checklist will display “Temp Basal rate”. Stop the temp basal rate delivery before the SmartGuard feature can be used or wait until the temp basal rate delivery is complete.

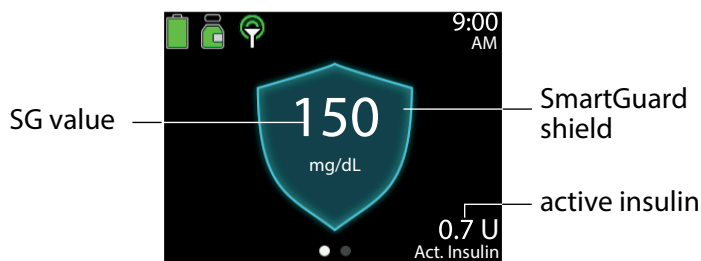
Line	Item	Instructions
8	SmartGuard updating 	If SmartGuard active insulin is updating, it will take up to five hours to complete. Wait for the update time to end before the SmartGuard feature can activate.
9	SmartGuard warming up 	Wait for the SmartGuard feature to gather insulin delivery history and determine the basal rate.

### To view the SmartGuard Checklist:

1. From the Home screen, press , and then select .
2. Select **SmartGuard Checklist**.

## Home screen with the SmartGuard feature


When the pump is using the SmartGuard feature, the **Home** screen displays a shield with the current SG, if available. When a BG is entered, it appears on the Home screen until either the next SG arrives or the BG is 12 minutes old.



## Using the SmartGuard feature
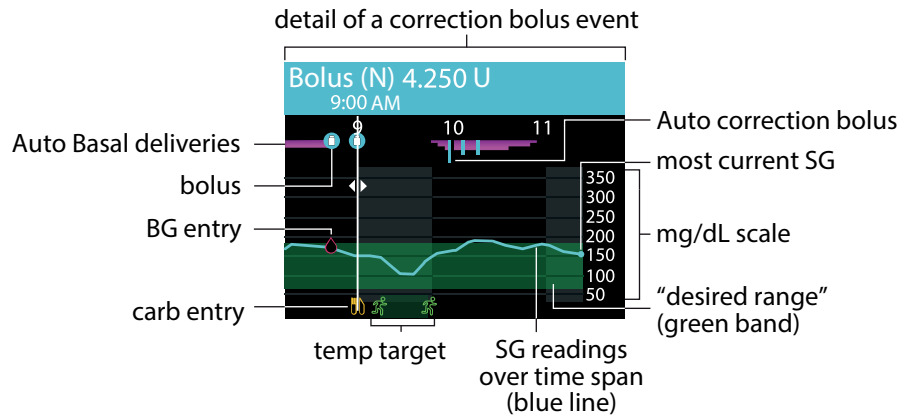
### The sensor graph with the SmartGuard feature



The sensor graph with the SmartGuard feature shows historical SG readings provided by the sensor.



The SmartGuard feature sensor graph includes the following information:

- When a location on the graph is selected, specific details of the SG or event appear, such as a correction bolus.
- Historical SG readings are displayed for the last 3-hour, 6-hour, 12-hour, or 24-hour periods. They appear as a blue line across the screen.
- Boluses are shown as white vials inside blue circles.
- Carb entries are shown as yellow knife and fork symbols. These represent any bolus amounts that include a carb entry.
- BG entries appear as red drop symbols.
- Magenta bands across the top represent Auto Basal deliveries provided by the SmartGuard feature.
- Blue vertical bars at the top represent Auto correction boluses delivered by the SmartGuard feature.
- A time change event appears as a white clock symbol.
- Temp target is shown as green runners.

### To view the sensor graph:

1. From the Home screen, press the  button to display the SG graph.  
A full-screen view of the 3-hour graph appears.
2. Press  to navigate to the 6-hour, 12-hour, and 24-hour graphs.



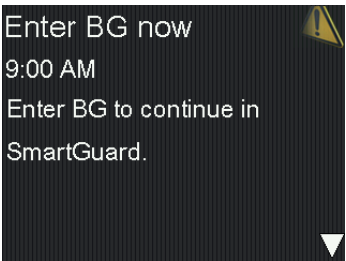
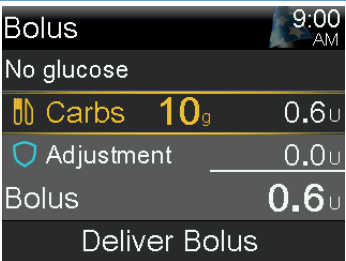
3. Press < to view SG readings and event details.
4. To exit the sensor graph, press ↶ or press the ⬢ button again.

## Entering a BG value in the SmartGuard feature

The pump may require a BG value to continue using the SmartGuard feature.

There are two ways to enter a BG value when using the SmartGuard feature. Manually enter a BG value or enter a BG value using the compatible Accu-Chek Guide Link meter. For more information on manually entering a BG, see *Entering a blood glucose (BG) meter reading*, page 96.

The following table shows when to use a blood glucose (BG) meter reading:

When to use a blood glucose (BG) meter reading	Examples
Anytime the system requests a blood glucose (BG) meter reading.	 A screenshot of the SmartGuard interface. It displays the text 'Enter BG now' at the top, followed by '9:00 AM' and 'Enter BG to continue in SmartGuard.' There is a yellow warning icon in the top right corner and a white downward arrow in the bottom right corner.
Anytime you deliver a bolus in SmartGuard when a sensor glucose (SG) value is not displayed on the bolus screen and you want to use a glucose for a correction.	 A screenshot of the SmartGuard bolus screen. It shows 'Bolus' at the top right with a clock icon and '9:00 AM'. Below that is 'No glucose'. Then, 'Carbs 10g' is shown with a carb icon and '0.6U'. Below that is 'Adjustment' with a blue circle icon and '0.0U'. At the bottom, 'Bolus' is shown with '0.6U' and a 'Deliver Bolus' button.
When using a medication that impacts glucose levels.	
When your sensor glucose (SG) values are different than the symptoms you are experiencing.	
The most recent sensor glucose (SG) reading is unavailable. Sensor glucose (SG) readings are unavailable in the following conditions: <ul style="list-style-type: none"><li>• A new sensor is started.</li><li>• A <b>Sensor Updating</b> notification appears.</li></ul>	

When to use a blood glucose (BG) meter reading	Examples
<ul style="list-style-type: none"><li>The sensor requires a new blood glucose (BG) meter reading to be entered because the system was unable to use the blood glucose (BG) meter reading that was entered to calibrate the sensor. All blood glucose (BG) meter readings that are entered are used to calibrate the sensor.</li></ul>	
There is doubt that sensor glucose (SG) values are correct.	



**WARNING:** Consult a healthcare professional if a medication that contains acetaminophen or paracetamol is taken while wearing the sensor. Medications that contain acetaminophen or paracetamol can falsely raise sensor glucose readings. The level of inaccuracy depends on the amount of acetaminophen or paracetamol active in the body and can differ for each person. Falsely elevated sensor readings can result in over-delivery of insulin, which can cause hypoglycemia. Medications that contain acetaminophen or paracetamol include, but are not limited to, cold medicines and fever reducers. Check the label of any medications being taken to see if acetaminophen or paracetamol is an active ingredient. Use additional blood glucose meter readings to confirm blood glucose levels.

While the SmartGuard feature is active, if acetaminophen or paracetamol is taken, program a temp target for up to eight hours, or the amount of time recommended by a healthcare provider. For more information, see *Setting a temp target, page 196*. Use blood glucose values instead of sensor glucose readings to calculate a meal bolus or correction bolus for up to eight hours, or the duration recommended by a healthcare provider, after taking acetaminophen or paracetamol.

### Delivering a bolus in the SmartGuard feature

A current BG or SG reading is used to determine the bolus amount. A carb amount can be entered for a food bolus.



**WARNING:** Do not use the SmartGuard feature for a period of time after giving a manual injection of insulin by syringe or pen. Manual injections are not accounted for in the active insulin amount. Using the SmartGuard feature after a manual injection may result in over-delivery of insulin. Too much insulin may cause hypoglycemia. Consult a healthcare professional for how long to wait after a manual injection before resuming the SmartGuard feature.

If the BG or SG is under 120 mg/dL, or if the bolus is zero after the pump accounts for active insulin, or if the SmartGuard feature estimates current basal delivery is sufficient, no correction is recommended.

A paired Accu-Chek Guide Link meter sends BG readings directly to the pump. Confirm the BG reading for use in the SmartGuard feature. If an Accu-Chek Guide Link meter is not used, the BG value must be manually entered on the BG screen or the Bolus screen while using the SmartGuard feature.



**Note:** Do not use a blood glucose (BG) meter reading if more than 12 minutes have passed since the last BG meter reading was taken. That BG meter reading and the calculated bolus amount may no longer be accurate.



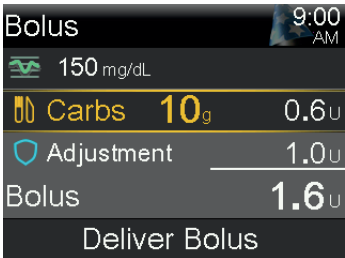

**WARNING:** SG readings are used to calculate meal boluses or correction boluses when delivering a bolus in the SmartGuard feature. SG is not the same as BG. Sensor performance may occasionally vary from sensor to sensor and in different situations for a sensor, such as on the first day of use.




When SG readings are used for meal boluses and for correction boluses, there is a risk of both hypoglycemia and hyperglycemia. If an SG reading is much lower than a BG reading would be at that time, there is a risk of hyperglycemia, because the amount of insulin delivered could be smaller. If an SG reading is much higher than a BG and there are symptoms of feeling low, but the SG reading is not low, and if there are symptoms of a severe hypoglycemic event, a severe hyperglycemic event, or diabetic ketoacidosis, a BG meter reading is needed.

This can also occur when SG readings are used when the Auto correction feature is turned on. For example, when an SG reading is much higher than a BG reading at that time, there is a risk of hypoglycemia, because the amount of insulin delivered could be larger.

If there are symptoms of feeling low, but the SG reading is not low, and if there are symptoms of a severe hyperglycemic event or diabetic ketoacidosis, a BG meter reading is needed.





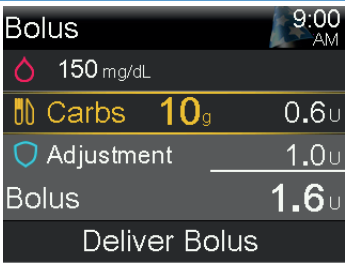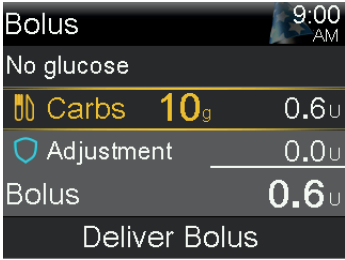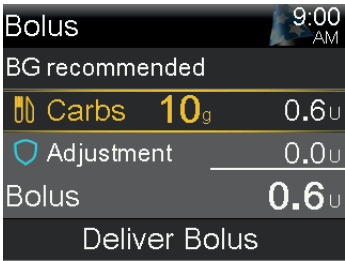The following table describes how glucose readings are shown on the SmartGuard bolus screen.

Bolus screen	Glucose reading information
	The  icon indicates the sensor glucose (SG) value when there is no recent blood glucose (BG) meter reading available during the last 12 minutes. A blood glucose (BG) meter reading can be entered to calculate a correction bolus. The

Bolus screen	Glucose reading information
	correction bolus is included in the Adjustment.
	A blood glucose (BG) meter reading is available to calculate a correction bolus. The correction bolus is included in the Adjustment.
	There are no blood glucose (BG) meter readings or sensor glucose (SG) values available. You can enter a carb amount for a food bolus or a blood glucose (BG) meter reading for a correction bolus.
	<p>The BG recommended message indicates that neither a blood glucose (BG) meter reading nor a sensor glucose (SG) reading is available to calculate a correction bolus.</p> <p><b>Note:</b> If a sensor glucose (SG) value shows on the Home screen, but does not show on the Bolus screen, the system determined that the sensor glucose (SG) value is not optimal to use to calculate a correction bolus. Enter a blood glucose (BG) meter reading if a correction bolus is desired.</p>

### Bolus adjustments in the SmartGuard feature

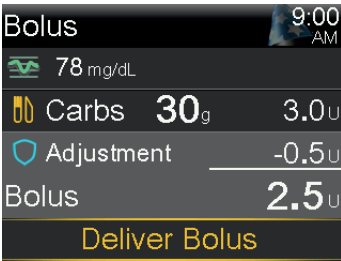
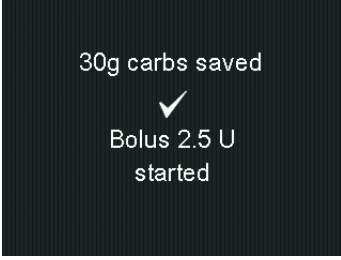
The SmartGuard feature calculates a bolus based on the current BG or SG reading and carbs, and may make an additional adjustment to the bolus.



Bolus adjustment	Example screens
<p>The bolus amount is adjusted down if the SmartGuard feature predicts a risk of hypoglycemia after the meal. Carbs are saved for use in future bolus adjustment calculations.</p>	 
<p>If the bolus amount is adjusted down to 0.0 for the bolus, no bolus is delivered. Carbs are saved for use in future bolus adjustment calculations.</p>	 

**Bolus adjustment**

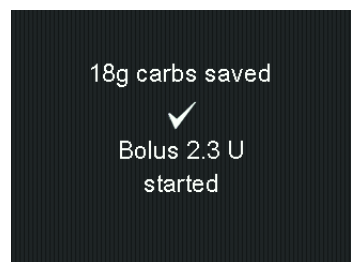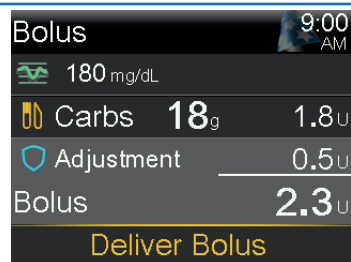
The bolus amount is adjusted up if a correction bolus is calculated based on high glucose and low active insulin.

Carbs are saved for use in future bolus adjustment calculations.

**Example screens****To deliver a bolus with the SmartGuard feature:**

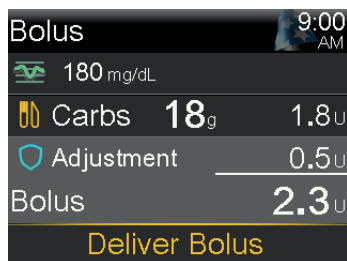
1. From the Home screen, press , and then select .
2. Select **Bolus**.
3. Enter a carb amount, if desired.



The screen indicates the amount of the calculated bolus.



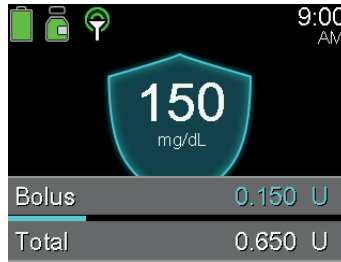4. Select **Deliver Bolus**.

A screen appears briefly to indicate the bolus delivery has started. The Home screen appears and shows the progress of the bolus delivery.



**Note:** To stop a bolus, press  from the Home screen, select  and then select **Stop Bolus**. Select **Yes** to confirm.



## Setting a temp target

A temporary target (temp target) of 150 mg/dL can be set for events such as exercise or other times when less insulin is needed. Consult a healthcare professional before using a temp target.



**Note:** The Auto correction feature is not active during an active temp target. It resumes after the temp target completes.

### To set a temp target:

1. From the Home screen, press , and then select .
2. Select **Temp Target** to turn the feature on or off.

Temp Target	
Duration	2:00 hr
Set target to 150 mg/dL	
<b>Start</b>	

3. Set the Duration, from 30 minutes to 24 hours, in 30-minute increments.
4. Select **Start**.

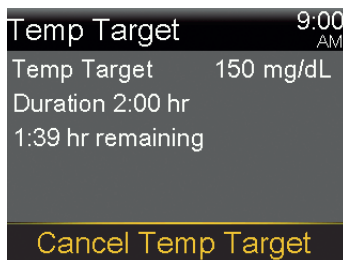


The screen shows a Temp Target Started message, and then changes to the Home screen, where a banner shows the remaining temp target time.



### To cancel a temp target:

1. From the Home screen, press , and then select .




2. Select **Cancel Temp Target**.

### Staying in the SmartGuard feature

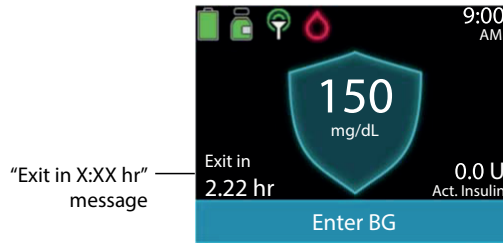When the pump requires an action to stay in the SmartGuard feature, it delivers insulin at a fixed basal rate for up to a maximum of 4 hours.

The message “Exit in X:XX hr” appears on the Home screen, showing the time remaining before the pump enters Manual mode. The basal rate delivered during this time is based on insulin delivery history and represents a delivery rate that minimizes the risk of hypoglycemia in situations when SG values are temporarily unavailable. The pump provides a notification of any required actions.



The pump resumes using SG readings for basal insulin delivery when certain conditions are met. The following table describes these conditions and the notification and required action to resume using SG readings for basal insulin delivery.

Condition	Notification and action
The SmartGuard feature has reached the time limit for minimum delivery. The minimum delivery time is 3 to 6 hours, depending on the reason.	A SmartGuard min delivery alert appears. Enter a BG.
The SmartGuard feature has been delivering basal insulin at its maximum limit for 7 hours.	A SmartGuard max delivery alert appears. Check the SmartGuard Checklist to determine the required steps. Enter a BG.
SG readings may be lower than actual glucose values.	An Enter BG now alert appears. Enter a BG.
No SG data has been received for more than 5 minutes.	<ul style="list-style-type: none"> <li>If SG data is not available, three dashes appear on the screen in place of the SG data. If the loss of SG data is intermittent, no action is required.</li> <li>If an action is required, an alert appears such as a Lost sensor signal alert or the Enter BG now alert. Follow the instructions on the screen.</li> </ul>



**Note:** To stay in the SmartGuard feature when changing the sensor, be sure the sensor warmup completes within four hours of the last available SG reading.

## Exiting the SmartGuard feature

The SmartGuard feature may stop functioning under the following conditions:

- The SmartGuard feature is turned off.
- The pump is delivering basal insulin based on insulin delivery history, and not SG readings, for four hours. See *Staying in the SmartGuard feature*, page 197.
- All insulin delivery has been manually suspended and has not resumed for four hours.
- The Sensor feature is turned off.

The SmartGuard feature can be turned off at any time. For more information, see *Setting up the SmartGuard feature*, page 183.

### Returning to the SmartGuard feature after an exit

The pump indicates any required actions on the Home screen, after an exit from the SmartGuard feature. In the example below, a BG entry is needed. Once the BG is entered, the pump resumes using the SmartGuard feature.



While in Manual mode, resume using the SmartGuard feature by meeting all requirements in the SmartGuard Checklist. For more information, see *SmartGuard Checklist*, page 184.

The SmartGuard feature can be resumed under the following conditions:

- The SmartGuard feature is turned on.
- The sensor is providing SG readings.
- A bolus is not in progress.
- A temp basal rate is not in progress.
- The 48-hour warm-up is complete.

- The SmartGuard feature is not in a 5-hour warm-up period.
- A new BG reading is entered.

If any of these conditions are not met, the SmartGuard feature cannot restart.

### **Using Block mode with the SmartGuard feature**

Block mode lets caregivers lock the pump to restrict access to critical pump features. While the pump is locked, Auto Basal delivery is active, and Auto correction boluses can occur if the feature is turned on. BG readings received from the Accu-Chek Guide Link meter can be confirmed. For more information on Block mode, see *Block mode*, page 204.

### **Alert silence feature**

The Alert silence feature silences certain sensor alerts for a set period of time. For more information, see *Silencing sensor alerts*, page 168.

9



# 9



## General settings

This chapter provides information about common tasks for various settings.

### Time and date

Confirm that the time and date are always set correctly on the MiniMed 780G insulin pump. Incorrect time and date settings can affect basal insulin delivery and the accuracy of pump history. Change the time or the date to match the time zone or daylight saving time. After the time and date are changed, the pump adjusts all settings automatically.



#### To change the time and the date:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Time & Date**.
3. Select and change the **Time**, **Time Format**, or **Date** as necessary. If a 12-hour clock is being used, specify AM or PM.
4. Select **Save**.

### Display options

The brightness of the pump screen can be controlled from the Display Options screen. The duration the backlight is on can also be adjusted.

## To adjust the display options:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Display**.
3. Select **Brightness** to adjust the brightness of the screen. A level from 1 to 5 can be set, or select **Auto** for the screen to automatically adjust to the current environment.
4. Select **Backlight** to adjust the timeout for the backlight on the pump screen. Select 15 seconds, 30 seconds, 1 minute, or 3 minutes.
5. Select **Save**.



**Note:** The brightness and backlight can affect the life of the battery. Use a lower brightness level setting, and set the backlight timeout to 15 or 30 seconds to help the battery last longer.



**CAUTION:** If you have not pressed Save after settings are entered and the screen goes dark, the entered settings will not be saved.

## Block mode

Block mode lets caregivers lock the pump to restrict access to critical pump features. While the pump is in Block mode, the pump automatically locks two minutes after the screen goes dark from inactivity.



**WARNING:** Always monitor the pump while it is locked. The pump can still be manually suspended while locked using the shortcut to the Status screen, which could result in hyperglycemia and ketoacidosis.

The following are examples of functions that are blocked while the pump is locked:

- Access the Menu screen.
- Deliver a bolus
- Start a new basal pattern





- Start a new temp basal delivery
- Change settings

The following are examples of important functions that remain available while the pump is locked:

- Previous bolus and basal deliveries continue normally
- Stop a bolus delivery using the shortcut to the Status screen
- Suspend and resume insulin delivery using the shortcut to the Status screen
- Receive sensor glucose (SG) values and receive and accept blood glucose (BG) meter readings
- Clear alarms and alerts

### To turn Block mode on or off:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Block Mode**.
3. Select **Block Mode** to turn the feature on or off.
4. Select **Save**.



The pump is in Block mode, but it is not yet locked.

### To lock the pump:

Press and hold  to manually enter Sleep mode.

The pump locks when it goes to sleep. While the pump is locked,  appears on the Home screen.

### To unlock the pump:

1. Press any button to wake up the pump.
2. Press .
3. Press and hold .

The Screen locked message appears.



**Note:** When the pump goes to sleep it will lock again.

## Self Test

The **Self Test** option can be used to confirm the pump is operating properly. Self test is additional to the routine tests that run independently while the pump operates.





**Note:** Insulin delivery is suspended for up to two minutes while the pump runs a self test.

The **Self Test** option includes the following tests. Observe the pump during these tests.

Test	Description
Display	The display turns on for up to 45 seconds.
Notification light	The notification light turns on for three seconds, and then it turns off.
Vibration	Two vibration tones are generated.
Tone	An alert tone, an Easy bolus step tone, and an alarm tone are generated.

### To run the self test:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Self Test**.

A message confirms self test is in progress.

Self test takes up to two minutes to complete. During that time, the display briefly turns white, the notification light blinks, the pump vibrates and then beeps.

If self test does not detect a problem, the Device Settings screen appears. If a problem is detected, a message appears with more information.

If an error message appears or the pump does not perform as indicated during the test, contact 24-Hour Technical Support.

## Manage Settings

The Manage Settings screen includes the following options:





- **Save Settings**
- **Restore Settings**
- **Clear All Settings**
- **Clear Active Insulin**
- **Settings History**
- **Max basal rate**

For information on how to use these options, see the procedures in this section.

### Saving the settings

The Save Settings option saves a record of the settings to restore the settings at a later date, if necessary.

#### To save the current settings:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Manage Settings**.
3. Simultaneously press and hold  and  until the Manage Settings screen appears.
4. Select **Save Settings**.





If these are the first settings saved, a message confirms that the settings are saved.

If the settings have been saved previously, a screen asks to replace the previous settings with the current settings. Select **Yes** to accept. Select **No** to cancel.

### Restoring the settings

The **Restore Settings** option replaces the current pump settings with the last settings that were saved. The **Restore Settings** option is available only if settings were previously saved.

### To restore the previous settings:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Manage Settings**.
3. Simultaneously press and hold  and  until the Manage Settings screen appears.
4. Select **Restore Settings**.  
A screen asks to confirm.
5. Select **Yes** to accept. Select **No** to cancel.

### Clearing the settings





The **Clear All Settings** option erases the current settings and returns them to the factory defaults. After the settings are cleared, the Startup Wizard appears and pump settings can be re-entered. The settings must be entered to continue using the pump.

The Clear All Settings option does not delete paired devices, such as the transmitter or meter.



**CAUTION:** Do not clear the pump settings unless directed by a healthcare professional. If pump settings are cleared they must be re-programmed as directed by a healthcare professional.

### To clear all settings:

1. Disconnect the pump from the body.
2. From the Home screen, press , and then select .
3. Select **Device Settings** > **Manage Settings**.
4. Simultaneously press and hold  and  until the Manage Settings screen appears.
5. Select **Clear All Settings**.  
A screen asks to confirm.
6. Select **Yes** to continue. Select **No** to cancel.

After the settings are cleared, the Startup Wizard appears. For more details on entering the startup settings, see *Startup settings*, page 69.





## Clearing the active insulin

Use the **Clear Active Insulin** option to use the pump with insulin for the first time. This option clears the SmartGuard therapy history and any active insulin values that the pump has tracked.

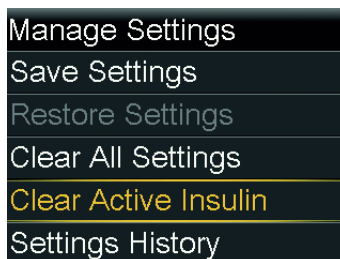
After the existing insulin values are cleared, it sets the active insulin value to zero. If bolus delivery was practiced with the pump prior to using the pump with insulin, the active insulin must be cleared. Clearing active insulin confirms that the Bolus Wizard feature has an accurate active insulin amount for bolus calculations.

Active insulin can be cleared only once. After the active insulin is cleared, this option is no longer available.

### To clear the active insulin:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Manage Settings**.
3. Simultaneously press and hold  and  until the Manage Settings screen appears.

The Manage Settings screen appears. If the active insulin has never been cleared, the **Clear Active Insulin** option appears.



**Note:** If the **Clear Active Insulin** option does not appear on the Manage Settings screen, the active insulin has already been cleared.

4. Select **Clear Active Insulin**.

A screen asks to confirm.





5. To clear the active insulin, select **Clear**. If the active insulin should not be cleared, select **Cancel**.

A message confirms that the active insulin is cleared.

## Viewing the pump setting history

The **Settings History** option shows a history of activities performed through the Manage Settings screen, such as when pump settings were saved, restored, or cleared.

### To view the pump setting history:





1. From the Home screen, press , and then select .
2. Select **Device Settings > Manage Settings**.
3. Simultaneously press and hold  and  until the Manage Settings screen appears.
4. Select **Settings History**.

## Max basal rate

If after consulting with your healthcare professional, you require a Manual mode Max basal rate higher than 10 units per hour, use **Manage Settings > Max Basal Rate**.

To set your Max basal rate, see *Max basal rate*, page 82.

### To increase the Max basal rate above 10 U/hr:

1. From the Home screen, press , and then select .
2. Select **Device Settings > Manage Settings**.
3. Simultaneously press and hold  and  until the Manage Settings screen appears.
4. Select **Max Basal Rate**.

A message appears indicating that changing this value will also change the Max basal rate setting in Delivery Settings. Consult a healthcare professional before changing this value.

Select **Continue** to change the Max basal rate setting or click **Cancel** to return to the previous screen.

5. Change the Max basal rate setting and select **Save**.



The Max basal rate setting applies to Manual mode only. SmartGuard insulin delivery is not affected by this value.

## Auto suspend

Auto suspend is a safety feature that stops all insulin delivery and sounds an alarm if a button is not pressed within a specified period of time. Consult a healthcare professional about how to best use this feature.

Auto suspend continues to work if the SmartGuard feature is active.



### To set up auto suspend:

1. From the Home screen, press , and then select .
2. Select **Device Settings** > **Auto Suspend**.
3. Select **Alarm**.
4. Select **Time** and enter the number of hours.
5. Select **Save**.

## Language

The language that the pump uses to show information can be updated after the startup.

### To change the language:

1. From the Home screen, press , and then select .  
A checkmark indicates which language is active.
2. Select **Device Settings** > **Language**.
3. Select a language.  
A screen asks to confirm.
4. Select **Yes** to accept. Select **No** to cancel.





# 10

History and graph



# 10

## History and graph

This chapter provides information about how to read historical data in the MiniMed 780G system.

### Introduction

The History screens provide details about personal therapy history in the MiniMed 780G insulin pump. The SG Review and Graph screens are available if the Sensor feature is turned on. The Time in Range screen shows the percent of time glucose levels are between 70 mg/dL and 180 mg/dL.

### History & Graph menu



The History & Graph menu provides information about insulin delivery, blood glucose (BG) meter readings, sensor glucose (SG) values, paired transmitter, and any alarms and alerts received.

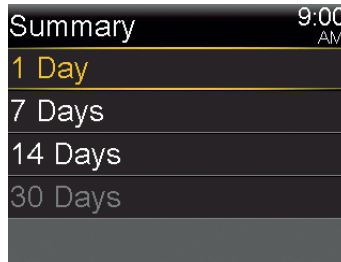
### History

#### Summary screen

The Summary screen displays information about past insulin deliveries, SG readings, and meter readings. Historical details can be viewed for a single day or for multiple days.

#### To view the Summary screen:

1. From the Home screen, press , and then select .
2. Select **History** > **Summary**.



3. Select the desired time period for the Summary screen.

The Summary screen appears and displays information for the number of days selected.



4. Scroll down to view the entire screen. In the **1 Day** view, use the < and > buttons on the pump to view the history of a specific day.

## Understanding the Summary screen

The Summary screen separates information into the following categories:

- Time in range information
- Insulin delivery overview
- Bolus Wizard
- Bolus in the SmartGuard feature
- BG
- Sensor
- Low management mode

## Summary screen: Time in SmartGuard and Time in range information

The following table describes the Time in SmartGuard, Time in Target Range, Time below range, and Time above range portions of the Summary screen.

Name	Description
Time in SmartGuard	number of hours / percent of time in the SmartGuard feature

Name	Description
Time in Target Range	number of hours / percent of time in target range (70 mg/dL to 180 mg/dL)
Time below range	number of hours / percent of time below target range (below 70 mg/dL)
Time above range	number of hours / percent of time above target range (above 180 mg/dL)

### Summary screen: insulin delivery overview

This summary screen appears when in Manual mode. If **1 Day** view is selected, the values are shown for that day. If multiple days are selected, the values shown are an average of the values for the selected number of days.

Name	Description
TDD	Total daily dose of insulin units.
Basal	<ul style="list-style-type: none"> <li>Insulin units devoted to basal delivery.</li> <li>Percentage of insulin devoted to basal delivery.</li> </ul>
Bolus	<ul style="list-style-type: none"> <li>Insulin units devoted to bolus delivery.</li> <li>Percentage of insulin devoted to bolus delivery.</li> </ul>
Total Carbs	Daily carbohydrate amount, in grams.

### Summary screen: Bolus Wizard

This summary screen appears when in Manual mode. If **1 Day** view is selected, the values are shown for that day. If multiple days are selected, the values shown are an average of the values for the selected number of days.

Name	Description
Carb bolus	<ul style="list-style-type: none"> <li>Total insulin units delivered using the Bolus Wizard feature in Manual mode with a food amount or with a food and glucose correction.</li> <li>Number of times the Bolus Wizard feature delivered a food bolus or a food plus correction bolus in Manual mode.</li> </ul>
Glucose correction only	<ul style="list-style-type: none"> <li>Total insulin units delivered using the Bolus Wizard feature in Manual mode or a bolus with BG correction amount only.</li> <li>Number of times the Bolus Wizard feature delivered a correction bolus in Manual mode.</li> </ul>

## Summary screen: SmartGuard

If **1 Day** view is selected, the values are shown for that day. If multiple days are selected, the values shown are an average of the values for the selected number of days.

Name	Description
Auto Correction	Total insulin units delivered by the Auto correction feature.
Bolus	<ul style="list-style-type: none"><li>Total insulin units delivered using the SmartGuard bolus feature.</li><li>Number of times the SmartGuard bolus feature was used.</li></ul>

## Summary screen: BG

The pump is only compatible with the Accu-Chek Guide Link meter.

Name	Description
BG	Total number of BG meter readings, including readings from an Accu-Chek Guide Link meter and BG meter readings entered manually.
Average BG	Average BG meter readings.
BG Std. Dev.	Standard deviation of BG meter readings.
Low BG	Lowest BG meter reading.
High BG	Highest BG meter reading.

## Summary screen: sensor

The sensor portion appears if a sensor has been used at least once.

Name	Description
SG Average	Average SG reading.
SG Std. Dev.	Standard deviation of the SG readings.

## Summary screen: low management mode

For information about the Suspend before low and Suspend on low features, see *Low SG settings, page 158*.



Name	Description
Suspend before low	The average number of Suspend before low events per day.
Suspend on low	The average number of Suspend on low events per day.
Time suspended by sensor	The average duration (amount of time) suspended as a result of Suspend before low or Suspend on low events per day.

## Daily History screen

Actions performed on the pump can be viewed on the Daily History screen for the selected day. The list shown on the screen provides further details and shows the most recent action first.

Daily History	9:00 AM
Temp Target Comp..	10:45 PM
Temp Target	10:40 PM
SmartGuard Active	10:35 PM
SmartGuard Exit	10:30 PM
◀ Thu, Jan 22 ▶	



### To view the Daily History screen:

1. From the Home screen, press , and then select .
2. Select **History** > **Daily History**.  
A list of dates appears.
3. Select a specific date. A list appears with any pump actions or events entered on the specified day.
4. Select any item in the list to open the Detail screen and view more information about the selected action or event.

## Alarm History screen

Select a specific day to view the history of alarms and alerts that occurred on the selected day. The list provides further details and shows the most recent alarm or alert first.

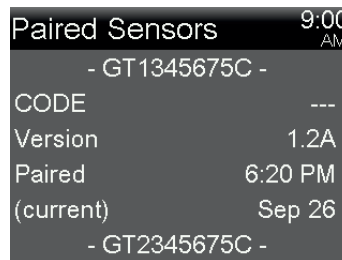
### To view the Alarm History screen:

1. From the Home screen, press , and then select .
2. Select **History** > **Alarm History**.  
A list of dates appears.



3. Select a specific date. A list appears showing any alarms or alerts that occurred on the specified day.
4. Select any alarm or alert in the list to open the Detail screen and view more information about the selected alarm or alert.

## Paired Sensors screen

The Paired Sensors screen displays the serial number, date, and time of the current sensor paired to the pump. The screen also provides a history of the sensors that were paired with and unpaired from the pump. The code is not applicable for the MiniMed 780G insulin pump with the Guardian 4 sensor.



### To view the Paired Sensors screen:

1. From the Home screen, press , and then select .
2. Select **History** > **Paired Sensors**.  
A list of transmitters appears.
3. Scroll down to view the entire screen.



## SG Review screen

Pair the pump with a sensor to view a graph of SG history based on high and low limits entered. Information can be viewed for 1 day or refer to an average of SG data over multiple days.

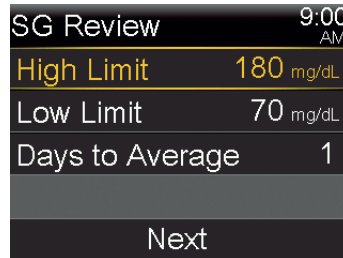
High and low limits set in the SG Review screen are only used to view SG data. These limits are not the same as the high and low glucose limits used for SG alerts. Changing the limits in the SG Review screen will not affect the high and low glucose limits used for the SG alerts.



### To review the SG history:

1. From the Home screen, press , and then select .
2. Select **Sensor Glucose Review**.

The SG Review screen appears. The high and low limits that appear are either the values entered for the last SG Review, or the default values of 180 mg/dL for the high limit and 70 mg/dL for the low limit.




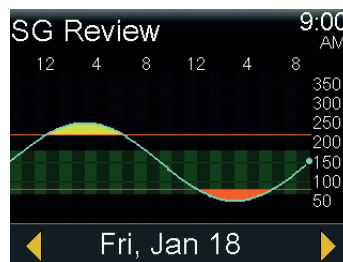
SG Review		9:00 AM
High Limit	180	mg/dL
Low Limit	70	mg/dL
Days to Average	1	
Next		

3. Enter the High Limit and Low Limit for the SG data review.

There must be a minimum of 20 mg/dL difference between the High Limit and the Low Limit.

4. Enter the number of days of SG history to average, and select **Next**.

If only one day is entered, the graph shows details about when the SG was above, below, or within the specified limits. Use the arrow keys to see the data for specific dates. Press  to see information about the time that SG was above, within, or below range. A message appears and states there is no data available if no data was saved.




If multiple days are entered, the pie chart shows the average percentage of time that the SG was above, below, or within the specific limits over an average of

multiple days. A message appears and states there is no data available if no data was saved.



## Graph screen

The graph shows information about the SG readings and trends, BG entries, auto correction bolus deliveries, and bolus entries. The below screen is an example of the graph screen using the SmartGuard feature.

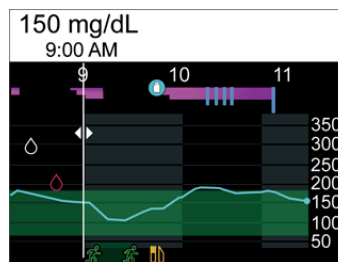



### To view the Graph screen:

- Press , or select **Graph** on the History & Graph screen.



## Time in Range screen

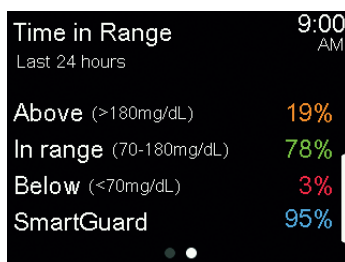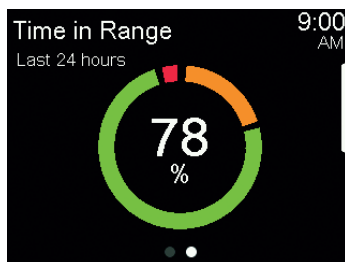Time in range is the percentage of time SG is between 70 mg/dL and 180 mg/dL. These values cannot be changed. Use the Time in Range screen to see how much time is spent below, above, and within range in the last 24 hours.

When using CGM, the following information can be viewed:



**To view the Time in Range screen:**

1. From the Home screen, press , and then select .
2. Select **Time in Range**.



11

Notifications and reminders



# 11

## Notifications and reminders

This chapter describes how to use reminders. It also covers the general behavior of the most common and the most serious notifications and how to resolve them.

### Notifications in the MiniMed Mobile app

If the MiniMed Mobile app is used, alarms, alerts, and messages can be viewed on the paired mobile device. For information about how to set the notification preferences in the app, see the MiniMed Mobile app user guide. For a table that describes the meaning, consequences, reasons, and resolutions for the most common or serious notifications, see *Pump alarms, alerts, and messages*, page 291.



**WARNING:** Do not rely on the MiniMed Mobile app to view all alerts. Alerts will not appear on the MiniMed Mobile app during reservoir set up. Some alerts may only appear on the pump. In some cases, alerts could be sent to the MiniMed Mobile app after they appear on the pump. Relying on the MiniMed Mobile app for all alerts could result in an alert being missed, which may lead to hypoglycemia or hyperglycemia.



### Reminders

There are several specific reminders that prompt a specific action. Personal reminders can be used for any purpose.

## Personal reminders

Up to 5 personal reminders can be set, along with the specific reminders for blood glucose (BG) meter readings and medication.

### To create a new Personal reminder:



1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Reminders** > **Personal**.
3. Select **Add New**.

The Select Name screen shows the available reminders.

4. Select a reminder.  
An edit screen appears for the selected reminder.
5. Enter the time the reminder should occur.
6. Select **Save**.

The Personal reminder occurs at the specified time each day unless it is edited or deleted.

### To edit, rename, or delete an existing Personal reminder:



1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Reminders** > **Personal**.
3. Select a reminder.
4. Do any of the following:
  - Select **Reminder** to turn the reminder on or off.
  - Select **Edit** to change the time of the reminder.
  - Select **Rename** to assign a different name to the reminder. When the Select Name screen appears, select any available name from the list.
  - Select **Delete** to delete the reminder.



## Bolus BG Check reminder

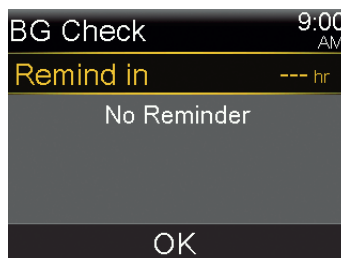
The Bolus BG Check reminder notifies when BG needs to be checked after a bolus delivery. After a bolus is started, the BG Check screen appears and the timer must be set for the reminder. The timer counts down from the time the bolus was started.

### To turn on or turn off Bolus BG Check reminders:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Reminders** > **Bolus BG Check**.
3. To turn the reminder on or off, select **Reminder**.
4. Select **Save**.

### To use a Bolus BG Check reminder if a bolus is being delivered:

1. If the Bolus BG Check reminder is on, the BG Check screen appears each time a bolus is started.





2. Enter a time between 30 minutes and 5 hours and select **OK**. If no reminder is necessary after the bolus delivery, select the dashes without adding a time, and select **OK**.



## Missed Meal Bolus reminder

Missed Meal Bolus reminders can be set up around typical meal times. Up to 8 reminders can be set.

### To create a new Missed Meal Bolus reminder:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Reminders** > **Missed Meal Bolus**.
3. Select **Add New**.
4. Select **Start Time** and enter a time.
5. Select **End Time** and enter a time.
6. Select **Save**.

### To turn on or off, edit, or delete existing Missed Meal Bolus reminders:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Reminders** > **Missed Meal Bolus**.
3. Select a reminder.
4. Change any of the following:
  - Select **Reminder** to turn this reminder on or off.
  - Select **Edit** to change the time of this reminder.
  - Select **Delete** to delete this reminder.

## Low Reservoir reminder

Set a Low Reservoir reminder to occur when the insulin level in the reservoir reaches a specified number of units and again when half of those units have been used.





**Note:** The number of units that remain in the reservoir can be found on the Pump status screen. For more information, see *Status screen*, page 76.



**WARNING:** Always check the amount of insulin left in the reservoir when the Low reservoir alert occurs. Confirm that the MiniMed 780G insulin pump has sufficient insulin. The insulin level in the reservoir can reach a low level during a bolus delivery or fill cannula delivery. If this occurs, the Low reservoir alert displays. If the pump does not have sufficient insulin, under-delivery of insulin can occur, which may cause hyperglycemia.



### To set up the Low Reservoir reminder:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Reminders** > **Low Reservoir**.
3. Select **Units** to enter the number of units. Set a value from 5 to 50 units.
4. Select **Save**.

## Set Change reminder

The Set Change reminder tracks the time between infusion set changes and provides a reminder to change the infusion set.



### To turn on or off, or edit the Set Change reminder:

1. From the Home screen, press , and then select .
2. Select **Alert Settings** > **Reminders** > **Set Change**.
3. Select **Reminder** to turn the reminder on or off.
4. Select **Time** and choose the number of days needed for the reminder.
5. Select **Save**.



**WARNING:** When changing the Set Change reminder, do not set a duration greater than what is indicated on the infusion set labeling. If the infusion set is labeled for three days then the reminder must only be set to two or three days.

## Alarms, alerts, and messages

The pump has a sophisticated safety network. If this safety network detects anything unusual, it communicates this information in the form of notifications. Notifications include alarms, alerts, and messages. When more than one notification is received, and there are multiple messages to view, a small white flap appears on the notification icon . When the first notification is cleared, the next notification becomes visible. A white triangle means that  must be pressed to continue.



**Note:** The notification light flashes when the pump has an alarm or alert.



**Note:** Promptly address all notifications and confirmations that appear on the pump screen. The notification will remain on the pump screen until it is cleared. When responding to a message, there may be times when another message appears.



**WARNING:** When the critical pump error occurs, the following screen appears and the pump siren goes off:

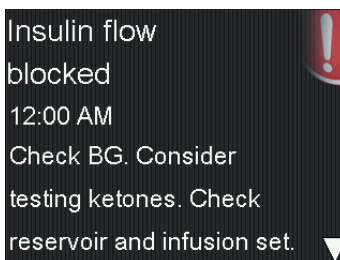


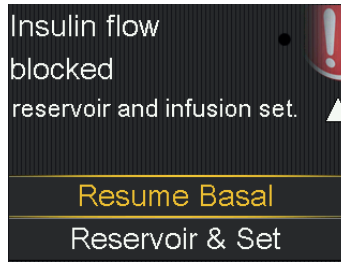
Immediately disconnect the pump and discontinue use. Contact 24-Hour Technical Support.

Insulin delivery is still required when the pump is removed. Consult a healthcare professional to determine an alternate method of insulin delivery while the pump is removed.

## Alarms

An alarm warns of a condition that requires immediate attention. Stopped insulin delivery and low glucose levels are the most common reasons for alarms.






**WARNING:** Always address alarms immediately when they occur. Ignoring an alarm can result in hyperglycemia or hypoglycemia.

When an alarm occurs:


**Display:** The pump displays a notification with a red icon and instructions.

**Notification light:** The red notification light blinks twice, followed by a pause, in a continuous repeating pattern.

**Audio:** Depending on the sound and vibration settings, the pump emits a series of alarm tones and vibrations.

The underlying problem that triggered the alarm must be resolved. In most cases, press  and then make a selection to clear the alarm. Sometimes the underlying problem is not resolved when the alarm is cleared. The alarm repeats until the underlying problem is fixed. If the alarm condition is not resolved after 10 minutes, the alarm tone escalates to a loud emergency siren.

## Alerts

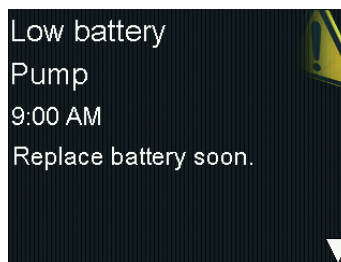


Alerts indicate that a situation may require attention. When an alert occurs, check the pump screen to see if any action is required.

When an alert occurs:

**Display:** The pump displays a notification with a yellow icon and instructions.

**Notification light:** The red notification light on the pump blinks once, followed by a pause, then blinks once again in a continuous repeating pattern.


**Audio:** Depending on the sound and vibration settings, the pump generates a series of tones and vibrations.

To clear an alert, press  and then make a selection. The pump beeps every 5 minutes or every 15 minutes, depending on the alert, until the alert is resolved. Some alerts will also escalate to a loud emergency siren after 10 minutes.



**Note:** If an alert occurs when the pump is on a screen other than the Home screen, the alert message may only appear after the pump returns to the Home screen.

## Messages

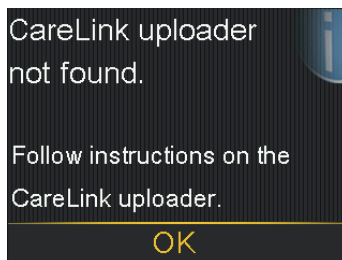

A message is a notification that shows the status of the pump or displays when a decision needs to be made.

When a message occurs:

**Display:** The pump displays a notification with a blue icon and instructions. Some messages show a yellow icon.

**Notification light:** The red notification light on the pump does not blink.

**Audio:** Depending on the sound and vibration settings, the pump emits a tone, a one-pulse-only vibration, or it emits a tone and a one-pulse-only vibration. To clear a message, press ▼ and then make a selection.

### **Pump alarms, alerts, and messages**

For a table that describes the meaning, consequences, reasons, and resolutions for the most common or serious notifications, see *Pump alarms, alerts, and messages, page 291*.



# 12

Additional basal features



# 12



## Additional basal features

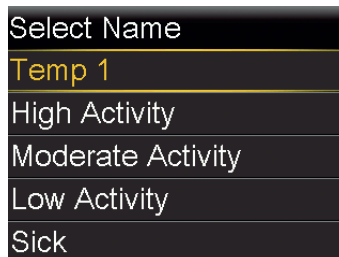
This chapter provides information about setting up additional features for basal insulin delivery.

### Preset temp basal rates

Set up preset temp basal rates for reoccurring short-term situations. Up to four preset temp basal rates can be set up for specific situations. There are also four additional preset temp rates available for use in other circumstances (Temp 1 through Temp 4).

#### To set up a preset temp basal rate:



1. From the Home screen, press , and then select .
2. Select **Delivery Settings** > **Preset Temp Setup**.
3. Select **Add New**.



4. Select a name for the preset temp basal rate.
5. Select **Type** to select Percent or Rate, and then enter the percentage or the rate in units per hour.

6. Set the **Duration** for the preset temp basal rate to be active.
7. Select **Save**.

**To edit, rename, or delete a preset temp basal rate:**

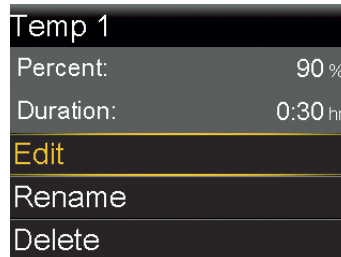
1. From the Home screen, press , and then select .
2. Select **Delivery Settings > Preset Temp Setup**.

The Preset Temp Setup screen appears and shows the settings for any existing preset temp basal rate.



3. Select a preset temp basal rate.

A screen appears that shows the preset temp basal rate information.



4. Do any of the following:
  - Select **Edit** to adjust the type (Percent or Rate), the percent or rate amount, and the duration.
  - Select **Rename** to assign a different name to the preset temp basal rate. When the Select Name screen appears, select any available name from the list.
  - Select **Delete** to delete the preset temp basal rate.

## Starting a preset temp basal delivery

Follow the steps to use the preset temp basal rate for basal insulin delivery. If a preset temp basal rate has not yet been set up, see *Preset temp basal rates*, page 239. After the preset temp basal delivery is completed or canceled, basal insulin delivery resumes using the programmed basal rate.

### To start a preset temp basal delivery:

1. From the Home screen, press , and then select .
2. Select **Basal > Preset Temp**.

The Preset Temp screen appears and shows the preset temp basal rates set up, along with their percentage or rate amounts.

Preset Temp	9:00 AM
Current rate:	0.025 U/hr
Temp 1	0.100 U/hr
High Activity	25 %
Moderate...	50 %



**Note:** If a percentage preset temp basal rate was set up so that it could exceed the current Max basal limit, that preset temp is grayed out in the list and cannot be selected.

3. Select a preset temp basal rate to start.
4. Select **Begin**.



The Temp Basal banner appears on the Home screen during delivery.

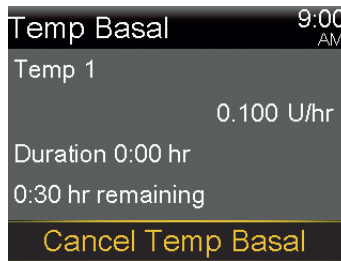


## Canceling a temp basal or preset temp basal

A temp basal rate or preset temp basal rate can be canceled at any time. After it is canceled, the scheduled basal pattern automatically resumes.

### To cancel a temp basal rate:

1. From the Home screen, press , and then select .
2. Select **Cancel Temp Basal**.  
The Temp Basal screen appears.



3. Select **Cancel Temp Basal**.

## Additional basal patterns



### Adding an additional basal pattern

This procedure shows how to add a new basal pattern after at least one basal pattern has been set. If this is the first time a basal pattern is being set, see *Setting up a basal pattern*, page 85.

The following basal patterns can be set up:



- Basal 1
- Basal 2
- Workday
- Day Off
- Sick Day

### To add an additional basal pattern:

1. From the Home screen, press , and then select .
2. Select **Basal > Basal Pattern Setup**.  
The Basal Pattern Setup screen appears.
3. To add a new basal pattern, select **Add New**.  
The Select Name screen appears.
4. Select a name for the basal pattern.
5. Set the basal rate.
6. Select **Review**.
7. Select **Save**.

### Editing, copying, or deleting a basal pattern

#### To edit, copy, or delete a basal pattern:

1. From the Home screen, press , and then select .
2. Select **Delivery Settings > Basal Pattern Setup**.  
The Basal Pattern Setup screen appears



3. Select a basal pattern.
4. Select **Options**.
5. Do any of the following:
  - Select **Edit** to adjust the end time or rate values.
  - Select **Copy** to copy the basal rate information from the selected basal pattern to a new basal pattern. When the Select Name screen appears, select any available name from the list.
  - Select **Delete** to delete the selected basal pattern. The active basal pattern cannot be deleted.

## Changing from one basal pattern to another

If more than one basal pattern has been set, the basal pattern can be changed. The MiniMed 780G insulin pump delivers basal insulin according to the selected basal pattern.

### To change to a different basal pattern:

1. From the Home screen, press , and then select .
2. Select **Basal > Basal Patterns**.

The Basal Patterns screen appears. A check mark displays next to the active basal pattern.

Basal Patterns		9:00 AM
Basal 1	1.125 U	✓
Basal 2	1.2 U	

3. Select a basal pattern.



Basal 2		9:00 AM
24 hr Total: 1.2 U		
Start	End	Rate (U/hr)
12:00 <sub>A</sub>	12:00 <sub>A</sub>	0.050
Begin		

4. Select **Begin**.



13



Additional bonus features



# 13

## Additional bolus features

This chapter provides information about additional features for bolus delivery. Square Wave, Dual Wave, Easy, Manual, and Preset bolus are only available in Manual mode. Since these bolus types are only available in Manual mode, remember that you must enter a blood glucose (BG) meter reading when setting up the bolus delivery. Do not use a sensor glucose (SG) value when delivering a bolus in Manual mode.

### Bolus types

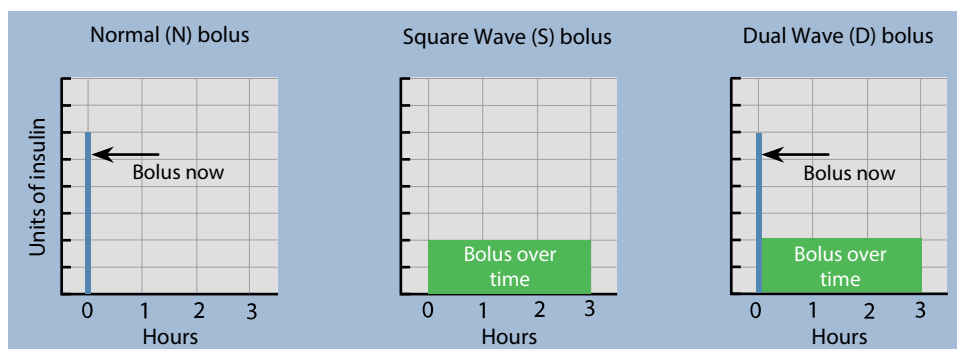
The following table provides general information about the available bolus types.

Bolus type	Description	Purpose
Normal	Normal bolus provides a single immediate dose of insulin.	<p>This is the typical bolus type used to cover food intake or to correct a high blood glucose (BG) meter reading.</p> <p>For details about delivering a normal bolus, see <i>Delivering a Normal bolus</i>, page 105.</p>
Square Wave bolus	Square Wave bolus delivers a single bolus evenly over an extended period of time from 30 minutes up to 8 hours.	<p>A Square Wave bolus can be used for the following reasons:</p> <ul style="list-style-type: none"> <li>• A delayed food digestion due to gastroparesis or meals high in fat.</li> <li>• Snacking over an extended period of time.</li> <li>• A normal bolus drops the BG too rapidly.</li> </ul> <p>For details about using the Square Wave bolus feature, see <i>Square Wave bolus</i>, page 254.</p>

Bolus type	Description	Purpose
Dual Wave bolus	Dual Wave bolus delivers a combination of an immediate normal bolus followed by a Square Wave bolus.	<p>A Dual Wave bolus can be used for the following reasons:</p> <ul style="list-style-type: none"> <li>When meals are high in carbs and fat, which may delay digestion.</li> <li>When a meal bolus is combined with a correction bolus for an elevated BG.</li> </ul> <p>For details about using a Dual Wave bolus, see <i>Dual Wave bolus</i>, page 258.</p>

## Bolus type example

The following example shows how the different bolus types work.





## Bolus settings

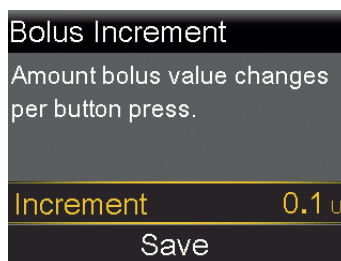
Additional settings are required to use the Bolus Wizard feature in Manual mode. These are described in the section, *Manual mode bolus delivery options*, page 97.

## Bolus increment

The Bolus increment is the number of units that are increased or decreased with each button press for the bolus delivery amount in the Bolus Wizard, Manual Bolus, and Preset Bolus screens. Depending on the typical bolus amount, the increment can be set to 0.1 units, 0.05 units, or 0.025 units.

### To set the bolus increment:

1. From the Home screen, press , and then select .
2. In Manual mode, select **Delivery Settings** > **Bolus Increment**.
3. Select **Increment** to set the desired increment value.






4. Select **Save**.

### Bolus speed

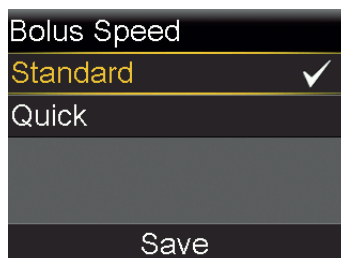The bolus speed sets the rate at which the pump delivers bolus insulin. Set a standard rate (1.5 units per minute), or a quick rate (15 units per minute).

### To set the bolus speed:

1. From the Home screen, press , and then select .
2. Select **Delivery Settings** > **Bolus Speed**.
3. Select **Standard** or **Quick**.



4. Select **Save**.




## Changing the Bolus Wizard settings in Manual mode

This section shows how to make changes to personal settings after the initial Bolus Wizard feature setup. Consult a healthcare professional before changes are made to the personal settings.

### Changing the carb ratio

The carb ratio can be set whether or not the Bolus Wizard feature is turned on.

#### To change the carb ratio:

1. From the Home screen, press , and then select .
2. In Manual mode, select **Delivery Settings > Bolus Wizard Setup > Carb Ratio**.
3. Select **Edit**.
4. Select the carb ratio. For one carb ratio, enter the g/U, and then press .

For more than one carb ratio, enter one carb ratio at a time to complete the full 24 hours, which ends at 12:00 A.





**Note:** For instructions on setting up more than one carb ratio over a 24-hour period, see *Settings covering a 24-hour period*, page 87.

5. Select **Save**.

### Changing the insulin sensitivity factor

The insulin sensitivity factor can be set only if the Bolus Wizard feature is turned on.

#### To change the insulin sensitivity factor:

1. From the Home screen, press , and then select .
2. In Manual mode, select **Delivery Settings > Bolus Wizard Setup > Insulin Sensitivity Factor**.
3. Select **Edit**.



4. Select the insulin sensitivity factor. For one insulin sensitivity factor, press  $\wedge$  and  $\vee$  to enter the mg/dL per U, and then press  $\odot$ .

For more than one insulin sensitivity factor, press  $\wedge$  or  $\vee$  to enter one insulin sensitivity factor at a time to complete the full 24 hours, which ends at 12:00 A.



**Note:** For instructions on setting up more than one insulin sensitivity factor over a 24-hour period, see *Settings covering a 24-hour period*, page 87.

5. Select **Save**.

## Changing the BG target

The BG target can be from 60 to 250 mg/dL. The BG target can be set only if the Bolus Wizard feature is turned on.

### To change the BG target:

1. From the Home screen, press  $\odot$ , and then select  $\text{⚙}$ .
2. In Manual mode, select **Delivery Settings > Bolus Wizard Setup > BG Target**.
3. Select **Edit**.
4. Select the BG target. For one BG target, enter the low BG limit and the high BG limit, and then press  $\odot$ .

For more than one BG target, enter one BG target at a time to complete the full 24 hours, which ends at 12:00 A.



**Note:** For instructions on setting up more than one BG target over a 24-hour period, see *Settings covering a 24-hour period*, page 87.

5. Select **Save**.



## Changing the active insulin time

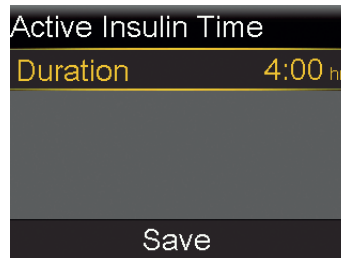
Active insulin is the bolus insulin that has been delivered by the pump and is still working to lower glucose levels. In the Bolus Wizard and SmartGuard Bolus feature, the

Active Insulin Time setting is used to calculate a correction bolus by subtracting the estimated active insulin from each bolus. In SmartGuard, auto correction boluses are delivered up to every 5 minutes. A shorter Active Insulin Time setting may result in more insulin being delivered in correction boluses.

A healthcare professional provides the personalized active insulin time based on historic glycemic control data for the individual user. When using SmartGuard, the recommended initial setting is an Active Insulin Time of 2-3 hours. The Active Insulin Time setting in the MiniMed 780G system is not necessarily reflective of the physiological insulin metabolism. Adjustments are not based on the pharmacokinetics and pharmacodynamics of the rapid-acting insulin. The current active insulin amount appears on the Home screen and includes only the bolus insulin received.

### To change the active insulin time:

1. From the Home screen, press , and then select .
2. In Manual mode, select **Delivery Settings** > **Bolus Wizard Setup** > **Active Insulin Time**.
3. Select **Duration**, and adjust the active insulin time in hours, using 15-minute increments.



4. Select **Save**.

## Square Wave bolus

A Square Wave bolus delivers a bolus evenly over a period of time from 30 minutes up to 8 hours.

When using the Bolus Wizard feature in Manual mode, a Square Wave bolus is available only when giving a food bolus without a correction for an elevated BG. A Square Wave

bolus is not available for a correction bolus alone or a correction bolus with food bolus. A normal bolus can be delivered while a Square Wave bolus is being delivered, as needed.

A Square Wave bolus can be useful in the following situations:



- Delayed food digestion due to gastroparesis or meals high in fat.
- When snacking over an extended period of time.
- A normal bolus drops BG too rapidly.

Since the Square Wave bolus extends delivery over a period of time, the insulin is more likely to be available as needed.

### Turning the Square Wave bolus feature on or off

A Square Wave bolus can be set up and delivered only after the Square Wave bolus feature is turned on.



#### To turn the Square Wave bolus feature on or off:

1. From the Home screen, press , and then select .
2. In Manual mode, select **Delivery Settings** > **Dual/Square Wave**.
3. Select **Square Wave** to turn the feature on or off.
4. Select **Save**.

### Delivering a Square Wave bolus using the Bolus Wizard feature

In Manual mode, the Bolus Wizard feature only delivers a Square Wave bolus if the Square Wave bolus feature is turned on and a carb value is entered. If a BG reading causes the Bolus Wizard feature to calculate that a correction bolus is necessary, then a Square Wave bolus cannot be delivered.

#### To deliver a Square Wave bolus using the Bolus Wizard feature:

1. From the Home screen, press , and then select .
2. In Manual mode, select **Bolus** > **Bolus Wizard**.

The Bolus Wizard screen appears.

Bolus Wizard		9:00 AM
BG	--- mg/dL	
Carbs	10 <sub>g</sub>	0.6 <sub>U</sub>
Adjustment	0.0 <sub>U</sub>	
Bolus	0.6 <sub>U</sub>	
Next		

- For a food bolus, select **Carbs** to enter the carb count of the meal.
- The calculated bolus appears in the Bolus field. To modify the bolus amount, select **Bolus**.
- Select **Next** to review the bolus information.

Bolus Wizard		9:00 AM
Bolus		0.8 <sub>U</sub>
Dual		Square
Deliver Bolus		

- Select **Square**.
- Select **Duration** to adjust the time period when the Square Wave bolus needs to be delivered.

Bolus Wizard		9:00 AM
Bolus		2.3 <sub>U</sub>
Duration		0:30 <sub>hr</sub>
Deliver Bolus		

- Select **Deliver Bolus** to start the bolus.





**Note:** To stop bolus delivery or to see details on the insulin that has been delivered, see *Stopping a Square Wave or Dual Wave bolus delivery*, page 268.

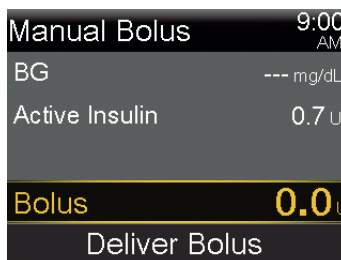
### Delivering a Square Wave bolus using the Manual bolus feature

The Square Wave bolus option is available in the Manual Bolus screen only after the Square Wave feature is turned on.


#### To deliver a Square Wave bolus using the Manual bolus feature:

1. From the Home screen, press , and then select .
2. In Manual mode, do one of the following:
  - Select **Bolus** if the Bolus Wizard feature is turned off.
  - Select **Bolus** > **Manual Bolus** if the Bolus Wizard feature is turned on.

The Manual Bolus screen appears.



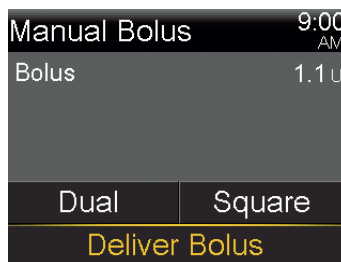3. Set the bolus delivery amount in units, and then select **Next**.



4. Select **Square**.
5. Select **Duration** to adjust the time period when the Square Wave bolus is to be delivered.
6. Select **Deliver Bolus** to start the bolus.



**Note:** To stop bolus delivery or to see details on the insulin that has been delivered, see *Stopping a Square Wave or Dual Wave bolus delivery*, page 268.

## Dual Wave bolus

The Dual Wave bolus feature meets both immediate and extended insulin needs by delivering a combination of an immediate normal bolus followed by a Square Wave bolus. A normal bolus can be delivered while the Square portion of a Dual Wave bolus is being delivered, as needed.



A Dual Wave bolus can be useful in these situations:

- When an elevated BG needs to be corrected before a meal, and a delayed bolus is needed for food that is absorbed slowly.
- When eating meals with mixed nutrients, such as carbs, fats and proteins, that are absorbed at different rates.

## Turning the Dual Wave bolus feature on or off

A Dual Wave bolus can be delivered only after the Dual Wave bolus feature is turned on.



### To turn the Dual Wave feature on or off:

1. From the Home screen, press , and then select .
2. In Manual mode, select **Delivery Settings** > **Dual/Square Wave**.
3. Select **Dual Wave** to turn the feature on or off.
4. Select **Save**.

## Delivering a Dual Wave bolus using the Bolus Wizard feature

In Manual mode, a Dual Wave bolus with the Bolus Wizard feature can be delivered only after the Dual Wave bolus feature is turned on.

### To deliver a Dual Wave bolus with the Bolus Wizard feature:

1. For a correction bolus or a food bolus with a correction, use a BG meter to check BG. For a food bolus only, go to step 2.
2. From the Home screen, press , and then select .
3. In Manual mode, select **Bolus** > **Bolus Wizard**.

The Bolus Wizard screen appears.





**Note:** For more information on how to manually enter the BG meter reading, see *Entering a blood glucose (BG) meter reading*, page 96.

4. For a food bolus, select **Carbs** to enter the carb count of the meal. For a correction bolus where no food was eaten, leave the carbs value as 0.

The calculated bolus appears in the Bolus field.

5. To modify the bolus amount, select **Bolus**.
6. Select **Next** to review the bolus information.

Bolus Wizard		9:00 AM
Bolus	0.8 U	
Dual	Square	
Deliver Bolus		

7. Select **Dual**.

The Bolus Wizard screen appears.

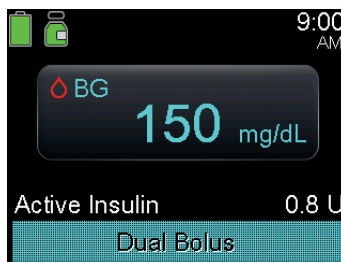
8. To change the amounts, select the area of the screen with the Now % and Square % values and adjust the **Now** % amount.

When adjusting the Now amount, the Square amount adjusts automatically.

Bolus Wizard		9:00 AM
Bolus	0.8 U	
Now	75 %	0.6 U
Square	25 %	0.2 U
Duration	0:30 hr	
Deliver Bolus		

9. Adjust the **Duration** of the square portion of the bolus to be delivered.
10. Select **Deliver Bolus** to start the bolus.







**Note:** To stop bolus delivery or to see details on the insulin that has been delivered, see *Stopping a Square Wave or Dual Wave bolus delivery*, page 268.

## Delivering a Dual Wave bolus using the Manual bolus feature

The Dual Wave bolus option is available in the Manual Bolus screen only after the Dual Wave feature is turned on.

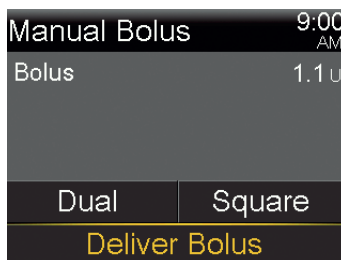
### To deliver a Dual Wave bolus using the Manual bolus feature:

1. From the Home screen, press , and then select .
2. In Manual mode, do one of the following:
  - Select **Bolus** if the Bolus Wizard feature is turned off.
  - Select **Bolus** > **Manual Bolus** if the Bolus Wizard feature is turned on.

The Manual Bolus screen appears.

3. Set the bolus delivery amount in units, and then select **Next**.

The Manual Bolus screen appears, with the option to select the bolus type.



4. Select **Dual**.

The Manual Bolus screen appears.

- To change the amounts, select the area of the screen with the Now % and Square % values and adjust the **Now** % value. When the Now amount is adjusted, the Square amount adjusts automatically.

Manual Bolus		9:00 AM
Bolus		0.8 U
Now	50 %	0.4 U
Square	50 %	0.4 U
Duration		0:30 hr
Deliver Bolus		

- Select **Duration** to adjust the time period when the Square Wave bolus is to be delivered.
- Select **Deliver Bolus** to start the bolus.

		9:00 AM
BG		150 mg/dL
Active Insulin		0.4 U
Dual Bolus		



**Note:** To stop bolus delivery or to see details on the insulin that has been delivered, see *Stopping a Square Wave or Dual Wave bolus delivery*, page 268.

## Easy bolus

The Easy bolus feature can be used to deliver a normal bolus using only the **^** button. The Easy bolus feature only works when the pump is in Sleep mode.

When the **^** button is pressed while the Easy bolus feature is used, the bolus amount increases by a certain amount. This amount, or step size, can be set from 0.1 to 2.0 units


of insulin. The pump makes a tone or vibration each time the  $\wedge$  button is pressed to help keep count of the steps.



**Note:** The step size cannot be greater than the Max bolus amount. The maximum number of steps is 20 for each bolus delivery.

## Setting up the Easy bolus feature

The following graph provides an example of setting up a bolus of 2.0 units of insulin using a step size of 0.5 units.

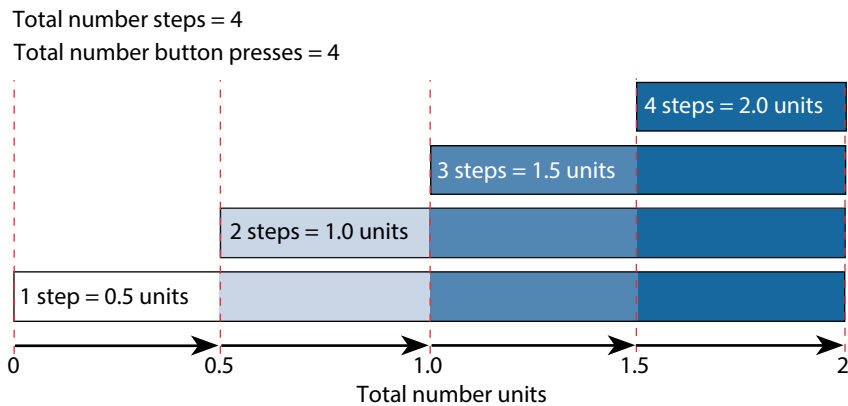

### To set up the Easy bolus feature:

1. From the Home screen, press  $\odot$ , and then select  $\text{⚙}$ .
2. Select **Device Settings** > **Easy Bolus**.
3. Select **Easy Bolus** to turn on the feature.
4. Set the **Step Size** amount in units.

Select a step size to a number that makes it easy to calculate the total bolus amount.



5. Select **Save**.

## Delivering a bolus using the Easy bolus feature



**WARNING:** Never rely on beeps or vibrations alone while using the Easy bolus feature. Always confirm the insulin delivery by looking at the pump screen. When using the Sound & Vibration options, it is possible that a sound or vibration notification may not occur as expected if the speaker or vibrator in the pump malfunctions. Relying on beeps or vibrations while using the Easy bolus feature may result in over-delivery of insulin.

### To deliver a bolus using the Easy bolus feature:

1. While the pump is in Sleep mode, press and hold  $\wedge$  for one second or until the pump beeps or vibrates. The bolus can now be set up.



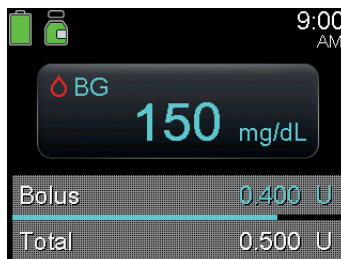
**Note:** If the pump does not respond when  $\wedge$  is pressed, it may not be in Sleep mode, even if the screen is dark. For more information, see *Sleep mode*, page 62.

2. Press  $\wedge$  the number of times needed to set the bolus amount. Count the tones or vibrations for each button press to confirm the total bolus amount.



**Note:** If  $\wedge$  is pressed too many times and the bolus amount is too high, press  $\vee$  to cancel the Easy bolus delivery and start at step 1 to set up a new bolus.

3. When the needed bolus amount is reached, press and hold  $\wedge$  to confirm the amount.
4. Press and hold  $\wedge$  for one second, or until the pump beeps or vibrates, to deliver the bolus.



**Note:** If the  $\wedge$  button is not pressed within 10 seconds after the bolus amount is confirmed, the bolus is canceled and a message appears that the bolus was not delivered.

## Preset bolus

The Preset Bolus feature allows frequently used bolus deliveries to be set up in advance. There are four preset bolus names that can be used to match a bolus to a meal that has a known carb content. Four additional preset bolus names can be set for other circumstances. These are numbered from Bolus 1 to Bolus 4.



**Note:** To set up a Preset bolus as a Dual Wave bolus or Square Wave bolus, the Dual Wave bolus feature or Square Wave bolus feature must be turned on.

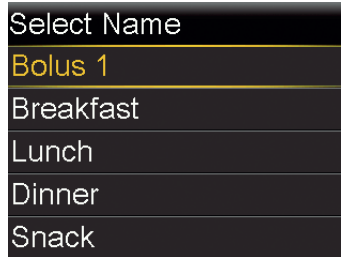
## Setting up and managing preset bolus deliveries

### To set up preset bolus amounts:


1. From the Home screen, press  $\odot$ , and then select  $\text{⚙}$ .
2. Select **Delivery Settings** > **Preset Bolus Setup**.



3. Select **Add New**.



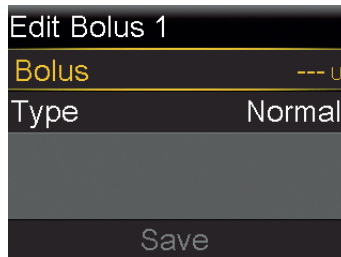4. Select a preset bolus.  
An edit screen appears.



5. Select **Bolus** to set the bolus amount.
6. Select **Type** to set this as a normal bolus, Square Wave bolus, or Dual Wave bolus.



**Note:** Square Wave and Dual Wave can be selected in the **Type** field only if the Square Wave bolus and Dual Wave bolus features are turned on.

If the type is set to Square or Dual, do the following:

- For a Square Wave bolus, set the **Duration** of time for the bolus delivery.
- For a Dual Wave bolus, adjust the **Now** % amount. When the **Now** amount is adjusted, the **Square** amount adjusts automatically. Then set the **Duration** of time for the Square portion of the bolus.



**Note:** If the Dual Wave bolus feature or Square Wave bolus feature is turned off, the existing Preset Bolus settings are still available for use.

7. Select **Save**.



## Editing, renaming, or deleting a preset bolus

Dual Wave Preset Boluses and Square Wave Preset Boluses can only be edited when the Dual Wave Bolus and Square Wave Bolus features are turned on.



**Note:** A preset bolus cannot be edited, renamed, or deleted during preset bolus delivery.


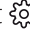
### To edit, rename, or delete a preset bolus:


1. From the Home screen, press , and then select .
2. Select **Delivery Settings > Preset Bolus Setup**.
3. Select a preset bolus.
4. Select **Options**.
5. Do any of the following:
  - Select **Edit** to adjust the bolus value and type, if applicable. If changing to a Square Wave bolus, enter the duration. If changing to a Dual Wave bolus, enter the Now and Square values and the Duration.
  - Select **Rename** to assign a different name to this preset bolus. When the Select Name screen appears, select any available name from the list.
  - Select **Delete** to delete this preset bolus.

## Delivering a preset bolus

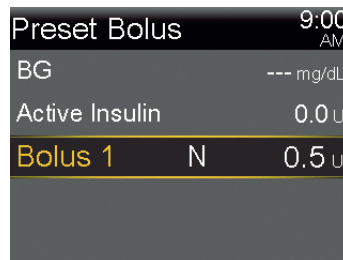
A preset bolus must be set before the Preset Bolus feature can be used. For more information, see *Setting up and managing preset bolus deliveries*, page 265.

### To deliver a preset bolus:



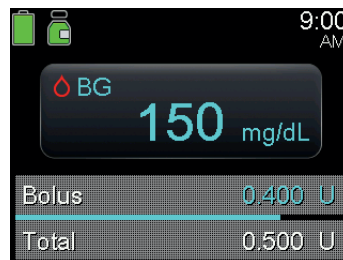1. From the Home screen, press , and then select .
2. Select **Bolus > Preset Bolus**.
3. Select the preset bolus to be delivered.



4. Review the bolus amount, and then select **Deliver Bolus** to start the bolus.



## Stopping a Square Wave or Dual Wave bolus delivery

This section describes how to stop a bolus in progress. It does not stop basal insulin delivery. To stop all insulin delivery, use the Suspend All Delivery feature (press , select , and select **Suspend All Delivery**).

This section describes how to stop the following bolus deliveries:




- A Dual Wave bolus during the Now portion delivery
- A Square Wave bolus delivery or a Dual Wave bolus during the Square portion delivery

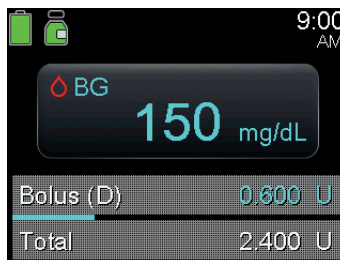
To stop a normal bolus delivery see *Stopping a bolus delivery*, page 108.




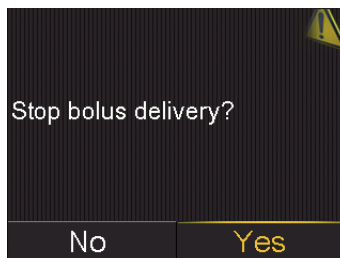
**Note:** When delivering a normal bolus and a Square Wave bolus at the same time, or a normal bolus and the Square portion of a Dual Wave bolus at the same time, both boluses are stopped.

### To stop a Dual Wave bolus delivery during the Now portion:

1. While the pump is delivering the Now portion of a Dual Wave bolus, press  from the Home screen.



2. Select .
3. Select **Stop Bolus**, then select **Yes** to confirm.



The Bolus Stopped screen appears and shows the amount of bolus delivered, and bolus amount that was originally set up.





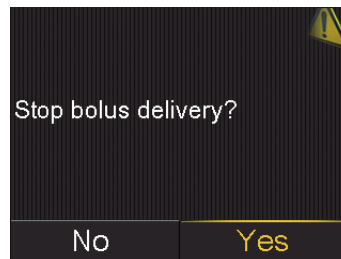
**Note:** When a Dual Wave bolus is stopped during the Now portion, the Now portion is stopped and the Square portion is canceled.



4. Select **Done**.


**To stop a Square Wave bolus delivery or the Square portion of a Dual Wave bolus delivery:**

1. While the pump is delivering a Square Wave bolus delivery or the Square portion of a Dual Wave bolus delivery, press  from the Home screen.
2. Select  and then select **Bolus**.
3. Select **Stop Bolus**, then select **Yes** to confirm.



The Bolus Stopped screen appears and shows the amount of bolus delivered, and the bolus amount that was originally set up.

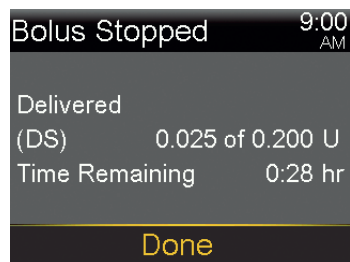4. Select **Done**.





# 14

Troubleshooting



# 14

## Troubleshooting

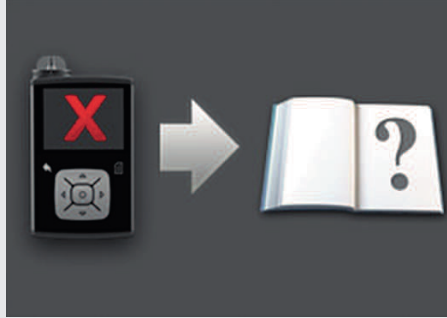
This chapter provides information about common MiniMed 780G insulin pump and sensor issues, as well as possible resolutions.

For a list of alarms, alerts, and messages, see *List of alarms, alerts, and messages*, page 291.

## Pump issues




**WARNING:** When the critical pump error occurs, the following screen appears and the pump siren goes off:



Immediately disconnect the pump and discontinue use. Contact 24-Hour Technical Support.




Insulin delivery is still required when the pump is removed. Consult a healthcare professional to determine an alternate method of insulin delivery while the pump is removed.

The following table provides troubleshooting information for the insulin pump:

Issue	Resolution
 appears on the Home screen or Bolus screens after an Active Insulin reset to zero alarm occurs.	<p>Select <b>OK</b> to clear the alarm.</p> <p>Contact 24-Hour Technical Support for assistance with the following steps:</p> <ol style="list-style-type: none"><li>1. Check the Daily History screen or the sensor graph for the recent bolus amounts, and when they were delivered, before giving any bolus.</li><li>2. Consult a healthcare professional for how long to wait after active insulin has been reset to zero before relying on the active insulin calculation of the Bolus Wizard feature. The active insulin tracked prior to the Active Insulin reset to zero alarm is not included in new Bolus Wizard calculations.</li><li>3. Check blood glucose (BG) using a blood glucose meter and treat as needed.</li></ol>




Issue	Resolution
	<div data-bbox="532 144 606 213"></div> <div data-bbox="629 144 1243 326"> <p><b>WARNING:</b> Do not rely on active insulin tracked in the pump when giving any bolus after active insulin has been reset to zero. Relying on the active insulin shown on the pump screen can result in the infusion of too much insulin, which can cause hypoglycemia.</p> </div>
<p>The pump buttons are stuck during airplane travel.</p>	<p>During atmospheric pressure changes, the pump buttons may not work for up to 45 minutes. For example, during airplane travel, pump buttons may get stuck and the pump will alarm. This is rare. If this occurs, either wait for the problem to correct itself, or confirm the AA battery connection:</p> <ol style="list-style-type: none"> <li>1. Remove the battery cap.</li> <li>2. Place the battery cap back onto the pump. The pump will check the AA battery power, and may require a new AA battery.</li> <li>3. If prompted, insert a new AA battery. For more information about changing the battery, see <i>Removing the battery</i>, page 68.</li> </ol> <p>If these steps do not correct the problem, contact 24-Hour Technical Support for assistance.</p>
<p>The pump was dropped or there are concerns that the pump may be damaged.</p>	<div data-bbox="532 885 606 954"></div> <div data-bbox="629 885 1243 1032"> <p><b>CAUTION:</b> Always inspect the pump for cracks before exposing the pump to water, especially if the pump was dropped or damaged. Water leakage can cause the pump to malfunction and result in injury.</p> </div> <ol style="list-style-type: none"> <li>1. Disconnect the pump from body. Confirm all infusion set and reservoir connections are secure.</li> <li>2. Disconnect the pump from body. Check the infusion set, including the tubing connector and tubing, for cracks or damage.</li> <li>3. Check the display, button area, and pump case for cracks or damage.</li> <li>4. Confirm the information on the Status screen.</li> <li>5. Confirm the settings for the basal rates and the pump are correct.</li> <li>6. Perform a self test. For more information, see <i>Self Test</i>, page 206.</li> <li>7. If necessary, contact 24-Hour Technical Support and check BG.</li> </ol> <p>For health-related questions or concerns, consult a healthcare professional.</p>

Issue	Resolution
The pump display times out too quickly.	In order to conserve battery, the pump display times out after 15 seconds. To increase the time, see <i>Display options, page 203</i> .
The pump displays a Check Settings alarm.	The pump has reset to factory settings. Review any settings that were not already set in the Startup Wizard and re-enter them, if necessary.
The pump settings have been cleared and need to be re-entered.	<p>Do not clear pump settings unless directed to do so by a healthcare professional. Certain pump errors may cause the pump to reset to factory default values, which clears the current pump settings. To restore saved pump settings, see <i>Restoring the settings, page 207</i>. Consult a healthcare professional to determine the necessary settings. Have the settings that need to be entered into the pump ready before starting the procedure below.</p> <p>Use the following procedure to re-enter personalized pump settings using the Startup Wizard:</p> <ol style="list-style-type: none"> <li>1. After the pump resets, the Startup Wizard appears. Select a language, and then press .</li> <li>2. Select a time format, and then press .</li> <li>3. Enter the current time, and then select <b>Next</b>.</li> <li>4. Enter the current date, and then select <b>Next</b>.</li> <li>5. Select the carb unit, and then press .</li> <li>6. When the Active Insulin Time screen appears, select Next. For more information, see <i>Bolus Wizard feature in Manual mode, page 99</i>.</li> <li>7. Enter the <b>Duration</b>, and then select <b>Next</b>.</li> <li>8. Enter the basal rates for the new basal pattern, and then select <b>Next</b>. For more information, see <i>Setting up a basal pattern, page 85</i>.</li> <li>9. Review the basal pattern information, and then select <b>Next</b>.</li> <li>10. On the Startup screen, a message displays to ask to set up Bolus Wizard now. Do one of the following: <ul style="list-style-type: none"> <li>• Select <b>Yes</b> to enter the Bolus Wizard settings in Manual mode. For more information, see <i>Bolus Wizard feature in Manual mode, page 99</i>.</li> <li>• Select <b>No</b> to skip the Bolus Wizard setup.</li> </ul> </li> </ol>

## Sensor issues

Issue	Resolution
The pump has lost connection with the sensor.	<p>After 30 minutes without a signal, the Lost sensor signal alert appears. Follow the steps on the pump screen or the steps below to try to resolve the issue.</p> <p><b>Note:</b> If alerts are silenced and a sensor alert occurs, the alert still appears on the screen.</p> <ol style="list-style-type: none"><li>1. Move the pump closer to the transmitter, and then select <b>OK</b>. It can take up to 15 minutes for the pump to find the transmitter signal. If the pump still cannot find the transmitter signal, the Possible signal interference alert appears.</li><li>2. Move away from electronic devices that may cause interference, and then select <b>OK</b>. Wait 15 minutes for the pump to locate the transmitter signal. If a signal is not found, the Check connection alert appears.</li><li>3. Confirm that the connection between the transmitter and sensor is secure, and then select <b>OK</b>. The "Check sensor insertion" message appears.</li><li>4. Do one of the following:<ul style="list-style-type: none"><li>• If the sensor connection is secure, select <b>Yes</b>. If the pump cannot find the sensor signal within 15 minutes or if the "Sensor signal not found - See User Guide" alert appears on the SG graph, contact 24-Hour Technical Support.</li><li>• If the sensor is not securely connected to the transmitter, select <b>No</b>. A Change sensor alert appears. Select <b>OK</b> and change the sensor.</li></ul></li></ol>
A calibration is not accepted.	<p>The sensor does not require calibration for use with the system. However, every BG reading entered into the pump is used to calibrate the sensor. A Calibration not accepted alert occurs in one of the following situations:</p> <ul style="list-style-type: none"><li>• The system cannot use the entered BG meter reading. Only a BG value between 50 mg/dL and 400 mg/dL can be used to calibrate the sensor. Wait at least 15 minutes, wash hands, and try again.</li><li>• The entered BG meter reading differs too greatly from the most recent SG reading. Check the accuracy of the BG meter reading and try again.</li></ul>

Issue	Resolution
	<ul style="list-style-type: none"> <li>The transmitter cannot receive the calibration BG meter readings from the pump due to a failed sensor signal. Troubleshoot the failed sensor signal.</li> </ul>
<p>The suspend by sensor icon appears with a red X. </p>	<p>The suspend by sensor icon appears with a red X when the Suspend before low or the Suspend on low feature is unavailable. This can occur in the following situations:</p> <ul style="list-style-type: none"> <li>A suspend event recently occurred. For information about the availability of the suspend functionality, see <i>The Suspend before low feature, page 160</i> or <i>The Suspend on low feature, page 162</i>.</li> <li>Sensor glucose (SG) values are unavailable.</li> </ul> <p>SG readings may be unavailable in the following situations:</p> <ul style="list-style-type: none"> <li>A BG meter reading is required.</li> <li>The pump has lost communication with the sensor. Restore pump communication with the sensor.</li> <li>The sensor is updating. Clear the alert and wait up to 3 hours for the SG values to resume.</li> </ul> <p>If necessary, insert a new sensor. If the issue continues after a new sensor is inserted, contact 24-Hour Technical Support.</p>

# 15



Pump Care



# 15

## Pump Care

The pump does not require preventative maintenance. This chapter provides information about caring for the components of the MiniMed 780G system.

### Clean, store, and dispose of the pump

#### Cleaning the pump

Prepare the following supplies to clean the pump:

- four small, clean, soft cloths
- mixture of water and mild detergent
- clean water
- 70% alcohol
- clean cotton swabs
- clean cotton balls



**CAUTION:** Never use organic solvents, such as lighter fluid, nail polish remover, or paint thinner to clean the MiniMed 780G insulin pump. Never use lubricants with the pump. When the pump is being cleaned, be sure to keep the reservoir compartment dry and away from moisture. If organic solvents are used to clean the pump, they can cause the pump to malfunction and result in minor injury.

### To clean the pump:

1. Dampen a cloth with water mixed with a mild detergent.
2. Use the cloth to wipe the outside of the pump while keeping the inside of the reservoir compartment dry.
3. Dampen a clean cloth with water and wipe to remove any detergent residue.
4. Dry with a clean cloth.
5. Wipe the pump with a 70% alcohol wipe.
6. Use a dry, clean cotton swab to remove any battery residue from the battery cap.
7. Use a dry, clean cotton swab to remove any battery residue from the battery compartment housing.

### Storing the pump

The pump can be stored when it is not in use.

If you place your pump in storage mode, it is important to insert a new AA battery for 8 to 12 hours every six months to ensure that the internal battery does not discharge to a deep discharge. A battery that is deeply discharged may experience decreased performance.



**WARNING:** After the pump is stored, do not rely on active insulin tracked in the pump when making new Bolus Wizard calculations. Storage mode clears active insulin. Inaccurate Bolus Wizard calculations may result in inaccurate insulin delivery and serious injury.

### To place the pump in storage mode:

1. Remove the AA battery from the pump. For details, see *Removing the battery*, page 68.





**Note:** When the battery is removed, the pump issues an Insert Battery alarm for 10 minutes or until the pump is in storage mode.

2. Press and hold  until the screen turns off.

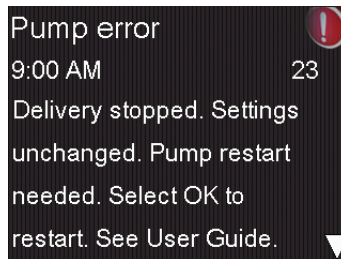


**CAUTION:** Never expose the pump to temperatures below -4 °F (-20 °C) or above 122 °F (50 °C). Storing the pump in temperatures outside of this range can damage the pump.

### To use the pump after it has been stored:


1. Insert a new AA battery into the pump. For details, see *Inserting the battery*, page 66.

A Pump Error alarm appears.




2. Select **OK**.

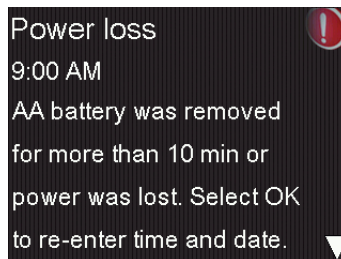The pump displays a Power Loss alarm.




3. Select **OK**.

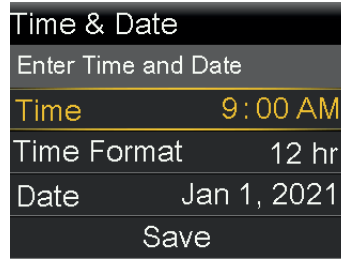The Time & Date screen appears.



4. Enter the current **Time**, **Time Format**, and **Date**.
5. Select **Save**.

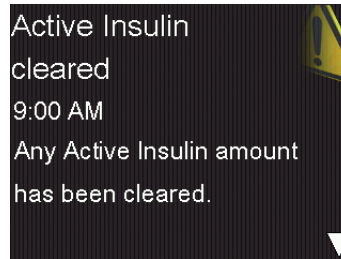The pump displays an Active Insulin Cleared alert.



6. Select **OK**.

Confirm that all the settings, such as basal rate, are set as desired. Use the Restore Settings option to reapply the last saved settings, if needed. For more information, see *Restoring the settings*, page 207.

7. Repeat the pairing process for the transmitter and meter. For transmitter details, see *Pairing the pump and transmitter*, page 137. For meter details, see *Pairing the pump and meter*, page 133.

## Pump disposal

Always follow local laws and regulations for the disposal of medical devices.

Do not dispose of the insulin pump in an unsorted municipal waste stream. The pump uses an AA battery and an internal battery.

Waste batteries, waste battery packs, electronics, and packaging are accepted by many recyclers. For information about recycling programs, contact local authorities. In addition, Medtronic and its distributors engage with many national recycling programs.

For information on the disposal of the components used with the MiniMed 780G system, refer to the corresponding user guide.



## Appendix A: List of alarms, alerts, and messages





# Appendix A: List of alarms, alerts, and messages

This appendix provides information about alarms, alerts, and messages that can occur in the MiniMed 780G system.

## Pump alarms, alerts, and messages

The following table lists the most common or serious alarms, alerts, and messages related to the MiniMed 780G insulin pump. The table also explains the meaning, consequences, and the reasons why these notifications appear, and provides steps for problem resolution.



**Note:** Use the MiniMed Mobile app to view the sensor graph on a mobile device. Always read and acknowledge all alarms and alerts on the pump. If the pump simultaneously generates more than one alarm or alert, only one of the alarms or alerts appears on the mobile device.

Title and text	Type	Explanation	Next steps
<b>Active Insulin cleared</b> Any Active Insulin amount has been cleared.	Alert	The pump shows the active insulin amount at 0 units. The pump shows this alert when the active insulin is cleared from the Clear Active Insulin option on the Manage Settings screen or if the pump	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>The active insulin tracked prior to pump restart is not included in new Bolus Wizard calculations. Consult a healthcare professional for how long to</li> </ul>

Title and text	Type	Explanation	Next steps
		has been shut down and is powered back on.	<p>wait after active insulin is cleared before relying on the active insulin calculation of the Bolus Wizard feature.</p> <ul style="list-style-type: none"> <li>• Check Daily History for the last bolus amount and when it was delivered.</li> </ul>
<b>Active Insulin reset to zero ?</b> Call local Medtronic support for assistance. See User Guide for phone numbers. Pump Active Insulin may be incorrect until XX:XX AM/PM due to a pump error. Monitor glucose.	Alarm	The pump shows the active insulin amount at 0 units. This occurs when a pump error clears active insulin in the pump. After the Active Insulin reset to zero alarm occurs, ? appears on the Home screen and Bolus screens until the time shown in the alarm.	<p>Select <b>OK</b> to clear the alarm. Contact 24-Hour Technical Support for assistance.</p>
<b>Active Insulin Reminder ?</b> Call local Medtronic support for assistance if needed. See User Guide for phone numbers. Pump Active Insulin was reset to zero at XX:XX AM/PM. Active insulin may be in-	Message	The Active Insulin Reminder message occurs when the Bolus Wizard screen or the Manual Bolus screen is accessed before the time shown in the message. After the Active Insulin reset to zero alarm occurs, ? appears on the Home screen and Bolus screens until the time shown in the Active Insulin reset to zero alarm or the Active Insulin Reminder message.	<p>Select <b>OK</b> to clear the message. Contact 24-Hour Technical Support for assistance with the following steps.</p> <ul style="list-style-type: none"> <li>• Check the Daily History screen or the sensor graph for the recent bolus amounts, and when they were delivered, before giving any bolus.</li> <li>• Consult a healthcare professional for how long to wait after active insulin</li> </ul>



Title and text	Type	Explanation	Next steps
correct until XX:XX AM/PM. Monitor glucose.			<p>has been reset to zero before relying on the active insulin calculation of the Bolus Wizard feature. The active insulin tracked prior to the Active Insulin reset to zero alarm is not included in new Bolus Wizard calculations.</p> <ul style="list-style-type: none"> <li>Check blood glucose (BG) using a blood glucose meter and treat as needed.</li> </ul>
<b>Auto Suspend</b> Insulin delivery suspended. No buttons pressed within time set in Auto Suspend.	Alarm	Insulin delivery is currently suspended by Auto Suspend. The Auto Suspend feature automatically suspends insulin delivery and triggers an alarm after no buttons are pressed for a specified period of time. Insulin delivery is suspended until the alarm is cleared and basal insulin delivery resumed.	<ul style="list-style-type: none"> <li>Select <b>Resume Basal</b> to clear the alarm and resume basal insulin delivery.</li> <li>Check BG and treat as needed.</li> </ul>
<b>Battery failed</b> Insert a new AA battery.	Alarm	The battery in the pump is low on power.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alarm.</li> <li>Remove the old battery and insert a new AA battery.</li> </ul>
<b>Battery not compatible.</b> See User Guide.	Alarm	The inserted battery is not compatible with the pump.	<ul style="list-style-type: none"> <li>Remove the incompatible battery to clear the alarm.</li> <li>Insert a new AA battery.</li> </ul>
<b>Bolus not delivered</b> Bolus entry timed out before deliv-	Alert	A bolus value was entered, but a bolus was not delivered within 30 seconds.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>If a bolus delivery was intended, check BG, re-en-</li> </ul>

Title and text	Type	Explanation	Next steps
ery. If bolus intended, enter values again.			ter bolus values and redeliver the bolus.
<b>Bolus stopped</b> Cannot resume bolus or cannula fill. XX.XXX of YY.YYY U delivered. ZZ.ZZZ U not delivered. If needed, enter values again.	Alarm	The battery power was exhausted while a bolus delivery or Fill Cannula procedure was in progress or the Resume bolus? message appeared and was not cleared.	<ul style="list-style-type: none"> <li>• Note the amount of insulin not delivered.</li> <li>• Replace the AA battery.</li> <li>• Select <b>OK</b> to clear the alarm.</li> <li>• Deliver the remaining bolus amount if needed.</li> </ul>
<b>Check settings</b> Startup Wizard settings complete. Check and set up your other settings.	Alert	Some settings have been cleared or reverted to factory default values.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to clear the alert.</li> <li>• Review any settings that have not already been set in Startup Wizard and re-enter the values if necessary.</li> </ul>
<b>Critical pump error</b> Delivery stopped. Pump not working properly. Stop using pump. Remove infusion set from body. Consider other insulin treatment. See User Guide.	Alarm	The pump has encountered an error that cannot be resolved. For example, the pump may have a mechanical problem.	The pump is not able to deliver insulin. Disconnect the infusion set and stop using the pump. <ul style="list-style-type: none"> <li>• Consider another form of insulin delivery.</li> <li>• Check BG, and treat as necessary.</li> <li>• Write down the error code that appears on the alarm screen.</li> <li>• Contact 24-Hour Technical Support for assistance with the pump.</li> </ul>
<b>Delivery limit exceeded</b> Delivery stopped. Check	Alarm	The pump has suspended insulin delivery because the hourly delivery limit was reached. This limit is based	<ul style="list-style-type: none"> <li>• Check BG.</li> <li>• Select <b>Resume Basal</b>.</li> </ul>

Title and text	Type	Explanation	Next steps
BG. See User Guide for more information.		on the maximum bolus and maximum basal setting. If this alarm occurs during a bolus, the bolus is canceled before it can complete.	<ul style="list-style-type: none"> <li>Check Bolus History and re-evaluate insulin needs.</li> <li>Continue to monitor BG.</li> </ul>
<b>Device Limit</b> You must delete an existing device (device type) before you can pair a new one (device type).	Message	The pump is already paired with the maximum number of devices for this type. The following list describes the maximum number of each <b>device type</b> to pair with the pump: <ul style="list-style-type: none"> <li>Meter—four Accu-Chek Guide Link meters</li> <li>CGM—one Guardian 4 transmitter</li> <li>Mobile Device—one compatible mobile device</li> </ul>	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the message.</li> <li>Go to the Paired Devices screen and select the device to unpair from the list of devices. Select <b>Unpair</b>, and then select <b>Yes</b> to confirm or <b>No</b> to cancel. Pair the pump and the desired device.</li> </ul>
<b>Device not compatible</b> Device cannot be used with this pump.	Alert	The pump cannot pair with the selected device.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Contact 24-Hour Technical Support for assistance.</li> </ul>
<b>Device not found</b> Make sure device is in range and in pairing mode.	Alert	The pump did not pair with the device.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Confirm that the device is not already paired with a pump.</li> <li>Confirm that the device is ready to pair with the pump.</li> <li>Make sure the pump is away from any electronic devices that might cause interference, such as cellular phones that</li> </ul>

Title and text	Type	Explanation	Next steps
			<p>are not paired with the MiniMed 780G system and other wireless devices.</p> <ul style="list-style-type: none"> <li>• Move the device closer to the pump.</li> <li>• Try to pair the pump with the device again.</li> </ul>
<b>Fill Cannula?</b> Select Fill to fill cannula or select Done if not needed.	Alarm	The Fill Cannula? screen has been active for 15 minutes.	<ul style="list-style-type: none"> <li>• To fill the cannula, select <b>Fill</b>.</li> <li>• If the cannula does not need to be filled, select <b>Done</b> to skip this process.</li> </ul>
<b>High BG XXX mg/dL</b> Check infusion set. Check ketones. Consider insulin injection. Monitor BG. Confirm BG?	Alert	<p>The BG meter reading is above 250 mg/dL.</p> <p>This alert appears in Manual mode. For High BG XXX mg/dL while the SmartGuard feature is on, see <i>SmartGuard feature alerts and messages, page 315</i>.</p>	<ul style="list-style-type: none"> <li>• Select <b>No</b> to prevent the remote BG from being used by the pump. Select <b>Yes</b> to confirm the BG reading.</li> <li>• Check BG and treat as necessary.</li> </ul>
<b>Insert battery</b> Delivery stopped. Insert a new battery now.	Alarm	<p>The battery was removed from the pump.</p> <p>If a bolus was in progress when the battery was removed, a Resume bolus? message appears and a tone sounds when a new battery is inserted. The message indicates how much of the bolus was delivered.</p>	<ul style="list-style-type: none"> <li>• Insert a new AA battery.</li> <li>• The alarm clears when a new battery is inserted.</li> <li>• The pump powers off after 10 minutes unless a new battery is inserted.</li> </ul>


Title and text	Type	Explanation	Next steps
<b>Insulin flow blocked</b> Check BG. Consider testing ketones. Check reservoir and infusion set.	Alarm	The pump has detected that the basal or bolus insulin flow was blocked.	<ul style="list-style-type: none"> <li>Check BG and ketones. Administer an insulin injection if necessary.</li> <li>Remove the infusion set and reservoir.</li> <li>Select <b>Reservoir &amp; Set</b> to start the process with a new infusion set and reservoir.</li> </ul> <p>If the alarm occurs during a bolus delivery:</p> <ul style="list-style-type: none"> <li>Check the Daily History screen for the amount of bolus already delivered before the pump alarmed.</li> <li>Consider delivering remaining bolus, if the bolus insulin was not included in an insulin injection.</li> </ul>



**WARNING:** Do not use the SmartGuard feature for a period of time after giving a manual injection of insulin by syringe or pen. Manual injections are not accounted for in the active insulin amount. Using the SmartGuard feature too soon after a manual injection may result in over-delivery of insulin and may cause hypoglycemia. Consult a healthcare professional for how long to wait after a manual injection before using the SmartGuard feature.

<b>Insulin flow blocked</b> Reservoir is empty. No insulin delivery. Change the reservoir and infusion set now. Check BG. Con-	Alarm	The pump has detected that the insulin flow is blocked and there is no insulin in the reservoir.	<ul style="list-style-type: none"> <li>Check BG and ketones. Administer an insulin injection if necessary.</li> <li>Remove the infusion set and reservoir.</li> <li>Select <b>Reservoir &amp; Set</b> to start the process with</li> </ul>
---	-------	--	--

Title and text	Type	Explanation	Next steps
sider testing ketones.			<p>a new infusion set and reservoir.</p> <p>If the alarm occurs during a bolus delivery:</p> <ul style="list-style-type: none"> <li>• Check the Daily History screen for the amount of bolus already delivered before the pump alarmed.</li> <li>• Consider delivering remaining bolus, if the bolus insulin was not included in an insulin injection.</li> </ul>
<b>Insulin flow blocked</b> Fill Cannula stopped. Restart the Reservoir & Set procedure.	Alarm	The pump has detected that the insulin flow is blocked while filling the cannula.	<ul style="list-style-type: none"> <li>• Check BG and ketones. Administer an insulin injection if necessary.</li> <li>• Remove the infusion set and reservoir.</li> <li>• Select <b>Reservoir &amp; Set</b> to start the process with a new infusion set and reservoir.</li> </ul>
<b>Insulin flow blocked</b> Fill Tubing stopped. Restart the Reservoir & Set procedure.	Alarm	The pump has detected that the insulin flow is blocked while filling the tubing. Possible connection issue between the tubing and reservoir.	<ul style="list-style-type: none"> <li>• Remove the reservoir and select <b>Reservoir &amp; Set</b> to restart the fill tubing process.</li> <li>• Disconnect the tubing from the reservoir.</li> <li>• Confirm that the tubing is not crimped or bent.</li> <li>• Continue to follow the steps displayed on the pump using the same infusion set and reservoir.</li> <li>• If this alarm occurs again, replace the infusion set.</li> </ul>

Title and text	Type	Explanation	Next steps
<b>Loading incomplete</b> Restart the Reservoir & Set procedure.	Alarm	 was pressed after loading began.	<ul style="list-style-type: none"> <li>Remove the reservoir to start again.</li> <li>Select <b>Reservoir &amp; Set</b> and follow the on-screen instructions.</li> </ul>
<b>Low battery Pump</b> Replace battery soon.	Alert	The battery in the pump is low on power. Remaining battery life is 10 hours or fewer.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Replace the AA battery as soon as possible. Otherwise, insulin delivery stops, and the Replace battery now alarm occurs.</li> <li>If the pump is delivering a bolus or filling the cannula, wait until delivery is complete to replace battery.</li> </ul>
<b>Low BG XX mg/dL</b> Treat Low BG. Do not bolus until BG is normal. Monitor BG. Confirm BG?	Alert	The BG meter reading is below 70 mg/dL.	<ul style="list-style-type: none"> <li>Select <b>No</b> to prevent the remote BG reading from being used by the pump. Select <b>Yes</b> to confirm the BG reading.</li> <li>Check BG and treat as necessary.</li> </ul>
<b>Low reservoir</b> XX units remaining. Change reservoir.	Alert	The reservoir is low on insulin, according to the number of units set in the Low Reservoir reminder.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Change the reservoir soon.</li> <li>If the reservoir is not changed after this alert is received, a second Low reservoir alert appears when the insulin level reaches half of the original alert amount.</li> </ul>

Title and text	Type	Explanation	Next steps
<b>Manage settings error</b> Delivery stopped. Backup settings cleared from Manage Settings. Current settings are working properly. Select OK to restart. See User Guide.	Alarm	A pump error occurred and the pump needs to be restarted. The backup settings have been lost, but the current settings are unchanged.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to restart the pump. The current settings are unchanged. Only the backup settings are lost.</li> <li>When the pump restarts, follow instructions on the pump display.</li> <li>If the pump was delivering a bolus or filling the cannula, check Daily History and evaluate if insulin is needed.</li> </ul>
<b>Max Fill reached</b> 3X.X U. Did you see drops at the end of tubing?	Alarm	The number of units expected to fill the tubing has been exceeded. By now, insulin should be visible at the end of the tubing.	<ul style="list-style-type: none"> <li>If there are drops of insulin at the end of the tubing, select <b>Yes</b>.</li> <li>If there are no drops of insulin at the end of the tubing, select <b>No</b>.</li> <li>Follow instructions displayed on the pump.</li> </ul>
<b>Max Fill reached</b> 4X.X U. Restart the Reservoir & Set procedure.	Alarm	The number of units expected to fill the tubing has been exceeded. By now, insulin should be visible at the end of the tubing.	<ul style="list-style-type: none"> <li>Remove the reservoir.</li> <li>Check if there is still insulin in the reservoir. If there is insulin in the reservoir the same reservoir can be used.</li> <li>Select <b>Reservoir &amp; Set</b> to restart the new reservoir procedure.</li> </ul>
<b>No reservoir detected</b> Restart the Reservoir & Set procedure.	Alarm	There is no reservoir in the pump or the reservoir is not properly locked into place.	<ul style="list-style-type: none"> <li>Select <b>Reservoir &amp; Set</b>.</li> <li>Confirm that the reservoir is filled with insulin.</li> <li>When prompted, confirm that the reservoir is insert-</li> </ul>



Title and text	Type	Explanation	Next steps
			ed and properly locked into place.
<b>Power error detected</b> Delivery stopped. Record your settings by uploading to CareLink or write your settings on paper. See User Guide.	Alarm	The internal power source in the pump is unable to charge. The pump is operating on the AA battery only.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to clear the alarm.</li> <li>• Check BG and treat as necessary.</li> <li>• Record the pump settings as soon as possible because the AA battery may not last long.</li> <li>• Contact 24-Hour Technical Support for assistance with the pump.</li> </ul>
<b>Power loss</b> AA battery was removed for more than 10 min or power was lost. Select OK to re-enter time and date.	Alarm	The battery has been out of the pump for more than ten minutes and the pump has lost power. The date and time must be reset.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to go to the Time &amp; Date screen.</li> <li>• Enter the current time, time format, and date.</li> </ul>
<b>Pump error</b> Delivery stopped. Current settings cleared. Pump restart needed. Select OK to restart and then re-enter your settings. See User Guide.	Alarm	The pump encountered an error and will restart. The pump settings will return to factory default values.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to restart the pump.</li> <li>• When the pump restarts, follow instructions on the pump display.</li> <li>• After the pump restarts, check settings and re-enter values as needed.</li> <li>• If the backup settings were recently saved in Manage Settings, use Restore Settings.</li> <li>• If the pump was delivering a bolus or filling the cannula, check Daily His-</li> </ul>

Title and text	Type	Explanation	Next steps
			<p>tory and re-evaluate if insulin is needed.</p> <ul style="list-style-type: none"> <li>• If this alarm recurs frequently, write down the error code on the alarm screen (it can also be found in the Alarm History) and contact 24-Hour Technical Support for assistance.</li> </ul>
<b>Pump error</b> Delivery stopped. Settings unchanged. Pump restart needed. Select OK to restart. See User Guide.	Alarm	A pump error has occurred, the pump needs to be restarted.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to restart the pump.</li> <li>• If the pump was delivering a bolus or filling the cannula, check Daily History and re-evaluate if insulin is needed.</li> <li>• If this alarm recurs frequently, write down the error code on the alarm screen (it can also be found in the Alarm History) and contact 24-Hour Technical Support for assistance.</li> </ul>
<b>Pump error</b> Delivery stopped. Settings unchanged. Select OK to continue. See User Guide.	Alarm	The pump encountered an error but a restart is not necessary. The issue is resolved. The settings are not changed.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to resume basal insulin delivery.</li> <li>• If the pump was delivering a bolus or filling the cannula, check Daily History and re-evaluate if insulin is needed.</li> <li>• If this alarm recurs frequently, write down the error code on the alarm screen (it can also be found in the Alarm History)</li> </ul>

Title and text	Type	Explanation	Next steps
			ry) and contact 24-Hour Technical Support for assistance.
<b>Pump restarted</b> Delivery stopped. Settings unchanged. Select OK to continue. See User Guide.	Alarm	The pump has encountered a problem and has restarted. The settings have not been changed.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to continue.</li> <li>If the pump was delivering a bolus or filling the cannula, check Daily History and re-evaluate if insulin is needed.</li> <li>If this alarm recurs frequently, write down the error code on the alarm screen (it can also be found in the Alarm History) and contact 24-Hour Technical Support for assistance.</li> </ul>
<b>Replace battery</b> Battery life less than 30 minutes. To ensure insulin delivery, replace battery now.	Alert	Battery power is low and will be exhausted within 30 minutes.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Replace the AA battery.</li> </ul>
<b>Replace battery now</b> Delivery stopped. Battery must be replaced to resume delivery.	Alarm	Insulin delivery has stopped due to low power. The battery was not replaced after the Low battery Pump alert.	Replace the battery immediately to resume insulin delivery.
<b>Reservoir estimate at 0 U</b> To ensure insulin delivery, change reservoir.	Alert	The reservoir level is estimated at 0 U.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Change the reservoir.</li> </ul>

Title and text	Type	Explanation	Next steps
<b>Resume bolus?</b> XXX of YYY U delivered. Resume delivery of ZZZ U?	Message	A normal bolus delivery has been interrupted because the pump battery was removed. If it is within ten minutes since this interruption, the bolus can be resumed.	<ul style="list-style-type: none"> <li>• Check the message to see how much of the bolus was delivered.</li> <li>• To cancel the remaining bolus delivery, select <b>Cancel</b>.</li> <li>• To resume the bolus delivery, select <b>Resume</b>.</li> </ul>
<b>Resume Dual bolus?</b> XX of YY U delivered. Resume delivery of ZZ U for XX:XX hr?	Message	The Square portion of Dual Bolus delivery has been interrupted. If it is within ten minutes since this interruption, the bolus can be resumed.	<ul style="list-style-type: none"> <li>• Check the message to see how much of the Dual Wave bolus was delivered.</li> <li>• To cancel the remaining bolus delivery, select <b>Cancel</b>.</li> <li>• To resume the bolus delivery, select <b>Resume</b>.</li> </ul>
<b>Resume Dual bolus?</b> XX of YY U delivered. Resume delivery of ZZ U now, and AA U Square for XX:XX hr?	Message	The Now portion of a Dual Wave bolus delivery has been interrupted because the pump battery was removed. If it is within ten minutes since this interruption, the bolus can be resumed.	<ul style="list-style-type: none"> <li>• Check the message to see how much of the Dual Wave bolus was delivered.</li> <li>• To cancel the remaining bolus delivery, select <b>Cancel</b>.</li> <li>• To resume the bolus delivery, select <b>Resume</b>.</li> </ul>
<b>Resume Square bolus?</b> XX of YY U delivered for XX:XX hr. Resume delivery of ZZ U for XX:XX hr?	Message	The Square Wave bolus delivery was interrupted. If it is within ten minutes since this interruption, the bolus can be resumed.	<ul style="list-style-type: none"> <li>• Check the message to see how much of the Square Wave bolus was delivered.</li> <li>• To cancel the remaining bolus delivery, select <b>Cancel</b>.</li> <li>• To resume the bolus delivery, select <b>Resume</b>.</li> </ul>

Title and text	Type	Explanation	Next steps
<b>Rewind required</b> Delivery stopped. Rewind was required due to pump error. Select OK to continue. See User Guide.	Alarm	The pump encountered an error.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alarm after the pump has completed rewinding.</li> <li>Select <b>Reservoir &amp; Set</b> from the Menu screen to start the new reservoir process with a new infusion set and reservoir. For details, see <i>Setting up the reservoir and infusion set</i>, page 114.</li> <li>If this alarm recurs frequently, contact 24-Hour Technical Support for assistance.</li> </ul>
<b>Stuck button</b> Button pressed for more than 3 minutes.	Alarm	The pump has detected that a button has been pressed for an unusually long time.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alarm.</li> <li>If this alarm occurs again, contact 24-Hour Technical Support for assistance with the pump.</li> </ul> <p>If the alarm cannot be cleared</p> <ul style="list-style-type: none"> <li>See <i>Pump issues</i>, page 276.</li> <li>Consider another form of insulin, because the pump is not delivering insulin.</li> <li>Check BG and treat as necessary.</li> <li>Contact 24-Hour Technical Support for assistance with the pump.</li> </ul>

## CGM (sensor) alarms, alerts, and messages

The following table lists the most common or serious alarms, alerts, and messages related to sensor glucose (SG) values, as well as the status of the transmitter and sensor. The table also explains the meaning, consequences, and the reasons why these notifications appear, and provides steps for problem resolution.

Title and text	Type	Explanation	Next steps
<b>Alert before high</b> Sensor glucose approaching High Limit. Check BG.	Alert	The SG reading is approaching the specified high limit.	<ul style="list-style-type: none"><li>• Select <b>OK</b> to clear the alert.</li><li>• Check BG.</li><li>• Follow instructions from a healthcare professional and continue to monitor BG.</li></ul>
<b>Alert before low</b> Sensor glucose approaching Low Limit. Check BG.	Alert	The SG reading is approaching the specified low limit.	<ul style="list-style-type: none"><li>• Select <b>OK</b> to clear the alert.</li><li>• Check BG.</li><li>• Follow instructions from a healthcare professional and continue to monitor BG.</li></ul>
<b>Alert on high XXX mg/dL</b> High sensor glucose. Check BG.	Alert	The SG reading is at or above the specified high limit.	<ul style="list-style-type: none"><li>• Select <b>OK</b> to clear the alert.</li><li>• Check BG.</li><li>• Follow instructions from a healthcare professional and continue to monitor BG.</li></ul>
<b>Alert on low XX mg/dL</b> Low sensor glucose. Check BG.	Alert	The SG reading is at or below the specified low limit.	<ul style="list-style-type: none"><li>• Select <b>OK</b> to clear the alert.</li><li>• Check BG.</li><li>• Follow instructions from a healthcare professional and continue to monitor BG.</li></ul>
<b>Alert on low XX mg/dL</b> Low sensor glucose. Insulin delivery suspended since XX:XX AM/PM. Check BG.	Alarm	The SG reading is at or below the specified low limit, and the pump has suspended insulin delivery due to a Suspend before low or Suspend on low.	<ul style="list-style-type: none"><li>• Select <b>OK</b> to clear the alarm.</li><li>• Check BG.</li><li>• Follow instructions from a healthcare professional and continue to monitor BG.</li></ul>

Title and text	Type	Explanation	Next steps
<b>Basal delivery resumed</b> Basal delivery resumed at XX:XX AM/PM after suspend by sensor. Check BG.	Message	The pump is resuming basal insulin delivery after a Suspend before low or Suspend on low event occurred.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the message.</li> <li>Check BG.</li> <li>Follow instructions from a healthcare professional and continue to monitor BG.</li> </ul>
<b>Basal delivery resumed</b> Basal delivery resumed at XX:XX AM/PM after suspend by sensor. Check BG.	Alert	The pump is resuming basal insulin delivery after a Suspend before low or a Suspend on low event occurred, because the Suspend before low or the Suspend on low feature was turned off.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Check BG.</li> <li>Follow instructions from a healthcare professional and continue to monitor BG.</li> </ul>
<b>Basal delivery resumed</b> Maximum 2 hour suspend time reached. Check BG.	Alert	The pump is resuming basal insulin delivery two hours after a Suspend before low or Suspend on low event occurred.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Check BG.</li> <li>Follow instructions from a healthcare professional and continue to monitor BG.</li> </ul>
<b>Basal delivery resumed</b> Maximum 2 hour suspend time reached. SG is still under Low limit. Check BG.	Alarm	The pump is resuming basal insulin delivery two hours after a Suspend before low or Suspend on low event occurred.	<ul style="list-style-type: none"> <li>The pump has resumed basal insulin delivery; however, the SG reading is still at or below the low limit.</li> <li>Select <b>OK</b> to clear the alarm.</li> <li>Check BG.</li> <li>Follow instructions from a healthcare professional and continue to monitor BG.</li> </ul>

Title and text	Type	Explanation	Next steps
<b>BG not received</b> Place pump close to transmitter. Select OK to resend BG to transmitter.	Alert	The transmitter was unable to receive the BG meter reading from the pump.	<ul style="list-style-type: none"> <li>• Move the pump and transmitter closer together.</li> <li>• Select <b>OK</b> to clear the alert, and then enter a new BG meter reading.</li> </ul>
<b>Calibration not accepted</b> Sensor information is unavailable for up to 2 hours. Entered BGs may not calibrate the sensor but can still be used for therapy.	Alert	The system was unable to use the BG meter readings entered to calibrate the sensor. This alert occurs on the first day only.	<ul style="list-style-type: none"> <li>• Wash and dry hands thoroughly.</li> <li>• Select <b>OK</b> to clear the alert.</li> <li>• Consider waiting up to two hours, and then enter a new BG meter reading.</li> <li>• Contact 24-Hour Technical Support for assistance, if needed.</li> </ul>
<b>Calibration not accepted</b> Wait at least 15 minutes. Wash hands, test BG again and calibrate.	Alert	The system was unable to use the BG meter readings entered to calibrate the sensor.	<ul style="list-style-type: none"> <li>• Wash and dry hands thoroughly.</li> <li>• Select <b>OK</b> to clear the alert.</li> <li>• After 15 minutes, enter a new BG meter reading. If a Calibration not accepted alert is received on the second calibration after 15 minutes, a Change sensor alert occurs.</li> <li>• Contact 24-Hour Technical Support for assistance, if needed.</li> </ul>
<b>Change sensor</b> Insert new sensor and Start New Sensor.	Alert	<b>No</b> was selected in the Check sensor insertion message, indicating that the sensor is not fully inserted.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to clear the alert.</li> <li>• Change the sensor. For details, see the sensor user guide.</li> <li>• After the sensor is changed, refer to <i>Starting the sensor</i>, page 172.</li> </ul>
<b>Change sensor</b> Second calibration not accepted. Insert new sensor.	Alert	A Calibration not accepted alert occurs if the entered BG meter reading differs too greatly from the most recent SG reading.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to clear the alert.</li> <li>• Change the sensor. For details, see the sensor user guide.</li> </ul>




Title and text	Type	Explanation	Next steps
		This alert occurs when two Calibration not accepted alerts are received in a row.	
<b>Change sensor</b> Sensor not working properly. Insert new sensor.	Alert	This alert occurs when the pump diagnoses a problem with the sensor that cannot be resolved.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Change the sensor. For details, see the sensor user guide.</li> </ul>
<b>Check connection</b> Ensure transmitter and sensor connection is secure, then select OK.	Alert	The pump fails to detect the transmitter and is unable to receive sensor signal.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>If the sensor is fully inserted, select <b>Yes</b>. If the sensor is not fully inserted, select <b>No</b>.</li> <li>If the sensor was not fully inserted, insert a new sensor.</li> <li>See <i>Pump issues</i>, page 276 for additional assistance, if needed.</li> </ul>
<b>Enter BG now</b> Enter BG to calibrate sensor. Sensor information is no longer available	Alert	A BG meter reading is required to calibrate the sensor. SG readings cannot be received until the sensor is calibrated.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert. If no BG meter reading is entered within 30 minutes, the Enter BG now alert occurs again.</li> <li>Select <b>Snooze</b>, enter the desired snooze time, and select <b>OK</b>. If no BG meter reading is entered before the Snooze time has ended, the Enter BG now alert occurs again.</li> <li>Enter a BG meter reading to calibrate the sensor.</li> </ul>
<b>High SG</b> Glucose was 250 mg/dL or higher for more than 3 hours. Check infusion set. Check ke-	Alert	SG was 250 mg/dL or higher for three hours.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Check BG and treat as necessary.</li> </ul>

Title and text	Type	Explanation	Next steps
tones. Monitor glucose.			
<b>Lost sensor signal</b> Move pump closer to sensor. May take 15 minutes to find signal.	Alert	A transmitter signal has not been received for 30 minutes during or after sensor initialization.	<ul style="list-style-type: none"> <li>• Move the pump closer to the transmitter. It can take up to 15 minutes for the pump to establish communication with the transmitter.</li> <li>• Select <b>OK</b> to clear the alert.</li> </ul>
<b>Low battery transmitter</b> Recharge transmitter within 24 hours.	Alert	The battery in the transmitter needs to be recharged within 24 hours.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to clear the alert.</li> <li>• Recharge the transmitter as soon as possible.</li> </ul>
<b>Low SG XX mg/dL</b> SG is under 64 mg/dL. Check BG and treat.	Alarm	The SG reading has fallen below 64 mg/dL. This alarm is factory set and cannot be changed or turned off. This alarm cannot be silenced and is always active, whether the pump is using the SmartGuard feature or Manual mode.	<ul style="list-style-type: none"> <li>• Select <b>OK</b> to clear the alarm.</li> <li>• Check BG and treat as necessary.</li> </ul>

**Note:** This alarm does not suspend insulin delivery.

**Note:** XX represents the current SG reading that appears on the pump. This alarm remains until the alarm is cleared, even if glucose values reach or rise above 64 mg/dL.

Title and text	Type	Explanation	Next steps
<div>  <b>WARNING:</b> For MiniMed 780G Users Ages 7-13: Do not rely solely on the use of a low sensor glucose (SG) value for “Alert on Low” or “Alert before Low” or the “Low SG” alarm. A low sensor glucose alert may not reflect the user’s true blood glucose at these levels, or may not alert. Do not ignore symptoms of low glucose. Always confirm sensor glucose readings with a blood glucose meter, and treat according to the recommendations of a healthcare professional. Solely relying on these sensor glucose alerts and readings for treatment decisions could result in missing severe hypoglycemia (low blood glucose) events. </div>			
<b>Max Basal</b>	Message	The Max basal rate setting is being changed.	<ul style="list-style-type: none"> <li>Select <b>Continue</b> or <b>Cancel</b> to clear the message.</li> </ul>
<b>Max Basal Error</b> Your max basal rate cannot be lower than an existing basal rate.	Alert	You are trying to lower the Max basal rate setting lower than existing basal rates.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Lower other basal rates before changing the Max basal rate.</li> </ul>
<b>Medical device</b> CALL FOR EMERGENCY ASSISTANCE. I have diabetes.	Alarm	The pump is suspended due to low SG and there has been no response to the alarm within 10 minutes.	<ul style="list-style-type: none"> <li>Select <b>Dismiss</b>.</li> <li>Immediately call for emergency assistance.</li> </ul>
<b>No calibration occurred</b> Confirm sensor signal. Calibrate by XX:XX AM/PM.	Alert	The sensor was unable to receive the calibration BG meter readings from the pump.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Check the status icons on the Home screen to confirm that the pump has a signal from the sensor. If there is no sensor signal, see <i>Sensor issues</i>, page 279.</li> <li>For SG readings to be monitored without interruption, enter or confirm a BG meter reading by the time displayed on the pump screen.</li> </ul>

Title and text	Type	Explanation	Next steps
<b>No calibration occurred</b> Confirm sensor signal. Check BG again to calibrate sensor.	Alert	The sensor was unable to receive the required calibration BG meter readings from the pump. Calibration is required by the system for SG readings to resume. "Calibration required" appears on the sensor graph.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Take another BG meter reading and calibrate again.</li> </ul>
<b>Not for SmartGuard</b> These settings apply to Manual mode only and do not affect how much insulin SmartGuard delivers.	Message	The Max basal or Max bolus rate settings or basal pattern settings are being accessed.	<ul style="list-style-type: none"> <li>Select <b>Continue</b> or <b>Cancel</b> to clear the message.</li> </ul>
<b>Outside typical range</b> Max Basal of 6.00 U/hr can allow basal delivery of up to 144 U/day in Manual mode. This amount of insulin may not be safe.	Alert	The Max basal rate setting is being increased past 6 U/hr.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Consult your healthcare professional before increasing the Max basal rate.</li> </ul>
<b>Possible signal interference</b> Move away from electronic devices. May take 15 minutes to find signal.	Alert	There may be interference from another electronic device that is affecting the communication be-	<ul style="list-style-type: none"> <li>Move away from other electronic devices. It can take up to 15 minutes for the pump to start communicating with the transmitter.</li> <li>Select <b>OK</b> to clear the alert.</li> </ul>

Title and text	Type	Explanation	Next steps
		tween the pump and the transmitter.	
<b>Rise Alert</b> Sensor glucose rising rapidly.	Alert	The SG reading has been rising as fast or faster than the preset Rise Alert limit.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Check BG using a meter.</li> <li>Follow instructions from a healthcare professional.</li> </ul>
<b>Sensor connected</b> If new sensor, select Start New. If not, select Reconnect.	Message	The transmitter has detected that a sensor is connected. The pump needs to know if this is a new sensor or if an old sensor has been reconnected.	<ul style="list-style-type: none"> <li>If a new sensor has been connected, select <b>Start New Sensor</b>.</li> <li>If a sensor that was already being used has been reconnected, select <b>Reconnect Sensor</b>.</li> <li>In either case, a "Sensor warm up X:XX hr" message appears for two hours. After the warm up, the pump starts receiving SG readings.</li> </ul>
<b>Sensor connected</b> Start new sensor.	Message	The pump has detected that this is a new sensor, which needs to be started and warmed-up.	<p>Select <b>Start New Sensor</b>.</p> <p>The alert closes and a "Warm-up" message appears on the sensor graph with a progress bar.</p>
<b>Sensor expired</b> Insert new sensor.	Alert	The sensor has reached the end of its useful life.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Change the sensor. For details, see the sensor user guide.</li> </ul>
<b>Sensor signal not found</b> See User Guide.	Alert	After multiple attempts, the pump failed to detect the sensor and is unable to receive sensor signal.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>If the pump still cannot find the sensor signal, contact 24-Hour Technical Support for assistance.</li> </ul>
<b>Sensor updating</b> Updating can take X. Monitor BG. Entered BGs will not calibrate the sen-	Alert	The SG reading is unavailable due to a temporary situation.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Follow the instructions on the pump screen. The sensor does not need to be changed.</li> </ul>

Title and text	Type	Explanation	Next steps
sor, but can still be used for therapy.			
<b>Suspend before low</b> Delivery stopped. Sensor glucose approaching Low Limit. Check BG.	Alert	The SG reading is falling. Insulin delivery is suspended according to the Suspend before low setting and the SG is approaching the specified low limit. The Suspend before low feature is not available with the SmartGuard feature.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Check BG. If necessary, treat BG as directed by a healthcare professional.</li> </ul>
<b>Suspend on low</b> Delivery stopped. Sensor glucose XX mg/dL. Check BG.	Alarm	The SG reading is at or below the specified low limit. The Suspend on low feature is not available with the SmartGuard feature.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alarm.</li> <li>Check BG. If necessary, treat BG as directed by a healthcare professional.</li> </ul>
<b>Transmitter battery depleted</b> Recharge transmitter now.	Alert	The battery in the transmitter needs to be recharged. SG readings cannot be recorded or transmitted until the transmitter is recharged.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Recharge the transmitter.</li> </ul>
<b>Very high Basal setting</b> Basal X will deliver YY.YY U per day when active in Manual mode. This basal pattern may	Message	When selecting, adding, or editing a basal pattern and the system calculates the pattern would deliver significantly more in-	<ul style="list-style-type: none"> <li>Select <b>Edit</b> or, if selecting a basal pattern, click <b>Back</b> to change the basal pattern.</li> <li>Select <b>Continue</b> to save or select the basal pattern.</li> </ul>

Title and text	Type	Explanation	Next steps
not be safe because it delivers significantly more insulin than you typically need. Consult healthcare professional before changing this value.		insulin than typically needed.	

### SmartGuard feature alerts and messages

The following table lists the most common or serious alerts and messages related to the SmartGuard feature. The table also explains the meaning, consequences, and the reasons why these notifications appear, and provides any necessary steps for problem resolution.

Title and text	Type	Explanation	Next steps
<b>Enter BG now</b> SmartGuard has been at maximum delivery rate for 7 hours. Enter BG to continue in SmartGuard.	Alert	SmartGuard has been delivering at the maximum SmartGuard basal delivery rate for seven hours. This rate is determined automatically by the system.	<ul style="list-style-type: none"><li>• Select <b>OK</b> to clear the alert.</li><li>• Enter a BG meter reading to return to Auto Basal.</li><li>• Follow instructions from a healthcare professional and continue to monitor BG.</li></ul>
<b>Enter BG now</b> SmartGuard has been at maximum delivery rate for 7 hours. Enter BG to continue in SmartGuard. This event occurred while pump was suspended, and ac-	Alert	The pump is suspended and the SmartGuard feature has been unable to lower the SG reading. SG is predicted to remain above the SmartGuard target.	<ul style="list-style-type: none"><li>• Select <b>OK</b> to clear the alert.</li><li>• Enter a BG meter reading.</li><li>• Follow instructions from a healthcare professional and continue to monitor BG.</li></ul>

Title and text	Type	Explanation	Next steps
tion is required to resume delivery.			
<b>Notes:</b>			
<ul style="list-style-type: none"> <li>The title of the alert appears the same as the previous SmartGuard max delivery alert in the table.</li> <li>If the pump is suspended, there will be no delivery. However, the alert may still occur.</li> </ul>			
<b>Enter BG now</b> SmartGuard has reached the time limit for minimum delivery rate. Enter BG to continue in SmartGuard.	Alert	The SmartGuard feature has reached the time limit for minimum delivery. The minimum delivery time is three to six hours, depending on the reason for the minimum delivery rate.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Enter a BG meter reading to return to Auto Basal.</li> <li>Follow instructions from a health-care professional and continue to monitor BG.</li> </ul>
<b>Enter BG now</b> SmartGuard has reached the time limit for minimum delivery rate. Enter BG to continue in SmartGuard. This event occurred while pump was suspended, and action is required to resume delivery.	Alert	SmartGuard has reached the time limit for minimum delivery. The minimum delivery time is three to six hours, depending on the reason for the minimum delivery rate.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Enter a BG meter reading.</li> <li>Follow instructions from a health-care professional and continue to monitor BG.</li> </ul>
<b>Notes:</b>			
<ul style="list-style-type: none"> <li>The title of the alert appears the same as the previous SmartGuard min delivery alert in the table.</li> <li>If the pump is suspended, there will be no delivery. However, the alert may still occur.</li> </ul>			
<b>Enter BG now</b> Enter BG to continue in SmartGuard.	Alert	The SmartGuard feature requires a BG reading to check the reliability of the sensor.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Enter a BG meter reading to return to Auto Basal, or to enter the</li> </ul>



Title and text	Type	Explanation	Next steps
			SmartGuard feature from Manual mode.
<b>Enter BG now</b> Glucose values may be lower than current SG. Enter BG to continue in SmartGuard.	Alert	The SmartGuard feature requires a BG reading to check the reliability of the sensor.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Enter a BG meter reading to return to Auto Basal, or to enter the SmartGuard feature from Manual mode.</li> </ul>
<b>High BG XXX mg/dL</b> Check infusion set. Check ketones. Monitor BG. Confirm BG?	Alert	The BG meter reading is above 250 mg/dL. This alert applies only to the SmartGuard feature. There is a similar alert for Manual mode when the SmartGuard feature is off. See <i>SmartGuard</i> , page 179.	<ul style="list-style-type: none"> <li>Select <b>No</b> to prevent the remote BG from being used by the pump.</li> <li>Select <b>Yes</b> to confirm the BG reading.</li> </ul>
<b>SmartGuard exit</b> Basal X started. Would you like to review the SmartGuard Checklist?	Alert	The pump has exited the SmartGuard feature because: <ul style="list-style-type: none"> <li>the sensor has been turned off</li> <li>the pump has been delivering basal insulin based on insulin delivery history, and not SG readings for the maximum of four hours.</li> </ul> This alert cannot be silenced, and is al-	<ul style="list-style-type: none"> <li>Select <b>No</b> to clear the alert. Select <b>Yes</b> to view the SmartGuard Checklist.</li> <li>Enter a BG meter reading.</li> <li>Follow instructions from a health-care professional and continue to monitor BG.</li> </ul> For details, see <i>Exiting the SmartGuard feature</i> , page 198 and <i>Returning to the SmartGuard feature after an exit</i> , page 199.

Title and text	Type	Explanation	Next steps
		ways active whenever the system is using the SmartGuard feature.	
<b>SmartGuard exit</b> Insulin delivery is still suspended.	Alert	<p>The pump has exited the SmartGuard feature because:</p> <ul style="list-style-type: none"> <li>the sensor has been turned off</li> <li>a suspend event message has not been cleared within four hours</li> <li>the pump has been delivering basal insulin based on insulin delivery history, and not SG readings for the maximum of four hours.</li> </ul> <p>This alert cannot be silenced, and is always active whenever the system is using the SmartGuard feature.</p>	<ul style="list-style-type: none"> <li>Enter a BG meter reading.</li> <li>Manually resume basal insulin delivery, when appropriate.</li> <li>Follow instructions from a health-care professional and continue to monitor BG.</li> </ul> <p>For details, see <i>Exiting the SmartGuard feature</i>, page 198 and <i>Returning to the SmartGuard feature after an exit</i>, page 199.</p>
<b>SmartGuard exit</b> Manual mode Basal X started. Check basal setting. Monitor glucose.	Alert	The pump has exited the SmartGuard feature and the system has detected a basal pattern that delivers significant-	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Check the active basal pattern.</li> </ul>

Title and text	Type	Explanation	Next steps
		ly more insulin than you typically need.	
<b>SmartGuard started</b> Current action canceled.	Alert	An operation that is not allowed while transitioning to the SmartGuard feature has been selected.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Allow the pump to complete its transition to the SmartGuard feature.</li> </ul>
<b>Very high Basal setting</b> SmartGuard is no longer active. Manual mode Basal X started and will deliver YY.YY U per day. Basal X may not be safe because it delivers significantly more insulin than you typically need. Consult healthcare professional for basal setting. Monitor glucose.	Alert	The SmartGuard feature has exited and the system has detected that the current basal pattern would deliver significantly more insulin than you typically need.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Check all basal patterns and consult with your healthcare professional.</li> </ul>
<b>Very high Basal setting</b> Basal X will deliver YY.YY U per day when active in Manual mode. Basal X may not be safe because it delivers significantly more insulin than you typically need. Consult healthcare professional for basal settings.	Alert	The system has determined that if SmartGuard exits, the active basal pattern will deliver significantly more insulin than you typically need.	<ul style="list-style-type: none"> <li>Select <b>Snooze</b> to postpone receiving this message for a set period of time.</li> <li>Select <b>OK</b> to clear the alert.</li> <li>Check all basal patterns and consult with your healthcare professional.</li> </ul>

Title and text	Type	Explanation	Next steps
<b>Very high Basal setting</b> Basal X is active and delivers YY.YY U per day. Basal X may not be safe because it delivers significantly more insulin than you typically need. Consult healthcare professional for basal settings. Monitor glucose.	Alert	The pump has detected the current Manual mode basal pattern that delivers significantly more insulin than you typically need.	<ul style="list-style-type: none"> <li>Select <b>Snooze</b> to postpone receiving this message for a set period of time.</li> <li>Select <b>OK</b> to clear the alert.</li> <li>Check all basal patterns and consult with your healthcare professional.</li> </ul>

## CareLink software alert and message

The following table lists the most common or serious alerts and messages related to CareLink software. The table also explains the meaning, consequences, and the reasons why these notifications appear, and provides steps for problem resolution. If an alarm, alert, or message occurs that is not listed, select **OK** to clear the notification and contact 24-Hour Technical Support.

Title and text	Type	Explanation	Next steps
<b>CareLink uploader not found.</b> Follow instructions on the CareLink uploader.	Message	The pump cannot find the CareLink uploader because the wrong pump code was entered, or the search timed out before the pump found the uploader.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the message.</li> <li>Follow the instructions on the CareLink uploader. For details, see <i>Uploading device data to CareLink software, page 144</i>.</li> </ul>
<b>Download slow</b> Insulin delivery not affected. CareLink download may take longer than usual.	Alert	The download of pump data is taking longer than expected. Data will not be affected.	<ul style="list-style-type: none"> <li>Select <b>OK</b> to clear the alert.</li> <li>Wait for the data to finish downloading.</li> <li>If problem still persists or if there is no progress in download, call</li> </ul>

Title and text	Type	Explanation	Next steps
Select OK to continue. See User Guide.			24-Hour Technical Support for assistance.









# Appendix B: Product specifications

This appendix provides detailed product specifications.

## Specifications and default settings

### Alarm and alert escalation

The following alerts may escalate to a siren if not cleared:

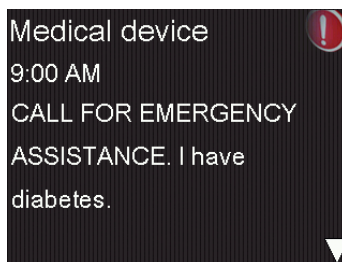
- Alert before high
  - Alert before low
  - Alert on high
  - Alert on low
  - Basal delivery resumed
  - Calibration not accepted
  - Change sensor
  - Enter BG now
  - Lost sensor signal
  - No calibration occurred
- Possible signal interference
  - High SG
  - Rise Alert
  - Sensor expired
  - Sensor signal not found
  - Low SG XX mg/dL (XX is a value below 64 mg/dL)
  - Sensor updating
  - Warm up not started
  - Very high basal setting

The MiniMed 780G insulin pump may generate a siren if the alert is not cleared within ten minutes. Before ten minutes, the pump beeps, vibrates, or both, depending on the sound and vibration settings.

Minutes	Sound	Vibration	Sound and vibration
0-5	Beep	Vibrate	Beep and vibrate
6-9	Beep and vibrate	Sound and vibrate	Beep and vibrate
10	Siren and vibrate	Siren and vibrate	Siren and vibrate



**Note:** The Medical device alarm plays a siren when this screen appears.



## Altitude range

- Operating range: 10.2 psiA (70.33 kPa) to 15.4 psiA (106.18 kPa)
- Storage range: 7.2 psiA (49.64 kPa) to 15.4 psiA (106.18 kPa)

## Backlight

Type	LED (Light-emitting Diode)
Time out	15 seconds (default), 30 seconds, one minute, three minutes
Time out when battery is low	15 seconds (default), 30 seconds

## Basal delivery

The basal delivery specifications do not affect the amount of insulin delivered by SmartGuard. These specifications relate to Manual mode only.

Delivery rate range	0 to 35 units per hour or the Max Basal Rate amount, whichever is lower. The default range is 0 to 10 U/hr.
Max Basal Rate default	2 units per hour
Basal patterns	Maximum of 8 patterns for Manual mode. Each pattern covers a 24-hour period and can have up to 48 rates. Rates are set in 30-minute increments.

Basal pattern names	Fixed names: Basal 1, Basal 2, Basal 3, Basal 4, Basal 5, Workday, Day Off, Sick Day
Increments	<ul style="list-style-type: none"><li>• 0.025 units per hour for basal amounts in the range 0 to 0.975 units</li><li>• 0.05 units per hour for basal amounts in the range 1 to 9.95 units</li><li>• 0.1 units per hour for basal amounts of 10 to 35 units</li></ul>

### BG meter reading

The BG meter reading refers to the most recent blood glucose (BG) meter reading received from the blood glucose meter after it is accepted by the pump user. When an Accu-Chek Guide Link meter is used, the reading appears on the Home screen when the Sensor feature is off. The reading also appears in the Bolus Wizard and SmartGuard Bolus screens when a bolus is programmed.

Expiration	12 minutes
Range	20 to 600 mg/dL

### Bolus delivery

Bolus Speed options	<ul style="list-style-type: none"><li>• Standard: 1.5 units/minute</li><li>• Quick: 15 units/minute</li></ul>
Bolus programming increments	<ul style="list-style-type: none"><li>• 0.025 units</li><li>• 0.05 units</li><li>• 0.1 units</li></ul>
Fluid delivered/stroke	<ul style="list-style-type: none"><li>• 0.25 µL (microliter) for 0.025 unit pump stroke</li><li>• 0.5 µL for 0.05 unit pump stroke</li><li>• 2.0 µL for 0.2 unit pump stroke</li></ul>

### Bolus Wizard feature in Manual mode default settings



**Note:** When using the SmartGuard feature, the Bolus Wizard feature is called the Bolus feature.

Item	Default	Limits	Maximum available segments	Increments
<b>Carb units</b>	grams	-	8	-
<b>Insulin to carb ratio</b>	None	1–200 g/U	8	0.1 g/U for 1–9.9 g/U; 1 g/U for ratios of 10 g/U to 200 g/U
<b>Insulin Sensitivity Factor*</b>	None	5–400 mg/dL	8	1 mg/dL
<b>BG Target*</b>	None	60–250 mg/dL	8	1 mg/dL
<b>Active Insulin Time</b>	4 hours	2 to 8 hours	1	15 minutes

\*Applies to Manual mode only.

## Bolus Wizard feature in Manual mode specifications

The Bolus Wizard feature uses four formulas to estimate a bolus, depending on the current BG reading. The following formulas apply only when the carb units are in grams.

1. If the current BG reading is higher than the High BG Target, the Bolus Wizard feature subtracts active insulin from the BG correction estimate, then adds this value to the food estimate to get the total bolus estimate. However, if the result of subtracting the active insulin amount from the BG correction estimate is a negative number (less than zero), the total bolus estimate is based only on the food estimate.

$$\text{total bolus estimate} = \frac{\text{(food estimate)}}{\frac{A}{B}} + \frac{\text{(correction estimate)}}{\frac{C - D}{E}} - \text{active insulin}$$

where: A = food (grams)  
 B = carb ratio  
 C = current BG  
 D = High BG Target  
 E = insulin sensitivity

Food estimate:

Carb grams ÷ Carb ratio = Units of insulin

Correction estimate:

(Current BG - High BG Target) ÷ Insulin sensitivity - Active insulin = Units of insulin

Total bolus estimate:

Food estimate + Correction estimate = Units of insulin

2. If the current BG is less than the Low BG Target, the Bolus Wizard feature adds the BG correction estimate to the food estimate to get the total bolus estimate.

(food estimate)

(correction estimate)

total bolus estimate =  $\frac{A}{B} + \frac{C - D}{E}$

where: A = food (grams)  
B = carb ratio  
C = current BG  
D = Low BG Target  
E = insulin sensitivity

Food estimate:

Carb grams ÷ Carb ratio = Units of insulin

Correction estimate:

(Current BG - Low BG Target) ÷ Insulin sensitivity = Units of insulin

Total bolus estimate:

Food estimate + Correction estimate = Units of insulin

3. If the current BG reading is within the High or Low BG Target, the total bolus estimate is based only on the food estimate.

(food estimate)

total bolus estimate =  $\frac{\text{food (grams)}}{\text{carb ratio}}$

Food estimate:

Carb grams ÷ Carb ratio = Units of insulin



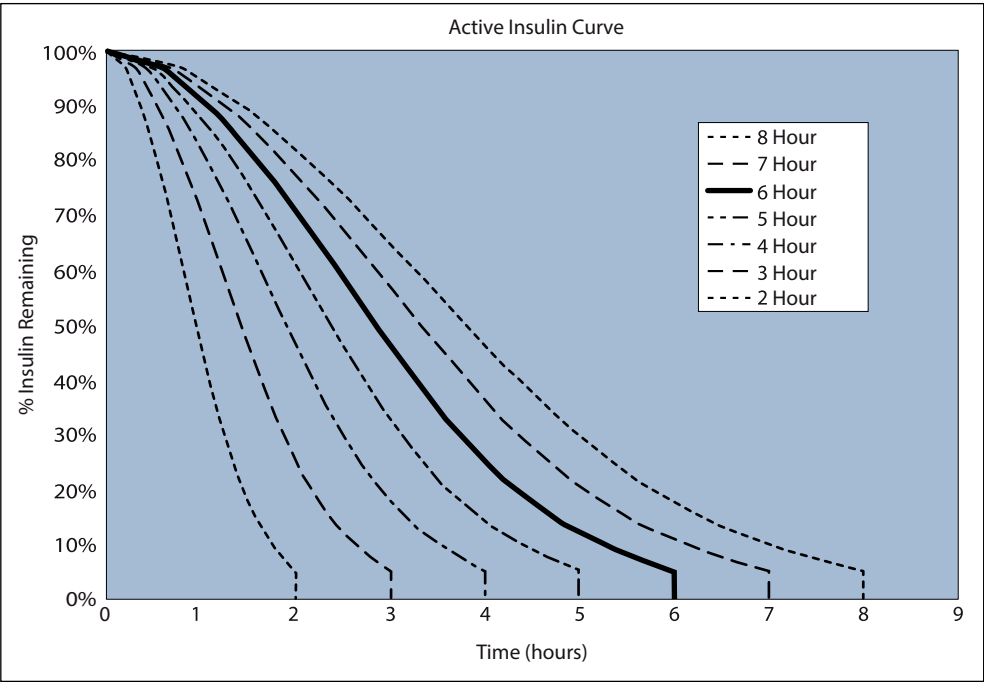
**Note:** When the current BG reading is below the Low BG Target, an active insulin amount is not considered in the Bolus Wizard feature calculations.

Total bolus estimate = Food estimate

4. If no BG reading is entered, the total bolus estimate is based only on the food estimate.

The following list includes additional conditions to consider when using the Bolus Wizard feature in Manual mode.

- If a Dual Wave bolus amount is less than the estimate due to the Max bolus limit or a change that is made, the Square portion of the bolus is reduced first.
- Active insulin is the bolus insulin that has been delivered by the pump and is still working to lower glucose levels. In the Bolus Wizard and SmartGuard Bolus feature, the Active Insulin Time setting is used to calculate a correction bolus by subtracting the estimated active insulin from each bolus. This is shown as Active Insulin, or Act. Insulin, on the Home screen, Bolus screen, Manual Bolus screen, Preset Bolus screen, and Daily History screen. This prevents over-infusion of insulin and reduces the risk of hypoglycemia.
- The Bolus Wizard feature may use the current BG reading, carb units, and active insulin to calculate the estimated bolus.
- The Active Insulin Curve graph shows how the Active Insulin Time setting affects the active insulin amount that is subtracted from correction boluses over time. The percentage of insulin remaining changes at varying rates depending on the Active Insulin Time setting.



Graph adapted from Mudaliar and colleagues, Diabetes Care, Volume 22, Number 9, Sept. 1999, page 1501.

Carb ratios

Maximum ratio settings	Range
8	1 to 200 g/U

Easy bolus feature

Use the Easy bolus feature to set up and deliver a normal bolus when the pump is in Sleep mode. This is done using ^ and with the help of sound and vibration cues.

Sound mode range	0 to 20 increments or Max bolus limit, whichever comes first
Vibrate mode range	0 to 20 increments or Max bolus limit, whichever comes first
Default step size	0.1 unit
Adjustable step size	0.1 to 2 units per increment up to Max bolus limit

## Environmental conditions

The MiniMed 780G system is designed to withstand most conditions encountered in daily life. For more details about environmental conditions, such as exposure to magnetic fields and radiation, waterproof capabilities, and extreme temperatures, see *User safety*, page 26.

- Pump storage and transport temperature range without a AA battery is from -4 °F (-20 °C) to 122 °F (50 °C).
- Pump operating temperature range is from 41 °F (5 °C) to 98.6 °F (37 °C).
- Operating air pressure range is from 10.2 psi (700 hPa) to 15.4 psi (1060 hPa).
- Storage and transport air pressure range is from 7.2 psi (496.4 hPa) to 15.4 psi (1060 hPa).
- Relative humidity (RH) range during operation is from 20% to 90%.
- RH range during storage and transport is from 5% to 95%.

## Essential performance

The pump will maintain the following functionalities to avoid under-infusion and over-infusion:

- Delivery accuracy
- Occlusion detection
- Empty reservoir detection
- Detection of power loss
- Pump therapy status—UI component: LCD
- Notification annunciation and display—UI components: piezo-electric speaker, LCD—applies to all features above

## Expected service life

The overall expected service life for the MiniMed 780G insulin pump is four years when used in accordance with this guide.

If there are concerns that the insulin pump may be damaged, contact 24-Hour Technical Support.



For additional information, see *Pump issues*, page 276.

For health-related questions or concerns, consult a healthcare professional.

## Filling the infusion set and cannula

- The cannula can be filled from 0.025 units to 5.1 units, in increments of 0.025 units.
- The standard fill rate is 1.5 units per minute.  
The quick fill rate is 15 units per minute.
- When filling the tubing, a warning occurs at 30 units. A second warning occurs at 40 units indicating that the pump must be rewound.
- Insulin used to fill the infusion set is recorded in the Daily History. This insulin is NOT included in the Total Daily Delivery (TDD) totals on the Summary screen.

## Insulin delivery default settings

### Bolus settings

Item	Default setting	Limits	Increments
Bolus Wizard feature:	Off	-	-
Easy bolus feature:	Off	-	-
Easy bolus step size:	0.1 U	0.1 U to 2 U	-
Bolus increment:	0.10 U	0.025 U 0.05 U 0.10 U	-
Dual/Square bolus:	Off	-	-
Max bolus:	10 U	0 to 25 U (per single bolus)	-
Bolus BG Check Reminder:	Off	0:30 to 5:00	0:30

### Basal settings for Manual mode

Item	Default setting	Limits	Increments
Max Basal Rate	2 U/hr	0–35 U/hr Default limit 10 U/hr	0.025 U for 0.025–0.975 U/hr 0.05 U for 1.00–9.95 U/hr

Item	Default setting	Limits	Increments
			0.1 U for rates of 10.0 U/hr or more
Basal Rate	0.000 U/hr	0.000 U/hr to Max Basal Rate	0.025 U for 0.025–0.975 U/hr 0.05 U for 1.00–9.95 U/hr 0.1 U for rates of 10.0 U/hr or more
Temp Basal Type	Percent	Percent, Rate	N/A
Temp Basal Percent	100%	0–200%	5%
Temp Basal Rate	Current basal rate	0.0 U/hr to Max Basal Rate	0.025 U for 0.025–0.975 U/hr 0.05 U for 1.00–9.95 U/hr 0.1 U for rates of 10.0 U/hr or more

## Low Reservoir reminder

The values are based on amount shown, not actual amount.

Alert range	Increment	Default value
The first reminder occurs at 5 to 50 units. The second reminder occurs at half of the remaining specified amount. The second reminder is automatic and cannot be changed.	1 unit	20 units

## Max bolus

The bolus delivery specifications do not affect the amount of insulin delivered by SmartGuard. These specifications relate to Manual mode only.

Range	0 to 25 units
Default	10 units

## Normal bolus

Range is 0.025 to 25 units of insulin, and limited by the Max bolus setting.

## Percent temp basal

The default value is 100 percent of basal programming. For example, if six units of basal insulin are delivered per day, the default temp basal amount will be six units per day.

Range	0 to 200%
Default	100% of basal programming
Increment	5%

## Program safety checks

A single fault condition causes the pump to suspend insulin delivery. Maximum infusion with a single fault condition is 0.2 units.

## Pump dimensions

The pump dimensions in inches are no greater than 4.0 length x 2.3 width x 1.1 depth.

The pump dimensions in centimeters are no greater than 10.2 length x 5.8 width x 2.8 depth.

## Pump memory

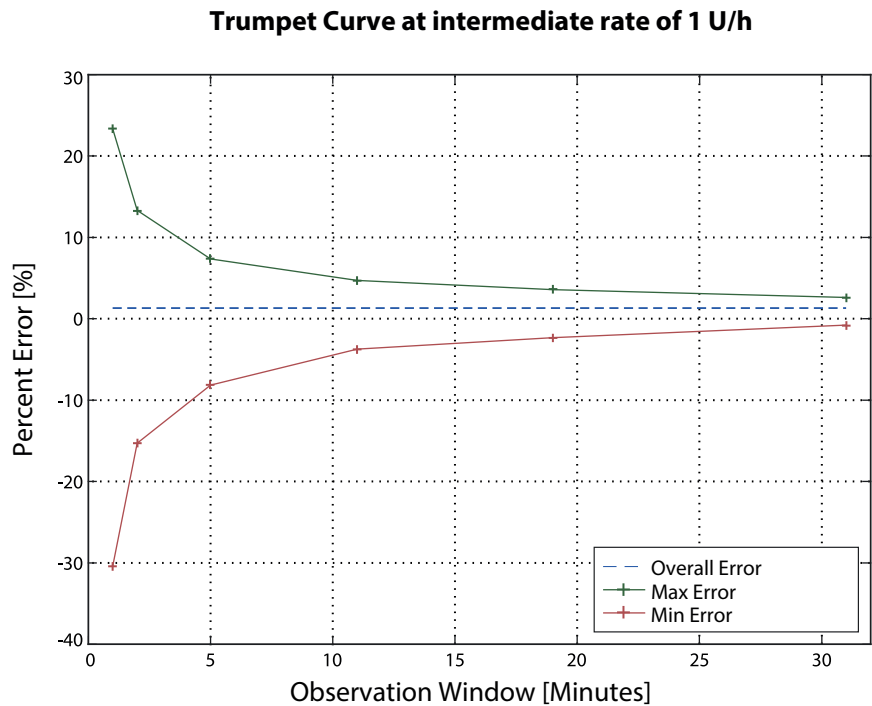
User settings and pump history are stored in pump memory. The pump keeps at least 35 days of history.

## Pump performance characteristics

### Delivery accuracy

- For a basal rate of 1.0 U/hr, the delivery accuracy is  $\pm 5\%$ .  
For a basal rate of 0.025 U/hr, the delivery accuracy is  $\pm 10\%$ .  
Delivery accuracy for bolus volumes  $< 0.1$  unit is  $\pm 20\%$  and delivery accuracy for bolus volumes  $\geq 0.1$  unit is  $\pm 5\%$ .
- All normal boluses are delivered within 16 minutes, 41 seconds  $\pm 3$  seconds at Standard rate (25 units, at 1.5 units per minute), and within 1 minute, 41 seconds  $\pm 3$  seconds at Quick rate (25 units, at 15 units per minute).
- During delivery, the maximum infusion pressure generated and the occlusion threshold pressure using a 3.0-mL reservoir does not exceed 13.15 psi (90.67 kPa). The average resulting bolus volume generated upon clearing the occlusion is 0.0112 mL (equivalent to 1.12 units of U-100 insulin).
- The following image is a representative delivery accuracy curve. The Trumpet Curve represents the maximum percentage change from the expected insulin dosage for a given time interval, known as the observation window, during the

infusion of insulin. The upper curve corresponds to positive changes, and the lower curve corresponds to negative changes.



Infusion pressure

The maximum infusion pressure and occlusion pressure during the fill tubing process are 25 psi (172.4 kPa).

Occlusion detection

When occlusion is detected, the Insulin flow blocked alarm occurs. The occlusion alarm is triggered by an average of 2.98 units of missed insulin (standard bolus) or 2.45 units of missed insulin (quick bolus). This table shows occlusion detection for four different situations when using U-100 insulin.

Rate	Minimum time before alarm	Average time before alarm	Maximum time before alarm
bolus delivery (10 units at standard speed)	77 seconds	124 seconds	173 seconds

Rate	Minimum time before alarm	Average time before alarm	Maximum time before alarm
bolus delivery (10 units at quick speed)	10 seconds	13 seconds	20 seconds
basal delivery (1.0 U/hr)	2 hours 27 minutes	3 hours 21 minutes	4 hours 21 minutes
basal delivery (0.025 U/hr)	131 hours 5 minutes	181 hours 16 minutes	211 hours 30 minutes



**Note:** Certain factors, such as ambient temperature changes or the presence of air in the infusion set or the reservoir, can delay an occlusion alarm.

## Sound frequency

The following table lists audible tones that the pump emits, and their corresponding frequencies:

Tone name	Frequency
Alarm	1655 Hz followed by 3310 Hz
Alternate Alarm	1850 Hz
Siren (escalated alarm)	1655 Hz, followed by 3310 Hz
Alert	934 Hz
High SG	1312 Hz, followed by 1410 Hz, 1500 Hz, 1619 Hz, 1722 Hz
Low SG	1722 Hz, 1619 Hz, 1500 Hz, 1410 Hz, 1312 Hz
Lost SG	1485 Hz, followed by 1395 Hz, 1320 Hz, 1395 Hz
Message tone	1655 Hz
Suspend message tone	2100 Hz, followed by 1800 Hz and 2100 Hz
Reminder tone	934 Hz
Fill tubing tone	1850 Hz
Bolus delivery cancellation tone	1485 Hz, followed by 1655 Hz and 1485 Hz
Loading complete tone	934 Hz
Reservoir loading in progress tone	1850 Hz
Easy bolus activation	1045 Hz
Easy bolus step 1 increment	1175 Hz
Easy bolus step 2 increment	1320 Hz
Easy bolus step 3 increment	1395 Hz

Tone name	Frequency
Easy bolus step 4 increment	1570 Hz
Easy bolus step 5 increment	1760 Hz

## Pump weight

The mass of the insulin pump without battery and consumables is less than 117 grams.

## Sensor default settings

High sensor settings			
Item	Default setting	Limits	Increments
High SG alert limit	250 mg/dL	100 to 400 mg/dL	5 mg/dL
High SG fixed alert	On (cannot be turned off)	250 mg/dL for 3 hours	-
Alert before high	Off	-	-
Alert on high	Off	-	-
Time before high	15 minutes	5 to 30 minutes	5 minutes
Rise Alert	Off	-	-
Rise Limit	Two up arrows	<ul style="list-style-type: none"> <li>1 up arrow (1 mg/dL/min)</li> <li>2 up arrows (2 mg/dL/min)</li> <li>3 up arrows (3 mg/dL/min)</li> <li>Custom limit (1.0 to 5.0 mg/dL/min)</li> </ul>	
High Snooze	1 hour	5 minutes to 3 hours	5 minutes
Low sensor settings			
Item	Default setting	Limits	Increments
Low SG alert limit	60 mg/dL	50 to 90 mg/dL	5 mg/dL
Low SG alarm	On (cannot be turned off)	64 mg/dL	-
Suspend before low	Off	-	-
Suspend on low	Off	-	-

Low sensor settings			
Item	Default setting	Limits	Increments
Alert before low	Off	-	-
Alert on low	Off	-	-
Low Snooze	20 minutes	5 minutes to 1 hour	5 minutes
Resume basal alert	Off	-	-

The SmartGuard feature settings			
Item	Default setting	Limits	Increments
SmartGuard	Off	-	-
Target	100 mg/dL	100 to 120 mg/dL	10 mg/dL
Auto Correction	On	120 mg/dL	-
Temp Target	Off	150 mg/dL	-
Temp Target Duration	2 hours	30 minutes to 24 hours	30 minutes

## IEC 60601-1-2

### IEC 60601-1-2, Special EMC Precautions for Medical Electrical Equipment

1. Special Precautions regarding Electromagnetic Compatibility (EMC): This body worn device is intended to be operated within a reasonable residential, domestic, public or work environment, where common levels of radiated “E” (V/m) or “H” fields (A/m) exist; such as cellular phones that are not paired with the MiniMed 780G system, Wi-Fi networks, Bluetooth wireless technology, electric can openers, microwave and induction ovens. This device generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the provided instructions, may cause harmful interference to radio communications.
2. Portable and mobile RF communications equipment can affect Medical Electrical Equipment as well. If RF interference from a mobile or stationary RF transmitter is encountered, move away from the RF transmitter that is causing the interference.

## IEC 60601-1-2

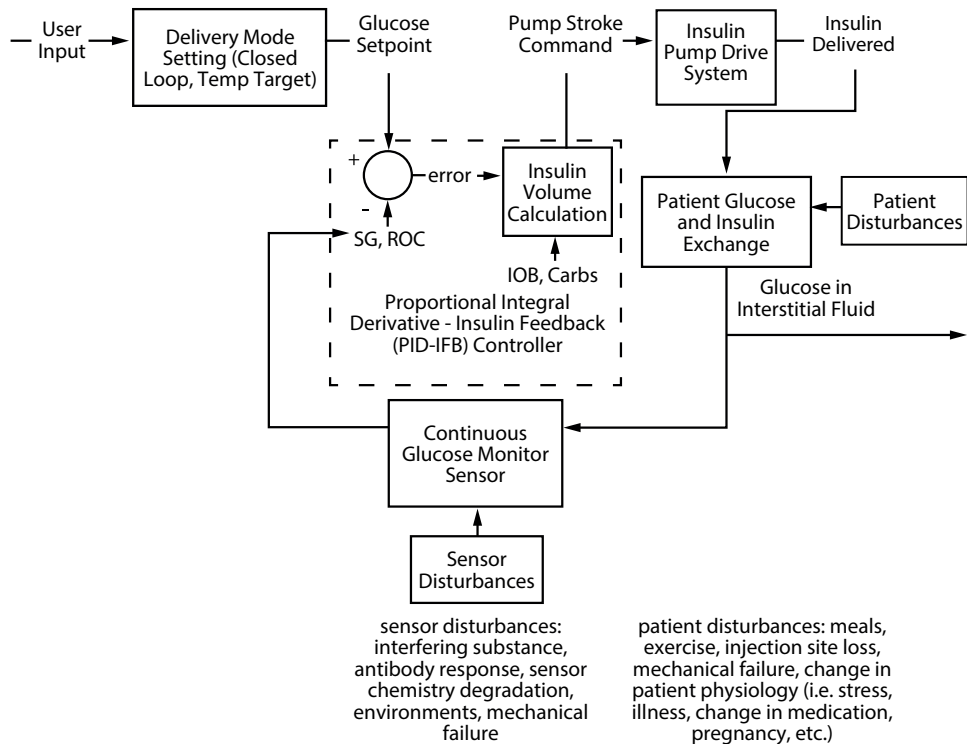
The MiniMed 780G system should not be used adjacent to other electrical equipment. If adjacent use becomes necessary, the MiniMed 780G system should be observed to confirm normal system operation.

## IEC 60601-1-10: PCLCS

The MiniMed 780G is a Physiological Closed-Loop Controlled system (PCLCS).

Auto Mode manages basal delivery using a closed loop control algorithm based on a Proportional Integral Derivative controller with insulin feedback (PID-IFB). The PID-IFB monitors the Rate Of Change (ROC) of sensor glucose (SG) and calculates the insulin volume using the Insulin On Board (IOB) and the reported Carbs. The closed loop controller uses continual feedback of SG values to calculate the insulin delivery rate for basal insulin control. The control algorithm is part of the pump application code. SG values are received by the pump via RF from the CGM sensor. This theory of operation is described in the following block diagram.





## Guidance and manufacturer's declaration


Guidance and Manufacturer's Declaration - Electromagnetic Emissions		
The MiniMed 780G insulin pump is intended for use in the electromagnetic environment specified below. Make sure that the MiniMed 780G insulin pump is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions Test: 47 CFR Part 15, Subpart C Section 15.247/FCC Part 15 Subpart B Section 15.109	<ul style="list-style-type: none"><li>6 dB and 99% Bandwidths: Complies</li><li>Maximum Output Power: Complies</li><li>TX Spurious Emissions: Complies</li></ul>	The MiniMed 780G insulin pump must emit electromagnetic energy in order to perform its intended function. Nearby electronic equipment may be affected.

Guidance and Manufacturer's Declaration - Electromagnetic Emissions		
	<ul style="list-style-type: none"> <li>Power Spectral Density: Complies</li> <li>Radiated Emissions at Band Edge: Complies</li> </ul>	
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Not applicable	
RF emissions CISPR 11	Complies Group 1 Class B	The MiniMed 780G insulin pump is suitable for use in aircraft and in all establishments, including domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RTCA DO 160G Radio Frequency Susceptibility (Radiated and Conducted) and Emission of Radio Frequency Energy	Complies	

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
<b>The MiniMed 780G insulin pump is intended for use in the electromagnetic environment specified below. Make sure that the MiniMed 780G insulin pump is used in such an environment.</b>			
Immunity Test	IEC 60601-1-2 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2, 4, 8, 15 kV air	±8 kV contact ±2, 4, 8, 15 kV air	For use in a typical domestic, commercial, or hospital environment.
Conducted disturbances induced by RF fields IEC 61000-4-6	3 V <sub>RMS</sub> 150 kHz to 80 MHz 6 V <sub>RMS</sub> ISM bands between 150 kHz to 80 MHz	Not applicable	Requirement does not apply to this battery powered device.
Electrical fast transient/burst IEC 61000-4-4	±2 kV 100 kHz repetition frequency	Not applicable	Requirement does not apply to this battery powered device.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
Surge IEC 61000-4-5	Line to Line: $\pm 0.5$ kV, $\pm 1$ kV Line to Ground: $\pm 0.5$ kV, $\pm 1$ kV, $\pm 2$ kV	Not applicable	Requirement does not apply to this battery powered device.
Voltage dips, short interruptions, and voltage variations on power supply lines IEC 61000-4-11	0% $U_T$ ; 0.5 cycle (at 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315°) 0% $U_T$ ; 1 cycle (at 0°) 70% for 25/30 cycles (at 0°) 0% for 250/300 cycles	Not applicable	Requirement does not apply to this battery powered device.
Power frequency (50/60 Hz) electromagnetic field IEC 61000-4-8	30 A/m (continuous field at 60 seconds)	30 A/m 400 A/m per IEC 60601-2-24	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Proximity fields from RF wireless communications equipment IEC 61000-4-39	IEC 60601-1-2	IEC 60601-1-2	For use in a typical domestic, commercial, or hospital environment.
<b>Note:</b> $U_T$ is the a.c. mains voltage prior to application of the test level.			

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
<b>The MiniMed 780G insulin pump is intended for use in the electromagnetic environment specified below. The customer or user of the MiniMed 780G insulin pump should assure that it is used in such an electromagnetic environment.</b>			
Immunity Test	IEC 60601-1-2 Test Level	Compliance Level	Electromagnetic Environment Guidance
Radiated RF IEC 61000-4-3 EN 301 489-17	10 V/m 80 MHz to 2.7 GHz	10 V/m 80 MHz to 2.7 GHz	Portable and mobile RF communications equipment should be used no closer to any part of the MiniMed 780G insulin pump,

Guidance and Manufacturer’s Declaration - Electromagnetic Immunity			
	80% AM at 1 kHz	80% AM at 1 kHz	<p>including cables, than the recommended separation distance of 12 in (30 cm). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol:</p> 

Wireless communication

The MiniMed 780G insulin pump communicates using smart device connectivity.

Operating frequency/Modulation type(s)	2.4 GHz band, GFSK
Effective radiated power (ERP)	1.48 mW (1.69 dBm)
Effective isotropic radiated power (EIRP)	2.42 mW (3.83 dBm)

FCC notice

This device complies with the United States Federal Communications Commission (FCC) and international standards for electromagnetic compatibility. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. These standards are designed to provide reasonable protection against excessive radio frequency interference, and to prevent undesirable operation of the devices from unwanted electromagnetic interference.



**Note:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Decrease the distance between the transmitter and the insulin pump to 6 feet (1.8 meters) or less.
- Decrease the distance between the meter and the insulin pump to 6 feet (1.8 meters) or less.
- Increase the separation between the transmitter and the device that is receiving/emitting interference.

**IMPORTANT:** Do not change or modify the internal RF transmitter or antenna unless expressly approved by Medtronic Diabetes. Doing so could interfere with your ability to operate the equipment.



**Note:** Harmful interference is defined by the FCC as follows. Any emission, radiation or induction that endangers the functioning of a radio navigation service or of other safety services or seriously degrades, obstructs or repeatedly interrupts a radio communications service operating in accordance with FCC rules.

## Open Source Software disclosure

This document identifies the Open Source Software that may be separately called, executed, linked, affiliated, or otherwise utilized by this product.

Such Open Source Software is licensed to users subject to the terms and conditions of the separate software license agreement for such Open Source Software.

Use of the Open Source Software shall be governed entirely by the terms and conditions of such license.

The source and object code, and applicable license for any Open Source Software can be obtained at the following site(s):

- LZ4-compression library (v1.9.1): <http://www.lz4.org>
- SWIG (v3.0.12): <http://www.swig.org>
- FNV-1 hash algorithm (v5.1): <http://www.isthe.com/chongo/tech/comp/fnv/> and <http://www.isthe.com/chongo/src/fnv/fnv64.c>
- CRC32 algorithm:  
<https://opensource.apple.com/source/xnu/xnu-792.13.8/bsd/libkern/crc32.c>





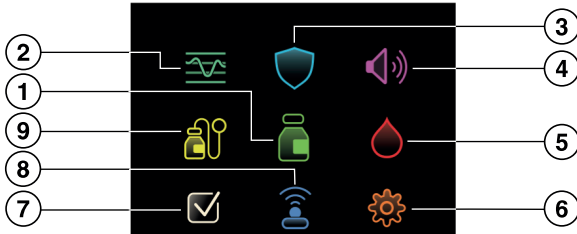


# Appendix C: Menu map

## Menu map

The following diagrams provide a map to the screens and features that are available from the Menu screen.

<div> <div>1</div> <div></div> </div> <p><b>Insulin</b></p> <p>Bolus/Stop Bolus*</p> <ul style="list-style-type: none"> <li>• Bolus Wizard*</li> <li>• Manual Bolus</li> <li>• Preset Bolus*</li> <li>• ➔ Delivery Settings</li> </ul> <p>Basal/Cancel Temp Basal*</p> <ul style="list-style-type: none"> <li>• Temp Basal</li> <li>• Preset Temp*</li> <li>• Basal Patterns</li> <li>• ➔ Basal Pattern Setup</li> </ul> <p>Suspend All Delivery/Resume Basal*</p> <p>➔ Delivery Settings</p>	<div> <div>2</div> <div></div> </div> <p><b>History &amp; Graph</b></p> <p>History</p> <ul style="list-style-type: none"> <li>• Summary</li> <li>• Daily History</li> <li>• Alarm History</li> <li>• Paired Sensors</li> </ul> <p>Sensor Glucose Review</p> <ul style="list-style-type: none"> <li>• High Limit</li> <li>• Low Limit</li> <li>• Days to Average</li> </ul> <p>Graph</p> <p>Time in Range</p>	<div> <div>3</div> <div></div> </div> <p><b>SmartGuard</b></p> <p>SmartGuard Checklist</p> <p>Temp Target/Cancel Temp Target*</p> <p>SmartGuard Settings</p> <ul style="list-style-type: none"> <li>• Target</li> <li>• Auto Correction</li> </ul> <p>SmartGuard (On and Off)</p> <div> <div>4</div> <div></div> </div> <p><b>Sound &amp; Vibration</b></p> <p>Silence Sensor Alerts</p> <p>Volume</p> <p>Sound</p> <p>Vibration</p> <p>➔ Alert Settings</p> <div> <div>5</div> <div></div> </div> <p><b>Blood Glucose</b></p> <p>BG</p>
<p>* indicates an item that is shown when applicable. For example, Stop Bolus is shown when a bolus is being delivered. Preset Bolus is shown when the Preset Bolus feature is set up.</p> <p>➔ This icon indicates a shortcut to the screen.</p>		

		
<div> <div>6</div> <div></div> </div> <p><b>Settings</b></p> <p>Alert Settings</p> <ul style="list-style-type: none"> <li>• High Alert</li> <li>• Low Alert</li> <li>• Snooze High &amp; Low</li> <li>• Reminders</li> </ul> <p>Delivery Settings</p> <ul style="list-style-type: none"> <li>• Bolus Wizard Setup</li> <li>• Basal Pattern Setup</li> <li>• Max Basal/Bolus</li> <li>• Dual/Square Wave</li> <li>• Bolus Increment</li> <li>• Bolus Speed</li> <li>• Preset Bolus Setup</li> <li>• Preset Temp Setup</li> </ul>	<div> <div>6</div> <div></div> </div> <p><b>Settings continued</b></p> <p>Device Settings</p> <ul style="list-style-type: none"> <li>• Sensor (On and Off)</li> <li>• Time &amp; Date</li> <li>• Device Info</li> <li>• Display</li> <li>• Block Mode</li> <li>• Self Test</li> <li>• Review Settings</li> <li>• Manage Settings</li> <li>• Easy Bolus</li> <li>• Auto Suspend</li> <li>• Language</li> </ul>	<div> <div>7</div> <div></div> </div> <p><b>Status</b></p> <p>Stop Bolus*</p> <p>Suspend All Delivery/Resume Basal*</p> <p>SmartGuard Checklist</p> <p>Pump</p> <p>Sensor</p> <div> <div>8</div> <div></div> </div> <p><b>Paired Devices</b></p> <p>Pair New Device</p> <p>Pair CareLink</p> <p>Mobile</p> <p>Meter</p> <p>CGM</p> <div> <div>9</div> <div></div> </div> <p><b>Reservoir &amp; Set</b></p> <p>New Reservoir &amp; Set</p> <p>New Reservoir Only</p> <p>New Set Only</p> <p>Fill Cannula</p>
<p>* indicates an item that is shown when applicable. For example, Resume Basal is shown when insulin delivery is suspended.</p>		







# Appendix D: Performance data

## MiniMed 780G system device performance

The MiniMed 780G system adjusts insulin delivery based on sensor glucose (SG) readings from CGM, while alleviating the complexity of trying to maintain glucose levels around meals. Clinical studies have shown that integrated insulin pump and CGM systems may provide better diabetes management, compared with multiple daily injections, or with a pump alone. Studies suggest that pump therapy, when regulated by sensor information, can improve HbA1C levels significantly without increasing the risk of hypoglycemia.<sup>1,2,3</sup>

The MiniMed 780G system continues to use the SmartGuard feature, which is designed to keep patient blood sugar levels in range by automatically adjusting basal insulin dosage every five minutes, delivering more or less insulin when it predicts that SG values are trending too high or too low. The SmartGuard feature has been updated to adjust how auto correction boluses and daily user adaptations are calculated. The system offers the following features:

---

<sup>1</sup> Bergenstal RM, Tamborlane WV, Ahmann A, et al. Effectiveness of sensor-augmented insulin-pump therapy in type 1 diabetes [STAR3 Study]. *N Engl J Med*.2010;363:311–320.

<sup>2</sup> Battelino T, Conget I, Olsen B, et al. The use and efficacy of continuous glucose monitoring in type 1 diabetes treated with insulin pump therapy [SWITCH study]. *Diabetologia*. 2012 Dec;55(12):3155-62. doi: 10.1007/s00125-012-2708-9. Epub 2012 Sept 11.

<sup>3</sup> Bergenstal RM, Klonoff DC, Bode BW, et al. Threshold-based insulin-pump interruption for reduction of hypoglycemia [ASPIRE in-home study]. *N Engl J Med*. 2013;369(3):224-232.

1. **Adjustable glycemic target settings.** With the help of a healthcare provider, patients can program the device to one of three setpoints to target their ideal SG value (100, 110, or 120 mg/dL). The device uses the programmed setpoint as a reference to adjust the rate of insulin delivered, which helps maintain control according to patient needs.
2. **Automatic correction boluses.** Mealtimes can be stressful and require that patients calculate boluses prior to and after meals to avoid hyperglycemia. The SmartGuard feature also includes an Auto correction feature that can calculate and deliver correction boluses every five minutes if the patient underestimates the amount of carbs in a meal or if they accidentally forget to deliver a meal bolus prior to eating.

The MiniMed 780G system retains the Suspend on low and Suspend before low features that were introduced in prior Medtronic insulin pumps. These features temporarily stop insulin delivery when SG values reach a preset low target (Suspend on low) or are predicted to reach the preset low target within 15 or 30 minutes (Suspend before low). Insulin delivery also resumes when SG values return to a safe range. These optional features are available when the pump is in Manual mode and function as a backup for the SmartGuard feature.

## The SmartGuard feature

### Clinical study overview

The SmartGuard feature with modifications (that controls insulin dosing in the MiniMed 780G system) was studied with subjects who wore the MiniMed 780G pump with the Simpler Sync sensor and the Extended infusion set and reservoir at home for 3 months.<sup>4</sup> The study did not include a control group. The study included subjects from different clinics around the US who were between 7 and 80 years old. Subjects had to have been diagnosed with type 1 diabetes mellitus for at least one year for subjects aged 7 to 13 years, and at least two years for subjects aged 14 to 80 years. All subjects in the study had to have used pump therapy for at least 6 months prior to screening and had an HbA1C value of less than 10.0% at the time of screening.

---

<sup>4</sup> Medtronic Inc., Clinical Study Report: CIP337 Safety and Effectiveness Evaluation of the MiniMed 780G System Used in Combination with the DS5 CGM.



This study started with a run-in (baseline) period. During run-in, subjects with no prior automated insulin delivery (AID) pump experience were instructed to use the MiniMed 780G system with only the sensor augmented pump (SAP) function activated (i.e., SmartGuard feature turned OFF). Subjects with AID pump experience were instructed to use the MiniMed 780G system with the SmartGuard feature turned ON and the Auto correction feature turned OFF. The intent of the run-in period was to allow subjects to become familiar with the new study devices while using their own insulin.

After the run-in period, subjects were instructed to use the study devices with both the SmartGuard feature and the Auto Correction feature turned ON during a study period comprising 3 stages. In the first two stages, subjects were instructed to use the study pump with the 120 mg/dL Auto Basal target setpoint and active insulin time set to 4 hours (stage 1), then to change the pump settings to the 100 mg/dL setpoint and active insulin time set to 2-3 hours (stage 2). In stage 3, subjects were instructed to use the study pump with the Auto Basal target setpoint and active insulin time set as considered best by the investigator for the individual subject. A total of 250 subjects were enrolled, and 212 subjects completed the study.

The updated MiniMed 780G system is compatible with the Guardian 4 Sensor and Guardian 4 Transmitter, but the clinical study tested a similar next-generation CGM known as the Simplerla Sync devices.

SmartGuard Performance: HbA1C and Time in Target Range

Table 1 shows the HbA1C from Baseline to the end of the 3-month study. This data helps explain how using the 780G SmartGuard feature with the Auto Correction feature enabled might affect a patient’s HbA1C.

Table 1. HbA1C from Baseline to End of 3-month Study Period

Category	Age 7-17 Years		Age 18-80 Years	
	Baseline	End of Study	Baseline	End of Study
HbA1C (%) Mean ± SD (Median) [N]	7.7 ± 1.0 (7.8) [112]	7.3 ± 0.8 (7.2) [111]	7.4 ± 0.9 (7.3) [110]	6.7 ± 0.5 (6.7) [106]

Table 2 reports the mean percentage of time spent in range (TIR, 70-180 mg/dL) in Stage 3 of the Study Period.

**Table 2.** Mean Percentage of Time Spent in Range (70-180 mg/dL) in Study Period Stage 3

Subject Age	Number of Subjects	Mean	95% Confidence Interval
7-17 Years	109	71.4%	(69.5%, 73.3%)
18-80 Years	107	80.2%	(78.7%, 81.8%)

Safety

Table 3 lists the device-related adverse events reported at screening, and during the run-in and study periods. Overall, 59 adverse events were reported. For subjects ages 7–17 years, there were no reports of device-related severe hypoglycemia, unanticipated serious or non-serious adverse device effects during the study. For adult subjects, there were no reports of unanticipated adverse device effects. Two severe hypoglycemia events and 1 diabetic ketoacidosis event were reported for adult subjects but none of these were device-related.

**Table 3.** Device Related Adverse Events

Adverse Events	Age 7–17 Years (N = 125)			Age 18–80 Years (N = 125)		
	Screening period	Run-in period	Study period	Screening period	Run-in period	Study period
Bleeding at infusion site	0	0	1	0	0	0
Bleeding at sensor site	0	2	2	0	1	1
Discomfort/irritation with infusion set	0	2	1	0	6	5
Discomfort/irritation with sensor	0	3	2	0	1	1
Infusion site infection	0	4	15	0	0	0
Mild ketonemia	0	0	0	0	1	0
Rash or contact dermatitis (infusion set related)	0	0	1	0	2	2
Rash or contact dermatitis (sensor or tape related)	0	0	1	0	0	0
Severe hyperglycemia	0	1	3	0	0	1
Total	0	12	26	0	11	10

Table 4 lists the study-period related adverse events. A total of 83 adverse events during the study period and one serious adverse event were reported from all investigational sites for 7–17-year-old study subjects enrolled in the study. There were 0 serious adverse events, no reports of severe hypoglycemia, 8 reports of severe hyperglycemia, no reports of diabetic ketoacidosis, and no reports of unanticipated adverse device effects (UADEs).

A total of 50 adverse events during the study period and three serious adverse events were reported from all investigational sites for 18–80-year-old study subjects enrolled in the study. Out of 50 events, there were 3 serious adverse events, 2 reports of

hypoglycemia and 1 report of a hyperglycemia event, 1 report of diabetic ketoacidosis events, and no reports of unanticipated adverse device effects (UADEs).

**Table 4.** Study Period Related Adverse Events

Category	Age 7–17 Years (N = 112)	Age 18–80 Years (N = 110)
	<b>Study Period</b>	<b>Study Period</b>
Total number of adverse events	83	50
Study Exit		
Led to study exit	0	1
Did not lead to study exit	83	49
Seriousness		
Serious adverse events (SAEs)	0	3
Death	0	0
Non-death	0	3
Non-serious adverse events	83	47
Diagnosis		
Severe hypoglycemia	0	2
Severe hyperglycemia	8	1
Diabetic ketoacidosis (DKA)	0	1
None of the above	75	46
Study procedure and device relatedness		
Related to study procedure only	0	0
Related to study device only	26	10
Unanticipated adverse device effects (UADE) / Unanticipated serious adverse device effects (USADE)	0	0
Unanticipated non-serious adverse device effects	0	0
Anticipated adverse device effects (ADEs)	26	10
Related to both study procedure and study device	0	0

**Table 4.** Study Period Related Adverse Events (continued)

Category	Age 7–17 Years (N = 112)	Age 18–80 Years (N = 110)
Unanticipated adverse device effects (UADE) / Unanticipated serious adverse device effects (USADE)	0	0
Unanticipated non-serious adverse device effects	0	0
Anticipated adverse device effects (ADEs)	0	0
Not related to study procedure or study device	57	40

**SmartGuard Use**

During the study period, subjects had access to the study device and were instructed to use SmartGuard with Auto Correction ON. *Table 5* presents the percentage of time that subjects spent using the sensor and the percentage of time spent using the SmartGuard feature with the Auto correction feature turned ON. This information shows that the SmartGuard feature was ON greater than 92% of the time.

**Table 5.** Sensor and Auto Mode Usage (Percentage of Time) During Study Period, Stage 3

Category	Age 7-17 Years (N = 111)	Age 18-80 Years (N = 110)
Time spent using sensor	92.9%	95.8%
Time spent not using sensor	7.1%	4.2%
Time spent in SmartGuard	93.5%	96.6%
Time spent in Manual mode	6.5%	3.4%

**SmartGuard Performance**

*Table 6* shows the mean percentage of SG values in specific glucose ranges during the run-in period and during stage 3 of the study period by all subjects using the 780G system with the Simpler Sync sensor. An international group of diabetes experts and the American Diabetes Association (ADA) consider patients in good control when

patients are in the target glucose range of 70–180 mg/dL for more than 70% of the day.<sup>5</sup>

The data in *Table 6* show that using the SmartGuard feature with the Auto correction feature kept SG values in range and reduced time above range. Specifically, adult subjects spent more time in range (70–180 mg/dL) and less time in hypoglycemia (<70 mg/dL) and hyperglycemia (>180 mg/dL) during stage 3 of the study period compared with the run-in period. Pediatric subjects spent more time in range (70–180 mg/dL) and less time in hyperglycemia (>180 mg/dL) without significantly increasing time in hypoglycemia (<70 mg/dL) during stage 3 of the study period compared with the run-in period.

**Table 6.** Percentage of SG values in Different Ranges during the Run-In Period and Study Period Stage 3

Category	SG Range (mg/dL)	Age 7-17 Years		Age 18-80 Years	
		Run-in Period (N = 112)	Study Period Stage 3 (N = 109)	Run-in Period (N = 110)	Study Period Stage 3 (N = 107)
Low SG Value	<54	0.3 ± 0.6 (0.2, 0.4)	0.4 ± 0.3 (0.3, 0.4)	0.3 ± 0.5 (0.2, 0.4)	0.2 ± 0.4 (0.1, 0.3)
	<70	1.6 ± 1.7 (1.3, 1.9)	1.9 ± 1.4 (1.7, 2.2)	1.7 ± 1.9 (1.4, 2.1)	1.5 ± 1.4 (1.3, 1.8)
Target SG Value	70 – 140	32.1 ± 14.1 (29.5, 34.7)	49.2 ± 9.7 (47.4, 51.0)	39.2 ± 13.0 (36.8, 41.7)	56.1 ± 10.5 (54.1, 58.1)
	70 – 180	54.4 ± 15.7 (51.5, 57.3)	71.4 ± 9.9 (69.5, 73.3)	66.5 ± 12.6 (64.1, 68.8)	80.2 ± 8.1 (78.7, 81.8)
High SG Value	>140	66.3 ± 14.7 (63.5, 69.0)	48.9 ± 10.0 (47.0, 50.8)	59.1 ± 13.9 (56.4, 61.7)	42.4 ± 11.0 (40.3, 44.5)
	>180	44.0 ± 16.1 (41.0, 47.0)	26.7 ± 10.1 (24.7, 28.6)	31.8 ± 13.1 (29.4, 34.3)	18.2 ± 8.4 (16.6, 19.9)
	>250	16.4 ± 11.1 (14.3, 18.5)	8.0 ± 6.6 (6.8, 9.3)	7.4 ± 6.1 (6.2, 8.5)	3.4 ± 3.0 (2.8, 4.0)
	>350	2.4 ± 3.5 (1.8, 3.1)	1.3 ± 2.2 (0.9, 1.8)	0.4 ± 0.7 (0.3, 0.5)	0.3 ± 0.5 (0.2, 0.4)

Note: Values are presented by Mean ± SD (95% CI) except Number of subjects.

*Table 7* shows the difference in mean sensor glucose from baseline to the end of the study period for all subjects using the 780G system with the Simplerla Sync sensor. The data shows that, compared to the run-in period, the subjects' mean glucose levels were closer to the center of the target euglycemic range.

<sup>5</sup> Battelino T, Danne T, Bergenstal RM, et al. Clinical Targets for Continuous Glucose Monitoring Data Interpretation: Recommendations From the International Consensus on Time in Range. *Diabetes Care*. 2019 Aug; 42(8): 1593-1603. doi: 10.2337/ dci19-0028. Epub 2019 Jun 8.

**Table 7.** Difference in Mean Sensor Glucose Values (mg/dL) between the Run-In Period and Study Period Stage 3

Category	Age 7-17 Years			Age 18-80 Years		
	Run-in Period (N = 112)	Study Period Stage 3 (N = 109)	Difference be- tween Run-in Pe- riod and Study Period Stage 3 (N = 109)	Run-In Period (N = 110)	Study Period Stage 3 (N = 107)	Difference be- tween Run-in Pe- riod and Study Period Stage 3 (N = 107)
Mean Glucose Value	180.4 ± 27.1 (175.3, 185.4)	154.4 ± 17.6 (151.0, 157.7)	-26.2 ± 22.2 (-30.4, -22.0)	161.0 ± 18.7 (157.5, 164.5)	142.2 ± 12.8 (139.7, 144.7)	-18.5 ± 14.0 (-21.2, -15.8)
Note: Values are presented by Mean ± SD (95% CI) except Number of subjects.						

During the study period, some subjects wore the study pump with the SmartGuard feature and the Auto correction feature turned ON, and with the target setpoint set to either 100 mg/dL, 110 mg/dL, 120 mg/dL, or 150 mg/dL (Temp Target) for at least an entire day. *Table 8* shows the mean sensor glucose (SG) value for each target setpoint option when that setpoint was used for the entire day during the overall study period. The data in *Table 8* shows that using the SmartGuard feature with the Auto correction feature and the 100 mg/dL target setpoint resulted in a lower mean SG value than when the features were used with the 120 mg/dL target setpoint.

**Table 8.** Mean Sensor Glucose Values (mg/dL) during SmartGuard Use Stratified by Target Glucose Setpoint during the Study Period

Category	Age 7– 17 Years					Age 18–80 Years				
	Overall (N = 112)	Target Glucose (mg/dL)				Overall (N = 109)	Target Glucose (mg/dL)			
		100 (N = 109)	110 (N = 12)	120 (N = 111)	150 (N = 52)		100 (N = 107)	110 (N = 5)	120 (N = 108)	150 (N = 48)
Mean Glucose Values During SmartGuard	153.6 ± 14.4 (150.9, 156.3)	151.9 ± 15.0 (149.1, 154.8)	149.5 ± 16.5 (139.0, 160.0)	157.8 ± 14.6 (155.1, 160.6)	157.3 ± 44.4 (145.0, 169.7)	143.8 ± 12.2 (141.4, 146.1)	141.0 ± 11.9 (138.7, 143.3)	139.8 ± 11.2 (125.9, 153.7)	150.5 ± 12.4 (148.1, 152.8)	137.5 ± 29.0 (129.1, 145.9)
Note 1: Values are presented by Mean ± SD (95% CI).										
Note 2: Analysis of data was only performed when SmartGuard Glucose target was used the entire day (e.g., 100 mg/dL set point used for entire day versus 110 mg/dL set point used for entire day versus 120 mg/dL set point used for entire day). Any day with partial usage was excluded from this analysis.										

*Figure 1* shows the percentage of subjects that had an HbA1C that was less than 7% during the run-in (baseline) and study periods. The ADA considers a HbA1C target of less than 7% appropriate for non-pregnant adults and many children.<sup>6,7</sup> *Figure 1* shows

that a greater percentage of subjects had an HbA1C that was less than 7% at the end of the study than at baseline.

**Figure 1.** Percentage of Patients with less than 7% HbA1C

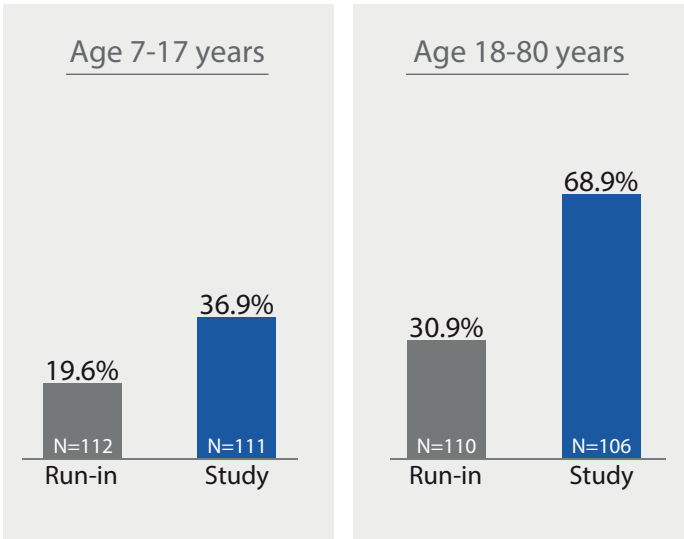


Table 9 shows the change in total daily dose of insulin (TDD) from the run-in period to study period stage 3, and the change in weight and BMI Z-score (for pediatrics) from baseline to the end of the study. Mean TDD increased for both pediatric and adult subjects. Mean weight increased slightly for pediatric subjects and remained unchanged for adults. In the pediatric population, 33 subjects gained more than 2.5 kg (5.5 lbs) in weight over the 3-month study period, and of these, 13 subjects gained 5 kg (11 lbs) or more. This data helps explain how using the SmartGuard feature with the

<sup>6</sup> American Diabetes Association Professional Practice Committee; 6. Glycemic Goals and Hypoglycemia: Standards of Care in Diabetes—2024. Diabetes Care 1 January 2024; 47 (Supplement\_1): S111–S125. <https://doi.org/10.2337/dc24-S006>

<sup>7</sup> American Diabetes Association Professional Practice Committee; 14. Children and Adolescents: Standards of Care in Diabetes—2024. Diabetes Care 1 January 2024; 47 (Supplement\_1): S258–S281. <https://doi.org/10.2337/dc24-S014>

Auto correction feature, childhood and pubertal growth, and elevated glucose levels may affect a patient’s TDD and weight.

**Table 9.** Changes in Mean TDD and Weight

Category	Age 7-13 Years		Age 14-17 Years		Age 18-80 Years	
TDD (U), Mean ± SD (Median)	Run-in Period (N = 57)	Study Period Stage 3 (N = 55)	Run-in Period (N = 55)	Study Period Stage 3 (N = 54)	Run-in Period (N = 110)	Study Period Stage 3 (N = 107)
	43.2 ± 24.0 (35.7)	50.3 ± 29.7 (40.6)	64.3 ± 23.7 (59.0)	75.0 ± 29.3 (71.9)	54.7 ± 27.1 (50.9)	57.8 ± 28.0 (50.0)
Weight* (kg), Mean ± SD (Median)	Baseline (N = 57)	End of Study (N = 56)	Baseline (N = 55)	End of Study (N = 55)	Baseline (N = 110)	End of Study (N = 108)
	47.1 ± 17.9 (44.1)	49.0 ± 19.2 (45.8)	68.7 ± 14.0 (66.7)	70.3 ± 14.6 (68.9)	84.8 ± 19.5 (82.3)	84.8 ± 19.3 (82.0)
BMI Z-score*, Mean ± SD (Median)	Baseline (N = 57)	End of Study (N = 56)	Baseline (N = 55)	End of Study (N = 55)	—	—
	0.6 ± 1.0 (0.6)	0.7 ± 1.1 (0.7)	0.5 ± 1.0 (0.6)	0.6 ± 1.0 (0.8)	—	—


\*Note: Weight and height were not collected in-clinic for some subjects.

Figure 2 shows the percentage of subjects that spent more than 70% of time in range (70-180 mg/dL), which is considered good glucose control by diabetes experts and the ADA, during the run-in (baseline) and the different study periods. The system offers three SG target setpoint options that allow users to customize insulin delivery. For the study period, percentages are shown for subjects that used the SmartGuard feature with the Auto correction feature.

In both pediatric and adult patients, the percentage of subjects who spent more than 70% of time in range when using the SmartGuard feature with the Auto correction feature increased in study period stage 3 from the run-in period. This data shows that more subjects using the SmartGuard feature with the Auto correction feature at the 100 mg/dL setpoint spent more time in range.



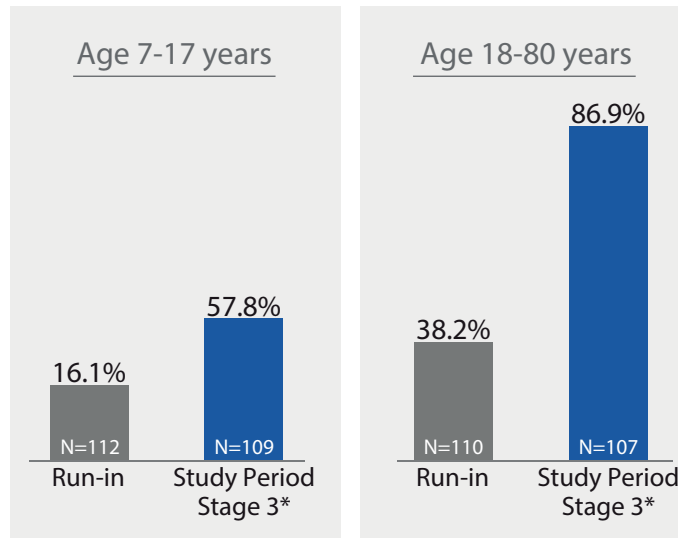**Figure 2.** Percentage of Subjects who spent More than 70% of Time in Range (70-180 mg/dL)



\*During this period, subjects were instructed to use the study device with the Auto basal target as well as Active insulin time (AIT) set to what is best for the individual subject, at the investigator's discretion. Note that 38% (41/109) of pediatric patients and 41% (44/107) of adult patients used the study device with the Auto Basal target with setpoint 100 mg/dL and AIT set to 2 hours during study period stage 3.

Overall, the clinical study suggested that the 780G system was safe, and subjects showed improvements in HbA1C (compared to baseline) and time in the target euglycemic range with use of the updated 780G pump and the Simpler Sync sensor. However, the study had the following limitation:

It did not compare subjects who were using the SmartGuard feature and the Auto correction feature to those who were not on a system (control group). Instead, the study compared how subjects did before using the Auto correction feature (run in period -2 weeks) against results while using the Auto correction feature (study period -3 months).

Due to this limitation, the results of the study should be interpreted with caution and you should understand that your individual results may vary.

## The Suspend before low feature

### Clinical study overview (Ages 14-75 Years)

The Suspend before low feature was evaluated for safety in a multi-center, single-arm, in-clinic study of the MiniMed 640G System.<sup>8</sup> This feature is the same in the MiniMed 780G system. Study subjects included persons aged 14 to 75 years diagnosed with type 1 diabetes mellitus who were on pump therapy at the time of screening.

A total of 71 subjects were subjected to hypoglycemic induction, followed by an observation period. For hypoglycemic induction, the target was set to 65 mg/dL, using the rate of change basal increase algorithm. The Suspend before low feature was activated with the Low Limit setting for the Suspend before low feature ON set to 65 mg/dL, and the subject was observed with frequent sample testing (FST, or frequent blood sampling for glucose measurements) for a maximum of 19 hours. The observation period included the suspension period, the insulin resumption period, and if applicable, an insulin resuspension after basal insulin delivery resumed.

### Feature performance and safety

Of the 71 subjects with induced hypoglycemia, 69 inductions were successful, 27 subjects experienced a hypoglycemic event and 42 subjects did not. At 120 minutes after the start of the pump suspension events, the mean reference glucose value (measured using a Yellow Springs Instrument [YSI]) was  $102 \pm 34.6$  mg/dL.

Five adverse events were reported during the study. Four adverse events were neither device nor procedure related. One adverse event was procedure related.

Data from this in-clinic study demonstrated that the Suspend before low feature is safe to use. Study success criteria, as defined in the protocol, were met (i.e. there were no device related serious adverse events, no diabetic ketoacidosis events related to the Suspend before low feature, and no unanticipated adverse device effects).

---

<sup>8</sup> Buckingham BA, Bailey TS, Christiansen M, et al. Evaluation of a Predictive Low-Glucose Management System In-Clinic. *Diabetes Technology and Therapeutics*. 2017;19(5):288-292.

## Clinical study overview (Ages 7-13 Years)

The Suspend before low feature was also evaluated in a study of the MiniMed 670G system that included subjects 7-13 years, diagnosed with type 1 diabetes mellitus.<sup>9</sup> This feature is the same in the MiniMed 780G system.

A total of 105 study subjects were observed overnight after exercise/activity while using the system with the Suspend before low feature activated. The Low Limit setting for the Suspend before low feature turned ON was set to 65 mg/dL and the subjects were observed with FST for a maximum of 12 hours.

## Feature performance and safety

In 79.7% of cases, after activation of the Suspend before low feature, the threshold of  $\leq 65$  mg/dL was avoided. Mean glucose levels up to six hours after the suspend feature was activated remained below the starting glucose levels.

Data from this in-clinic evaluation demonstrated that the Suspend before low feature is safe to use in a pediatric population.

## Guardian 4 System Performance



**Note:** You should review the information in this section with your healthcare professional to understand the performance of the Guardian 4 system.

## Clinical study overview

The performance of the Guardian 4 system was evaluated using data collected during a multi-center prospective clinical study.<sup>10</sup> The study included participants 7 to 80 years old. Within the 7 to 80 years age range, the study enrolled a total of 308 subjects previously diagnosed with type 1 or 2 diabetes and 267 of these subjects completed the study. Subjects ages 18 and older were instructed to wear a total of three sensors and transmitters in the abdomen and arm. Subjects ages 7 to 17 years old were instructed to wear a total of three sensors and transmitters in the arm and buttock. For

<sup>9</sup> Buckingham BA, Bailey TS, Christiansen M, et al. Evaluation of a Predictive Low-Glucose Management System In-Clinic. Diabetes Technology and Therapeutics. 2017;19(5):288-292

all subjects, the transmitters were used to record raw sensor signals during the study and there was no real-time calculation of sensor glucose values.

Frequent sample testing (FST) was performed on four occasions for subjects 14 and older and on two occasions for subjects 7 to 13 years of age.

Reference blood (plasma) glucose values were obtained with a YSI Glucose Analyzer every 5-15 minutes for subjects 7 years and older. During each FST, subjects with an established insulin sensitivity ratio and insulin carbohydrate ratio underwent a hypoglycemic challenge and a hyperglycemic challenge.

Data collected during the study was post-processed after the study using the Guardian 4 system sensor algorithm to convert the raw sensor information to sensor glucose values every five minutes. For the accuracy information presented in the following sections, YSI reference values were paired with the closest sensor glucose reading within five minutes of the time of the reference value measurement.

Sensor accuracy

Sensor accuracy was calculated for sensors compared to a YSI reference for subjects 7 years and older, in the arm insertion site. Do not insert the sensor into any other location.

Table 10. Overall Accuracy Compared to YSI

Patient Population	Insertion Site	Number of Subjects	Number of paired SG-YSI	Percent of SG within 20/20% of YSI	Mean Absolute Relative Difference (%)
Adults (18+)	Arm	153	20612	88.3 (88.0)	10.6
Pediatrics (7-17)	Arm	107	7702	85.6 (85.0)	11.6

\*CGM readings are within 50-400 mg/dL, inclusive, which are slightly truncated from GS3 (40-400 mg/dL).  
\*For 20% agreement, 20 mg/dL used when YSI < 70 mg/dL.

In Table 11 through Table 12, the agreement of the SG values to paired YSI values was assessed by calculating the percentage of SG values that were within 15%, 20%, and 40% of the paired YSI values. For SG readings less than 70 mg/dL, the absolute difference in mg/dL between the SG and paired YSI values was calculated.

<sup>10</sup> Medtronic Inc., Performance Evaluation of an Advanced Algorithm with CGM in Adults, Adolescents and Pediatrics 10838519DOC. September 2019.

**Table 11.** Overall accuracy of SG-YSI paired points within SG ranges; Adults, Arm

CGM Glucose Range (mg/dL)	Number of Subjects	Number of paired CGM-YSI	Percent within 15 mg/dL YSI	Percent within 20 mg/dL YSI	Percent within 40 mg/dL YSI	Percent within 15% YSI	Percent within 20% YSI	Percent within 40% YSI	Bias (mg/dL)	MARD (%)
A) < 54	46	252	77.4	87.7	98.0				-9.3	14.9
B) 54-69	99	2204	91.0	95.9	99.7				-2.3	10.2
C) 70-180	153	12893				72.7	84.1	98.6	-4.9	11.3
D) 181-250	146	3857				81.6	91.1	99.7	-11.4	9.3
E) > 250	97	1406				86.5	94.5	99.9	-7.3	8.2

\* For reference range < 70 mg/dL, agreement was based on 15/20/40 mg/dL. CGM readings are within 50-400 mg/dL, inclusive.

**Table 12.** Overall accuracy of SG-YSI paired points within SG ranges; Pediatrics, Arm

CGM Glucose Range (mg/dL)	Number of Subjects	Number of paired CGM-YSI	Percent within 15 mg/dL YSI	Percent within 20 mg/dL YSI	Percent within 40 mg/dL YSI	Percent within 15% YSI	Percent within 20% YSI	Percent within 40% YSI	Bias (mg/dL)	MARD (%)
A) < 54	28	103	62.1	78.6	98.1				-13.0	19.4
B) 54-69	53	562	79.2	89.1	98.9				-6.3	12.9
C) 70-180	106	3967				64.6	79.2	98.4	-9.6	13.1
D) 181-250	103	1992				78.6	90.2	99.9	-14.1	9.9
E) > 250	77	1078				86.9	94.4	100.0	-8.0	8.3

\* For reference range < 70 mg/dL, agreement was based on 15/20/40 mg/dL. CGM readings are within 50-400 mg/dL, inclusive.  
\* Includes pediatric subjects 7-17 years of age.

### Agreement when CGM reads “Below 50 mg/dL” or “Above 400 mg/dL”

The real-time CGM systems display glucose values between 50 mg/dL and 400 mg/dL. It displays “Below 50 mg/dL” when the SG value detected is below 50 mg/dL. It displays “Above 400 mg/dL” when the SG value detected is above 400 mg/dL. *Table 13* and *Table 14* illustrate the number and percentage of the paired YSI values in different blood glucose levels when the CGM system displays “Below 50 mg/dL” (LOW) or “Above 400 mg/dL” (HIGH).

**Table 13.** The number and percentage of YSI values collected when CGM displays “Below 50 mg/dL” (LOW)

CGM Display	Population	Insertion Site	CGM-YSI pairs	YSI (mg/dL)					Total
				<55	<60	<70	<80	≥80	
LOW	Adult	Arm	Cumulative, n	128	263	445	481	9	490
			Cumulative %	26%	54%	91%	98%	2%	
	Pediatrics	Arm	Cumulative, n	51	87	168	194	6	200
			Cumulative %	26%	44%	84%	97%	3%	

\*Includes pediatric subjects 7-17 years of age.

**Table 14.** The number and percentage of YSI values collected when CGM displays "Above 400 mg/dL" (HIGH)

CGM Display	Population	Insertion Site	CGM-YSI pairs	YSI (mg/dL)					
				>340	>320	>280	>240	≤240	Total
HIGH	Adult	Arm	Cumulative, n	20	21	21	21	0	21
			Cumulative %	95%	100%	100%	100%	0%	
	Pediatrics	Arm	Cumulative, n	32	32	32	32	0	32
			Cumulative %	100%	100%	100%	100%	0%	

\*Includes pediatric subjects 7-17 years of age.

## Concurrence of SG and YSI values

Table 15 and Table 16 show, for each SG range, the percentage of concurring data points where the paired YSI values were in different blood glucose ranges.

**Table 15.** Overall concurrence of YSI values and SG readings using SG ranges; Adults, Arm

Percent of matched pairs in each SMBG glucose range for each SG range (mg/dL)												
SG ranges (mg/dL)	Number of paired SG-YSI	YSI glucose ranges (mg/dL)										
		<50	≥50–60	>60–80	>80–120	>120–160	>160–200	>200–250	>250–300	>300–350	>350–400	>400
A) <50	490	6.7% (33/490)	47.6% (233/490)	43.9% (215/490)	1.8% (9/490)	0.0% (0/490)	0.0% (0/490)	0.0% (0/490)	0.0% (0/490)	0.0% (0/490)	0.0% (0/490)	0.0% (0/490)
B) ≥50–60	1036	7.0% (73/1036)	38.4% (398/1036)	50.8% (526/1036)	3.5% (36/1036)	0.1% (1/1036)	0.2% (2/1036)	0.0% (0/1036)	0.0% (0/1036)	0.0% (0/1036)	0.0% (0/1036)	0.0% (0/1036)
C) >60–80	2899	0.5% (15/2899)	17.8% (517/2899)	69.7% (2022/2899)	11.6% (335/2899)	0.2% (6/2899)	0.1% (4/2899)	0.0% (0/2899)	0.0% (0/2899)	0.0% (0/2899)	0.0% (0/2899)	0.0% (0/2899)
D) >80–120	4334	0.0% (2/4334)	0.7% (31/4334)	12.7% (552/4334)	63.0% (2730/4334)	22.5% (973/4334)	1.0% (43/4334)	0.1% (3/4334)	0.0% (0/4334)	0.0% (0/4334)	0.0% (0/4334)	0.0% (0/4334)
E) >120–160	5123	0.0% (0/5123)	0.0% (0/5123)	0.1% (6/5123)	11.8% (604/5123)	63.8% (3271/5123)	22.0% (1129/5123)	2.0% (105/5123)	0.2% (8/5123)	0.0% (0/5123)	0.0% (0/5123)	0.0% (0/5123)
F) >160–200	3337	0.0% (0/3337)	0.0% (0/3337)	0.0% (0/3337)	0.3% (10/3337)	13.4% (447/3337)	57.5% (1920/3337)	25.7% (856/3337)	2.8% (95/3337)	0.3% (9/3337)	0.0% (0/3337)	0.0% (0/3337)
G) >200–250	2477	0.0% (0/2477)	0.0% (0/2477)	0.0% (0/2477)	0.0% (0/2477)	0.3% (7/2477)	11.7% (291/2477)	61.1% (1514/2477)	23.8% (589/2477)	3.0% (74/2477)	0.1% (2/2477)	0.0% (0/2477)
H) >250–300	895	0.0% (0/895)	0.0% (0/895)	0.0% (0/895)	0.0% (0/895)	0.0% (0/895)	0.0% (0/895)	13.7% (123/895)	56.3% (504/895)	27.7% (248/895)	2.2% (20/895)	0.0% (0/895)
I) >300–350	390	0.0% (0/390)	0.0% (0/390)	0.0% (0/390)	0.0% (0/390)	0.0% (0/390)	0.0% (0/390)	1.0% (4/390)	25.4% (99/390)	55.4% (216/390)	16.7% (65/390)	1.5% (6/390)
J) >350–400	121	0.0% (0/121)	0.0% (0/121)	0.0% (0/121)	0.0% (0/121)	0.0% (0/121)	0.0% (0/121)	0.0% (0/121)	1.7% (2/121)	36.4% (44/121)	51.2% (62/121)	10.7% (13/121)
K) >400	21	0.0% (0/21)	0.0% (0/21)	0.0% (0/21)	0.0% (0/21)	0.0% (0/21)	0.0% (0/21)	0.0% (0/21)	0.0% (0/21)	4.8% (1/21)	52.4% (11/21)	42.9% (9/21)

**Table 16.** Overall concurrence of YSI values and SG readings using SG ranges; Pediatrics, Arm

Percent of matched pairs in each YSI glucose range for each SG range (mg/dL)											
SG Ranges (mg/dL)	Number of paired SG-YSI	YSI glucose ranges (mg/dL)									
		<50	≥50–60	>60–80	>80–120	>120–160	>160–200	>200–250	>250–300	>300–350	>350–400
A) <50	200	11.5% (23/200)	32.5% (65/200)	53.0% (106/200)	2.5% (5/200)	0.5% (1/200)	0.0% (0/200)	0.0% (0/200)	0.0% (0/200)	0.0% (0/200)	0.0% (0/200)
B) ≥50–60	346	4.0% (14/346)	26.0% (90/346)	64.2% (222/346)	5.2% (18/346)	0.6% (2/346)	0.0% (0/346)	0.0% (0/346)	0.0% (0/346)	0.0% (0/346)	0.0% (0/346)
C) >60–80	692	0.7% (5/692)	13.4% (93/692)	61.6% (426/692)	23.6% (163/692)	0.7% (5/692)	0.0% (0/692)	0.0% (0/692)	0.0% (0/692)	0.0% (0/692)	0.0% (0/692)
D) >80–120	1348	0.1% (1/1348)	1.2% (16/1348)	12.5% (168/1348)	57.0% (768/1348)	27.8% (375/1348)	1.4% (19/1348)	0.1% (1/1348)	0.0% (0/1348)	0.0% (0/1348)	0.0% (0/1348)
E) >120–160	1529	0.0% (0/1529)	0.0% (0/1529)	0.1% (1/1529)	7.5% (114/1529)	58.3% (891/1529)	30.0% (458/1529)	3.4% (52/1529)	0.8% (12/1529)	0.1% (1/1529)	0.0% (0/1529)
F) >160–200	1439	0.0% (0/1439)	0.0% (0/1439)	0.0% (0/1439)	0.1% (1/1439)	9.0% (129/1439)	51.4% (739/1439)	35.9% (516/1439)	3.3% (48/1439)	0.3% (5/1439)	0.0% (0/1439)
G) >200–250	1270	0.0% (0/1270)	0.0% (0/1270)	0.0% (0/1270)	0.0% (0/1270)	0.2% (3/1270)	9.9% (126/1270)	59.4% (754/1270)	28.7% (364/1270)	1.8% (23/1270)	0.0% (0/1270)
H) >250–300	695	0.0% (0/695)	0.0% (0/695)	0.0% (0/695)	0.0% (0/695)	0.0% (0/695)	0.1% (1/695)	13.1% (91/695)	58.6% (407/695)	25.8% (179/695)	2.4% (17/695)
I) >300–350	296	0.0% (0/296)	0.0% (0/296)	0.0% (0/296)	0.0% (0/296)	0.0% (0/296)	0.0% (0/296)	1.0% (3/296)	18.9% (56/296)	61.1% (181/296)	18.9% (56/296)
J) >350–400	87	0.0% (0/87)	0.0% (0/87)	0.0% (0/87)	0.0% (0/87)	0.0% (0/87)	0.0% (0/87)	0.0% (0/87)	5.7% (5/87)	24.1% (21/87)	60.9% (53/87)
K) >400	32	0.0% (0/32)	0.0% (0/32)	0.0% (0/32)	0.0% (0/32)	0.0% (0/32)	0.0% (0/32)	0.0% (0/32)	0.0% (0/32)	0.0% (0/32)	31.3% (10/32)
											68.8% (22/32)

\* Includes pediatric subjects 7-17 years of age.

## Trend accuracy

**Table 17.** Trend accuracy compared to YSI over time; Adults, Arm

SG Rate Ranges (mg/dL/min)	No. of Paired SG-YSI	YSI (mg/dL/min)					
		<-2	[-2,-1]	[-1, 0]	[0, 1]	(1, 2]	>2
<-2	179	53.6% (96/179)	38.0% (68/179)	7.3% (13/179)	1.1% (2/179)	0.0% (0/179)	0.0% (0/179)
[-2,-1]	1036	5.2% (54/1036)	49.9% (517/1036)	42.3% (438/1036)	2.5% (26/1036)	0.1% (1/1036)	0.0% (0/1036)
[-1, 0]	10059	0.2% (17/10059)	4.0% (401/10059)	79.1% (7958/10059)	16.3% (1640/10059)	0.4% (37/10059)	0.1% (6/10059)
[0, 1]	7342	0.0% (3/7342)	0.5% (38/7342)	22.6% (1656/7342)	69.9% (5129/7342)	6.6% (487/7342)	0.4% (29/7342)
(1, 2]	1513	0.0% (0/1513)	0.3% (5/1513)	2.0% (31/1513)	29.5% (446/1513)	58.6% (886/1513)	9.6% (145/1513)
>2	461	0.0% (0/461)	0.0% (0/461)	0.4% (2/461)	4.3% (20/461)	37.5% (173/461)	57.7% (266/461)

CGM readings are within 50-400 mg/dL, inclusive.

**Table 18.** Trend accuracy compared to YSI over time; Pediatrics, Arm

SG Rate Ranges (mg/dL/min)	YSI (mg/dL/min)						
	No. of Paired SG-YSI	<-2	[-2,-1]	[-1, 0]	[0, 1]	(1, 2]	>2
<-2	196	50.5% (99/196)	41.3% (81/196)	7.1% (14/196)	1.0% (2/196)	0.0% (0/196)	0.0% (0/196)
[-2,-1]	742	6.2% (46/742)	53.6% (398/742)	36.7% (272/742)	2.7% (20/742)	0.7% (5/742)	0.1% (1/742)
[-1, 0]	3103	0.3% (9/3103)	6.5% (201/3103)	76.2% (2363/3103)	15.8% (490/3103)	1.0% (31/3103)	0.3% (9/3103)
[0, 1]	2450	0.0% (0/2450)	0.9% (23/2450)	21.1% (517/2450)	68.6% (1680/2450)	8.9% (218/2450)	0.4% (11/2450)
(1, 2]	851	0.0% (0/851)	0.1% (1/851)	3.1% (26/851)	32.2% (274/851)	55.9% (476/851)	8.7% (74/851)
>2	354	0.0% (0/354)	0.0% (0/354)	0.3% (1/354)	4.2% (15/354)	28.0% (99/354)	67.5% (239/354)

CGM readings are within 50-400 mg/dL, inclusive.  
 \* Includes pediatric subjects 7-17 years of age.

## Accuracy over time

**Table 19.** Sensor Accuracy Compared to YSI Over Time; Adults, Arm

Wear Period	Number of paired SG-YSI	Percent of SG within 15/15% of YSI (%)	Percent of SG within 20/20% of YSI (%)	Percent of SG within 40/40% of YSI (%)	Mean Absolute Relative Difference (%)
Beginning	5678	65.8	78.3	98.3	13.9
Early Middle	5504	80.9	91	99.7	10.0
Late Middle	5142	81	91.5	99.7	9.7
End	4288	87.5	94.4	99.7	8.3

CGM readings are within 50-400 mg/dL, inclusive.  
 \* For reference range < 70 mg/dL, agreement was based on 15/20/40 mg/dL.

**Table 20.** Sensor Accuracy Compared to YSI Over Time; Pediatrics, Arm

Wear Period	Number of paired SG-YSI	Percent of SG within 15/15% of YSI (%)	Percent of SG within 20/20% of YSI (%)	Percent of SG within 40/40% of YSI (%)	Mean Absolute Relative Difference (%)
Beginning	3127	68.5	84.7	99.5	12.8
Early Middle	2546	74.4	85.5	98.9	11.5
Late Middle	1145	74.8	84.3	99.8	10.9
End	884	81.7	91	99.8	9.1

\* Includes pediatric subjects 7-17 years of age. CGM readings are within 50-400 mg/dL, inclusive.  
 \* For reference range < 70 mg/dL, agreement was based on 15/20/40 mg/dL.

## Precision

Precision of the system was evaluated by comparing the results from two separate sensors worn in the location on the same subject at the same time.

	Number of paired points	Percent Absolute Relative Difference (PARD)	Coefficient of variation (%CV)
7-17 YO Arm	10386	9.48	6.5
18+ YO Arm	23549	9.08	6.5

## Sensor life

### Adults (18-75 Years)



77.8% of sensors worn in the arm functioned more than six days and up to the full seven days of wear (144 to 168 hours). The median functional sensor life for sensors worn in the arm insertion site over the course of the study was 167.9 hours, with a mean functional life of 147.9 hours.

## Pediatrics

61.7% of sensors worn in the arm functioned more than six days and up to the full seven days of wear (144 to 168 hours). The median functional sensor life for sensors worn in the arm insertion site over the course of the study was 163.7 hours, with a mean functional life of 138.9 hours.

## Safety

There were no device-related or procedure-related serious adverse events, or unanticipated adverse device effects after seven days of use.

## Alert performance

CGM enables a device to display SG readings, glucose trend arrows, glucose trend graphs, and SG alerts, for example, High and Low Sensor Glucose alerts, High and Low Predicted alerts, and Rise and Fall alerts for rate-of-change.

The high and low SG alerts (**Threshold alerts**) let the user know when the SG is at or above the high limit or at or below the low limit. Using only a high or low Threshold alert may reduce the number of false alerts, but does not provide a warning before reaching a high or low limit. The default alert thresholds are highlighted in gray below.

**Predictive alerts** notify users that their SG level may soon reach a high or low limit setting. Users may select how early they would like to be notified before their SG level reaches a high or low limit. The earliest warning is 60 minutes before reaching a high or low limit, but users can reduce the amount of warning down to 10 minutes. Users receive a Predictive alert when their SG level is predicted to reach their high or low limit in the Time Before High or Time Before Low setting they select. In general, the earlier the warning, the more time a user has to react to a potential high or low, but this also increases the potential for false alerts.

A predictive alert is simply an estimation of a future SG level compared to the high or low limit setting. If the predicted future SG value is at or above the high limit or is at or

below the low limit, then a predictive alert is sounded even though the current SG level has not crossed the high or low limit. The predicted SG level is calculated using the current SG level, the derivative of current and previous SG readings (the trend or slope of the SG readings) and the Time Before High or Time Before Low duration the user selects.

The device always alerts the user with a Low SG alarm when the CGM reads that the user is at or below 64 mg/dL, regardless of the high/low threshold and/or predictive alerts that the user sets.

## Glucose TRUE Alert Rate

The glucose true alert rate is the rate at which the blood glucose (BG) confirmed that the CGM alert was triggered correctly. For example:

- **True Threshold Hypoglycemic alert rate** is a measure of how often the CGM read that the user was at or below the low threshold and the user's BG was actually at or below that low threshold.
- **True Threshold Hyperglycemic alert rate** is a measure of how often the CGM read that the user was at or above the high threshold and the user's BG was actually at or above that high threshold.
- **True Predictive Hypoglycemic alert rate** is a measure of how often the CGM predicted that the user would reach or go below the low threshold and the user's BG was actually at or below that low threshold within 15 or 30 minutes.
- **True Predictive Hyperglycemic alert rate** is a measure of how often the CGM predicted that the user would reach or go above the high threshold and the user's BG was actually at or above that high threshold within 15 or 30 minutes.

The true alert rate is important because it is necessary that users be notified when their blood glucose is low (or high) so that they can correct the low (or high) blood glucose. A high true alert rate indicates that when the CGM says that their glucose values are, or will reach a specified threshold, the user's blood glucose is likely to be at or approaching that threshold.

For example, per the following table, the low glucose alerts would have correctly indicated that the user was below (i.e. threshold only), or predicted to reach below the threshold (i.e. predictive only) or both (predictive and threshold) 85.9%, 64.1%, or 72.4%

of the time within 30 minutes (or 85.9%, 60.3% or 70.1% of the time within 15 minutes) when the user had BG values lower than 70 mg/dL for a sensor inserted in the arm.

**Table 21.** Glucose TRUE Alert Performance, Adults

Glucose TRUE Alert Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
50	Arm	27.3%	25.8%	28.3%	19.9%	28.0%	21.6%
60*	Arm	73.3%	71.9%	57.4%	53.0%	63.2%	59.9%
64	Arm	83.0%	82.4%	-	-	-	-
70	Arm	85.9%	85.9%	64.1%	60.3%	72.4%	70.1%
80	Arm	81.9%	80.8%	62.2%	58.9%	70.0%	67.6%
90	Arm	77.7%	77.4%	59.8%	57.0%	66.9%	65.1%
180	Arm	87.6%	87.2%	66.4%	64.0%	74.3%	72.7%
220	Arm	86.8%	86.8%	64.2%	62.1%	72.2%	70.9%
250	Arm	87.9%	87.0%	67.3%	62.0%	74.5%	70.8%
300	Arm	88.9%	88.9%	69.5%	65.5%	76.1%	73.5%

\*The default alert threshold is highlighted in gray.

**Table 22.** Glucose TRUE Alert Performance, Pediatrics (7-17 Years of Age)

Glucose TRUE Alert Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
50	Arm	23.8%	23.8%	13.3%	9.2%	16.4%	13.6%
60*	Arm	52.8%	52.8%	34.3%	32.9%	40.6%	39.6%
64	Arm	59.2%	59.2%	-	-	-	-
70	Arm	63.0%	63.0%	44.6%	40.8%	51.1%	48.6%
80	Arm	69.1%	69.1%	50.0%	44.5%	57.1%	53.6%
90	Arm	74.4%	72.2%	56.3%	51.7%	63.6%	60.0%
180	Arm	92.2%	92.2%	82.7%	78.3%	86.8%	84.3%
220	Arm	91.3%	90.2%	77.3%	74.0%	83.0%	80.5%
250	Arm	88.4%	88.4%	70.7%	68.4%	77.4%	75.9%
300	Arm	85.2%	85.2%	60.5%	59.9%	68.9%	68.5%

\* The default alert threshold is highlighted in gray.  
\* Includes pediatric subjects 7-17 years of age.

## Glucose FALSE Alert Rate

The glucose false alert rate is the rate at which the BG did not confirm that the CGM alert was triggered correctly. For example:

- **False Threshold Hypoglycemic alert rate** is a measure of how often the CGM read that the user was at or below the low threshold, but the user's BG was actually above that low threshold.
- **False Threshold Hyperglycemic alert rate** is a measure of how often the CGM read that the user was at or above the high threshold, but the user's BG was actually below that high threshold.

- **False Predictive Hypoglycemic alert rate** is a measure of how often the CGM predicted that the user would be at or below the low threshold, but the user's BG was actually above that low threshold within 15 or 30 minutes.
- **False Predictive Hyperglycemic alert rate** is a measure of how often the CGM predicted that the user would be at or above the high threshold, but the user's BG was actually below the high threshold within 15 or 30 minutes.

The false alert rate is important because it is necessary that users be correctly notified when their glucose is low or high so that they can correct the low or high glucose. A low false alert rate indicates that when the CGM says that their glucose values are, or will reach a specified threshold, the user's glucose is likely to be at or approaching that threshold.

For example, per the following table, the high glucose threshold alerts would have incorrectly indicated that the user was above (i.e. threshold only), or predicted to reach above the threshold (i.e. predictive only), or both (threshold and predictive) for adult 12.4%, 33.6% or 25.7% of the time within 30 minutes (or 12.8%, 36.0%, or 27.3% of the time within 15 minutes) when the user had BG less than 180 mg/dL for a sensor inserted in the arm.

**Table 23.** Glucose FALSE Alert Performance, Adults

Glucose FALSE Alert Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
50	Arm	72.7%	74.2%	71.7%	80.1%	72.0%	78.4%
60*	Arm	26.7%	28.1%	42.6%	47.0%	36.8%	40.1%
64	Arm	17.0%	17.6%	-	-	-	-
70	Arm	14.1%	14.1%	35.9%	39.7%	27.6%	29.9%
80	Arm	18.1%	19.2%	37.8%	41.1%	30.0%	32.4%
90	Arm	22.3%	22.6%	40.2%	43.0%	33.1%	34.9%
180	Arm	12.4%	12.8%	33.6%	36.0%	25.7%	27.3%
220	Arm	13.2%	13.2%	35.8%	37.9%	27.8%	29.1%
250	Arm	12.1%	13.0%	32.7%	38.0%	25.5%	29.2%
300	Arm	11.1%	11.1%	30.5%	34.5%	23.9%	26.5%
*The default alert threshold is highlighted in gray.							

**Table 24.** Glucose FALSE Alert Performance, Pediatrics (7-17 Years of Age)

Glucose FALSE Alert Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
50	Arm	76.2%	76.2%	86.7%	90.8%	83.6%	86.4%
60*	Arm	47.2%	47.2%	65.7%	67.1%	59.4%	60.4%
64	Arm	40.8%	40.8%	-	-	-	-
70	Arm	37.0%	37.0%	55.4%	59.2%	48.9%	51.4%

**Table 24.** Glucose FALSE Alert Performance, Pediatrics (7-17 Years of Age) (continued)

Glucose FALSE Alert Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
80	Arm	30.9%	30.9%	50.0%	55.5%	42.9%	46.4%
90	Arm	25.6%	27.8%	43.7%	48.3%	36.4%	40.0%
180	Arm	7.8%	7.8%	17.3%	21.7%	13.2%	15.7%
220	Arm	8.7%	9.8%	22.7%	26.0%	17.0%	19.5%
250	Arm	11.6%	11.6%	29.3%	31.6%	22.6%	24.1%
300	Arm	14.8%	14.8%	39.5%	40.1%	31.1%	31.5%
* The default alert threshold is highlighted in gray.							
* Includes pediatric subjects 7-17 years of age.							

## Glucose Correct Detection Rate

Glucose Correct Detection Rate is the rate that the device alerted when it should have alerted. For example, the BG was at or below the hypoglycemic threshold, or at or above the hyperglycemic threshold, and the device sounded a threshold or predictive alert.

The correct detection rates are important because it is necessary that users be notified when their glucose is low or high so that they can correct the low or high glucose. A high glucose correct detection rate indicates that users can have confidence that they will be notified by the device if their glucose is low or high.

For example, per the following table, the threshold alert, the predictive alert, or both (threshold and predictive) for adults notified the user 83.0%, 94.2% or 94.4% of the time within 30 minutes (or 81.8%, 90.4% or 91.3% within 15 minutes) when the user had BG less than 180 mg/dL in a sensor inserted in the arm.

**Table 25.** Glucose CORRECT DETECTION Alert Performance, Adults

Glucose CORRECT DETECTION Alert Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
50	Arm	37.5%	35.4%	89.6%	72.9%	89.6%	79.2%
60*	Arm	65.1%	64.0%	86.0%	83.7%	86.6%	84.3%
64	Arm	75.5%	75.0%	-	-	-	-
70	Arm	87.3%	86.8%	94.8%	93.4%	95.3%	94.3%
80	Arm	84.2%	83.1%	91.4%	86.8%	91.4%	88.7%
90	Arm	87.3%	86.3%	93.3%	90.2%	93.7%	91.7%
180	Arm	83.0%	81.8%	94.2%	90.4%	94.4%	91.3%
220	Arm	79.4%	78.3%	92.4%	90.6%	92.7%	90.9%
250	Arm	72.5%	71.1%	88.3%	85.3%	88.3%	86.1%
300	Arm	62.3%	60.1%	83.3%	80.4%	83.3%	80.4%
*The default alert threshold is highlighted in gray.							

**Table 26.** Glucose CORRECT DETECTION Alert Performance, Pediatrics (7-17 Years of Age)

Glucose CORRECT DETECTION Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
50	Arm	66.7%	66.7%	86.7%	60.0%	86.7%	73.3%
60*	Arm	64.6%	64.6%	80.0%	75.4%	80.0%	76.9%
64	Arm	74.6%	74.6%	-	-	-	-
70	Arm	83.3%	81.0%	92.9%	86.9%	92.9%	91.7%
80	Arm	86.6%	85.7%	97.3%	89.3%	97.3%	92.9%
90	Arm	89.9%	87.9%	97.3%	89.9%	98.0%	93.3%
180	Arm	90.3%	89.0%	96.2%	91.3%	96.9%	94.1%
220	Arm	88.0%	84.9%	95.9%	92.1%	96.2%	94.5%
250	Arm	77.1%	74.8%	91.6%	85.0%	91.6%	87.4%
300	Arm	72.2%	72.2%	93.8%	88.7%	93.8%	89.7%

\* The default alert threshold is highlighted in gray.  
\* Includes pediatric subjects 7-17 years of age.

Glucose Missed Detection Rate

The Missed Detection Rate is the rate that the device did not alert when it should have. For example, the BG was at or below the hypoglycemic threshold, or at or above the hyperglycemic threshold, and the device did not sound a threshold or predictive alert.

Missed detection rates are important because it is necessary that users be notified when their glucose is low or high, so that they can correct the low or high glucose. A low missed detection rate indicates that users can have confidence that they will be notified by the device if their glucose is low or high.

For example, per the following table, the threshold alert, predictive alert, or both alerts (threshold and predictive) for adults did not sound 12.7%, 5.2% or 4.7% of the time within 30 minutes (or 13.2%, 6.6% or 5.7% within 15 minutes) when the user had BG less than 70 mg/dL in a sensor inserted in the arm.

**Table 27.** Glucose MISSED Detection Performance, Adults

Glucose Missed Detection Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
50	Arm	62.5%	64.6%	10.4%	27.1%	10.4%	20.8%
60*	Arm	34.9%	36.0%	14.0%	16.3%	13.4%	15.7%
64	Arm	24.5%	25.0%	-	-	-	-
70	Arm	12.7%	13.2%	5.2%	6.6%	4.7%	5.7%
80	Arm	15.8%	16.9%	8.6%	13.2%	8.6%	11.3%
90	Arm	12.7%	13.7%	6.7%	9.8%	6.3%	8.3%
180	Arm	17.0%	18.2%	5.8%	9.6%	5.6%	8.7%
220	Arm	20.6%	21.7%	7.6%	9.4%	7.3%	9.1%
250	Arm	27.5%	28.9%	11.7%	14.7%	11.7%	13.9%
300	Arm	37.7%	39.9%	16.7%	19.6%	16.7%	19.6%

**Table 27.** Glucose MISSED Detection Performance, Adults (continued)

Glucose Missed Detection Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
*The default alert threshold is highlighted in gray.							

**Table 28.** Glucose MISSED Detection Performance, Pediatrics (7-17 Years of Age)

Glucose MISSED Detection Rate							
mg/dL	Insertion Site	Threshold Only		Predictive Only		Threshold and Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
50	Arm	33.3%	33.3%	13.3%	40.0%	13.3%	26.7%
60*	Arm	35.4%	35.4%	20.0%	24.6%	20.0%	23.1%
64	Arm	25.4%	25.4%	-	-	-	-
70	Arm	16.7%	19.0%	7.1%	13.1%	7.1%	8.3%
80	Arm	13.4%	14.3%	2.7%	10.7%	2.7%	7.1%
90	Arm	10.1%	12.1%	2.7%	10.1%	2.0%	6.7%
180	Arm	9.7%	11.0%	3.8%	8.7%	3.1%	5.9%
220	Arm	12.0%	15.1%	4.1%	7.9%	3.8%	5.5%
250	Arm	22.9%	25.2%	8.4%	15.0%	8.4%	12.6%
300	Arm	27.8%	27.8%	6.2%	11.3%	6.2%	10.3%
* The default alert threshold is highlighted in gray.							
* Includes pediatric subjects 7-17 years of age.							









# Glossary

<b>active insulin</b>	Active insulin is bolus insulin delivered by the insulin pump that continues to lower BG levels. Active insulin is not necessarily reflective of the pharmacokinetics and pharmacodynamics of rapid acting insulins.
<b>active insulin time</b>	Affects the length of time that bolus insulin is tracked as active insulin.
<b>activity guard</b>	An attachment that secures the reservoir during activity or when the insulin pump is worn by a child.
<b>alarm</b>	An audible beep or vibration with a message that requires immediate attention.
<b>alarm history</b>	A feature that stores information about recent alarms and alerts.
<b>alert</b>	An audible beep or vibration with a message to inform of a situation that may require attention.
<b>alert before low</b>	An alert that occurs when the low SG reading is being approached.
<b>alert limits</b>	The settings that determine when low and high SG alerts are triggered.
<b>alert on low</b>	An alert that occurs when the SG reading reaches or falls below the low limit.
<b>auto basal</b>	The automatically adjusted basal insulin delivered by the SmartGuard feature based on the current SG readings.
<b>Auto correction</b>	A correction bolus automatically delivered by the MiniMed 780G system to maximize time in range. Auto correction only occurs when using the SmartGuard feature.

<b>auto suspend</b>	A feature that suspends insulin delivery and triggers an alarm if no buttons are pressed for the specified period of time. Insulin delivery resumes when the alarm is cleared.
<b>awake mode</b>	A state in which the pump screen is on. The Home screen appears unless another screen is being used.
<b>basal insulin</b>	Insulin that is continuously delivered by the insulin pump to meet insulin needs between meals and during sleep.
<b>basal pattern</b>	A set of one or more basal rates that covers a 24-hour period.
<b>basal rate</b>	The setting for the amount of continuous basal insulin to be delivered per hour.
<b>BG</b>	The acronym for blood glucose. For more information, see <b>blood glucose (BG)</b> .
<b>BG meter</b>	A device that measures glucose levels in the blood.
<b>BG targets</b>	The high and low BG readings used for BG correction when using the Bolus Wizard feature in Manual mode.
<b>Block mode</b>	A feature that restricts the ability to change all settings. Certain functions can still be performed, such as suspend insulin delivery or clear alarms and alerts.
<b>blood glucose (BG)</b>	Glucose that is present in the blood, commonly measured by a BG meter.
<b>bolus BG check reminder</b>	A reminder for a BG check after programming a bolus. The reminder appears when the specified time period has passed.
<b>bolus insulin</b>	Insulin used to cover an expected rise in BG levels due to carbohydrates, or to lower a high BG reading down to the BG target range.
<b>bolus speed</b>	The delivery speed for bolus insulin.

<b>Bolus Wizard feature</b>	In Manual mode, a feature that uses individual Bolus Wizard settings to calculate an estimated bolus amount based on the BG value and the entered carbs. These settings include carb ratio, insulin sensitivity factor, BG target range, and active insulin time.
<b>calibrate</b>	The process of using a blood glucose (BG) meter reading to help the sensor glucose (SG) readings more closely match the glucose measured in your blood.
<b>cannula</b>	Short, thin, and flexible tube placed in the tissue below the skin. Insulin is delivered through the cannula into the body.
<b>carb ratio</b>	The number of grams of carbohydrates covered by one unit of insulin. The carb ratio is used to calculate bolus amounts.
<b>CGM</b>	The acronym for continuous glucose monitoring. For more information, see <b>continuous glucose monitoring (CGM)</b> .
<b>continuous glucose monitoring (CGM)</b>	A monitoring tool that uses a glucose sensor placed below the skin to continuously measure the amount of glucose in the interstitial fluid.
<b>correction bolus</b>	Insulin used to lower a high BG or SG reading down to a target value.
<b>CT scan</b>	The acronym for computed tomography scan.
<b>daily history</b>	Details of the events entered or actions performed using the insulin pump.
<b>diabetic ketoacidosis</b>	A serious condition that occurs when insulin levels are low, BG levels are elevated, and the body uses fat for energy. This process produces ketones, which upset the acid-base balance in the body, leading to a potentially life-threatening situation.
<b>Dual Wave bolus</b>	A type of bolus that provides a dose of insulin delivered as a combination of a normal bolus followed by a Square Wave bolus.
<b>Easy bolus</b>	A feature that delivers a normal bolus in preset increments using sound or vibrate confirmation.

<b>EMC</b>	The acronym for electromagnetic compatibility.
<b>ESD</b>	The acronym for electrostatic discharge.
<b>food bolus</b>	A dose of insulin given to cover an expected rise in glucose levels from carbohydrates.
<b>GPS</b>	The acronym for global positioning system.
<b>high limit</b>	The setting the insulin pump uses to determine when to alert for a high SG condition.
<b>infusion set</b>	Tubing that connects to the reservoir on one end, and has a needle or cannula on the other end, that is inserted into the body. Insulin travels from the insulin pump through the infusion set into the body.
<b>infusion site</b>	The location on the body where the infusion set is inserted.
<b>insulin sensitivity factor</b>	The amount that BG is reduced by one unit of insulin. The insulin sensitivity factor is used to calculate correction bolus amounts.
<b>insulin stacking</b>	Occurs when a bolus is delivered while active insulin from a previous bolus is still lowering glucose levels. Insulin stacking can result in hypoglycemia.
<b>interstitial fluid</b>	The fluid that surrounds the cells in the body.
<b>IV</b>	The acronym for intravenous.
<b>lock</b>	A feature that prevents accidental button presses.
<b>low limit</b>	The setting the insulin pump uses to determine when to alert for a low SG condition and suspend insulin delivery.
<b>Manual bolus</b>	A feature to manually enter and deliver a dose of insulin.
<b>Manual mode</b>	Manual mode refers to system functions that are used when the SmartGuard feature is not active.
<b>Max basal rate</b>	The maximum amount of basal insulin that can be programmed by the user per hour in Manual mode.
<b>Max bolus</b>	The maximum bolus amount that can be programmed by the user in one dose.

<b>meter</b>	A term for any BG meter.
<b>missed meal bolus reminder</b>	A reminder when a bolus is not delivered during the specified time period, which is often around meal times.
<b>MRI</b>	The acronym for magnetic resonance imaging.
<b>NiMH</b>	The acronym for nickel-metal hydride.
<b>normal bolus</b>	A type of bolus that provides an entire dose of insulin immediately.
<b>notifications</b>	All notifications are designed to get attention and convey different types of information. They include alarms, alerts, reminders, and messages.
<b>occlusion</b>	A blockage or crimp of the cannula or tubing that prevents proper insulin flow.
<b>piston</b>	The part of the insulin pump that engages the reservoir and moves insulin through the tubing.
<b>power save mode</b>	A state in which the insulin pump is fully functional, but the screen goes dark to save power.
<b>preset bolus</b>	A feature to set up and save a bolus for specific meals or snacks that are frequently consumed.
<b>preset temp basal</b>	A feature to set up and save temporary basal rates for repeated use.
<b>reminder</b>	A type of notification to help remember an action.
<b>reservoir</b>	The small container that is filled with insulin and inserted into the insulin pump.
<b>Resume basal alert</b>	An alert that occurs when the insulin pump has automatically resumed basal insulin delivery after a Suspend before low or Suspend on low event because the SG readings have met the necessary criteria. This alert always occurs if basal insulin delivery has resumed because the two-hour maximum suspend time has elapsed.

<b>rewind</b>	A feature that returns the piston to its start position to place a new reservoir into the insulin pump.
<b>RF</b>	The acronym for radio frequency.
<b>rise alert</b>	An alert that occurs if the SG reading is rising rapidly.
<b>sensitivity</b>	For more information, see <b>insulin sensitivity factor</b> .
<b>sensor (glucose sensor)</b>	The small part of the CGM system that is inserted just below the skin to measure glucose levels in the interstitial fluid and collect the sensor data. The sensor wirelessly sends the collected sensor data to the pump or other compatible mobile device.
<b>sensor glucose (SG)</b>	Glucose that is present in the interstitial fluid and is measured by a glucose sensor.
<b>set change reminder</b>	A reminder to change the infusion set.
<b>SG</b>	The acronym for sensor glucose. For more information, see <b>sensor glucose (SG)</b> .
<b>Sleep mode</b>	A state in which the insulin pump is fully functional, but the screen is dark. The insulin pump automatically enters Sleep mode when no buttons are pressed for about two minutes.
<b>SmartGuard bolus feature</b>	A feature that assists to calculate a recommended bolus amount based on optional carbohydrate intake and optional BG or SG measurement. One or both of the two optional values may be entered.
<b>SmartGuard feature</b>	An insulin delivery feature that automatically controls basal insulin delivery to regulate BG levels to a target SG value.
<b>SN</b>	The acronym for serial number.
<b>Square Wave bolus</b>	A bolus delivered evenly over the specified time period.
<b>suspend</b>	Suspend features include the Suspend before low feature and the Suspend on low feature.
<b>Suspend before low</b>	A feature that suspends insulin delivery when the sensor predicts the SG reading is approaching the low limit.



<b>suspend delivery</b>	A feature that stops all insulin delivery until it is resumed. Only the basal insulin restarts when delivery is resumed.
<b>Suspend on low</b>	A feature that suspends insulin delivery when the SG reading reaches or falls below the low limit.
<b>TDD</b>	The acronym for total daily dose.
<b>temp basal rate (temporary basal rate)</b>	A feature that temporarily increases or decreases the current basal rate for the specified duration of time.
<b>transfer guard</b>	The plastic piece that comes attached to the reservoir. It is used to connect the reservoir to the insulin vial while the reservoir fills with insulin.







# Index

## A

accuracy over time

    performance data 372

acetaminophen 31, 151, 190

active insulin

    about 100

    clearing 209

    viewing amount 71

Active Insulin cleared alert 291

Active Insulin reminder alert 292

active insulin reset to zero

    icon 72

Active Insulin reset to zero 292

Active Insulin Time

    about 100

    changing 253

adjustments

    bolus 193

Alarm History 219

alarms

    about 232, 233

    Active Insulin reset to zero 292

    audio options 234

    Auto Suspend 293

    Battery failed 293

    Battery not compatible 293

    Bolus stopped 294

Critical pump error 294

Delivery limit exceeded 294

Fill Cannula? 296

Insert battery 296

Insulin flow blocked 297, 298

Loading incomplete 299

Manage Settings error 300

Max Fill reached 300

Medical device 311

No reservoir detected 300

notification light 234

Power error detected 301

Power loss 301

pump 291

Pump error 301, 302

Pump restarted 303

red icon 234

Replace battery now 303

Rewind required 305

Stuck button 305

alert

    Active Insulin reminder 292

    max basal error 311

    performance data 373

- Alert before high
  - about 154
  - description 306
  - setting up 155
- Alert before low
  - about 161
  - description 339
  - setting up 166
- Alert on high
  - about 154
  - setting up 155
- Alert on low
  - about 164
  - setting up 166
- alerts
  - about 232, 234
  - Active Insulin cleared 291
  - Alert before high 306
  - Alert before low 306
  - audio options 235
  - Basal delivery resumed 307
  - BG not received 308
  - Bolus not delivered 293
  - Calibration not accepted 308
  - Change sensor 308, 309
  - Check connection 309
  - Check settings 294
  - Device not compatible 295
  - Device not found 295
  - Download slow 320
  - Enter BG now 309, 315, 316, 317
  - High BG 317
  - High SG 309
  - Lost sensor signal 310
  - Low battery Pump 299
  - Low battery transmitter 310
  - Low reservoir 299
  - No calibration occurred 311, 312
  - notification light 235
  - outside typical range 312
  - overview 234
  - Possible signal interference 312
  - pump 291, 320
  - Replace battery 303
  - Reservoir estimate 303
  - Rise alert 313
  - Sensor expired 313
  - Sensor signal not found 313
  - Sensor updating 313
  - silencing 168
  - SmartGuard 315
  - SmartGuard exit 317, 318
  - SmartGuard started 319
  - Suspend before low 314
  - Suspend on low 314

- Transmitter battery
  - depleted 314
  - Very high basal setting 319, 320
  - yellow icon 235
- Alert silence
  - in SmartGuard feature 200
- Alert Silence
  - about 168
  - canceling 170
- areas, for inserting infusion set 125
- arrows
  - trend arrows 176
- audio
  - settings 75
  - silencing 168
- Auto Basal 181
- Auto Correction
  - about 181
  - setting up 183
- Auto suspend 211
- Auto Suspend alarm 293
- B**
- backlight
  - battery use 66
  - setting up 203
- basal
  - about 81
  - daily delivery amount 217
  - history 217
  - Home screen 71
  - Max basal rate 82
  - patterns 84
  - preset temp basal rate 239
  - rate 81
  - summary 215
  - Temp Basal 94
- Basal delivery resumed
  - alarm 307
  - alert 307
  - message 307
- basal patterns
  - about 84
  - changing 244
  - copying 243
  - deleting 243
  - editing 243
  - setup 85
- Basal Patterns screen
  - about 91
- basal rate
  - about 81
- Basal screen
  - about 91
- battery
  - about 66
  - alarm 293, 296, 303
  - alert 299, 303
  - compartment location 63

- disposal 68
- icon 72
- inserting 66
- pump clip 68
- removing 68
- replacing 66
- Battery failed alarm 293
- Battery not compatible alarm 293
- BG
  - confirming 96
  - Entering a BG value in Manual mode 96
  - Entering a BG value in the SmartGuard feature 189
  - Home screen 71
  - Summary screen 218
- BG not received alert 308
- BG target
  - about 100
  - changing 253
- Block mode
  - about 204
  - icon 72
  - in SmartGuard 200
- blood glucose
  - confirming 96
- bolus
  - about 97
  - Bolus Wizard 97, 99
  - daily delivery amount 217
  - Dual Wave 258
  - Easy bolus 262
  - history 217, 218
  - increment 250
  - Manual bolus 97, 107
  - Max bolus setting 97
  - normal bolus 105, 107
  - options in Manual mode 97
  - preset 265
  - settings 250
  - shortcut 72
  - SmartGuard 190
  - speed 251
  - Square Wave 254
  - stopping 108, 268
  - summary 215, 217, 218
  - types 249
- Bolus BG Check reminder 229
- Bolus not delivered alert 293
- Bolus stopped alarm 294
- Bolus Wizard
  - about 97
  - Active Insulin Time 100
  - BG target 100
  - changing 252
  - Dual Wave 259
  - Insulin Sensitivity Factor 100
  - normal bolus 105



- settings 99
  - Square Wave 255
  - turning off 105
- buttons, pump 62
- C**
- calibrate
  - error 309
- calibration
  - calibrating the sensor 153
- Calibration not accepted alert 308
- cannula
  - about 65
  - filling 126
  - stop filling 128
- cannula alarm 296
- carb ratio
  - about 100
  - changing 252
  - setting up 101
- CareLink software 144
- CareLink uploader not found.
  - message 320
- CGM agreement
  - performance data 369
- Change sensor alert 308, 309
- Check connection alert 309
- Check settings alert 294
- cleaning
  - pump 283
  - transmitter 140
- clearing
  - active insulin 209
  - pump settings 208
- clinical study overview
  - performance data 367
- concurrence of SG and YSI values
  - performance data 370
- consumables
  - infusion set 47
  - reservoir 47
- continuous glucose monitoring (CGM)
  - high SG settings 153
  - home screen 174
  - low SG settings 158
  - pairing pump, transmitter 137
  - Sensor feature, turning on 173
  - sensor glucose reading 175
  - sensor graph 175
  - settings 153
- conventions, user guide 24
- Critical pump error alarm 294
- D**
- Daily History 219

- date
  - changing 203
  - Startup Wizard 69
- delivery 75
- Delivery limit exceeded alarm 294
- Delivery settings
  - SmartGuard 57
- Device Limit message 295
- Device not compatible alert 295
- Device not found alert 295
- display options 203
- Download slow alert 320
- Dual Wave
  - about 258
  - Bolus Wizard 259
  - example 250
  - Manual bolus 261
- Dual Wave bolus
  - stopping 268

**E**

- Easy bolus
  - about 262
  - delivering 264
  - setting up 263
- electromagnetic emissions
  - about 341
- electromagnetic immunity
  - about 342
- emergency kit 25

- Enter BG
  - in Manual mode 96
  - in SmartGuard 189
- Enter BG now alert 309, 315, 316, 317

## F

- Fill Cannula? alarm 296
- food estimate 329

## G

- glucose alert rate, false
  - performance data 375
- glucose alert rate, true
  - performance data 374
- glucose correct detection rate
  - performance data 377
- glucose missed detection rate
  - performance data 378
- graph
  - CGM 175
  - history 221
  - SmartGuard 187

## H

- High BG
  - alert 317
- High Limit
  - about 154
- High SG alert
  - about 154, 309

- high SG settings
  - about 153
  - setting up 155
- High Snooze 158
- history
  - Alarm History 219
  - Daily History 219
  - Graph 222
  - Paired Sensors 220
  - pump settings, viewing 210
  - SG Review 220
  - summary 215
  - Time in Range 222
- History & Graph 215
  - introduction 215
- Home screen
  - CGM 174
  - Manual mode 71
  - SmartGuard 187
- hydroxycarbamide 31
- hydroxyurea 31
- I
- icons
  - menu 74
  - SmartGuard graph 187
  - status 72
  - suspend 163
- infusion set
  - about 114
  - fill cannula alarm 296
  - fill tubing 121
  - inserting 125
  - overview 64
  - reservoir alarm 299
  - setting up 114
  - type 47
- Insert battery alarm 296
- inserting
  - infusion set 125
  - reservoir 121
  - sensor 170
- insulin
  - basal 81
  - basal patterns 84
  - bolus 97
  - flow blocked alarm 297, 298
  - icon 73
  - Max Fill reached alarm 300
  - setting up delivery 81
- insulin delivery
  - Block mode 204
  - resuming 91, 164, 168
  - stopping 91, 108, 268
  - suspending 91, 160, 162, 211
- Insulin flow blocked alarm 297, 298

Insulin Sensitivity Factor

- about 100
- changing 252
- setting up 102

## L

language

- changing 211
- setting up 69

light, notification 234

Loading incomplete alarm 299

Lost sensor signal alert 310

Low battery Pump alert 299

Low battery transmitter alert 310

Low limit

- about 159

low management mode

- summary 218

Low Reservoir

- alert 230, 299
- reminder 230

Low SG alarm

- about 159

low SG settings

- about 158
- changing 167
- setting up 165

Low Snooze 167

## M

Manage Settings 207

Manage Settings error alarm 300

Manual bolus

- about 97
- Dual Wave 261
- Easy bolus feature 264
- normal bolus 107
- Square Wave 257

Max basal

- about 82
- alarm 294

max basal error alert 311

max basal message 311

max basal rate

- manage settings 210

Max bolus

- about 98
- alarm 294

Max Fill reached alarm 300

meal bolus

- in SmartGuard 195
- reminder 229

Medical device alarm 311

medications

- acetaminophen 31
- hydroxycarbamide 31
- hydroxyurea 31
- paracetamol 31

## menu

- accessing 74
- Blood Glucose screen 75
- History & Graph screen 75
- icons 75
- Insulin screen 75
- map 75, 349
- Paired Devices screen 75
- Reservoir & Set screen 75
- Settings screen 75
- SmartGuard screen 75
- Sound & Vibration screen 75
- Status screen 75

## messages

- about 232, 235
- Basal delivery resumed 307
- CareLink uploader not found. 320
- CGM 306
- Device Limit 295
- max basal 311
- not for smartguard 312
- overview 235
- pump 291, 320
- Resume bolus 304
- Resume Dual bolus 304
- Resume Square bolus 304
- Sensor connected 313
- very high basal setting 314

## meter

- deleting 136
- pairing 133
- unpairing 135
- MiniMed 780G system
  - performance data 355
- Missed Meal Bolus reminder 229
- mobile device

- uploading 141

## modes

- Block 204
- Manual 71
- Sleep 62

**N**

- No calibration occurred
- alert 311, 312
- No reservoir detected alarm 300
- normal bolus
  - about 105
  - Bolus Wizard 105
  - delivery 105
  - example 250
  - Manual bolus 107
- not for smartguard message 312
- notification light 235
- notifications
  - about 227, 232
  - light 234

## O

### Occlusion detection

- alarm 337

- optional items 48

### options

- basal patterns 84

- bolus speed 327

- display 203

- Manage Settings 207

- Manual mode bolus delivery 97

- Self Test 206

- Sound & Vibration 75

- Status screen 76

- outside typical range alert 312

## P

- pair device 137, 141

- meter 133

- mobile device 141

- pairing transmitter 137

- unpairing transmitter 139

- Paired Sensors 220

- paracetamol 31, 151, 190

- patterns, basal

- about 84

- adding 242

- changing 244

- copying 243, 244

- deleting 243

- editing 243

- setup 85

- performance data 355

- accuracy over time 372

- alert 373

- CGM agreement 369

- clinical study overview 367

- concurrence of SG and YSI

- values 370

- glucose alert rate, false 375

- glucose alert rate, true 374

- glucose correct detection

- rate 377

- glucose missed detection

- rate 378

- MiniMed 780G system 355

- precision 372

- safety 373

- sensor 367

- sensor accuracy 368

- sensor life 372

- SmartGuard 356

- suspend before low 366

- trend accuracy 371

- personal reminders 228

- Possible signal interference

- alert 312

- Power error detected alarm 301

- Power loss alarm 301
- precision
  - performance data 372
- preset bolus
  - about 265
  - changing 267
  - delivering 268
  - setting up 265
- preset temp basal rate
  - about 239
  - canceling 242
  - setting up 239
  - starting 241
- Preset Temp basal rate
  - editing 240
- pump
  - alarms 291
  - alerts 291
  - buttons 61
  - cleaning 283
  - connection status 73
  - deleting meter 136
  - disposal 286
  - messages 291
  - overview 63
  - pairing, meter 133
  - pairing, mobile device 141
  - pairing, transmitter 137
  - parts illustration 63
  - rewinding 114
  - storing 284
  - unpairing, meter 135
  - unpairing, transmitter 139
  - updating the software 141
- pump care
  - removing the battery 68
- pump clip
  - removing battery cap with 68
  - tightening battery cap with 67
- Pump error alarm 301, 302
- Pump restarted alarm 303
- pump security 44
- pump settings
  - clearing 208
  - managing 207
  - restoring 207
  - saving 207
  - viewing history 210
- R**
- rate
  - basal 81
  - Max basal 82
  - preset temp basal 239
  - temp basal 94
- red light 235
- reminders
  - about 227
  - Bolus BG Check 229

- Low Reservoir 230
- Missed Meal Bolus 229
- Personal 228
- Set Change 231
- Replace battery alert 303
- Replace battery now alarm 303
- reservoir
  - about 114
  - connecting 117
  - diagram 63
  - filling 117
  - icon 73
  - inserting 121
  - Loading incomplete alarm 299
  - Low Reservoir alert 299
  - Max Fill reached alarm 300
  - models 47
  - No reservoir detected
    - alarm 300
  - overview 65
  - plunger 118
  - removing 114
  - setting up 114
  - transfer guard 118
- Reservoir estimate alert 303
- restoring pump settings 207
- Resume basal alert 164
- Resume bolus message 304
- Resume Dual bolus message 304
- Resume Square bolus message 304
- rewinding 121
- Rewind required alarm 305
- Rise Alert
  - about 154
  - description 313
  - setting up 156
- Rise Limit
  - about 154
  - setting up 156
- S**
- safety
  - contraindications 28
  - information 24
  - insulin guidelines 46
  - performance data 373
  - risks 29
  - warnings 33
- saving pump settings 207
- screen
  - Block mode 72, 200, 204, 205
  - Home screen 71
  - Language 69
  - Select Time Format 69
  - Summary 215, 219, 220, 221, 223
  - Time and Date 203
  - Welcome screen 69
- self test 206



- sensor
  - about 27
  - BG not received alert 308
  - calibration alert 308, 309
  - communication 73
  - components 170
  - expired alert 313
  - feature, deactivating 173
  - feature, turning on 173
  - Home screen 174
  - inserting 170
  - life icon 73
  - pairing the transmitter 137
  - performance data 367
  - removing 140, 174
  - starting 172
  - Summary screen 218
- sensor accuracy
  - performance data 368
- sensor alerts
  - canceling silence 170
  - silencing 168, 170
- Sensor connected message 313
- Sensor expired alert 313
- sensor glucose (SG)
  - graph 175
  - high SG settings 153
  - history 175
  - low SG settings 158
- sensor graph
  - about 174
  - CGM 175
- sensor life
  - performance data 372
- sensor performance
  - precision 372
- Sensor signal not found alert 313
- Sensor updating alert 313
- Set Change reminder 231
- settings
  - 24-hour period 87
  - bolus 250
  - Bolus Wizard 99
  - CGM 153
  - low SG 158
  - managing 207
  - startup 69
- SG Review 220
- shortcuts 72
- silencing alerts
  - about 168
- Sleep mode 62
- SmartGuard
  - about 179
  - alerts 315
  - Block mode in 200
  - bolus 190
  - bolus adjustments 193

- Checklist 184
- conditions to activate 183
- delivery settings 57
- entering a BG value 189
- exiting 198
- Home screen 187
- introduction 179
- messages 315
- performance data 356
- preparing 181
- returning to 199
- sensor graph 187
- Setting Temp Target 196
- setting up 183
- staying in 197
- suspend features 184
- Temp Target 196
- using 187
- SmartGuard exit alert 317, 318
- SmartGuard started alert 319
- software
  - updating the pump
  - software 141
- Sound & Vibration
  - adjusting 76
  - screen 75
- Square Wave
  - about 254
  - Bolus Wizard 255
  - example 250
  - Manual bolus 257
  - setting up 255
- Square Wave bolus
  - stopping 268
- starting
  - sensor 172
- Startup Wizard
  - about 69
- status icons
  - about 72
  - Block mode 72
  - CGM 174
  - Connection icon 73
  - insulin amount 73
  - sensor life 73
  - time 71
- Status screen
  - about 76
  - pump 76
  - sensor 76
  - shortcut 72
  - SmartGuard Checklist 76
- stopping
  - bolus 108, 268
  - insulin delivery 91
- storage mode 284
- Stuck button alarm 305

- Summary screen
  - about 215
  - BG 218
  - bolus 217, 218
  - Bolus Wizard 217
  - insulin delivery overview 217
  - low management mode 218
  - overview 216
  - sensor 218
  - SmartGuard 218
  - time in range 216
- Suspend before low
  - about 160
  - alert 314
  - in SmartGuard 184
  - performance data 366
  - setting up 166
  - unavailable 163
- Suspended delivery
  - about
    - resuming 92
  - after a suspend event
    - resuming 164
  - manually
    - resuming 168
- suspending insulin delivery
  - Auto suspend 211
  - in SmartGuard 184
  - manually 91
  - resuming automatically 164
  - resuming manually 168
  - Suspend before low 160
  - Suspend on low 162
- Suspend on low
  - about 162
  - alarm 314
  - in SmartGuard 184
  - setting up 166
  - unavailable 163
- system overview 63
- system security 44
- T**
  - target
    - BG 100
    - SmartGuard 180
  - TDD 217
  - temp basal rate
    - about 94
    - starting 95
    - types 94
  - Temp Target
    - canceling 197
    - setting up 196
  - time
    - changing 203
    - Startup Wizard 203

- Time before high
  - about 154
  - setting up 156
- time in range
  - about 222
  - shortcuts 72
- transfer guard 38, 117
- transmitter
  - charging 137
  - connecting 172
  - disconnect from the sensor 140
  - pairing 137
  - unpairing 139
- Transmitter battery depleted alert 314
- trend accuracy
  - performance data 371
- trend arrows 176
- troubleshooting
  - pump issues 276
  - sensor issues 279

- tubing
  - filling 121
  - image 63
  - Max Fill reached alarm 300
- tubing connector 33, 38, 63, 117

## U

- Updating
  - pump software 141

## V

- Very high basal setting alert 319, 320
- very high basal setting message 314
- vibration settings 75
- viewing
  - basal delivery information 91
  - pump settings history 210

## W

- wake the pump 62



# Medtronic



## Medtronic MiniMed

18000 Devonshire Street

Northridge, CA 91325

USA

1 800 646 4633

+1 818 576 5555

[www.medtronicdiabetes.com](http://www.medtronicdiabetes.com)

RF: M994838A001 or M994838A002 (Refer to device label)

**R<sub>x</sub>** Only

© 2025 Medtronic  
M056377C001\_3  
2025-04-22

**REF**

MMT-1884, MMT-1894



\*M056377C001\*

MiniMed™ 780G